**User Manual** 

# Tektronix

# RSA3303B & RSA3308B 3 GHz & 8 GHz Real-Time Spectrum Analyzers 071-2363-00

This document applies to firmware version 4.0 and above.

www.tektronix.com

Copyright © Tektronix. All rights reserved. Licensed software products are owned by Tektronix or its subsidiaries or suppliers, and are protected by national copyright laws and international treaty provisions.

Tektronix products are covered by U.S. and foreign patents, issued and pending. Information in this publication supercedes that in all previously published material. Specifications and price change privileges reserved.

TEKTRONIX and TEK are registered trademarks of Tektronix, Inc.

### **Contacting Tektronix**

Tektronix, Inc. 14200 SW Karl Braun Drive P.O. Box 500 Beaverton, OR 97077 USA

For product information, sales, service, and technical support:

- In North America, call 1-800-833-9200.
- Worldwide, visit www.tektronix.com to find contacts in your area.

### End User License Agreement for Microsoft Software

You have acquired a device ("DEVICE") that includes software licensed by Tektronix Japan, Ltd. ("TJ") from Microsoft Licensing Inc. or its affiliates ("MS"). Those installed software products of MS origin, as well as associated media, printed materials, and "online" or electronic documentation ("SOFTWARE") are protected by international intellectual property laws and treaties. The SOFTWARE is licensed, not sold. All rights reserved.

IF YOU DO NOT AGREE TO THIS END USER LICENSE AGREEMENT ("EULA"), DO NOT USE THE DEVICE OR COPY THE SOFTWARE. INSTEAD, PROMPTLY CONTACT TJ FOR INSTRUCTIONS ON RETURN OF THE UNUSED DEVICE(S) FOR A REFUND. **ANY USE OF THE SOFTWARE, INCLUDING BUT NOT LIMITED TO USE ON THE DEVICE, WILL CONSTITUTE YOUR AGREEMENT TO THIS EULA (OR RATIFICATION OF ANY PREVIOUS CONSENT).** 

GRANT OF SOFTWARE LICENSE. This EULA grants you the following license:

You may use the SOFTWARE only on the DEVICE.

#### NOT FAULT TOLERANT.

THE SOFTWARE IS NOT FAULT TOLERANT. TJ HAS INDEPENDENTLY DETERMINED HOW TO USE THE SOFTWARE IN THE DEVICE, AND MS HAS RELIED UPON TJ TO CONDUCT SUFFICIENT TESTING TO DETERMINE THAT THE SOFTWARE IS SUITABLE FOR SUCH USE.

#### NO WARRANTIES FOR THE SOFTWARE.

THE SOFTWARE IS PROVIDED "AS IS" AND WITH ALL FAULTS. THE ENTIRE RISK AS TO SATISFACTORY QUALITY, PERFORMANCE, ACCURACY, AND EFFORT (INCLUDING LACK OF NEGLIGENCE) IS WITH YOU. ALSO, THERE IS NO WARRANTY AGAINST INTERFERENCE WITH YOUR ENJOYMENT OF THE SOFTWARE OR AGAINST INFRINGEMENT. IF YOU HAVE RECEIVED ANY WARRANTIES REGARDING THE DEVICE OR THE SOFTWARE, THOSE WARRANTIES DO NOT ORIGINATE FROM, AND ARE NOT BINDING ON, MS.

### NOTE ON JAVA SUPPORT.

The SOFTWARE may contain support for programs written in Java. Java technology is not fault tolerant and is not designed, manufactured, or intended for use or resale as online control equipment in hazardous environments requiring fail-safe performance, such as in the operation of nuclear facilities, aircraft navigation or communication systems, air traffic control, direct life support machines, or weapons systems, in which the failure of Java technology could lead directly to death, personal injury, or severe physical or environmental damage. Sun Microsystems, Inc. has contractually obligated MS to make this disclaimer.

### NO LIABILITY FOR CERTAIN DAMAGES.

EXCEPT AS PROHIBITED BY LAW, MS SHALL HAVE NO LIABILITY FOR ANY INDIRECT, SPECIAL, CONSEQUENTIAL OR INCIDENTAL DAMAGES ARISING FROM OR IN CONNECTION WITH THE USE OR PERFORMANCE OF THE SOFTWARE. THIS LIMITATION SHALL APPLY EVEN IF ANY REMEDY FAILS OF ITS ESSENTIAL PURPOSE. IN NO EVENT SHALL MS BE LIABLE FOR ANY AMOUNT IN EXCESS OF U.S. TWO HUNDRED FIFTY DOLLARS (U.S.\$250.00).

#### LIMITATIONS ON REVERSE ENGINEERING, DECOMPILATION, AND DISASSEMBLY.

You may not reverse engineer, decompile, or disassemble the SOFTWARE, except and only to the extent that such activity is expressly permitted by applicable law notwithstanding this limitation.

#### SOFTWARE TRANSFER ALLOWED BUT WITH RESTRICTIONS.

You may permanently transfer rights under this EULA only as part of a permanent sale or transfer of the Device, and only if the recipient agrees to this EULA. If the SOFTWARE is an upgrade, any transfer must also include all prior versions of the SOFTWARE.

#### EXPORT RESTRICTIONS.

You acknowledge that SOFTWARE is of US-origin. You agree to comply with all applicable international and national laws that apply to the SOFTWARE, including the U.S. Export Administration Regulations, as well as end-user, end-use and country destination restrictions issued by U.S. and other governments. For additional information on exporting the SOFTWARE, see http://www.microsoft.com/exporting/.

#### LIMITATION ON SOFTWARE PROGRAMS USED ON THE DEVICE.

The combination of software programs you use on the DEVICE shall address not more than two (2) general office automation or consumer computing functions. Such functions include, but are not limited to: email, word processing, spreadsheets, database, network browsing, scheduling, and personal finance.

#### STORAGE/NETWORK USE.

The SOFTWARE may not be installed, accessed, displayed, run, shared or used concurrently on or from different computers, including a workstation, terminal or other digital electronic device ("Computing System"). Notwithstanding the foregoing and except as otherwise provided below, any number of Computing Systems, may access or otherwise utilize the file and print services and internet information services of the SOFTWARE, if included.

You may use the SOFTWARE on a single DEVICE as interactive workstation software, but not as server software. However, you may permit a maximum of ten (10) Computing Systems to connect to the DEVICE to access and use services of the SOFTWARE, such as file and print services and internet information services. The ten-connection maximum includes any indirect connections made through other software or hardware which pools or aggregates connections.

#### Warranty 2

Tektronix warrants that this product will be free from defects in materials and workmanship for a period of one (1) year from the date of shipment. If any such product proves defective during this warranty period, Tektronix, at its option, either will repair the defective product without charge for parts and labor, or will provide a replacement in exchange for the defective product. Parts, modules and replacement products used by Tektronix for warranty work may be new or reconditioned to like new performance. All replaced parts, modules and products become the property of Tektronix.

In order to obtain service under this warranty, Customer must notify Tektronix of the defect before the expiration of the warranty period and make suitable arrangements for the performance of service. Customer shall be responsible for packaging and shipping the defective product to the service center designated by Tektronix, with shipping charges prepaid. Tektronix shall pay for the return of the product to Customer if the shipment is to a location within the country in which the Tektronix service center is located. Customer shall be responsible for paying all shipping charges, duties, taxes, and any other charges for products returned to any other locations.

This warranty shall not apply to any defect, failure or damage caused by improper use or improper or inadequate maintenance and care. Tektronix shall not be obligated to furnish service under this warranty a) to repair damage resulting from attempts by personnel other than Tektronix representatives to install, repair or service the product; b) to repair damage resulting from improper use or connection to incompatible equipment; c) to repair any damage or malfunction caused by the use of non-Tektronix supplies; or d) to service a product that has been modified or integrated with other products when the effect of such modification or integration increases the time or difficulty of servicing the product.

THIS WARRANTY IS GIVEN BY TEKTRONIX WITH RESPECT TO THE PRODUCT IN LIEU OF ANY OTHER WARRANTIES, EXPRESS OR IMPLIED. TEKTRONIX AND ITS VENDORS DISCLAIM ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. TEKTRONIX' RESPONSIBILITY TO REPAIR OR REPLACE DEFECTIVE PRODUCTS IS THE SOLE AND EXCLUSIVE REMEDY PROVIDED TO THE CUSTOMER FOR BREACH OF THIS WARRANTY. TEKTRONIX AND ITS VENDORS WILL NOT BE LIABLE FOR ANY INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES IRRESPECTIVE OF WHETHER TEKTRONIX OR THE VENDOR HAS ADVANCE NOTICE OF THE POSSIBILITY OF SUCH DAMAGES.

# **Table of Contents**

General Safety Summary	
Preface	xix
About This Manual	xix
Related Documents	XX
Conventions	xxi

### **Getting Started**

Product Overview         Features         Application         Difference Between RSA3303B and RSA3308B         Real-Time Analysis         Architecture	<b>1-1</b> 1-1 1-2 1-2 1-3 1-7
Installation	1-11
Unpacking to Check Contents Applying Power Setting Up the Stand Functional Check Powering Off the Analyzer Restart Backing Up User Files About Installation of Other Applications	1-11 1-12 1-15 1-16 1-20 1-21 1-22 1-22
Calibration         Cal Menu         Calibrating Gain         Calibrating Center Offset         Calibrating DC Offset         Adjusting Display Brightness         Confirming Performance	<b>1-23</b> 1-24 1-25 1-26 1-27 1-28 1-28

## **Operating Basics**

Interface Maps		
Controls and Connectors	2-2	
Display Screen	2-8	
Functional Overview	2-13	
Menu Keys	2-13	
Menu Operations	2-16	
Measurement Basics	2-21	
Setting System Parameters	2-23	
Using Windows XP	2-26	

Tutorial
Preparations
Displaying Spectrum
Using Markers and Peak Search
Using Averaging and Comparison Displays
Displaying DPX Spectrum
Displaying a Spectrogram
Spectrum Analysis
Modulation Analysis
Turning Off the Power

### Reference

Spectrum Analysis (S/A Mode)	3-1
Measurement Screen Layout	3-2
Spectrum Analysis	3-4
DPX Spectrum Display	3-14
Spectrogram Display	3-17
Real-Time Analysis	3-19
Zoom Function	3-23
Modulation Analysis (Demod Mode)	
Measurement Screen Layout	3-29
Analog Modulation Analysis	3-37
Digital Modulation Analysis (Option 21 Only)	3-44
RFID Analysis (Option 21 Only)	3-61
Time Analysis (Time Mode)	3-83
Measurement Screen Layout	3-83
Time Variation Measurement	3-84
CCDF Measurement	3-88
Pulse Measurements	3-93
Signal Source Analysis (Option 21 Only)	3-103
Setting Frequency and Span	3-119
Frequency and Span Setting Menu	3-120
Using the Channel Table	3-122
Using the Marker and Peak Search	3-123
Frequency Setting Range	3-124
Vector Span	3-126
Setting Amplitude	3-127
Amplitude Menu	3-127
Over-Voltage Input	3-131
Amplitude Correction	3-132
Setting Acquisition and Analysis Parameters	3-141
Acquisition/Analysis Menu	3-141
Seamless Acquisition	3-145
Trigger	3-147
Trigger Menu	3-14/
Creating a Trigger Mask (Option 02 Only)	3-140
Trigger Point Indicator	3-155
Synchronous Operation with External Instruments	3-162

FFT and RBW	3-163
RBW/FFT Menu	3-164
FFT Points	3-168
FFT Window	3-169
Trace Comparison Display and Average Function	3-175
Trace/Avg Menu	3-176
Displaying Trace 1 and 2	3-179
Averaging the Waveform	3-180
Saving/Loading Waveform Data	3-183
Trace Compression	3-183
Setting Views	3-187
View Menu	3-188
Basic Procedure	3-189
Spectrum View Setting	3-190
DPX Spectrum View Setting	3-191
Spectrogram View Setting	3-193
Time Domain View Setting	3-195
CCDF View Setting	3-197
Constellation View Setting (Option 21 Only)	3-198
EVM View Setting (Option 21 Only)	3-199
Symbol Table Setting (Option 21 Only)	3-201
Eye Diagram Setting (Option 21 Only)	3-202
AM/AM View Setting (Option 21 Only)	3-204 3-205
AM/PM View Setting (Option 21 Only)	3-205
PDF View Setting (Option 21 Only) Noisogram View Setting (Option 21 Only)	3-200
Display Line	3-209
Display Line Operation (Other than Real Time S/A)	3-210
Multi Display Lines (Real Time S/A Only)	3-213
Marker Operation and Peak Search	3-217
Markers Menu	3-218
Using Markers	3-220
Peak Search	3-225
Using the Online Help	3-227
Displaying the Online Help	3-227
Using a Mouse and Keyboard	3-230
	3-231
Selecting Input Source	3-231
Input Menu	
File Operations	3-233
File Type	3-233
Load/Save Menu	3-234
Saving and Loading Files	3-236
Entering a File Name	3-242
Deleting a File	3-244
Making/Deleting a Directory	3-244

	File Format3-Data File Format3Trace File Format3Bitmap Trace File Format3CSV File Format3MAT File Format3
	Screen Copy3.Print Menu3Printing a Screen Copy3Creating a File3
Appendices	
	Appendix A: SpecificationsElectrical CharacteristicsPhysical CharacteristicsEnvironmental Characteristics
	Appendix B: Menu Tree
	Appendix C: Digital Demodulation Symbol Mapping
	Appendix D: Restoring the Instrument Operating System Restoring the Operating System from the Instrument Hard Disk
	Appendix E: Inspection and CleaningInspecting the ExteriorCleaning the ExteriorLubricationCleaning the Interior
	Appendix F: Options and Accessories         Options         Standard Accessories         Optional Accessories
Glossary	

Index

# **List of Figures**

Figure 1-1: Concept of the swept spectrum analyzer	1-3
Figure 1-2: Sweep of resolution filter	1-3
Figure 1-3: Concept of the real-time spectrum analyzer	1-4
Figure 1-4: Concurrent acquisition	1-4
Figure 1-5: Frame acquisition	1-5
Figure 1-6: DPX processing	1-5
Figure 1-7: DPX spectrum display	1-5
Figure 1-8: Block diagram	1-0
Figure 1-9: AC Input (rear panel)	1-12
Figure 1-9: AC Input (rear panel) Figure 1-10: Principal power switch (rear panel)	1-12
	1-13
Figure 1-11: Front panel power switch (ON/STANDBY switch)	
Figure 1-12: Initial screen	1-14
Figure 1-13: RF Input connector	1-14
Figure 1-14: Setting up the stand	1-15
Figure 1-15: Spectrum of the calibration signal (50 MHz, about	1 10
-20 dBm)	1-16
Figure 1-16: Setup display	1-17
Figure 1-17: Reference level setting and overrange indicator	1-18
Figure 1-18: Spectrogram display	1-19
Figure 1-19: Calibration menu structure	1-24
Figure 1-20: UNCAL display	1-25
Figure 1-21: Center offset	1-26
Figure 1-22: DC offset	1-27
Figure 1-23: System menu	1-28
Figure 2-1: Front panel	2-2
Figure 2-2: Rear panel	2-4
Figure 2-3: Side panel	2-5
Figure 2-4: Connecting USB devices	2-6
Figure 2-5: Operation with mouse and keyboard	2-7
Figure 2-6: Display screen configuration	2-8
Figure 2-7: Status display	2-10
Figure 2-8: Pre- and post-trigger regions	2-10 2-10
Figure 2-9: Key lock display	2-10 2-11
	2-11 2-12
Figure 2-10: Setup display	2-12

Figure 2-11: Menu keys	2-14
Figure 2-12: Example menu item display	2-16
Figure 2-13: Menu item types	2-17
Figure 2-14: Numeric setting menu	2-17
Figure 2-15: Changing value with the knob	2-18
Figure 2-16: Changing value with the keypad	2-19
Figure 2-17: Numeric keypad	2-19
Figure 2-18: Step size for the center frequency	2-20
Figure 2-19: Changing the step size to the predetermined values	2-20
Figure 2-20: Selecting the measurement mode	2-21
Figure 2-21: Starting/Stopping data acquisition	2-22
Figure 2-22: Restoring default settings	2-22
Figure 2-23: System menu	2-23
Figure 2-24: Displaying versions and options	2-24
Figure 2-25: Displaying the Windows XP accessory menu	2-27
Figure 2-26: Cable connection	2-30
Figure 2-27: Principal power switch (rear panel)	2-31
Figure 2-28: Power switch (On/Standby)	2-31
Figure 2-29: Initial screen	2-32
Figure 2-30: Frequency, span, and amplitude settings	2-33
Figure 2-31: Menu items for numeric value input	2-34
Figure 2-32: Numeric value entry keypad	2-34
Figure 2-33: Center frequency of 100 MHz, span of 15 MHz	2-35
Figure 2-34: Center frequency of 100 MHz, span of 20 kHz	2-36
Figure 2-35: Setting amplitude	2-37
Figure 2-36: Reference level of 10 dBm	2-38
Figure 2-37: Status indicator	2-39
Figure 2-38: Measurement with a single marker	2-40
Figure 2-39: Measurement with the delta marker	2-41
Figure 2-40: Searching for the peak	2-43
Figure 2-41: Averaging a waveform	2-44
Figure 2-42: Comparison display with averaged waveform	2-45
Figure 2-43: DPX spectrum display	2-46
Figure 2-44: Persistence display	2-47
Figure 2-45: Concurrent display of spectrum and spectrogram	2-48
Figure 2-46: Tall display of spectrum and spectrogram	2-49
Figure 2-47: Spectrogram view	2-50
Figure 2-48: Channel power measurement	2-51

= 40 kHz)2-5Figure 2-50: Carrier frequency measurement2-5Figure 2-51: AM signal measurement2-5Figure 2-52: Analysis range settings2-5Figure 2-53: Frame and block2-5Figure 2-54: Changing the acquisition length2-5Figure 2-55: Single view display2-5
Figure 2-51: AM signal measurement2-5Figure 2-52: Analysis range settings2-5Figure 2-53: Frame and block2-5Figure 2-54: Changing the acquisition length2-5Figure 2-55: Single view display2-5
Figure 2-52: Analysis range settings2-5Figure 2-53: Frame and block2-5Figure 2-54: Changing the acquisition length2-5Figure 2-55: Single view display2-5
Figure 2-53: Frame and block2-5Figure 2-54: Changing the acquisition length2-5Figure 2-55: Single view display2-5
Figure 2-54: Changing the acquisition length2-5Figure 2-55: Single view display2-5
Figure 2-55: Single view display2-5
Figure 2-56: Scale settings2-5
Figure 3-1: S/A menu structure    3-
Figure 3-2: S/A measurement screen   3-
Figure 3-3: "Flex" grid style   3-
Figure 3-4: Channel power measurement
Figure 3-5: ACPR measurement band power markers
Figure 3-6: Example of ACPR measurement
Figure 3-7: C/N measurement band power markers
Figure 3-8: Example of C/N measurement   3-
Figure 3-9: OBW measurement band power marker
Figure 3-10: OBW measurement   3-
Figure 3-11: Carrier frequency measurement       3-1
Figure 3-12: EBW measurement band power markers
Figure 3-13: EBW measurement
Figure 3-14: Setting up spurious signal measurement       3-1
Figure 3-15: Example of spurious signal measurement
Figure 3-16: DPX process flow
Figure 3-17: DPX spectrum display   3-1
Figure 3-18: Persistence display     3-1
Figure 3-19: Spectrum and spectrogram concurrent display 3-1
Figure 3-20: View orientation   3-1
Figure 3-21: Single display   3-1
Figure 3-22: Difference between normal and real-time modes 3-1
Figure 3-23: Real-time mode   3-2
Figure 3-24: Zoom function   3-2
Figure 3-25: Zoom domain settings using the Acquisition/Analysis
menu
Figure 3-26: Zoomed area settings using the marker       3-2
Figure 3-27: Demod menu structure3-2
Figure 3-28: Demod mode screen    3-2

Figure 3-29: Analysis range setting in the overview	3-30
Figure 3-30: Markers keys	3-31
Figure 3-31: Specifying the range with the marker and the	
reference cursor	3-32
Figure 3-32: Specifying the analysis range using the mouse	3-33
Figure 3-33: FFT processing range setting on the overview	3-34
Figure 3-34: Changing the overview and subview	3-35
Figure 3-35: One-view display	3-36
Figure 3-36: AM demod measurement	3-39
Figure 3-37: FM demod measurement	3-40
Figure 3-38: PM demod measurement	3-41
Figure 3-39: IQ versus time measurement	3-42
Figure 3-40: Pulse spectrum measurement	3-43
Figure 3-41: Process flow for digitally-modulated signals	3-48
Figure 3-42: Constellation analysis	3-50
Figure 3-43: EVM measurement	3-51
Figure 3-44: IQ/Frequency versus time measurement	3-52
Figure 3-45: Power variance measurement	3-53
Figure 3-46: Symbol table analysis	3-54
Figure 3-47: Eye diagram analysis	3-55
Figure 3-48: Setting Linear Signal Region	3-56
Figure 3-49: AM/AM measurement	3-57
Figure 3-50: AM/PM measurement	3-58
Figure 3-51: CCDF measurement	3-59
Figure 3-52: PDF measurement	3-60
Figure 3-53: RFID waveform and measurement items	3-62
Figure 3-54: Carrier measurement	3-65
Figure 3-55: Spurious measurement	3-66
Figure 3-56: ACPR measurements	3-68
Figure 3-57: Modified Miller code	3-71
Figure 3-58: Transmission power on/down measurement	
parameters	3-72
Figure 3-59: Power on/down measurements	3-75
Figure 3-60: RF envelope / FSK pulse measurement (main view)	3-77
Figure 3-61: RF envelope / FSK pulse measurement parameters	3-78
Figure 3-62: Constellation, eye diagram, and symbol table	
(main view)	3-80
Figure 3-63: Definition of the modulation depth and index	3-82
Figure 3-64: Time menu structure	3-83

Figure 3-65: IQ versus Time measurement	3-85
Figure 3-66: Power variance measurement	3-86
Figure 3-67: Frequency measurement	3-87
Figure 3-68: CCDF calculation process	3-88
Figure 3-69: CCDF single view	3-91
Figure 3-70: CCDF multi-view	3-92
Figure 3-71: Definition of pulse characteristics	3-94
Figure 3-72: Definition of pulse-on/off	3-94
Figure 3-72: Setting for View Results For	3-98
Figure 3-74: Setting the View Define menu	3-99
Figure 3-75: Waveform display in the subview	3-102
Figure 3-76: Phase noise measurement parameters	3-102
Figure 3-77: Phase noise measurement	3-100
Figure 3-77: Thase holse measurement	3-107
Figure 3-79: Real-time phase noise measurement parameters	3-110
Figure 3-80: Real-time phase noise measurement	3-112
Figure 3-81: Settling time displayed in the subview	3-113
Figure 3-82: Real-time spurious measurement	3-114
Figure 3-83: Frequency versus Time measurement parameters	3-110
Figure 3-84: Frequency versus Time measurement parameters	3-118
Figure 3-85: Setting frequency and span	3-110
Figure 3-86: Setting the center frequency using Marker ⇒	3-123
Figure 3-87: Relationship between the frequency and span settings	3-125
Figure 3-88: Vector mode and scalar mode	3-125
Figure 3-89: Amplitude menu structure	3-120
Figure 3-90: Setting the amplitude	3-127
Figure 3-91: Overrange indicator	3-120
Figure 3-92: The amplitude correction concept	3-131
Figure 3-93: Amplitude correction example	3-134
Figure 3-94: Amplitude correction data input	3-137
Figure 3-95: Amplitude offset	3-139
Figure 3-96: Frequency offset	3-139
Figure 3-97: Amplitude correction setup display	3-140
Figure 3-98: Acquisition/Analysis menu structure	3-141
Figure 3-99: Acquisition and analysis parameters	3-143
Figure 3-100: Parameters in the zoom mode	3-143
Figure 3-101: Frame cycle	3-145
Figure 3-101: Frame cycle      Figure 3-102: Seamless acquisition	3-145
11gur - 102. Scamicos acquisitivii	5-145

Fi	gure 3-103: Trigger menu structure	3-148
Fi	gure 3-104: Trigger position	3-150
Fi	gure 3-105: Power trigger detection range	3-151
Fi	gure 3-106: Trigger level and slope	3-152
Fi	gure 3-107: Acquiring and displaying data by Trigger and	
	Repeat modes	3-153
Fi	gure 3-108: Trigger mask	3-155
Fi	gure 3-109: Filling operation for creating a mask	3-156
Fi	igure 3-110: Example mask	3-157
Fi	igure 3-111: Default mask	3-158
Fi	gure 3-112: Changing the position of Point A	3-158
Fi	gure 3-113: Changing the position of Point B	3-159
Fi	gure 3-114: Adding Point C	3-159
Fi	gure 3-115: Trigger point display	3-161
Fi	gure 3-116: Synchronous operation of two analyzers	3-162
Fi	gure 3-117: FFT and RBW process	3-163
Fi	gure 3-118: RBW/FFT menu tree	3-164
Fi	gure 3-119: Process flow when RBW/FFT = FFT	3-166
Fi	gure 3-120: FFT overlap	3-167
Fi	gure 3-121: Generating spurious signals by increasing the	
	number of FFT points	3-168
Fi	gure 3-122: Windowing process of time domain data	3-170
175	gure 3-123: Comparison display of Trace 1 and 2	
F	Baree 1200 comparison anspiraj or frace france 2	3-175
	gure 3-124: Trace/Avg menu structure	3-175 3-176
Fi		
Fi Fi	gure 3-124: Trace/Avg menu structure	3-176
Fi Fi Fi	gure 3-124: Trace/Avg menu structure	3-176 3-181
Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure igure 3-125: Displaying an averaged waveform igure 3-126: Compared display	3-176 3-181 3-182
Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure igure 3-125: Displaying an averaged waveform igure 3-126: Compared display igure 3-127: Relationships between frame, bin, and pixel	3-176 3-181 3-182 3-184
Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structureigure 3-125: Displaying an averaged waveformigure 3-126: Compared displayigure 3-127: Relationships between frame, bin, and pixeligure 3-128: Compression method for displaying the waveform	3-176 3-181 3-182 3-184 3-185
Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure igure 3-125: Displaying an averaged waveform igure 3-126: Compared display igure 3-127: Relationships between frame, bin, and pixel igure 3-128: Compression method for displaying the waveform igure 3-129: Setting the scale in spectrum view	3-176 3-181 3-182 3-184 3-185 3-190
Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure	3-176 3-181 3-182 3-184 3-185 3-190 3-192
Fi Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure igure 3-125: Displaying an averaged waveform igure 3-126: Compared display igure 3-127: Relationships between frame, bin, and pixel igure 3-128: Compression method for displaying the waveform igure 3-129: Setting the scale in spectrum view igure 3-130: Setting the scale in spectrum view igure 3-131: Scale and format settings in spectrogram view	3-176 3-181 3-182 3-184 3-185 3-190 3-192 3-194
Fi Fi Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure	3-176 3-181 3-182 3-184 3-185 3-190 3-192 3-194 3-196
Fi Fi Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure igure 3-125: Displaying an averaged waveform	3-176 3-181 3-182 3-184 3-185 3-190 3-192 3-194 3-196 3-197
Fi Fi Fi Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure	3-176 3-181 3-182 3-184 3-185 3-190 3-192 3-194 3-196 3-197 3-198
Fi Fi Fi Fi Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure	3-176 3-181 3-182 3-184 3-185 3-190 3-192 3-194 3-196 3-197 3-198 3-199 3-200
Fi Fi Fi Fi Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure	3-176 3-181 3-182 3-184 3-185 3-190 3-192 3-194 3-196 3-197 3-198 3-199 3-200 3-201
Fi Fi Fi Fi Fi Fi Fi Fi Fi Fi Fi	igure 3-124: Trace/Avg menu structure	3-176 3-181 3-182 3-184 3-185 3-190 3-192 3-194 3-196 3-197 3-198 3-199 3-200

Figure 3, 140. Eye diagram view	2 202
Figure 3-140: Eye diagram view	3-203
Figure 3-141: Vector and dot displays for AM/AM view	3-204 3-205
Figure 3-142: Vector and dot displays for AM/PM view	
Figure 3-143: Setting the scale in PDF view	3-206
Figure 3-144: Setting the scale in noisogram view	3-208
Figure 3-145: Display line	3-209
Figure 3-146: Lines menu structure (other than Real Time S/A)	3-210
Figure 3-147: Two horizontal lines	3-211
Figure 3-148: Two horizontal and two vertical lines	3-212
Figure 3-149: Multi display lines (Real Time S/A mode)         Diagram (A) = 100 Million (A)	3-213
Figure 3-150: Lines menu structure (Real Time S/A)         Sile	3-214
Figure 3-151: Marker display	3-217
Figure 3-152: Markers menu structure	3-218
Figure 3-153: Measurement with a single marker	3-220
Figure 3-154: Using the delta marker to take measurements	3-221
Figure 3-155: Measurement with the reference cursor	3-222
Figure 3-156: Changing the Trace	3-223
Figure 3-157: Interlocked markers	3-224
Figure 3-158: Marker movement on the Bitmap trace	3-224
Figure 3-159: Functions of the peak search keys	3-225
Figure 3-160: Setting the minimum frequency jump	3-226
Figure 3-161: Online help for the front panel key	3-228
Figure 3-162: Online user manual	3-229
Figure 3-163: Word search using the keyboard	3-230
Figure 3-164: Input menu tree	3-231
Figure 3-165: Load/Save menu structure	3-234
Figure 3-166: Save to the preset file (lower right of the screen)	3-237
Figure 3-167: Selecting the folder	3-238
Figure 3-168: Save to File menu	3-239
Figure 3-169: Load from the preset file (lower right of the screen).	3-240
Figure 3-170: Load from File menu	3-241
Figure 3-171: Displaying the file name on the Load Data side key .	3-241
Figure 3-172: Alphanumeric keypad	3-242
Figure 3-173: Data file structure	3-246
Figure 3-174: Adding dummy frames	3-248
Figure 3-175: Data block	3-250
Figure 3-176: Trace file (.TRC) structure	3-254
Figure 3-177: Bitmap trace file (.DPT) structure	3-256
B	

Figure 3-178: Data sequence for writing to the Bitmap trace file	3-257
Figure 3-179: Print menu structure	3-259
Figure B-1: Menu keys	B-1
Figure C-1: Symbol mapping: GFSK/BPSK, QPSK, 8PSK,	
16QAM, 32QAM, and 64QAM	C-1
Figure C-2: Symbol mapping: 256QAM	C-2

# **List of Tables**

Table 1-1: Span and RBW	1-17
Table 2-1: Key functions of the keyboard	2-7
Table 2-2: Status display	2-10
Table 2-3: Setup display	2-12
Table 2-4: Menu key summary	2-14
Table 2-5: Measurement modes	2-21
Table 3-1: Standard menu in the S/A mode	3-1
Table 3-2: Measurement items in the S/A mode	3-4
Table 3-3: Features of the real-time mode	3-20
Table 3-4: Span setting range	3-20
Table 3-5: Standard menu in the Demod mode	3-28
Table 3-6: Communication standard and parameters	3-45
Table 3-7: Modulation type and available measurement item	3-45
Table 3-8: Compatible RFID standards	3-61
Table 3-9: Measurement items for RFID analysis	3-62
Table 3-10: Standard settings	3-69
Table 3-11: Details of the RF envelope / FSK pulse measurement	
table	3-77
Table 3-12: Measurement result display items	3-81
Table 3-13: Symbol value definition	3-82
Table 3-14: Signal source measurement items	3-103
Table 3-15: Phase noise measurement frequency band	3-108
Table 3-16: Bin width for each decade	3-108
Table 3-17: The number of waveform points	3-109
Table 3-18: Frequency and span setting range	3-124
Table 3-19: Reference level setting range	3-128
Table 3-20: RF attenuation level setting	3-129
Table 3-21: Mixer level setting	3-129
Table 3-22: Vertical scale setting range	3-130
Table 3-23: How to acquire a waveform	3-149
Table 3-24: Trigger level setting range	3-151
Table 3-25: Characteristics and usage of FFT windows	3-169
Table 3-26: FFT window and bandpass filter	3-171
Table 3-27: Trace type	3-177
Table 3-28: Averaging method	3-177

Table 3-29: Number of bins (FFT points: 1024)	3-183
Table 3-30: Bitmap color scheme	3-191
Table 3-31: Front panel keys for online help	3-229
Table 3-32: File save operation	3-236
Table 3-33: Preset file names	3-237
Table 3-34: File load operation	3-239
Table 3-35: Combinations of valid A, P, I and Q	3-251
Table A-1: Frequency	A-2
Table A-2: Spectrum purity	A-3
Table A-3: Noise sideband	A-3
Table A-4: Input	A-4
Table A-5: Amplitude	A-5
Table A-6: Dynamic range	A-5
Table A-7: Spurious response	A-6
Table A-8: Acquisition	A-7
Table A-9: Sampling rate	A-7
Table A-10: Frame time	A-8
Table A-11: Analog demodulation	A-8
Table A-12: RBW (Resolution Bandwidth)	A-9
Table A-13: Trigger	A-10
Table A-14: Measurement function	A-10
Table A-15: Digital demodulation (Option 21 only)	A-11
Table A-16: ACLR measurement (Option 30 only)	A-13
Table A-17: Pulse measurement	A-13
Table A-18: Digital phosphor spectrum processing (DPX)	A-13
Table A-19: Display   Image: Control of the second secon	A-14
Table A-20: Marker, trace, and display line	A-14
Table A-21: Controller and interface	A-14
Table A-22: Power requirements	A-15
Table A-23: Power connector	A-15
Table A-24: Physical characteristics	A-16
Table A-25: Environmental characteristics	A-16
Table A-26: Certifications and compliances	A-17
-	
Table C-1: π/4 shift DQPSK	C-3
Table C-2: GMSK	C-3
	<b>n</b> 4
Table E-1: External inspection checklist	E-1

Table F-1: Power cord identification	F-2
Table F-2: Standard accessories	<b>F-3</b>
Table F-3: Optional accessories	<b>F-3</b>
Table F-4: Accessories for optional software	F-4

Table of Contents

# **General Safety Summary**

Review the following safety precautions to avoid injury and prevent damage to this product or any products connected to it. To avoid potential hazards, use this product only as specified.

Only qualified personnel should perform service procedures.

To Avoid Fire or<br/>Personal InjuryUse Proper Power Cord. Use only the power cord specified for this product and<br/>certified for the country of use.

**Connect and Disconnect Properly.** Do not connect or disconnect probes or test leads while they are connected to a voltage source.

**Ground the Product.** This product is grounded through the grounding conductor of the power cord. To avoid electric shock, the grounding conductor must be connected to earth ground. Before making connections to the input or output terminals of the product, ensure that the product is properly grounded.

**Observe All Terminal Ratings.** To avoid fire or shock hazard, observe all ratings and markings on the product. Consult the product manual for further ratings information before making connections to the product.

**Power Disconnect.** The power switch disconnects the product from the power source. See instructions for the location. Do not block the power switch; it must remain accessible to the user at all times.

**Do Not Operate Without Covers.** Do not operate this product with covers or panels removed.

**Avoid Exposed Circuitry.** Do not touch exposed connections and components when power is present.

**Do Not Operate With Suspected Failures.** If you suspect there is damage to this product, have it inspected by qualified service personnel.

Do Not Operate in Wet/Damp Conditions.

Do Not Operate in an Explosive Atmosphere.

Keep Product Surfaces Clean and Dry.

**Provide Proper Ventilation.** Refer to the manual's installation instructions for details on installing the product so it has proper ventilation.

### Symbols and Terms



**WARNING.** Warning statements identify conditions or practices that could result in injury or loss of life.



**CAUTION.** Caution statements identify conditions or practices that could result in damage to this product or other property.

Terms on the Product. These terms may appear on the product:

**Terms in this Manual.** These terms may appear in this manual:

DANGER indicates an injury hazard immediately accessible as you read the marking.

WARNING indicates an injury hazard not immediately accessible as you read the marking.

CAUTION indicates a hazard to property including the product.

Symbols on the Product. The following symbols may appear on the product:





Ο



Protective Ground (Earth) Terminal

CAUTION Refer to Manual

Mains Disconnected OFF (Power)

Mains Connected ON (Power)

# Preface

This is the user manual for the RSA3303B and RSA3308B Real-Time Spectrum Analyzers. It covers the following information:

- Product features
- Installation
- Operation
- Menu structure
- Accessories

### **About This Manual**

This manual is composed of the following chapters:

- *Getting Started* describes the functions, architecture, installation, and calibration of the analyzer.
- Operating Basics explains the functions of the front, rear, and side panels and the menus. This section also provides tutorials for beginners. It gives step-by-step procedures for measurement using a signal generator.
- *Reference* explains the basic concepts of measurement processes and application-specific operations. It also describes the front panel keys and menus.
- *Appendices* provide additional information including the menu tree, accessories, and cleaning information.

First time users should install the analyzer in *Getting Started*, and then go to *Operating Basics* and perform the tutorials beginning on Page 2–29.

The analyzer uses Microsoft Windows XP as the operating system. This manual does not describe common usage of Windows XP. Refer to your Windows manuals as necessary.

### **Related Documents**

In addition to this user manual, the following documentation is available for your analyzer. Tektronix part numbers are supplied in *Appendix F: Options and Accessories*.

- RSA3000B Series Programmer Manual (standard accessory, PDF) contains an alphabetical listing of the programming commands and other information related to controlling the analyzer over the GPIB interface.
- RSA3303B & RSA3308B Service Manual (optional accessory) describes how to verify the characteristics of, adjust, disassemble, assemble, and troubleshoot the analyzer, and contains the information required for repair, including module replacement, and calibration.

**PDF Manual** The Programmer Manual described above (PDF only) is included in the Documents CD (Tektronix part number 063-4089-XX).

# Optional Software When using the optional software, refer to the following manuals: User Manuals

- RSA3303B, RSA3308B, & RSA3408B Option 24 User Manual provides operating basics for the GSM/EDGE analysis.
- RSA3303B, RSA3308B, & RSA3408B Option 25 User Manual provides operating basics for the cdma2000 forward and reverse link analysis.
- RSA3303B, RSA3308B, & RSA3408B Option 26 User Manual provides operating basics for the cdma2000 1xEV-DO forward and reverse link analysis.
- RSA3303B, RSA3308B, & RSA3408B Option 28 User Manual provides operating basics for the TD-SCDMA analysis.
- RSA3303B, RSA3308B, & RSA3408B Option 30 User Manual provides operating basics for the 3GPP Release 99 and Release 5 uplink and downlink analysis.
- RSA3303B, RSA3308B, & RSA3408B Option 40 User Manual provides operating basics for the 3GPP Release 6 analysis.

### **Conventions**

This manual uses the following conventions:

Front-panel key and control labels are printed in the manual in upper case text. For example, Span, Peak, Print. If it is part of a procedure, the key or control label is printed in boldface. For example:

Press Span.

• To identify keys on the front panel, the area name label precedes the key. For example:

Press the Mode: **Demod** key.

Menu and on-screen form titles are printed in the manual in the same case (initial capitals) as they appear on the analyzer screen, such as Span, Source, and Channel Power. If it is part of a procedure, the menu title is shown in boldface. For example:

Press the **Source** side key.

A series of keys, controls, and/or menu items separated by an arrow symbol
 (→) indicates the order in which to perform the listed tasks. For example:

Select **RBW/FFT**  $\rightarrow$  **Filter Shape...** $\rightarrow$  **Gaussian**.

Preface

# **Getting Started**

# **Product Overview**

The RSA3303B and RSA3308B are portable real-time spectrum analyzers having RF measurement capability up to 3 GHz and 8 GHz, respectively. The real-time spectrum analyzer has a different architecture from traditional tools, and is uniquely capable of capturing continuous, intermittent, or random signals with equal ease. The data acquired seamlessly is analyzed in time, amplitude, phase, and frequency domains correlating with time. Furthermore, DPX spectrum processing provides an intuitive understanding of time-varying RF signals with color graded displays based on frequency of occurrence.

### **Features**

- DC to 3 GHz (RSA3303B) / 8 GHz (RSA3308B) measurement frequency range
- 100 Hz to 3 GHz measurement span and 15 MHz vector span
- Real-time analysis for seamless capture of time-varying RF signals
- DPX spectrum display facilitates observation of intermittent signals
- Spectrogram provides a 3-D representation of time-varying spectrum
- Spectrum analysis of power, ACPR, C/N, OBW, EBW, and spurious signals
- Analog modulation analysis of AM, PM, FM, ASK and FSK signals
- Digital modulation analysis ranging from BPSK to 256QAM (Option 21)
  - Constellation analysis
  - EVM analysis
  - AM/AM and AM/PM distortion analyses
  - RFID analysis
- Time characteristic analysis
  - Pulse measurements
  - Signal source analysis
- CCDF analysis
- 8.4 inch TFT color display and sturdy cabinet
- USB, LAN, and GPIB interfaces

### Application

The RSA3303B and the RSA3308B can perform real-time analysis for the following purposes:

- Signal quality analysis of analog and digital modulation
- Understanding frequency and spectral occupancy behavior over time
- Capture and characterization of undesired, unknown, or interfering signals
- Device/system design or operational diagnostic measurement
- Finding answers to elusive EMI problems
- VCO/synthesizer design
- RFID device characterization
- General purpose digital modulation vector signal analysis (Option)
- Spectrum monitoring
- Radar measurements

### **Difference Between RSA3303B and RSA3308B**

The RSA3303B and the RSA3308B have the same functions except for their measurement frequency ranges:

RSA3303B	 DC to 3 GHz
RSA3308B	 DC to 8 GHz


Descriptions in this manual apply to both the RSA3303B and the RSA3308B unless otherwise noted.

### **Real-Time Analysis**

This section explains real-time analysis with the comparison of a conventional swept spectrum analyzer and a real-time spectrum analyzer.

Conventional Swept Spectrum Analyzer

Figure 1-1 is a block diagram of the conventional swept spectrum analyzer. There are two RF input signals in this example. The RF signal is converted to IF (intermediate frequency) by the swept local oscillator. IF output goes through a bandpass filter, where resolution of the spectrum analyzer is defined.

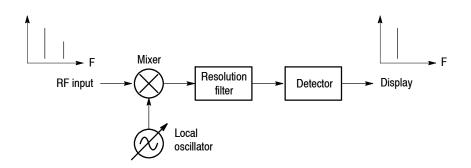

#### Figure 1-1: Concept of the swept spectrum analyzer

The filter is swept from  $F_{start}$  to  $F_{stop}$ . See Figure 1-2. Only signals within the filter bandwidth are observed at one point in time. Signal A is detected first and then signal B is detected and displayed.

**NOTE**. An intermittent signal, such as a burst phenomenon, will not be detected unless it is present at the exact moment that the filter is being swept past it.

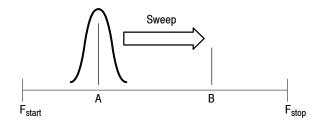


Figure 1-2: Sweep of resolution filter

### Real-Time Spectrum Analyzer

The real-time spectrum analyzer is conceptually equipped with a series of bandpass filters as shown in Figure 1–3. Signals passing through those filters are concurrently observed and recorded continuously. Signals A and B are concurrently acquired and displayed as shown in Figure 1–4.

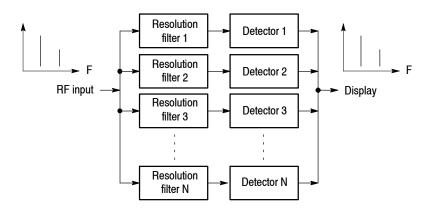


Figure 1-3: Concept of the real-time spectrum analyzer

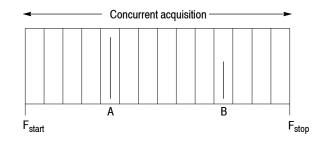


Figure 1-4: Concurrent acquisition

The method used to concurrently acquire signals in a certain frequency range is not bandpass filters, but FFT (Fast Fourier Transformation). The RSA3303B or RSA3308B first acquires a series of frames of data in the time domain, as shown in Figure 1-5, and then performs the FFT process for each frame. This method enables continuous analysis of spectra and ensures the capture of real-time phenomenon such as burst signals in digital mobile communication. The RSA3303B and the RSA3308B are equipped with a 51.2 MHz A/D converter to analyze the spectrum by single scanning for spans up to 15 MHz.

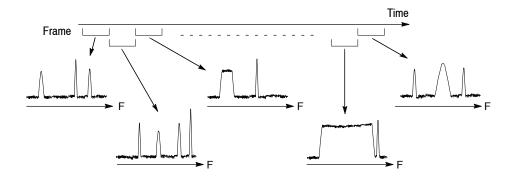


Figure 1-5: Frame acquisition

**DPX Spectrum Display** Along with the real-time processing capability, the analyzer has the proprietary DPX (Digital Phosphor) waveform processing technology that enhances visibility of transient phenomena and intermittent events. The DPX process emulates traditional variable persistence CRTs that store waveforms on the phosphor coating by employing the proprietary DPX engine in the real-time spectrum analyzer. The DPX engine produces more than 50,000 FFT spectrums per second and records the data occurrences for each pixel to make a histogram, which is displayed on the color-graded bitmap. Figure 1-6 illustrates the concept of DPX processing and Figure 1-7 shows an example of the DPX spectrum.


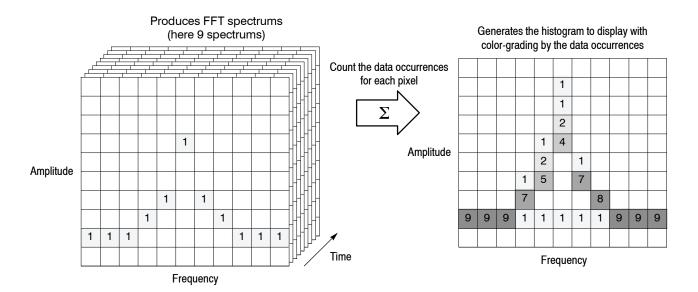


Figure 1-6: DPX processing

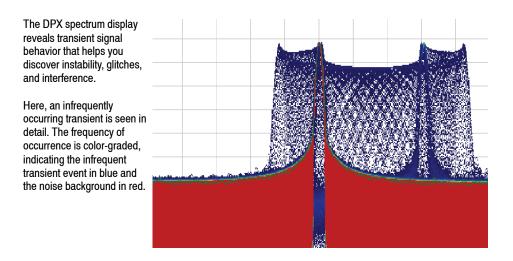

### Figure 1-7: DPX spectrum display

The screen display is generally updated at certain intervals to rewrite waveforms. But the persistence mode in the DPX spectrum display enables a data point once displayed to remain visible for a decay period you specify. You can adjust the persistence time to optimize display characteristics for varying signal conditions, from a live RF view of dynamic signals to the discovery of single occurrences. Persistence is what enables the spectrum analyzer to display events that occur only rarely.

### Architecture

Figure 1-8 shows the block diagram of the signal processing system.

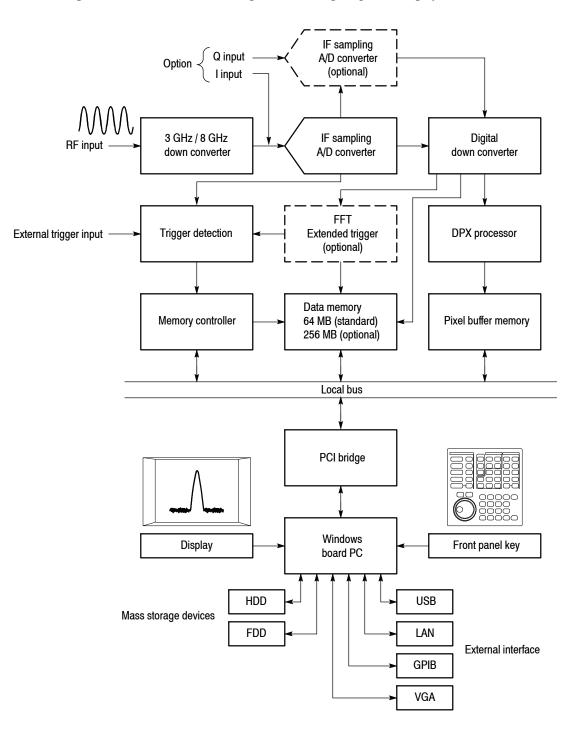


Figure 1-8: Block diagram

3 GHz / 8 GHz Down Converter	Converts the RF signal, input through the front-panel INPUT connector, into a 20 MHz IF signal. This converter processes signals up to 3 GHz for RSA3303B and 8 GHz for RSA3308B by IF conversions in three stages. A 10 MHz oscillator provides a high-accuracy reference for all stages of IF conversion. The converter also adjusts signals for A/D conversion using a low-noise amplifier, fine-tuning attenuator, and anti-alias filter. The output signal from the down converter is sent to the IF sampling A/D converter.
IF Sampling A/D Converter	The IF analog output from the down converter enters the A/D converter through the fine-tuning attenuator, low-noise amplifier, and anti-alias filter, to be converted into a digital signal. The sampling rate of the A/D converter is 51.2 MHz with a resolution of 14 bits.
Digital Down Converter	The digital down converter splits the real signal from the A/D converter into complex (I and Q) components and limits the frequency span of the resulting signal.
	With Option 03, you can enter I and Q signals at this point in the block diagram.
	The down converter performs span and center frequency fine-tuning. This converter consists of two primary stages. In the baseband, the first stage converts a 0 to 15 MHz real signal into a complex signal of $\pm 7.5$ MHz. The second stage converts frequency to set any center frequency.
	A decimating filter between stages changes the span by effectively reducing the sampling rate. A 503 tap FIR filter and four stage comb filter allows extremely accurate filtering with minimal spurious emissions.
	The data streams from the digital down converter are then divided into frames and saved in data memory.
FFT/Extended Trigger (Optional)	Option 02 provides a real-time digital trigger function that monitors the frequency spectrum for the occurrence of specific events. A trigger mask is used to set conditions for trigger.
	The FFT processor performs high-speed calculations to create extended trigger signals. The FFT processor performs 1024-point complex FFT at high speed to create extended trigger signals. The FFT processor consists of an input buffer, FFT calculation DSP, output buffer, and timing control circuit. Performing a 1024-point complex FFT at 12,500 times per second allows the trigger to operate in real time in spans up to 15 MHz.
	Because the trigger comparator is continuously operating at the maximum rate, no events will be missed. The pre-trigger and post-trigger positions can be set as needed; events before and after the trigger event can be measured.

Data Memory	A high-speed 64 MB standard SDRAM stores spectrum data. It is expanded to 256 MB with Option 02. For each data point, I and Q data use 2 bytes. This memory can store 16,000 frames for 1 waveform = 1024 points of analysis, and 64,000 frames for Option 02. For example, W-CDMA communication data can be captured for up to 2.5 seconds on standard and 10 seconds on Option 02. The memory is accessed from the system controller via the ISA/PCI bridge.
DPX Processor/ Pixel Buffer Memory	In the DPX spectrum mode, time-domain data from the digital down converter goes to the proprietary DPX processor that performs real-time FFT over 50,000 times per second. The spectrum traces are accumulated in the pixel buffer memory, where a counter is incremented each time a trace writes to a point on the display. A color is assigned to each display point based on the value of its counter. Thus, as acquisitions occur over time, a color-graded waveform, the Bitmap, develops on the display that shows the frequency of data occurrence (signal density).
Windows Board PC	The system controller board has an Intel Pentium III CPU. It runs on Win- dows XP, and controls menu operation from the front panel keys. It is equipped with a 20 GB hard disk and 3.5 inch disk drive to store data and settings. Waveforms, menus, and measurement results are shown on the color display, which uses an 8.4 inch XGA TFT-LCD module.

The standard analyzer has the following external interfaces:

- USB (for mouse, keyboard, and/or printer)
- LAN (Ethernet 10/100BASE-T)
- GPIB
- VGA (for an external monitor)

Product Overview

# Installation

This section describes how to install this instrument. The topics are organized as follows:

- Unpacking to check contents
- Applying the power
- Setting up the stand
- Functional check
- Powering off the analyzer
- Restart
- Backing up the user file

Before starting installation, you should become familiar with the *General Safety Summary* on page xvii.

### **Unpacking to Check Contents**

- **1.** This product is packed in a cardboard box for delivery. Before opening the box, make sure that there is no damage on the surface.
- 2. Open the box, check that the product has no damage and that all the standard accessories are found inside. For a list of accessories, refer to *Standard Accessories* on page F-3. If you find any damaged or missing components, contact your local Tektronix representative.
- **3.** It is recommended to keep the box and packing materials. You may need them to send this product to Tektronix for calibration or repair.



**CAUTION.** The analyzer has exhaust fans on the side panel. Leave a space of at least 5 cm (2 in) on both sides for proper air circulation.

### **Applying Power**

Power on the analyzer using the following procedure:

**AC Power Requirements** The analyzer operates from an AC line frequency of 47-63 Hz, over the range of 90-250 Volts, without the need for configuration, except the power cord. Refer to page F-2 for the power cord options.


The maximum power consumption is 350 W. Refer to *Appendix A, Specifications,* for additional information on power and environmental requirements.



**CAUTION.** Use only power cords that are approved for the country of use. Using non-approved power cords could result in fire or shock hazard.

### Connecting the Power Cord

1. Plug the power cord into the AC input on the rear panel.

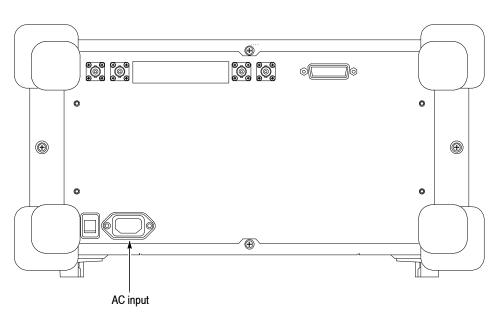

#### Figure 1-9: AC Input (rear panel)

2. Connect the plug of the power cord to a properly grounded outlet.

#### Turning on the Analyzer

1. Turn on the principal power switch on the rear panel.

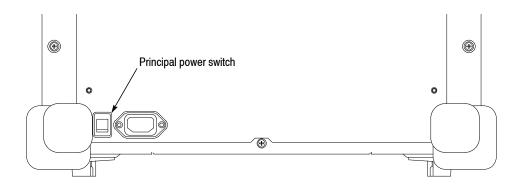


Figure 1-10: Principal power switch (rear panel)

When you turn on the principal power switch, a voltage is applied to the analyzer standby circuit. Make sure the LED next to the power switch on the front panel lights up in orange.

**2.** Turn on the power switch (ON/STANDBY) on the lower left of the front panel. The LED next to the power switch changes to green.

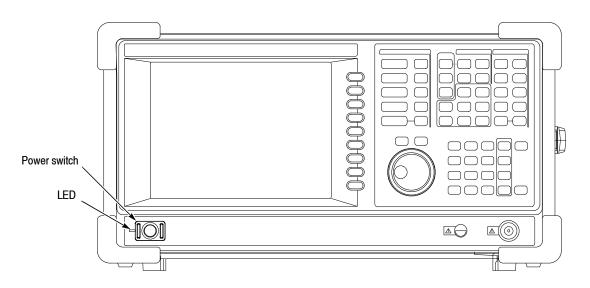


Figure 1-11: Front panel power switch (ON/STANDBY switch)

When you turn on the analyzer, Windows XP boots up. After several minutes, the analyzer application starts up.

The initial screen appears as shown in Figure 1-12. The displayed spectrum represents the noise floor of the analyzer.

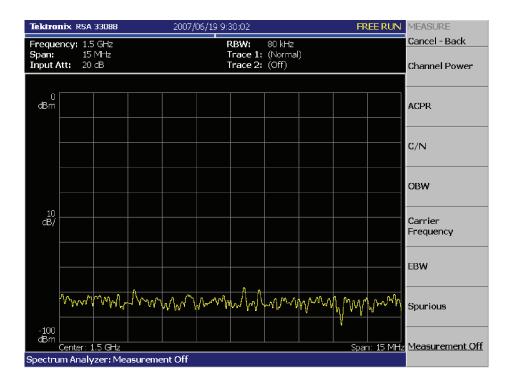


Figure 1-12: Initial screen

If "UNCAL" is displayed on top of the screen, run the gain calibration routine (refer to *Calibrating Gain* on page 1-25).



**CAUTION.** Never apply signals with a combined amplitude greater than +30 dBm to the RF INPUT connector. If you exceed this input rating, you can permanently damage the analyzer. (The RF Input connector is shown in Figure 1–13.)

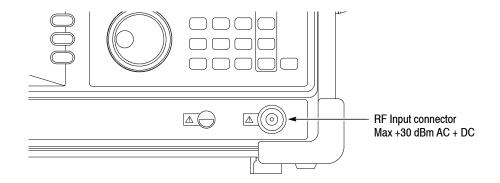


Figure 1-13: RF Input connector

### Setting Up the Stand

To set up the stand, place the analyzer on a table. Lift the front of the analyzer and pull out the stand until it is perpendicular to the analyzer.

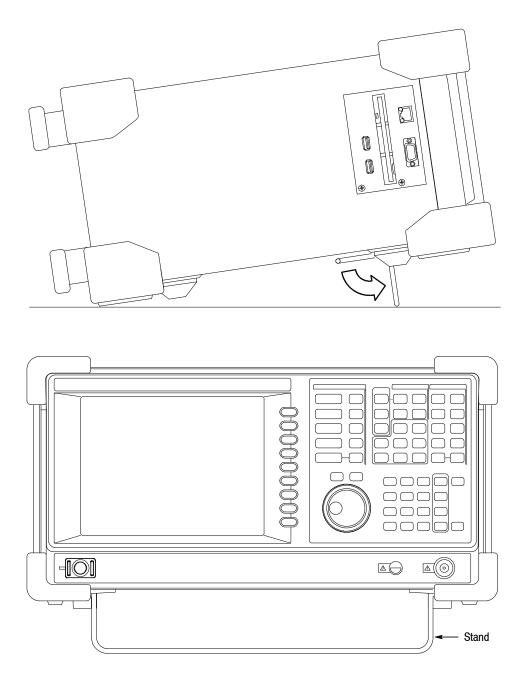


Figure 1-14: Setting up the stand

### **Functional Check**

The analyzer has a built-in calibration signal source with frequency of 50 MHz and amplitude of approximately -20 dBm. Using this source, perform this quick functional check to verify that your instrument is operating correctly.

- 1. Turn on the analyzer.
- 2. Display the spectrum of the calibration signal:
  - **a.** Press the **S**/**A** key on the front panel and then press the **Spectrum Analyzer** side key.
  - **b.** Press the **Preset** key on the front panel to reset the analyzer.
  - c. Press the Input key on the front panel.
  - d. Press the Signal Input Port... side key to select Cal.

The spectrum of the calibration signal appears.

e. Check that "INPUT: CAL" and "FREE RUN" are displayed in the status indicator at the upper right of the screen (see Figure 1-15).


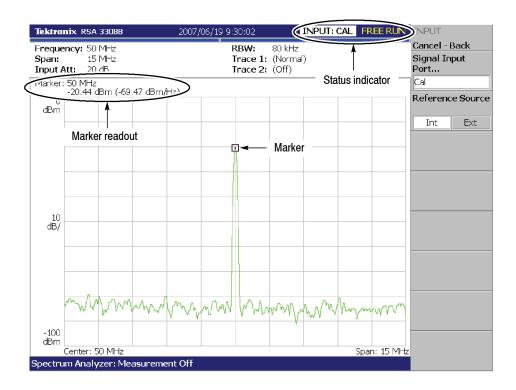


Figure 1-15: Spectrum of the calibration signal (50 MHz, about - 20 dBm)

- 3. Check the center frequency and peak amplitude using the marker:
  - **a.** Press the **Peak** key on the front panel to place the marker on the peak (see Figure 1-15).
  - **b.** Check the marker readouts on screen. The frequency should be 50 MHz and the amplitude should be approximately -20 dBm.
  - c. Press the Marker Setup key on the front panel and then the Markers side key to select Off. Check that the marker disappears.
- 4. Check the RBW (Resolution Bandwidth) while changing the span setting.
  - **a.** Press the **Span** key on the front panel.
  - **b.** Confirm that the span is 15 MHz and the RBW is 80 kHz in the setup display on the upper part of the screen (see Figure 1-16).

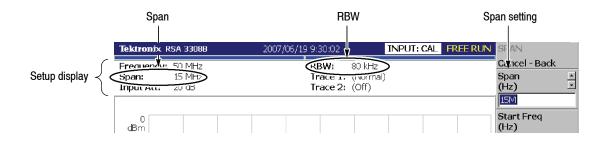

#### Figure 1-16: Setup display

**c.** Using the general purpose knob, change the span setting as listed in Table 1-1 and check that the RBW is displayed correctly.

#### Table 1-1: Span and RBW

Span	RBW
15 MHz	80 kHz
5 MHz	20 kHz
100 kHz	500 Hz
1 kHz	20 Hz

**d.** Using the numeric keypad, set the span back to 15 MHz. (Press  $1 \rightarrow 5 \rightarrow MHz$ , in that order, on the keypad.)

- **5.** Check the reference level:
  - a. Press the Amplitude key on the front panel.
  - **b.** Make sure that the reference level is set to 0 dBm with the **Ref Level** side key. Check that 0 dBm is displayed on the upper left side of the graticule (see Figure 1-17).
  - c. Use the general purpose knob to set the reference level to -30 dBm.
  - **d.** Confirm that "Overrange increase RefLev or Atten" is indicated in the red box at the top center of the screen. Make sure that -30 dBm is displayed on the upper left side of the graticule and that the spectrum waveform is distorted as shown in Figure 1-17.

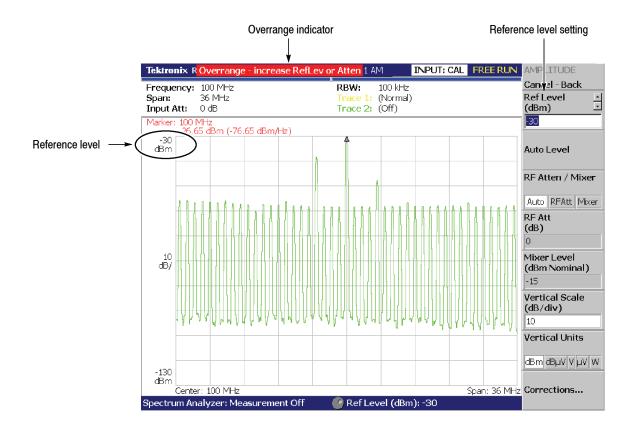


Figure 1-17: Reference level setting and overrange indicator

e. Using the numeric keypad, set the reference level back to 0 dBm. (Press  $0 \rightarrow Enter$ , in that order, on the keypad.)

- **6.** Check the spectrogram display:
  - **a.** Press the **S**/**A** key on the front panel.
  - **b.** Press the **S/A with Spectrogram** side key. Check that the spectrogram is displayed on the lower side of the screen (see Figure 1-18).
  - **c.** Press the **Run/Stop** key on the front panel to stop data acquisition. Confirm that the trace display freezes and PAUSE is displayed in the status indicator at the top right of the screen.

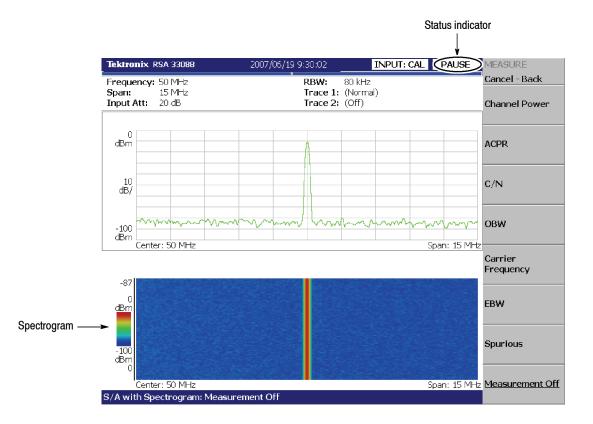


Figure 1-18: Spectrogram display

### **Powering Off the Analyzer**

Turn off the power switch on the front panel.



**CAUTION.** When you power on or off the analyzer, you must use the front-panel power switch. Failure to do so may cause the operating system to shut down improperly.

When you power on the analyzer again, wait at least 10 seconds after the last power off.

When you press the front-panel ON/STANDBY switch, the analyzer starts a shutdown process (including a Windows shutdown) to preserve settings and then power off. The LED next to the power switch changes to orange. Avoid using the rear-panel power switch or disconnecting the power cord to power off the analyzer.

To completely remove power to the analyzer, perform the shutdown just described, and then set the power switch on the rear panel to off.

**NOTE**. Turning off the front-panel power switch does not shut down the principal power supply completely. To turn off the principal power supply, press the principal power switch on the rear panel. Turning off the principal power switch turns off the front-panel LED. When you do not use the analyzer for a long time or in case of emergency, you should unplug the power cord.

### Restart

When the analyzer operates abnormally, use the following procedure to turn the analyzer off and on again.

**NOTE**. When the analyzer operates abnormally, it will not be shut down by turning off the front-panel power switch alone.

1.	Make sure that the front-	-panel power switch	is in the turned-off position.
	mane sare mat me nom	panel ponel on ten	is in the turned on position.

- 2. Turn off the principal power switch on the rear panel.
- 3. Wait at least 10 seconds and then turn on the principal power switch again.
- 4. Turn on the front-panel power switch.

**When Scan Disk Appears** If the analyzer was not shut down properly, Windows Scan Disk may run when you turn on the analyzer. When the Scan Disk screen appears, wait until the Scan Disk is completed. If an error is detected, refer to the Windows manual for treatment. For more information on accessing Windows on the analyzer, refer to *Using Windows XP* on page 2–26.

When the Display<br/>Brightness Is Not EvenIt is a characteristic of the LCD (Liquid Crystal Display) panel to sometimes<br/>have uneven brightness, dead pixels (dots that never turn on) or stuck pixels<br/>(dots that always stay on). This is neither a malfunction nor a defect, and not a<br/>cause for repair or exchange.

**Restoring the Instrument** Refer to Appendix D: Restoring the Instrument Operating System. **Operating System** 

### **Backing Up User Files**

You should back up your user files on a regular basis in case of system failures. The Back Up tool is located in the System Tool folder in the Windowss Accessory folder. Start this tool to select the files and folders to back up. For more information, use Windows online help. For more information about access to Windows on the analyzer, refer to *Using Windows XP* on page 2–26.

The following files should be backed up frequently:

- Status files (\*.sta)
- Data files (\*. iqt)
- Trace files (\*.trc)
- Correction files (\*.cor)

Refer to page 3-233 for details on file operations.

**Using LAN** The analyzer is equipped with LAN Ethernet interface as standard, allowing you to save data in peripheral devices such as other PCs, hard disks, and MO via network. Refer to page 2–5 for details on connecting to LAN.

### About Installation of Other Applications

The analyzer incorporates Windows XP as the operating system. Some combinations of internal measurement applications and external applications may cause deterioration in the basic performance or conflicts between these applications.

It is not recommended that you install other applications, including Microsoft Word, Excel, and Outlook, on the analyzer. If you install an external application, you do so at your own risk, keeping in mind that it may lower the performance of the analyzer.

# Calibration

Perform these routines to optimize the analyzer performance:

- Gain calibration
- Center offset calibration
- DC offset calibration
- Display brightness adjustment

Each item is explained in this section.

### **Cal Menu**

Use the Cal key to calibrate the analyzer.

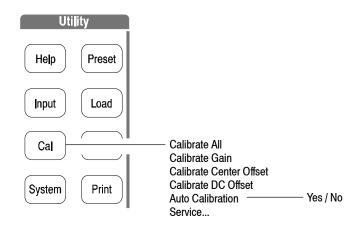


Figure 1-19: Calibration menu structure

The Cal menu contains the following controls:

Calibrate All. Performs all possible calibration operations.

**Calibrate Gain.** Calibrates the internal gain steps. Refer to page 1-25 for the gain calibration.

**Calibrate Center Offset.** This calibration cancels a center offset. Refer to page 1–26 for the center offset calibration.

**Calibrate DC Offset.** This calibration cancels DC offset in baseband. Refer to page 1-27 for the DC offset calibration.

**Auto Calibration.** Determines whether to automatically perform all possible calibration operations. The default setting is On.

**Service...** This menu item is only to be used by qualified personnel for repair and calibration. Refer to the *RSA3303B and RSA3308B Service Manual* (optional accessory) for details.

If you want to perform all possible calibration operations together, press the **Calibrate All** side key. When **Auto Calibration** is set to Yes, they will be performed automatically any time the analyzer is in an uncal state.

**NOTE**. When you run the calibration during signal acquisition, the calibration starts after the acquisition is completed.

### **Calibrating Gain**

The gain calibration calibrates the analyzer's amplifier gain using the internal signal generator. Run this internal calibration routine as required when you boot the analyzer or when UNCAL (uncalibrated) is displayed during operation.

Allow the analyzer to warm up for 20 minutes before you begin the calibration. The warm-up period allows electrical performance of the analyzer to stabilize.

During normal operation, when the ambient temperature changes by more than  $\pm 5$  °C from the temperature at the last calibration, UNCAL is displayed in the yellow box at the top of the screen (see Figure 1-20). If this happens, run the gain calibration.

	When UNCAL is displayed, run the gain calibration 		
	↓		
	UNCAL	FREE RUN	
Frequency Span Input Att			

#### Figure 1-20: UNCAL display

To run the gain calibration, do the following:


- 1. Press the Cal key on the front panel.
- 2. Press the Calibrate Gain side key.

The calibration runs. It takes several seconds to complete the process.

### **Calibrating Center Offset**

When you display a spectrum and no input signal is present, a spurious emission at the center frequency may appear regardless of frequency settings. The center offset calibration cancels those spurious emissions. If the spurious emission is too obvious when you narrow the span, run the calibration.

**NOTE**. Option 03 only. When you input I and Q signals from the rear panel connectors, set the IQ input signal level to zero externally.

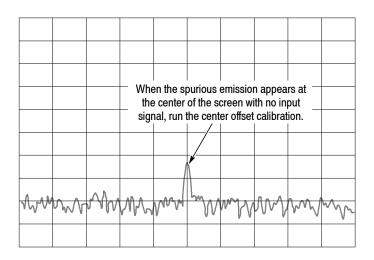

#### Figure 1-21: Center offset

To run the center offset calibration, do the following:

- 1. Press the Cal key on the front panel.
- 2. Press the Calibrate Center Offset side key.

The calibration runs. It takes several seconds to complete the process.

### **Calibrating DC Offset**

The DC offset calibration cancels DC offset that appears at 0 Hz in the baseband (DC to 20 MHz). When you change the amplitude setting and the DC offset is too obvious, run the DC offset calibration.

	When DC offset appears at 0 Hz in								
-									
WW	$\mathcal{W}$	M		MwM	~ <u>~</u> //~	W	$\mathcal{M}$	m	ms.

### Figure 1-22: DC offset

To run the DC offset calibration, do the following:

- 1. Press the Cal key on the front panel.
- 2. Press the Calibrate DC Offset side key.

The calibration runs. It takes several seconds to complete the process.

### **Adjusting Display Brightness**

Adjust the brightness of the display according to your environment.

1. Press the System key on the front panel (see Figure 1-23).

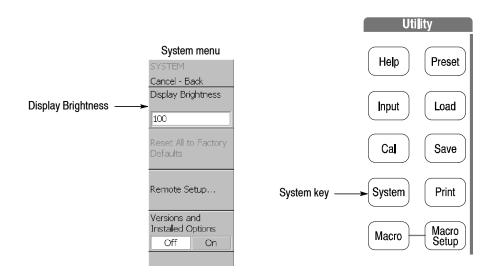


Figure 1-23: System menu

- 2. Press the Display Brightness side key.
- **3.** Turn the general purpose knob to adjust the brightness. The setting range is 0 to 100.

### **Confirming Performance**

The electrical characteristics described in *Appendix A: Specifications* can be checked only by our service personnel. If you need any service, contact your local Tektronix representative.

# **Operating Basics**

# **Interface Maps**

This section describes the controls, connectors, and display:

- Controls and connectors
  - Front panel
  - Rear panel
  - Side panel
  - Using a mouse and keyboard
- Display screen
  - Elements of the display
  - Status display
  - Front panel key lock
  - Setup display

For the analyzer specifications, refer to Appendix A: Specifications.

### **Controls and Connectors**

Figures 2-1 through 2-3 on the following pages show the controls and connectors on the front, the side, and the rear panels.

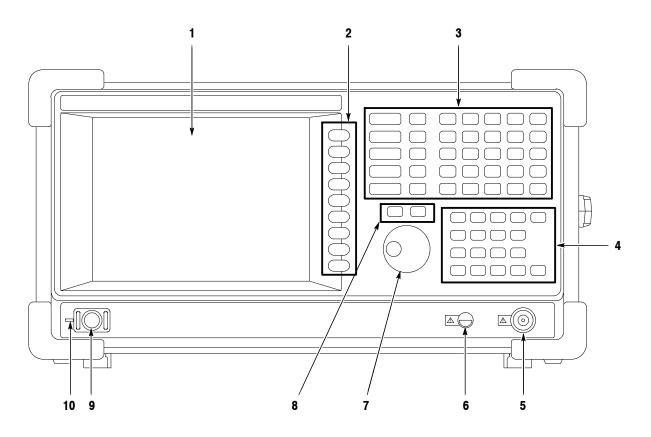


Figure 2-1: Front panel

#### **Front-Panel Interface**

 Display. LCD (Liquid Crystal Display). Size: 21.3 cm (8.4 inch) Resolution: 800 × 600 dots Color: 256 colors maximum

- 2. Side Keys. Select menu items associated with menu keys.
- 3. Menu keys. Select menus. Refer to page 2-13 for details.

Refer to page 2-16 for menu operations.

**4. Keypad.** Enters alphanumeric characters. Refer to page 2-17 for entering a numeric value and page 3-242 for entering a file name.

#### **Finput connector**. Connects an input signal. Type: N-type connector Input impedance: 50 Ω Maximum capacity of non-breakdown input: 30 dBm




**CAUTION.** Do not apply a signal of more than +30dBm; to do so may damage the instrument.

- **6. Preamplifier power source.** Provides power source for a preamplifier (optional accessory). Refer to Table A-23 on page A-15 for specifications.
- 7. General purpose knob. Changes a setting.
- 8. Up/Down keys. Increases or decreases a value.

Refer to page 2-18 for changing a setting using the general purpose knob or the up/down keys.

- 9. Power switch (ON/STANDBY). Refer to page 1-12 for applying power.
- 10. LED. Green on operating, orange on standby.

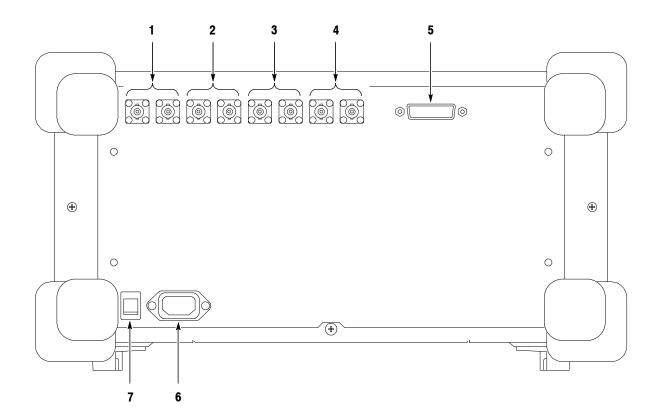

#### Figure 2-2: Rear panel

#### Rear-Panel Interface

- 1. **REF IN/OUT connector.** A 50  $\Omega$  BNC connector for reference signal input/output. Refer to page A-15 for specifications.
- 2. I+/I- INPUT connector (Option 03 only). A 50  $\Omega$  BNC connector for I signal differential input. When using one connector as a single-ended input, terminate the other end in 50  $\Omega$ .
- 3. Q+/Q- INPUT connector (Option 03 only). A 50  $\Omega$  BNC connector for Q signal differential input. When using one connector as a single-ended input, terminate the other end in 50  $\Omega$ .
- 4. TRIG IN/OUT connector. A 50  $\Omega$  BNC connector for trigger signal input/output. Refer to page A-10 for specifications.
- **5. GPIB connector.** Used to control the analyzer from an external controller. Refer to the RSA3303B and RSA3308B *Programmer Manual* for GPIB control.
- 6. AC input. Connect AC power cable.
- 7. **PRINCIPAL POWER SWITCH.** When this switch is on, the internal standby circuit is energized. Refer to page 1-12 for a power on procedure.

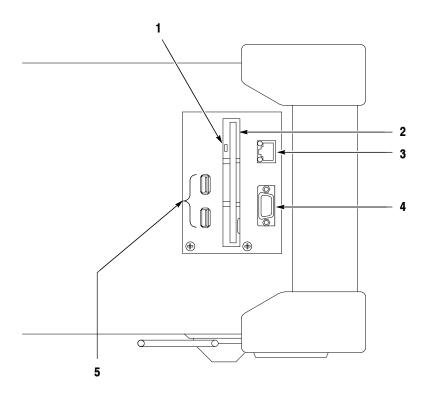


Figure 2-3: Side panel

### Side-Panel Interface

1. Indicator. Indicates whether the disk drive is active.



**CAUTION.** When this indicator is lit, you must not remove the disk from the disk drive. If you do so, the data stored on the medium may be destroyed or an error may occur.

- **2.** Floppy disk drive. A 3.5-inch 2HD (1.44MB) or 2DD (720KB) disk formatted for MS-DOS can be used to save and load data and settings.
- **3.** LAN Ethernet connector. 10/100BASE-T connector. When you reboot the analyzer after connecting to LAN, the analyzer recognizes the network speed automatically and sets it to 10 Mbps or 100 Mbps. You can share resources such as files or disks on the network.


**NOTE**. The Windows XP default network settings are configured when the instrument is shipped from the factory. Contact your system administrator for information about setting network parameters in your LAN.

For using Windows XP on the analyzer, refer to page 2-26.

**4.** VGA output connector. 15 pin D-sub connector (female) to send the display of this instrument to another monitor.

**NOTE**. When you use an external monitor, connect a VGA cable to the analyzer and to the monitor, and then turn on the monitor before you turn on the analyzer.

5. USB connectors (two ports). Connect USB devices such as a mouse, keyboard, and printer. You can connect the devices any time, with no setup (plug & play), to either or both ports. You can also connect a mouse to a USB port on a keyboard (see Figure 2-4).

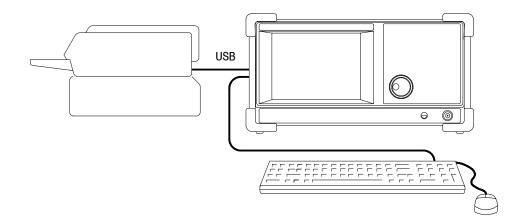



### Using a Mouse and Keyboard

You can operate the analyzer using the standard accessory mouse and keyboard instead of the side keys and the front panel keypad.

The mouse and keyboard operations are as follows:

- Click the menu item instead of pressing the side keys.
- If the menu item has arrow buttons, click them to select the value.



#### Figure 2-5: Operation with mouse and keyboard

Key	Purpose	Function
Numeric keys	Numeric input	Enter a numeric value in a numeric input field.
Left/Right arrow keys	Caret move	Moves the caret in an alpha or numeric input field.
Home	Caret move	Moves the caret to the beginning of an input field.
End	Caret move	Moves the caret to the end of an input field.
Backspace	Alphanumeric input	Deletes the character before the caret.
Delete	Alphanumeric input	Deletes the character after the caret.
ESC	Alphanumeric input	Aborts a numeric entry and restores the original value.
ENTER	Alphanumeric input	Accepts a value in the input field.
K or k key	Alphanumeric input	Kilo (10 <sup>3</sup> ). Press Enter to complete entry of the value.
M key	Alphanumeric input	Mega (10 <sup>6</sup> ). Press Enter to complete entry of the value.
G or g key	Alphanumeric input	Giga (10 <sup>9</sup> ). Press Enter to complete entry of the value.
m key	Alphanumeric input	milli (10 <sup>-3</sup> ). Press Enter to complete entry of the value.
U or u key	Alphanumeric input	micro (10 <sup>-6</sup> ). Press Enter to complete entry of the value.
N or n key	Alphanumeric input	nano (10 <sup>-9</sup> ). Press Enter to complete entry of the value.

#### Table 2-1: Key functions of the keyboard

### **Display Screen**

Figure 2-6 shows the elements of the display screen.

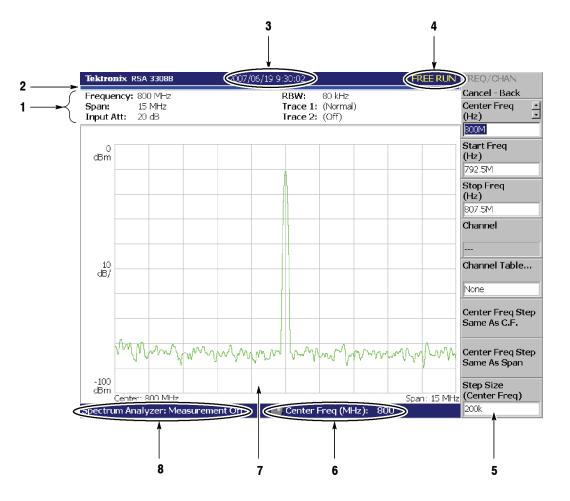


Figure 2-6: Display screen configuration

#### **Elements of the Display**

- 1. Setup display area. Displays the current hardware value. Refer to *Setup Display on* page 2-12.
- 2. **Progress bar.** Indicates the progress of the acquisition cycle on the left bar and the measurement cycle on the right bar. The progress fills up in blue from left to right.
- 3. Date/Time display area. Shows the current date and time.
- **4. Status display area.** Shows the trigger status. Refer to *Status Display* on page 2-10.
- **5.** Side menu display area. When you press a menu key on the front panel, the menu associated with that key is displayed. Refer to page 2-13 for details of menu items.
- 6. Menu setting display area. Displays the last setting of the menu item that can be set with the general purpose knob.
- 7. View. The View window displays the waveform or the measurement results. Multiple views can be displayed on one display screen, depending on the measurement mode.

For more information, refer to:	
Spectrum Analysis (S/A mode)	page 3-1
Modulation Analysis (Demod mode)	page 3-27
Time Analysis (Time mode)	page 3-83
Setting Views	page 3-187

**8. Measurement function display area.** Displays the measurement function currently in use (the settings of the Mode and Measure menus).

#### **Status Display**

The status display area in the upper right side of the screen (see Figure 2-7) shows the instrument status as listed in Table 2-2.


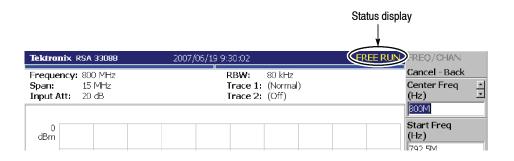



Figure 2-7: Status display

#### Table 2-2: Status display


ltem	Description
ARM	The pre-trigger portion of the acquisition record is filling. A trigger event occurring during this state will not be recognized.
READY	Pre-trigger data has been acquired, and the instrument is waiting for a trigger event.
TRIG'D	Pre-trigger data has been acquired, and a trigger event has been detected. The instrument is now acquiring post-trigger data.
FREE RUN	The instrument acquires and measures without waiting for a trigger event.
PAUSE	The user has temporarily stopped acquisition/measurement cycling.

The acquired data is stored in the data memory from address zero in order of acquisition. When you set a trigger condition, the acquired data is stored in the pre-trigger region until the trigger event occurs. Thereafter, it is stored in the post-trigger region (see Figure 2–8).

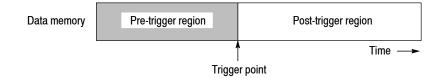

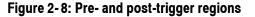

**Front Panel Key Lock** When controlling this instrument through GPIB, you can disable all the keys on the front panel except the power switch using the :SYSTEM:KLOCK command. At this time, the message "PANEL LOCK" is displayed on the top side key (see Figure 2-9).

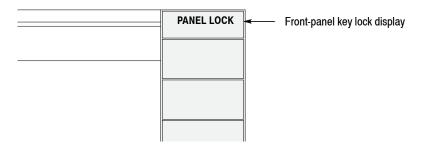

#### Figure 2-9: Key lock display

To cancel the key lock, the following two methods can be used:

- Use the :SYSTEM:KLOCK command to cancel.
- Turn off the power and then on.

Refer to the *RSA3300B Series & RSA3408A/B Programmer Manual* for information about the GPIB commands.

**Setup Display** The setup display area in the upper part of the screen shows the analyzer hardware settings (see Figure 2-10). The contents differ depending on the measurement mode: spectrum analysis (S/A), modulation analysis (Demod), or time analysis (Time), as shown in Table 2-3.

	Tektronix RSA 33088		2007/0	06/19 9	9:30:02	
Setup display {	Frequency: Span: Input Att:	15 MHz			RBW: Trace 1: Trace 2:	)
	0 dBm					

Figure 2-10: Setup display

#### Table 2-3: Setup display

ltem	Description	Mode
Frequency	Indicates the center frequency.	All
Span	Indicates the span.	
Input Att.	Indicates the attenuation of an input signal before it enters the internal mixer.	
RBW	Indicates RBW (Resolution Bandwidth) for compatibility with swept spectrum analyzers. Refer to <i>FFT and RBW</i> on page 3-163.	S/A
NBW	Indicates NBW (Noise Bandwidth) instead of RBW when FFT-processed data does not go through RBW process.	
Trace 1 and 2	Indicates the Trace 1 and 2 trace type.	
Spectrum Length	Indicates time length of a 1024-point FFT frame. It is determined by the span.	Real Time S/A
Spectrum Interval	Indicates time interval between FFT frames. Refer to <i>FFT Start Point</i> on page 3-167.	
Acquisition Length	Indicates time to acquire a block of data. It can be set in the Acquisition/Analysis menu.	Demod and Time

# **Functional Overview**


This section provides the operating fundamentals:

- Menu Keys
- Menu Operations
- Measurement Basics
- Setting System Parameters
- Using Windows XP

## Menu Keys

Menu keys on the front panel (shown in Figure 2-11 on page 2-14) are divided into the following three functional groups:

- Settings. Sets frequency, amplitude, and time parameters for specified measurement and controls data acquisition.
- Mode. Selects a measurement mode.
- View. Selects a view content and scales a measurement graph.
- Markers. Controls markers.
- Utility. Provides system initialization, waveform storage facilities, instrument calibration, screen hard copy, and miscellaneous other functions.

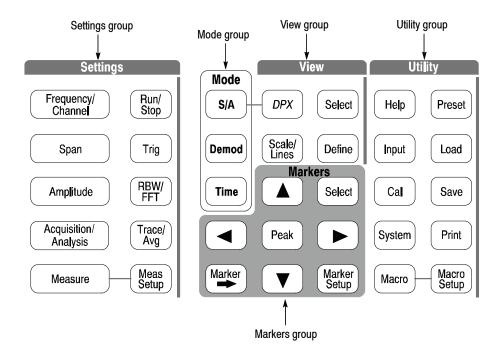

#### Figure 2-11: Menu keys

Table 2-4 summarizes the front-panel menu key functions. For details, refer to each section listed in the Reference column.

Menu group	Menu key	Function	Reference	
Settings	Frequency /Channel	Sets frequency or channel.	Setting Frequency and Span (page 3-119)	
	Span	Sets span.		
	Amplitude	Sets amplitude.	Setting Amplitude (page 3-127)	
	Acquisition/Analysis	Sets acquisition and analysis parameters.	Setting Acquisition and Analysis Parameters (page 3-141)	
	Run/Stop	Start or stop data acquisition.	Starting/Stopping Data Acquisition (page 2-22)	
	Trig	Controls trigger.	Trigger (page 3-147)	
	RBW/FFT	Sets RBW and FFT parameters.	FFT and RBW (page 3-163)	
	Trace/Avg	Controls trace display and average function.	Trace Comparison Display and Average Function (page 3-175)	

Table 2-4: Menu key summary

Menu group	Menu key	Function	Reference		
Settings	Measure	Selects a measurement item.	Spectrum Analysis (S/A Mode)		
	Meas Setup	Sets measurement parameters.			
Mode	S/A	Selects spectrum analysis.	(page 3-27) Time Analysis (Time Mode)		
	Demod	Selects modulation analysis.			
	Time	Selects time analysis.	(page 3-83)		
View	DPX	Displays the DPX spectrum.	DPX Spectrum Display (page 3-14)		
	Select	Selects a view on screen.	Spectrum Analysis (S/A Mode) (page 3-1)         Modulation Analysis (Demod Mode) (page 3-27)         Time Analysis (Time Mode) (page 3-83)         DPX Spectrum Display (page 3- Setting Views (page 3-187)         Display Line (page 3-209)         Marker Operation and Peak Sear (page 3-217)         Using Online Help (page 3-227)         Selecting Input Source (page 3-2         Calibration (page 1-23)         Setting System Parameters (page 2-23)         Restoring Default Settings (page 3-233)         Screen Copy (page 3-259)         For installing macros, contact you		
	Define	Sets view format.	Display Line (page 3-209)		
	Scale/Lines	Scales a graph and controls display lines.			
Markers	Peak	Searches for peak on waveform.	Marker Operation and Peak Search (page 3-217)		
	▲▼◀►	Move the marker to the next peak.			
	Select	Selects a marker.			
	Marker →	Sets parameters based on the marker position.			
	Marker Setup	Sets marker parameters.			
Utility	Help	Displays online help.	Using Online Help (page 3-227)		
	Input	Selects the input source.	Selecting Input Source (page 3-231)		
	Cal	Calibrates the analyzer.	Calibration (page 1-23)		
	System	Controls system parameters.			
	Preset	Returns settings to the factory defaults.			
	Load	Load data from a file.	File Operations		
	Save	Save data to a file.			
	Print	Prints the screen image.	Screen Copy (page 3-259)		
	Macro	Shows a menu listing all macros.	For installing macros, contact your local Tektronix sales office.		
	Macro Setup	Shows a menu to configure macros.			

Table 2-4: Menu key summary (Cont.)

## **Menu Operations**

This section describes basic operations of the analyzer menus and how to select menu items and input numeric values.

Menu Item Information Up to nine soft keys can be displayed down the right side of the screen (see Figure 2-12). Cancel - Back is always displayed at the top, and the other eight keys select menu items.

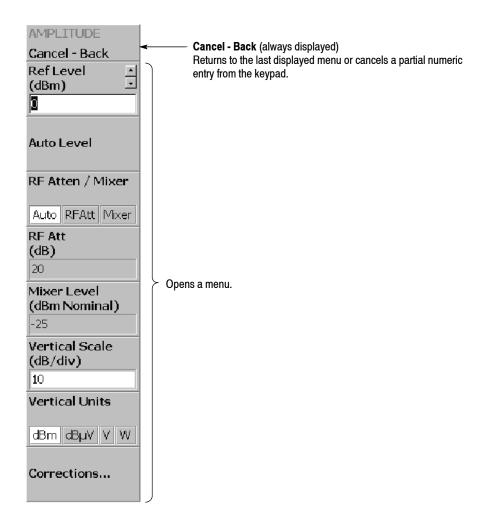


Figure 2-12: Example menu item display

**NOTE**. When the setting is prohibited or is not available, the item is displayed in gray.


#### **Menu Item Types**

The different types of menu items are shown in Figure 2-13.

Ref Level (dBm)	<b>Numeric entry</b> The current value of parameter is displayed. To change the value, press the associated side key and use the general purpose knob, up/down keys, or the keypad.
RF Atten / Mixer	<b>Toggle</b> You can switch selection items by pressing the associated side key.
Channel Power	Function execution The function indicated on the label is executed by pressing the associated side key. In this example, the "Channel Power" measurement is performed.
Corrections	Move to sub-menu If the label is followed by "", you can move to a lower-level menu by pressing the associated side key.
Go to page 2 (of 2)	Move between pages When the number of menu items is greater than seven, this menu appears. To move to the next page, press the associated side key. If you are on the last page, you will return to the first page.
	Invalid Menu items without labels do not cause any action.

#### Figure 2-13: Menu item types

**Numeric Input** An example numeric input field is shown in Figure 2-14. In this type of field, you can change the numeric value by turning the general purpose knob, by pressing the up/down ( $\Delta \nabla$ ) keys, or by entering a value using the keypad.

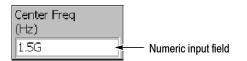



#### Changing a Setting Using the General Purpose Knob or the Up/Down Keys.

1. Press the side key for setting a numeric value. For example, press Frequency/Channel → Center Freq to set a center frequency.

The menu item changes to the display shown in Figure 2-15.


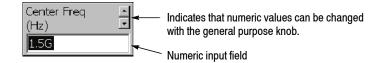


Figure 2-15: Changing value with the knob

2. Turn the general purpose knob to increase or decrease the value.

You can also use the up and down  $(\blacktriangle \nabla)$  keys to increase or decrease a setting value, respectively.

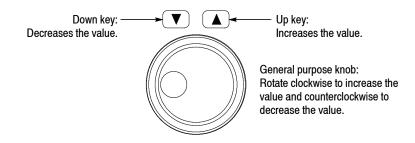

The up and down keys have the same functions as the general purpose knob except the step size (the amount per click by which the general purpose knob changes a setting value or the amount per press for the up and down keys) as follows:


- For the general purpose knob, the step size is determined internally. You cannot change the step size.
- For the up and down keys, the step size is set with the **Step Size** side key. Refer to *Changing the Step Size* on page 2-20 for the detail.

The changed value is immediately reflected on the analyzer settings and displays.

**Entering a Value Using the Keypad.** You can enter values using the front panel keypad, shown in Figure 2–17.

1. Press the side key for setting a numeric value. For example, press Frequency/Channel → Center Freq to set a center frequency.

The menu item changes to the display shown in Figure 2-16.

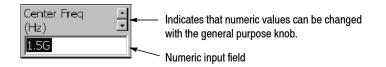

#### Figure 2-16: Changing value with the keypad

2. Press the keys required to enter the desired numeric value. For example, to enter the frequency 123.45 MHz, press 1 2 3 . 4 5 MHz.

ABC DEF Unlock BKSP 7 9 8 GHz Deletes a character before the cursor. GH JKL MNO ms 4 5 6 MHz PQRS WXYZ μs TUV 2 1 3 kHz #%&/ Caps Lock ns/dBm Next Enter +∕\_ 0 Accepts the input. Ηz • Numeric value keys The unit keys also function as the Enter key. Input is immediately accepted when one of those keys is pressed.

To delete an entered number, press the  $BK\ SP\ (Backspace)$  key.

#### Figure 2-17: Numeric keypad

**3.** Confirm the input by pressing the unit key or **Enter** key. The confirmed value is immediately reflected to the analyzer settings and display.

Press the **Cancel - Back** side key to cancel the change.

**Changing the Step Size** When using the up and down  $(\blacktriangle \nabla)$  keys to increase or decrease a setting value, you can change the step size (the amount per press by which the up or down key changes the setting value) with the **Step Size** side key. (The step size cannot be changed for menu items that do not display **Step Size**.)

In the example shown in Figure 2-18, the step size for the center frequency is set to 100 kHz; the displayed frequency set value changes by 100 kHz step for each press of the up or down key.



Figure 2-18: Step size for the center frequency

**Step Size for Center Frequency.** The step size is set with the **Step Size** side key. The center frequency step size can be also set with two side keys in the Frequency/Channel menu (see Figure 2–19 below).

- Center Freq Step Same As C.F. Useful for quickly locating harmonics of a signal seen at the center frequency.
- Center Freq Step Same As Span. Useful for quickly analyzing a larger frequency area without overlapping span windows.


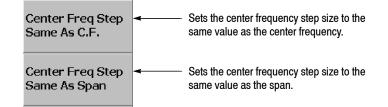



Figure 2-19: Changing the step size to the predetermined values

# **Measurement Basics**

#### Selecting the Measurement Mode

When you start a measurement, first select the measurement mode. For details on each mode, refer to the section shown in Table 2–5.

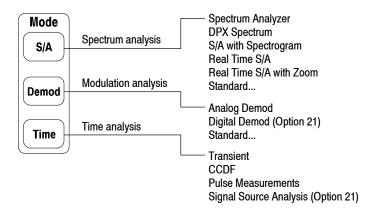

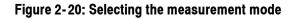


Menu key	Measurement mode	Description	Reference	
Spectrum and	alysis			
S/A	Spectrum Analyzer	General spectrum analysis	Spectrum Analysis	
	DPX Spectrum	DPX (Digital Phosphor) Spectrum analysis	<i>(S/A Mode)</i> – page 3-1 –	
	S/A with Spectrogram	Spectrum analysis with spectrogram		
	Real Time S/A	Real-time spectrum analysis		
	Real Time S/A with Zoom	Real-time spectrum analysis with zoom function		
	Standard <sup>1</sup> Communication standard specific analysis		1	
Modulation a	nalysis		·	
Demod	Analog Demod	Analog modulation analysis	Modulation Analysis	
	Digital Demod	Digital modulation analysis (Option 21)	(Demod Mode) page 3-27	
	Standard <sup>1</sup>	Communication standard specific analysis	- page 0-21	
Time analysis	3			
Time	Transient	Time characteristics analysis	Time Analysis	
	CCDF	CCDF analysis	<i>(Time Mode)</i> page 3-83	
	Pulse Measurements	Pulse characteristics analysis	- puge 0-00	
	Signal Source Analysis	Signal source analysis (Option 21)	]	

#### Table 2-5: Measurement modes

<sup>1</sup> Contains optional analysis functions. Refer to Appendix F: Options and Accessories for optional software.

#### Starting/Stopping Data Acquisition

Use the Run/Stop key on the front panel to start or stop data acquisition.

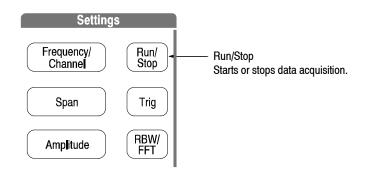




If acquisition and measurements are waiting for a trigger, or are paused/stopped, pressing this key will start the acquisition. If acquisition and measurements are running, pressing this key will stop acquisition and measurements, aborting the current acquisition.

#### **Restoring Default Settings**

The analyzer saves the settings when it is turned off. When you turn on the analyzer, it starts with the settings as they were when it was shut down.

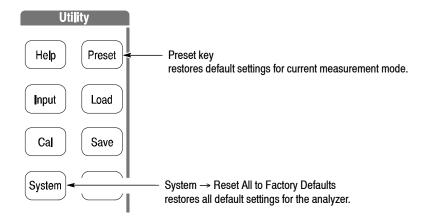


#### Figure 2-22: Restoring default settings

- Pressing the Preset key returns the instrument settings to the factory defaults for the current measurement mode.
- Pressing System → Reset All to Factory Defaults sets all parameters for all measurements and modes to the default values.

# **Setting System Parameters**

The System menu sets system-wide parameters such as the display brightness, factory defaults, and GPIB address.

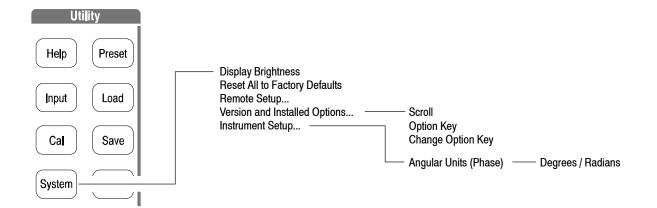

#### Figure 2-23: System menu

**System Menu** The System menu contains the following controls:

Display Brightness. Adjusts the brightness of display. Setting range: 0 to 100.

**Reset All to Factory Defaults.** Sets all parameters for all measurements and modes to the default values.

**Remote Setup...** Sets GPIB parameters. Refer to the *RSA3303B Programmer Manual* for details.

**Versions and Installed Options.** Displays the current versions of all standard and optional software, and any third party software licenses. Refer to page 2-24 for information about displaying versions and installed options.


Instrument Setup... Sets the fundamental parameters for measurements.

• Angular Units. Selects degree (default) or radian for the angular unit.

#### Displaying Versions and Installed Options

You can see the system versions and installed options. For information on options for the analyzer, refer to *Appendix F: Options and Accessories*.

- 1. Press the System key on the front panel.
- **2.** Press the **Versions and Installed Options...** side key. The screen appears as shown in Figure 2–24.

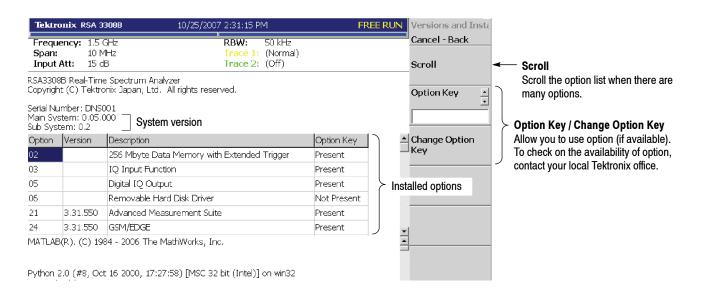

#### Figure 2-24: Displaying versions and options

The following information is shown on the screen:

- Version
  - Main System: Basic application software version
  - Sub System: Firmware version
- Option

The table shows version and name of the option installed. The "Option Key" field tells you the following information:

- Present: Indicates that you can use the option.
- Not Present: Indicates that you cannot use the option. For information about how to remove the protection and use the option, refer to *Enabling the Option*.
- 3. To exit the version display, press any front-panel key.

#### Enabling the Option.

- **1.** Contact your local Tektronix office to order the option. You will receive an option key (alphanumeric code).
- 2. Press the **Option Key** side key (see Figure 2-24) and enter the option key using the alphanumeric keypad on the front panel.

**NOTE**. To input "-" (hyphen), press the "." (period) key.

3. Press the Change Option Key side key to accept the input.

Now you can use the option on the analyzer.

# **Using Windows XP**

This analyzer operates under Windows XP. You can switch to a Windows XP
desktop screen or execute a Windows XP application program.

Displaying the<br/>Windows DesktopTo display the Windows desktop screen, follow these steps:Undesktop1. Connect a USB mouse and keyboard to the analyzer.

- For the USB connectors on the side panel, see Figure 2-3 on page 2-5.
- **2.** Use a mouse to move the pointer to the bottom of screen (see Figure 2-25). The task bar appears.
- **3.** Locate the pointer on the **RSA3303B** or **RSA3308B** icon in the task bar and right-click. A menu appears.
- **4.** Select **Close** from the menu. The analyzer system program terminates and the Windows XP desktop screen appears.

Switching from Windows Desktop to the Analyzer View. To switch from the Windows desktop screen to the analyzer view display, from the task bar, select Start  $\rightarrow$  Program  $\rightarrow$  TEK RTSA. The analyzer program starts.

- Starting Windows XP<br/>ApplicationsWhen you use a mouse to move the pointer to the bottom of screen, a task bar<br/>appears (see Figure 2-25). The task bar contains Start and the analyzer applica-<br/>tion icons. Follow the Windows XP operating procedure and access Windows<br/>applications using the Start menu.
- **Setting Date and Time** The screen displays the date and time managed by the Windows XP operating system. You can use the Windows time setting program to set the date and time.

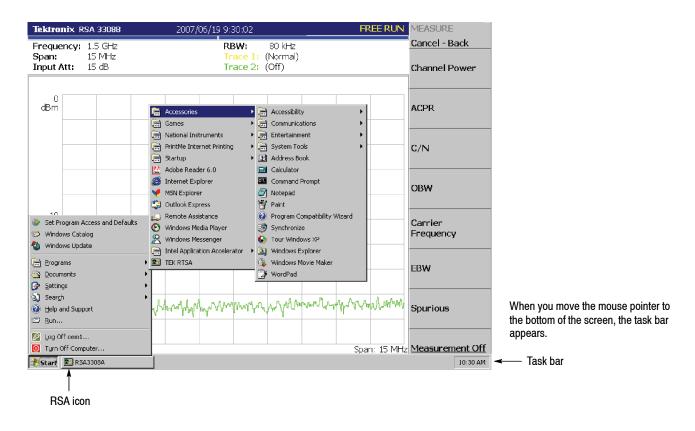


Figure 2-25: Displaying the Windows XP accessory menu

Functional Overview

# Tutorial

This section describes how to perform basic procedures; it contains practical examples of applying power, displaying the results of measurements, and shutting off the analyzer. This section uses default settings in most examples.

- Preparations: Connecting devices and applying the power
- Displaying spectrum
- Using markers and peak search
- Using averaging and comparison displays
- Displaying spectrogram
- Spectrum analysis
- Modulation analysis
- Shutting off the power

You must complete the installation procedures (starting on page 1-11) before performing any of the following steps.


# **Preparations**

This tutorial uses a digitally modulated signal. The following equipment is used as a signal source.

- Analog signal generator (example: HP8657B)
- One 50  $\Omega$  coaxial cable

#### Connecting the Signal Generator

1. Connect the output of the signal generator to the RF INPUT connector on the analyzer front panel using a coaxial cable (see Figure 2-26).

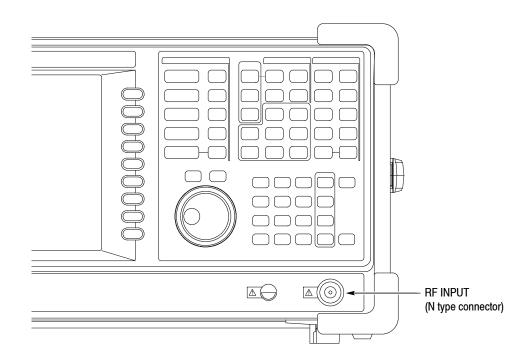

#### Figure 2-26: Cable connection

2. Set the signal generator as follows:

Center frequency	100 MHz
Output level	-10 dBm
Modulation type	AM
Modulation source	Internal 10 kHz
Modulation depth	50%

#### **Applying the Power**

- **1.** Apply power to the signal generator.
- 2. Turn on the principal power switch (PRINCIPAL POWER SWITCH) on the rear panel, as shown in Figure 2-27. The orange LED on the front panel lights.

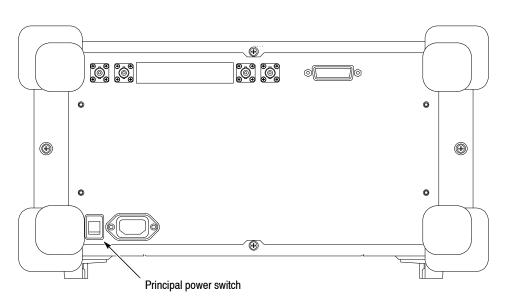


Figure 2-27: Principal power switch (rear panel)

**3.** Turn on the power switch (On/Standby) on the front panel, as shown in Figure 2-28.


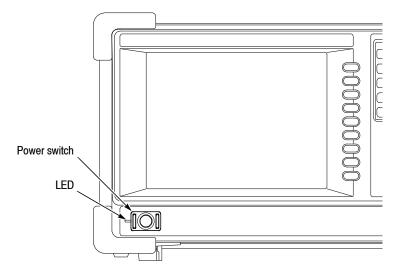


Figure 2-28: Power switch (On/Standby)

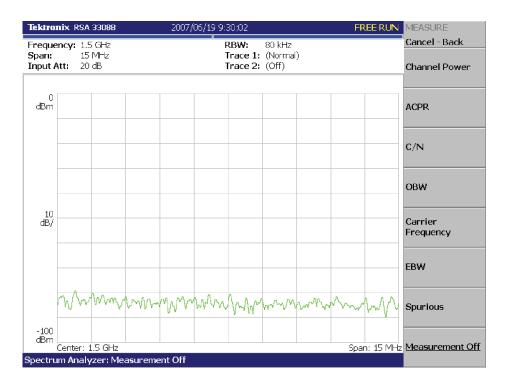

The initial screen shown in Figure 2-29 appears after Windows XP starts up (the background color is in white on this manual to make figures more viewable).


#### Figure 2-29: Initial screen

#### **Restoring Default Settings**

The analyzer saves the settings when it is turned off. When you turn on the analyzer, it starts with the settings as they were when it was shut down.

This tutorial starts with the factory default settings. Perform the following steps to restore the factory default settings:

**1.** Press the **System** key.

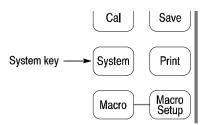

2. Press the Reset All to Factory Defaults side key.

The instrument is now ready to take measurements.

# **Displaying Spectrum**

This section describes how to set the frequency, span, and amplitude first, and then appropriately display the spectrum. Figure 2–30 shows the settings.


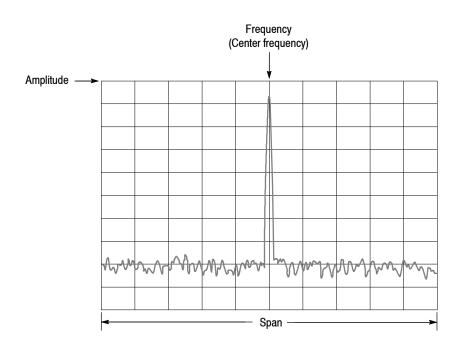


Figure 2-30: Frequency, span, and amplitude settings


#### Setting Center Frequency and Span

The center frequency is set to 1.5 GHz, and span is set to 15 MHz when you power on the analyzer. Change the settings of center frequency and span to display the waveform around 100 MHz.

1. Press the Frequency/Channel key on the front panel.

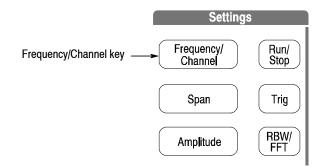

The Frequency/Channel menu, shown in Figure 2-31, is displayed on the right side of the screen. Notice that the Frequency menu item is available for entering the numeric value for center frequency.

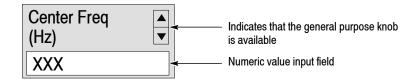

#### Figure 2-31: Menu items for numeric value input

You can change the value using the general purpose knob or enter the value using the numeric value entry keypad, shown in Figure 2-32.

2. Enter the new center frequency at 100 MHz. As the interval between 1.5 GHz (current setting) and 100 MHz is wide, it is more convenient to use the keypad here.

Press 1 0 0 MHz, in order, on the keypad.

The **GHz**, **MHz**, **kHz**, and **Hz** keys also function as the enter key. When you press any of these keys, the numeric value you enter will be set immediately.

If you enter incorrect values, clear them using the **BKSP** (back space) key and then enter the correct digit.

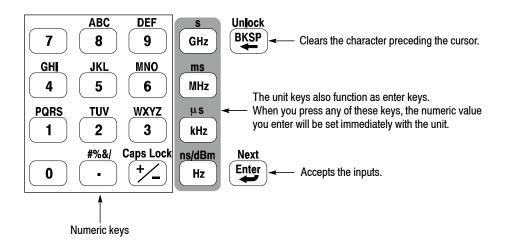


Figure 2-32: Numeric value entry keypad

The spectrum waveform is displayed on the screen as shown in Figure 2-33. The current setting is displayed on the bottom of the screen ("Center Freq (MHz): 100" in this example).


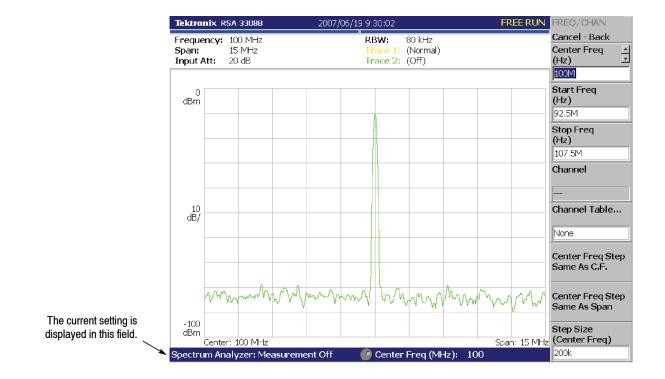



Figure 2-33: Center frequency of 100 MHz, span of 15 MHz

The next step is to set the span. It is currently set to 15 MHz as the default value. Change the setting to 20 kHz.

3. Press the Span key on the front panel.

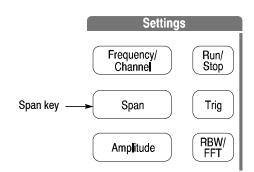

The Span menu item is selected.

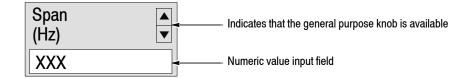

**4.** Turn the general purpose knob to the left to select **20k**. The hardware is set with the selected value immediately.

The spectrum waveform is displayed on the screen as shown in Figure 2-34.


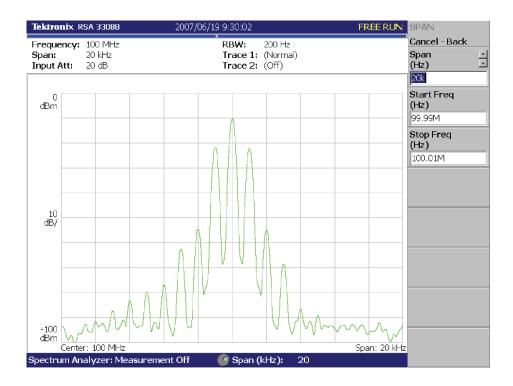



Figure 2-34: Center frequency of 100 MHz, span of 20 kHz

# **Setting Up Amplitude** In Figure 2-35, the vertical scale shows the spectrum view, set up for 10 dB per division. Ref level is the maximum value of the vertical axis, and it is set to 0 dBm when the analyzer is powered on. Observe the change in the waveform as you change this setting in the following procedure:




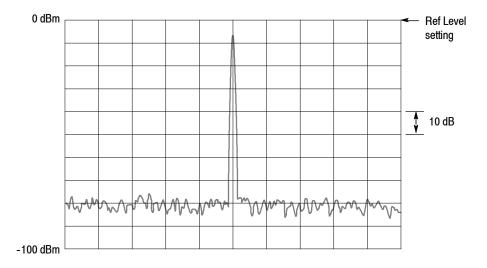

1. Press the Amplitude key on the front panel.

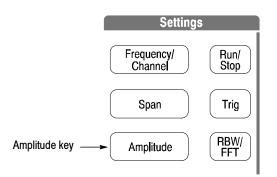

The Amplitude menu is displayed on the right side of the screen. Notice that the Ref Level item is available for entering amplitude.

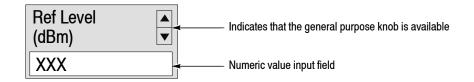

- 2. Observe the changes in the waveform as you turn the general purpose knob.
  - As you turn the knob to the right, the amplitude increases and waveform shifts to a relatively lower position.
  - As you turn the knob to the left, the amplitude decreases and waveform shifts to a relatively higher position.

In Figure 2–36, the amplitude is set to 10 dBm. When the amplitude is set to 0 dBm or higher, a blue standard line is displayed indicating 0 dBm.

3. Set the Ref Level back to 0 dBm after confirming the amplitude operation.


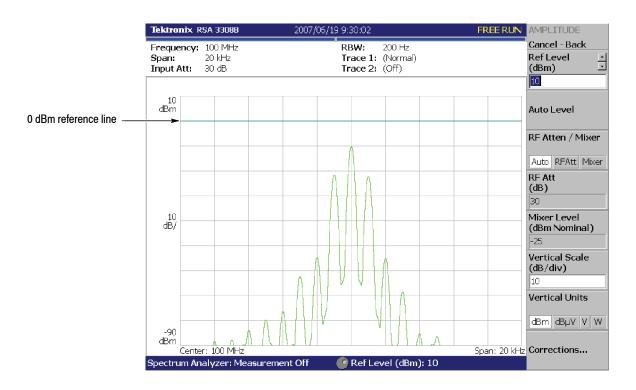


Figure 2-36: Reference level of 10 dBm

#### Starting and Stopping Data Acquisition

Use the **Run/Stop** key to start or stop data acquisition. There are two acquisition modes: the continuous mode, in which data is acquired repeatedly and the single mode, in which one waveform is acquired. You can select the mode using the Trig menu.

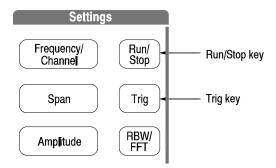

By default, the analyzer is now acquiring data in the continuous mode.

1. Press the **Run/Stop** key to stop data acquisition.

When acquisition stops, "PAUSE" is shown in the status indicator on screen (see Figure 2-37).

				Status indica	ator	
Tektronix R	SA 3308B	2007/06/19 9:30:0	2	PAUSE	MPLITUDE	
	100 MHz 20 kHz 20 dB		200 Hz 1: (Normal) 2: (Off)		Cancel - Back Ref Level (dBm)	× •

#### Figure 2-37: Status indicator

- 2. Acquire data in the single mode using the following steps:
  - **a.** Press the **Trig** key on the front panel.
  - b. Press the Repeat... side key and select Single.
  - **c.** Press the **Run/Stop** key to acquire data. One waveform is acquired and displayed every time you press the key.
- **3.** Press the **Repeat...** side key again and select Continuous to return to the continuous mode.

When acquisition starts, "READY" is displayed in the status indicator, indicating that the analyzer is ready for trigger.


# **Using Markers and Peak Search**

The markers are used to measure amplitude or frequency (and also to find the peak signal).

One or two markers can be displayed as Marker 1 and 2. To measure an absolute value, only Marker 1 is used, and that is called "single marker mode." To measure a relative value, Markers 1 and 2 are used, and that is called "delta marker mode." To position the markers, use the general purpose knob or the numeric input keypad.

#### 1. Press the Marker Setup key on the front panel.


Measuring with a Single Marker

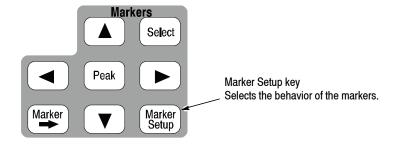

2. Press the Markers side key to select Single.

The marker  $(\Box)$  appears at center on the waveform.

**3.** The **Marker X Position** menu item is selected by default. Turn the general purpose knob to move the marker to a measurement point (see Figure 2-38).

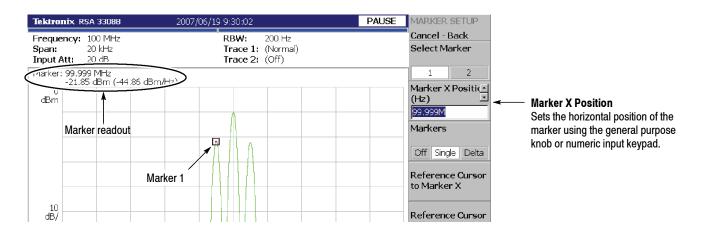




#### Measuring Difference with Delta Marker

Turn on Marker 1 and 2 to measure the difference in amplitude and frequency. On the screen, the symbol " $\Box$ " represents the active marker and " $\diamondsuit$ " represents the fixed marker. You can operate only the active marker.

- 1. Press the Marker Setup key on the front panel.
- 2. Press the Markers side key to select Delta.

The fixed marker ( $\diamondsuit$ ) appears at the active marker position.

- **3.** Make sure that **1** (Marker 1) is selected in the **Select Marker** menu item. This means that Marker 1 is the active marker.
- 4. The Marker X Position menu item is selected. Using the general purpose knob or the numeric input keypad, move the marker to a reference point (see Figure 2-39).

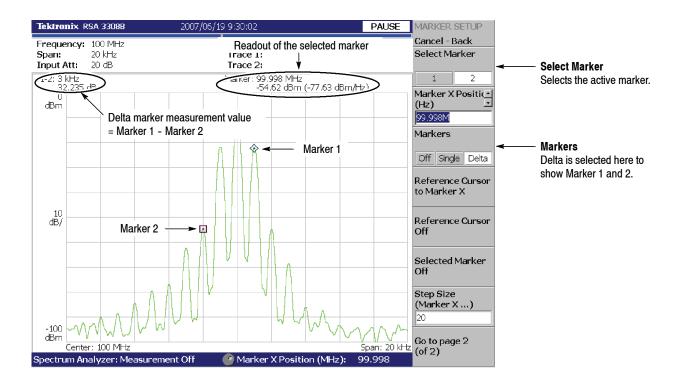

#### Figure 2-39: Measurement with the delta marker

5. Press the Select Marker side key (top) to select 2, making Marker 2 active.

You can also select the marker using the Markers: **Select** key on the front panel. The Marker: **Select** key and the **Select Marker** side key have the same function.

6. The Marker X Position menu item is selected. Using the general purpose knob or the numeric input keypad, move the marker to a measurement point (see Figure 2-39).

The difference between the two marker positions is indicated in the upper left part of the screen.


7. Press the Markers side key to select Single.

The analyzer returns to the single marker mode.

**Searching for the Peak** Measure the frequency interval between the maximum intensity spectrum and the peak to its left by using the delta marker and peak search functions together.

1. Press the **Peak** key on the front panel.

Marker 1 moves to the maximum intensity spectrum.

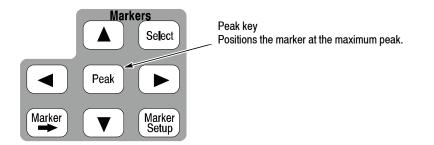

2. Press the Markers side key to select Delta.

The fixed marker ( $\diamondsuit$ ) appears at the active marker position.

- 3. Press the Select Marker side key to select 2, making Marker 2 active.
  - a. Press the marker right key (▶) to move the marker to the next signal peak to the right. Try it several times.
  - **b.** Press the marker left key (◀) to move the marker to the next signal peak to the left. Try it several times.
  - **c.** Using the marker left or right key, put the marker to the measurement peak.

The difference between the two marker positions is shown in the upper left part of the screen (see Figure 2-40).

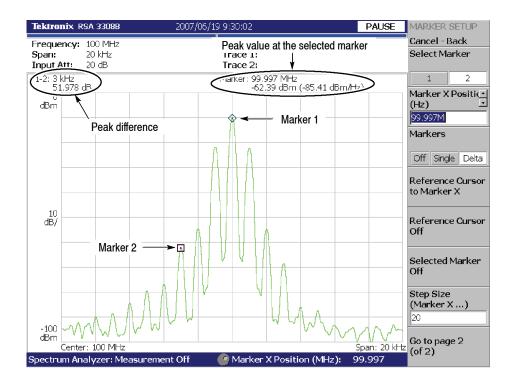



Figure 2-40: Searching for the peak

4. Press the Markers side key to select Off. Both markers disappear.

# **Using Averaging and Comparison Displays**

This section describes how to use the averaging function, which displays the waveform with reduced noise. The averaged waveform can be displayed with the original waveform.

- Averaging There are several averaging methods; select RMS (Root-Mean-Square) in this example.
  - 1. Press the Trace/Avg key.

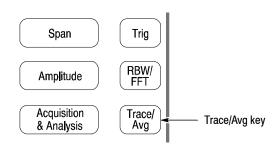

- 2. Press the Trace Type side key to select Average.
- **3.** Press the **Number of Averages** side key to specify how many traces are accumulated to create the averaged waveform.

In this example, enter 64 with the numeric keypad. Press **6 4 Enter** in order.

4. Press the Run/Stop key on the front panel to acquire the waveform.

The averaged waveform is displayed on the screen. The average count is shown in the upper right part of the screen (see Figure 2-41). When the waveform is acquired in the free run mode, the averaging is done with exponential RMS (root-mean-square). This method continues the average with an exponential weighting applied to old values, using the number of averages (64 in this example) as the weighting factor.

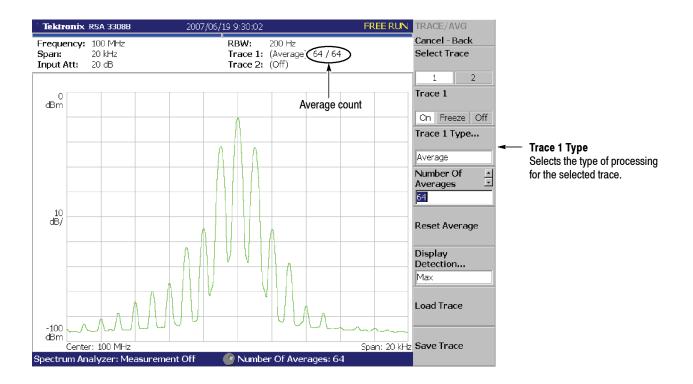


Figure 2-41: Averaging a waveform

5. Press the **Reset Average** side key to restart averaging.

**Comparison Display** You can display two different types of traces concurrently on the screen. In this procedure, you will display the currently acquired waveform and an averaged waveform for comparison.

- 1. Press the Trace/Avg key on the front panel.
- 2. Make sure that 1 (Trace 1) is selected in the Select Trace menu item.
- **3.** Press the **Trace Type...** side key and select **Normal** to define Trace 1 as the currently acquired waveform.
- 4. Press the Select Trace side key to select 2 (Trace 2).
- **5.** Press the **Trace 1 Type...** side key and select **Average** to define Trace 2 as the averaged waveform.
- 6. Press the Run/Stop key on the front panel to acquire waveform.

The currently acquired waveform (Trace 1 in yellow) is displayed with the averaged waveform (Trace 2 in green). See the display in Figure 2-42.

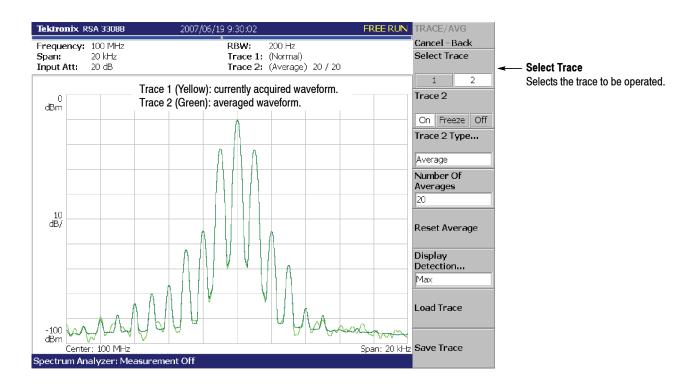



Figure 2-42: Comparison display with averaged waveform

7. Press the Trace 2 Type... side key again and select Off to remove Trace 2.


# **Displaying DPX Spectrum**


The DPX spectrum measurement records data occurrences for each pixel on screen to show the bitmap trace with color-grading according to the density.

1. Press the DPX key on the front panel to display the DPX spectrum.

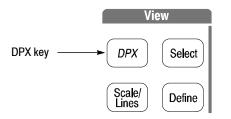

Two traces are displayed by default with Trace 1 as Bitmap and Trace 2 as +Peak. The Bitmap trace indicates high density in red and low in blue. (See Figure 2-43.)

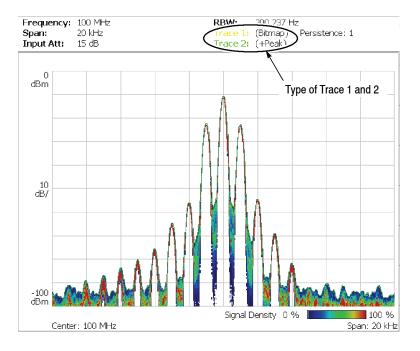

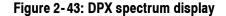

#### Changing the Trace Type

Change the Trace 1 type.

- 1. Press the Trace/Avg key on the front panel.
- 2. Press the Select Trace side key to select 1 (Trace 1).

**3.** Press the **Trace 1 Type...** side key to select the trace type. For example, select **MaxHold** and observe the waveform. The MaxHold trace holds the maximum amplitude for each frequency.

**Persistence Display** When persistence is enabled, signal events fade gradually from the display. You can limit how long a point will be displayed by adjusting the Variable Persistence control, selecting Infinite Persistence, or turning Persistence off. The persistence controls are valid for the Bitmap trace.

- 1. Press the Trace/Avg key on the front panel.
- 2. Press the Trace 1 Type... side key to select Bitmap.
- **3.** Press the **Dot Persistence** side key to select **Variable** (variable persistence display).
- **4.** Press the **Persistence Value** side key and set the value ranging from 1 to 1000 (no unit). For example, set 20.
- 5. Slowly change the amplitude or frequency of the signal generator to vary the output waveform and observe the transient on the DPX display.

Change the persistence value several times and observe how the waveform fades out. Use the **Reset Bitmap** side key to restart the display.

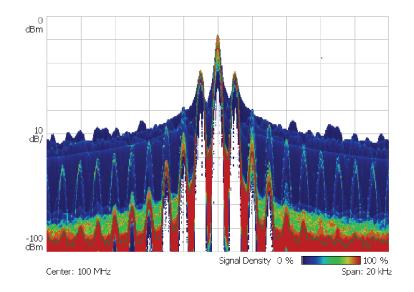


Figure 2-44: Persistence display


6. Reset the signal generator to the original setting for the next step. (Amplitude: -10 dBm, Frequency 100 MHz)

## **Displaying a Spectrogram**

The spectrogram view is a useful tool to observe spectrum changes in time series in a three-dimensional view. The horizontal and vertical axes indicate frequency and frame number, respectively, and the color axis represents amplitude.

Display the spectrogram using the following steps:

1. Press the Mode: S/A key.

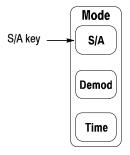

- 2. Press the S/A with Spectrogram side key.
- 3. If no waveform is displayed, press the **Run/Stop** key to acquire data.

The spectrum and spectrogram are displayed concurrently (see Figure 2-45).


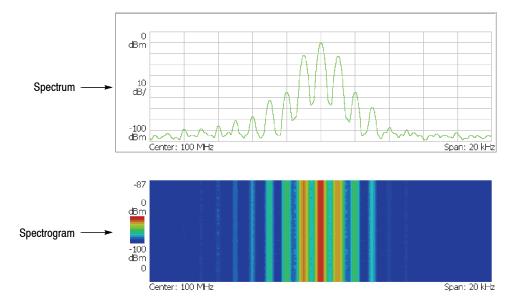


Figure 2-45: Concurrent display of spectrum and spectrogram

- 4. Display the spectrum and spectrogram side-by-side:
  - **a.** Press the View: **Define** key.

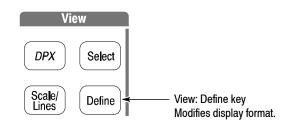

b. Press the View Orientation side key to select Tall.


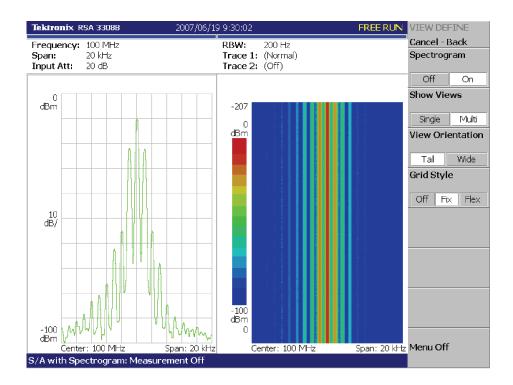


Figure 2-46: Tall display of spectrum and spectrogram


c. Press the View Orientation side key to select Wide.

- 5. Display spectrogram only:
  - **a.** Press the View: **Select** key and then select the spectrogram on the display.

The selected view is framed in white on screen.

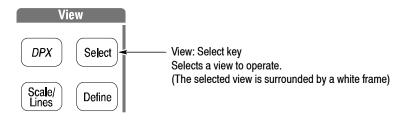

**b.** Press the View: **Define** key and then press the **Show Views** side key to select **Single**.

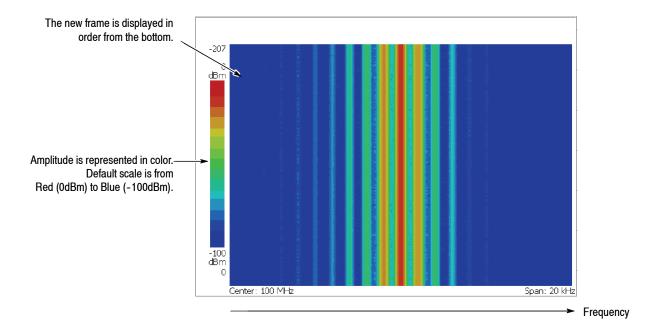

Only the spectrogram is displayed (see Figure 2-47).

Figure 2-47: Spectrogram view


c. Press the Show Views side key again to return to Multi.

# **Spectrum Analysis**

There are several measurement items in the spectrum analysis, such as ACPR (Adjacent Channel Leak Power Ratio), C/N (Carrier vs. Noise Power Ratio), and OBW (Occupied Band Width). These items enable you to perform measurements using simple key operations. Here measure channel power and carrier frequency as examples.

# Measuring Channel Power

1. Press the Measure key on the front panel.

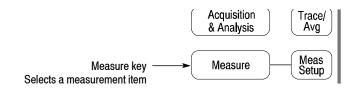

Measurement items are displayed in the menu on the right side of the screen.

2. Press the Channel Power side key.

A band power marker, indicating measurement range, is displayed on the spectrum waveform. The measurement result is displayed below the waveform (see Figure 2-48).

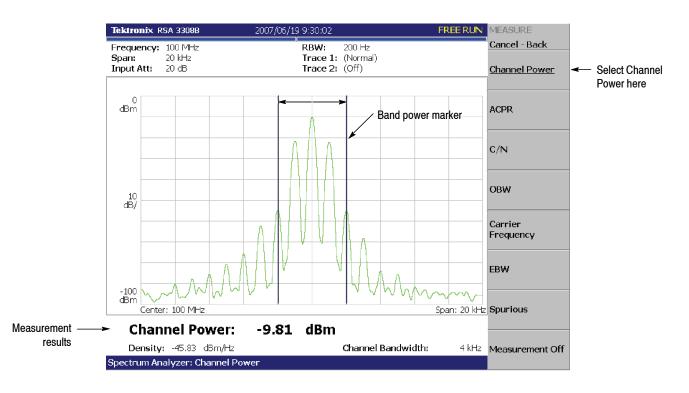


Figure 2-48: Channel power measurement

# Changing a Measurement Parameter

Change a measurement parameter:

- 1. Press the Meas Setup key on the front panel.
- 2. The Channel Bandwidth menu item is selected. Using the general purpose knob, set the measurement range to 40 kHz for example. See the display in Figure 2-49.

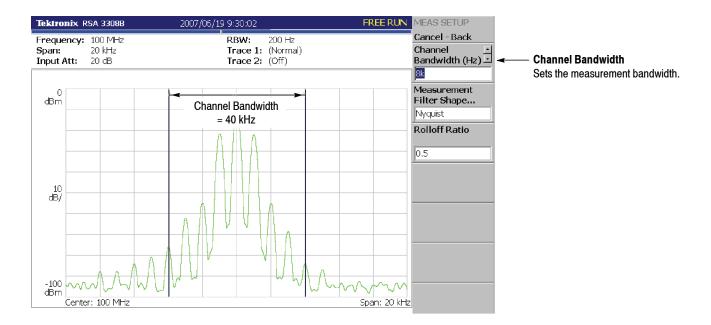



Figure 2-49: Channel power measurement (Channel Bandwidth = 40 kHz)

Measuring Carrier Frequency

Carrier frequency can be measured accurately, using the counter function:

- 1. Press the **Measure** key on the front panel.
- 2. Press the Carrier Frequency side key.

The measurement result is displayed at the bottom of the screen (see Figure 2-50).

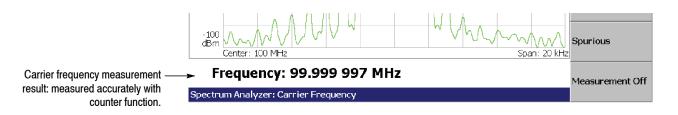




# **Modulation Analysis**

This section describes how to analyze a modulated signal. Settings are the same as those in the previous section; center frequency at 100 MHz, span of 20 kHz, and amplitude of 0 dBm.

**Selecting Analysis Mode** 

The analyzer functions are classified into three functional groups, selected using the Mode keys:

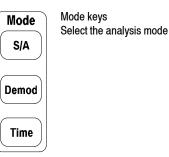

Spectrum analysis Mode: S/A
 Performs general spectrum analysis. All operations in the tutorial up to this point have been performed in this mode.

- Modulation analysis Mode: Demod Performs analog and digital (Option 21 only) modulation analyses.
- Time analysis Mode: Time Performs time characteristic analysis, including CCDF measurement.


The measurement procedures for modulation analysis and the time analysis are similar. Select the analog modulation analysis here:

- 1. Press the **Demod** key on the front panel.
- 2. Press the Analog Demod side key.

## Selecting a Measurement Item

Observe an AM signal change in time series as an example.

- 1. Press the AM Demod side key.
- 2. Press the Run/Stop key to acquire the signal (see Figure 2-51).

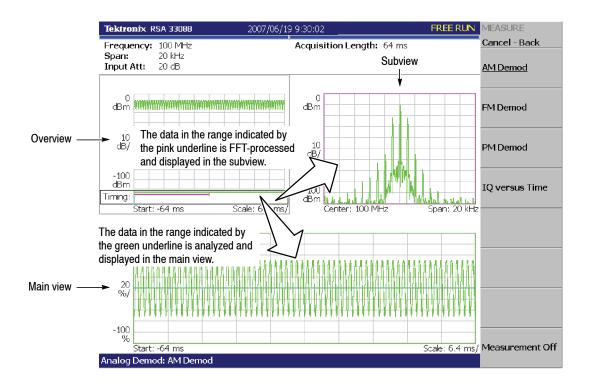

#### Figure 2-51: AM signal measurement


Three views are displayed on the screen in the Demod mode:

- Overview: All data of the specified block are displayed in the time domain. The measurement range indicated by the green underline is specified in this view.
- Main view: Displays the measurement result and waveform of the range specified in the overview (they may be displayed in separate views). In this case, level changes of I and Q signals are displayed as the yellow and green traces, respectively, in the main view.
- *Subview:* Spectrum is displayed as an auxiliary view by default. You can specify the range in the overview to create the spectrum for the subview.

# **Setting Analysis Range**

Set an analysis range, which will be displayed with a green underline in the overview.

1. Press the Acquisition/Analysis key on the front panel.

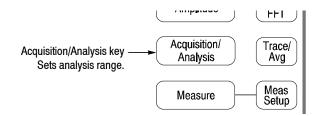

Two vertical green lines appear, indicating the analysis range, shown in Figure 2-52.

- **2.** Press the **Acquisition History** side key to specify the number of the block to be analyzed. Leave it at "0" (the latest block), its default value.
- **3.** Press the **Analysis Length** side key and specify the time length of the analysis range. For example, enter **32 ms** using the numeric keypad.
- 4. Press the **Analysis Offset** side key and specify the starting point of the range. For example, enter **24 ms** using the numeric keypad.


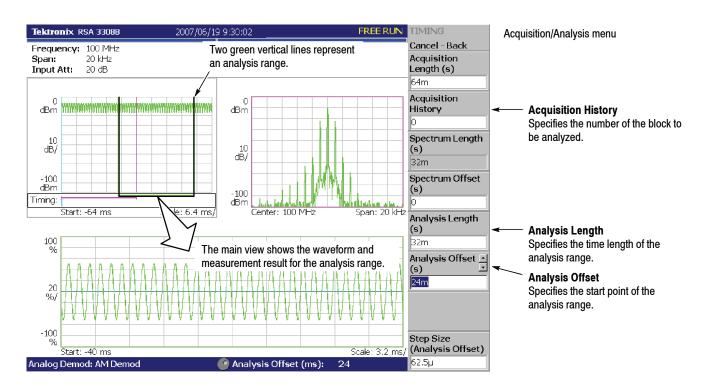


Figure 2-52: Analysis range settings

## Setting the Acquisition Length

One frame consists of 1024 data points and one block consists of several frames. Data are acquired in block units in the Demod mode. The number of frames in one block is referred to as block size. Frame and block size are shown in Figure 2-53.

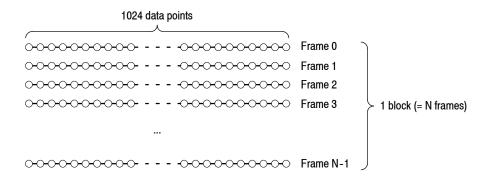

#### Figure 2-53: Frame and block

Suppose that one block contains N frames. The acquisition length for a block is calculated using this equation:

(One block acquisition length) =  $N \times$  (One frame acquisition length)

One block acquisition length is set using **Acquisition Length** in the Acquisition/ Analysis menu. One frame acquisition length is set internally, depending on span, and displayed in **Spectrum Length** in the Acquisition/Analysis menu.

To set the acquisition length, perform the following steps:

1. Press the Acquisition/Analysis key on the front panel.

By default, Acquisition Length is set to 64 ms and Spectrum Length 32 ms. The number of frames in a block is  $64 \div 32 = 2$ .

**2.** Press the **Acquisition Length** side key to change the value. For example, set 256 ms by turning the knob (see Figure 2-54).

In this case, Acquisition Length is 256 ms and Spectrum Length 32 ms, so the data is displayed in the overview for  $256 \div 32 = 8$  frames (8192 points). The green underline indicating the analysis range gets narrower than the previous setting: one block = 2 frames.

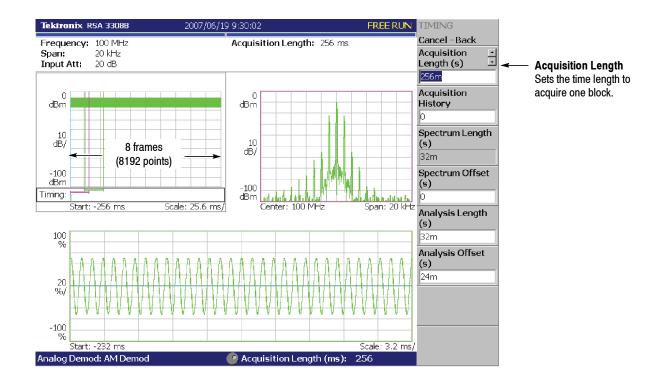


Figure 2-54: Changing the acquisition length

# Displaying Single View and Changing Scale


Although three views are normally displayed in the Demod mode, you can select one view to display full screen. This section describes how to display one view and change the vertical and horizontal scales using the View keys.

- 1. Press the Run/Stop key to stop data acquision and observe the waveform.
- **2.** Press the View: **Select** key to select the main view. The selected view is surrounded with a white frame.

Vi	ew	View k Contro
DPX	Select	
Scale/ Lines	Define	

View keys Control the scale and format of views **3.** Press the **Define** key, and then the **Show Views** side key, shown in Figure 2-55, to select **Single**.

Only the main view is displayed, with the size enlarged on the screen.

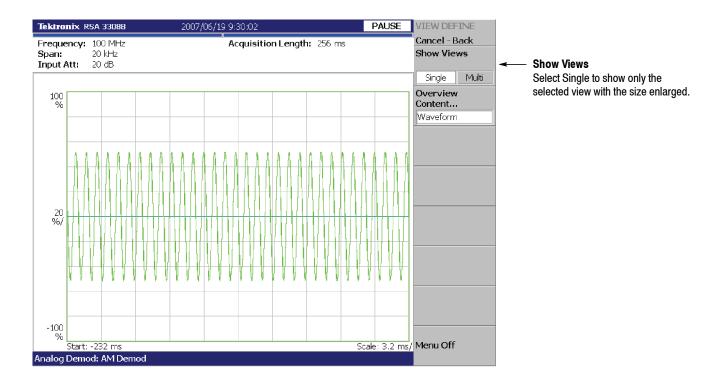

#### Figure 2-55: Single view display

- **4.** Change the scale:
  - a. Press the Scale/Lines key in the View area.
  - **b.** Press the **View Scale...** side key.
  - **c.** Press the **Horizontal Scale** side key and change the scale of the horizontal axis (see Figure 2-56).

Try several settings by turning the general purpose knob; observe the changes in the display.

d. Press the Vertical Scale side key and change the vertical axis scale.

Try several settings by turning the general purpose knob; observe the changes in the display.

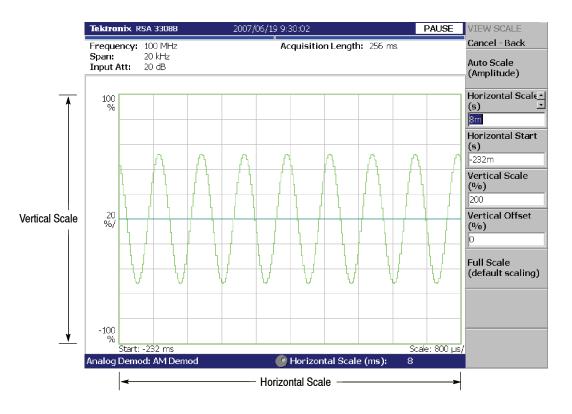


Figure 2-56: Scale settings

# Completing the Measurement

Complete the measurement with the following procedure:

- 1. Press the **Measure** key on the front panel.
- 2. Press the Measurement Off side key.

The display returns to the spectrum view. However, the analysis mode is still in the Demod mode.

# **Turning Off the Power**

When you have completed the measurements, turn off the power:

1. Press the power switch (On/Standby) at the lower left of the front panel.

The shutdown process of Windows XP runs and the power source goes to standby state, with the orange LED on.


**2.** Turn off the signal generator.

You have now completed the tutorial.

# Reference

# **Spectrum Analysis (S/A Mode)**

This section describes how to perform measurements in the spectrum analysis (S/A) mode.

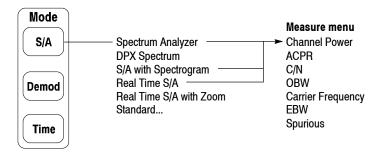

#### Figure 3-1: S/A menu structure

There are four items in the S/A menu:


- Spectrum Analyzer. Performs general spectrum analysis. Refer to page 3-4 for details.
- DPX Spectrum. Performs analysis displaying the DPX (Digital Phosphor) spectrum view. Refer to page 3-14 for details.
- S/A with Spectrogram. Performs spectrum analysis with spectrogram. Refer to page 3-17 for details.
- Real Time S/A. Performs real-time spectrum analysis with spectrogram. Refer to page 3-19 for details.
- **Real Time S/A with Zoom.** Performs real-time spectrum analysis with zoom function. Refer to page 3-23 for details.
- Standard.... Performs spectrum analysis according to a communication standard. For details, refer to the user manual included in each option as shown in Table 3-1:

Table 3-1:	Standard	menu in	the S/A	mode
------------	----------	---------	---------	------

Menu item	Option	Description
W-CDMA-UL	Option 30	W-CDMA uplink analysis
3GPP-R5-DL	Option 30	3GPP Release 5 downlink analysis
3GPP-R5-UL	Option 30	3GPP Release 5 uplink analysis

# **Measurement Screen Layout**

Figure 3-2 shows the basic screen layout in the spectrum analysis (S/A mode). The spectrum waveform and measurement results are displayed. You can display the spectrum and the spectrogram at the same time. Refer to *Spectrogram Display* on page 3-17.

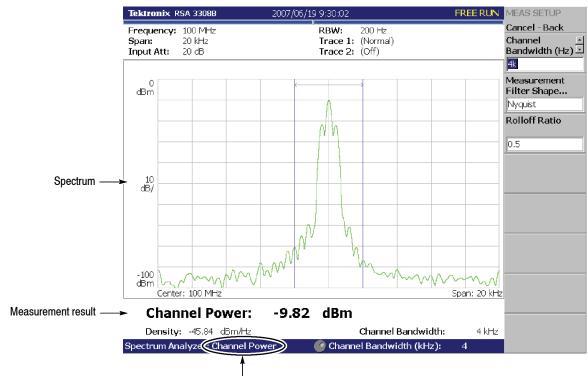

Measurement item

Figure 3-2: S/A measurement screen

# **Changing the Grid Style**


The graticule is a  $10 \times 10$  grid by default. You can select a different grid style by pressing View: **Define**  $\rightarrow$  **Grid Style** as follows:

**NOTE**. The Grid Style menu item is available only in the S/A mode except Real Time S/A.

**Off.** Turns off the graticule.

**Fix.** Always displays a  $10 \times 10$  grid.

**Flex.** Displays the graticule with the horizontal scale (per division) in a 1-2-5 sequence.

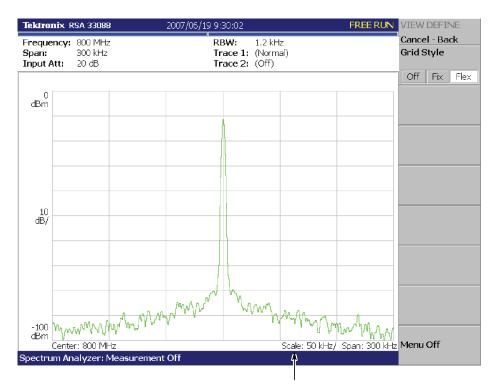

The horizontal scale is set to 50 kHz/div for a span of 300 kHz in the "Flex" grid style.

Figure 3-3: "Flex" grid style

# **Spectrum Analysis**

You can select following measurement items with the **Measure** key when you select **S**/**A** and then the **Spectrum Analyzer** side key.

Refer to:
age 3-5
age 3-6
age 3-7
age 3-8
age 3-10
age 3-10
age 3-12
08

Table 3-2: Measurement items in the S/A mode

Pressing the **Measurement Off** side key stops the measurement and returns to the original spectral display.

- **Basic Procedure** 1. Press the S/A key and then the Spectrum Analyzer side key.
  - 2. Select the measurement item (refer to Table 3-2).
  - 3. Display the spectrum waveform of the measurement signal:
    - a. Press the Run/Stop key on the front panel to start data acquisition.

**NOTE**. For details on setting frequency, span, and amplitude, refer to Setting Frequency and Span on page 3-119 and Setting Amplitude on page 3-127.

- **b.** Set the frequency by pressing the **Frequency/Channel** key on the front panel.
- c. Set the span by pressing the **Span** key on the front panel.
- d. Set the amplitude by pressing the Amplitude key on the front panel.
- 4. Set the following Measurement Setup controls by pressing the **Meas Setup** key on the front panel.

**Channel Power Measurement Measurement Measurement This measurement calculates power in the range specified using a band power marker (shown in Figure 3-4).** 

**Meas Setup Menu.** The Meas Setup menu for the channel power measurement contains the following controls:

**Integration Bandwidth.** Sets frequency range for power measurement (see Figure 3-4).

Measurement Filter Shape... Selects a filter shape from these types:

- Rect (Rectangular)
- Gaussian
- Nyquist
- Root Nyquist

**Rolloff Ratio.** Enters rolloff ratio when Nyquist or Root Nyquist filter is selected. Range: 0.0001 to 1 (default value: 0.5).

An example of the channel power measurement is shown in Figure 3-4.

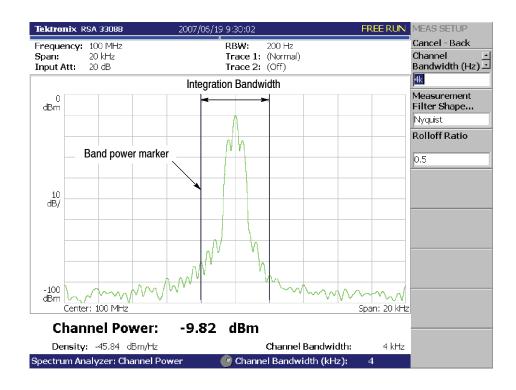


Figure 3-4: Channel power measurement

**ACPR Measurement** The power ratio of the carrier signal to the signal appearing in its adjacent frequency band (leak signal) is measured in the ACPR (Adjacent Channel Power Ratio) measurement. Frequency ranges are set using three band power markers.

**Meas Setup Menu.** The Meas Setup menu for the ACPR measurement contains the following controls (see Figure 3-5):

Main Channel Bandwidth. Sets frequency range of main channel.


Adjacent Channel Bandwidth. Sets frequency range of adjacent channel.

Chan Spacing. Sets a frequency interval between two adjacent channels.

Measurement Filter Shape... Selects a filter shape from these types:

- Rect (Rectangular)
- Gaussian
- Nyquist
- Root Nyquist

**Rolloff Ratio.** Enters rolloff ratio when Nyquist or Root Nyquist filter is selected. Range: 0.0001 to 1 (default value: 0.5).

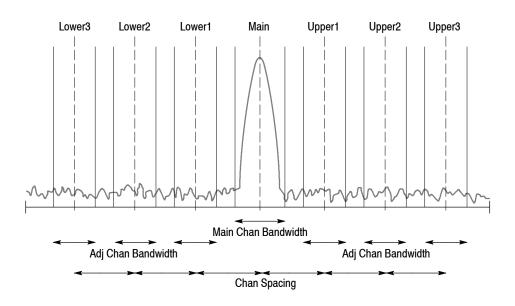

#### Figure 3-5: ACPR measurement band power markers

Figure 3-6 shows an example of the ACPR measurement.

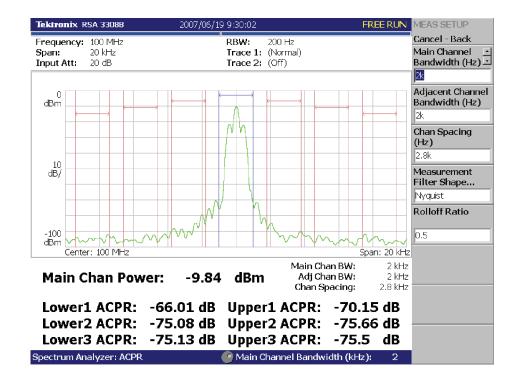


Figure 3-6: Example of ACPR measurement

**C/N Measurement** Measures Carrier to Noise ratio (C/N).

**Meas Setup Menu.** The Meas Setup menu for the C/N measurement contains the following controls (see Figure 3-7).

**Offset Frequency.** Sets the offset frequency from carrier to noise. Range: -Span/2 to +Span/2.


Noise Bandwidth. Sets the noise bandwidth.

Carrier Bandwidth. Sets the carrier bandwidth.

Measurement Filter Shape... Selects a filter shape from these types:


- Rect (Rectangular)
- Gaussian
- Nyquist
- Root Nyquist


**Rolloff Ratio.** Enters rolloff ratio when Nyquist or Root Nyquist filter is selected. Range: 0.0001 to 1 (default value: 0.5).

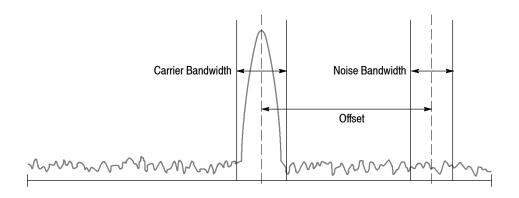

## Figure 3-7: C/N measurement band power markers

Figure 3-8 shows an example of the C/N measurement.

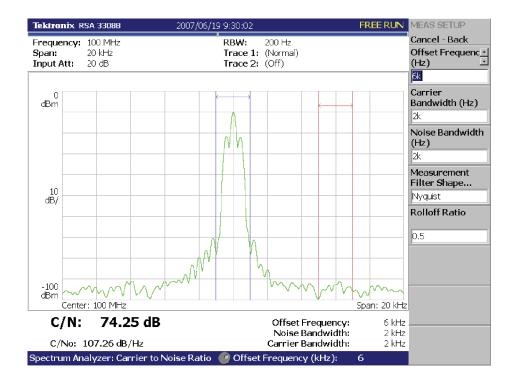


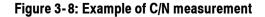

**OBW Measurement** The OBW (Occupied Bandwidth) measures the frequency bandwidth, using the ratio you specify for carrier signal power / power within the span setting.

**Meas Setup Menu.** The Meas Setup menu for the OBW measurement contains the following control.

**Power Ratio.** Specifies the power ratio of the carrier and span regions for calculating OBW (see Figure 3-9). The default setting is 99% as defined in T-53 or IS-95 standard. Range: 80 to 99.99%.

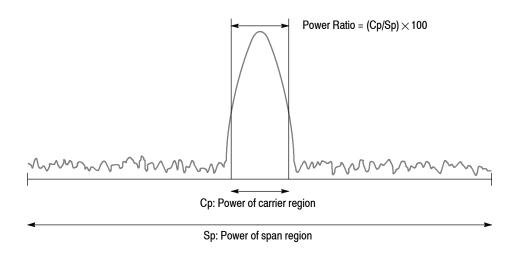

#### Figure 3-9: OBW measurement band power marker

Figure 3-10 shows an example of the OBW measurement.

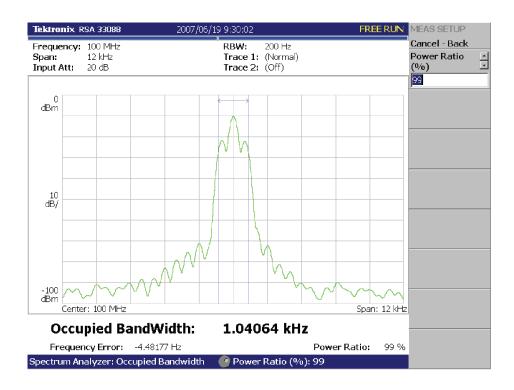


Figure 3-10: OBW measurement

# Carrier Frequency Measurement

Carrier frequency is accurately measured using the counter function.

**NOTE**. You must set the frequency and the span to display only the spectrum with the carrier you want to measure. Measurement error will occur if other frequency components are displayed together.

**Meas Setup Menu.** The Meas Setup menu for the Carrier Frequency measurement contains the following control:

**Counter Resolution.** Sets the counter resolution. The measurement result is displayed on the bottom of the view with this resolution.

Range: 1 mHz to 1 MHz in a tenfold sequence (default: 1 Hz).

Figure 3-11 shows an example of the carrier frequency measurement.

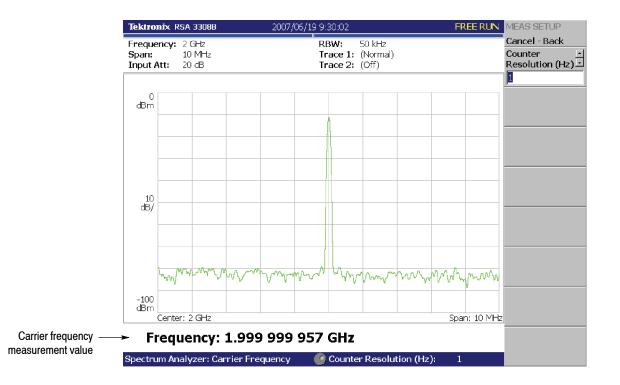



Figure 3-11: Carrier frequency measurement

## EBW Measurement

The frequency bandwidth between the maximum peak of the spectrum and the specified dB value is measured in EBW (Emission Bandwidth).

**Meas Setup Menu.** The Meas Setup menu for the EBW measurement contains the following control:

**Measurement Level.** Specifies how far down from the peak level the bandwidth is measured (see Figure 3-12). Range: -100 to -1 dB (default: -30 dB).

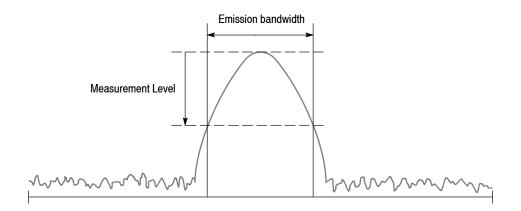

## Figure 3-12: EBW measurement band power markers

Figure 3-13 shows an example of the EBW measurement.

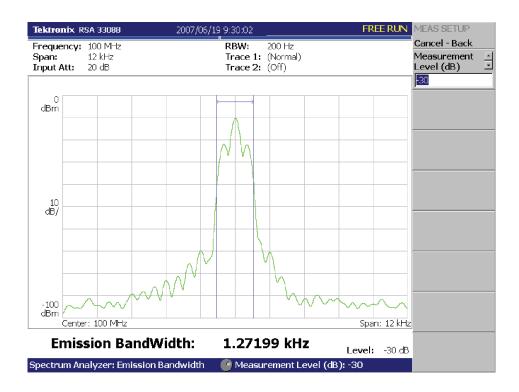


Figure 3-13: EBW measurement

# Spurious Signal Measurement

The spurious signal measurement detects up to 20 spurious signals and calculates the frequency difference and the amplitude ratio relative to the normal signal by setting the level conditions.

**NOTE**. For detecting stationary spurious signals, waveform averaging is effective for reducing noise. Refer to page 3-175 for averaging.

**Meas Setup Menu.** The Meas Setup menu for the spurious signal measurement contains the following controls:


**Signal Threshold.** Sets the threshold for detecting the standard signal (see Figure 3-14). A signal whose amplitude is larger than this threshold is regarded as a standard signal. Range: -100 to +30 dBm.

**Ignore Region.** Sets the frequency range with the carrier (standard signal) peak as the center in which spurious signals are not to be detected, for avoiding false recognition of spurious signals (see Figure 3-14). Range: 0 to Span/2 Hz (default: 0 Hz)

**Spurious Threshold.** Sets the threshold for detecting spurious signals (see Figure 3-14). Enter the value relative to the peak of standard signal. Range: -90 to -30 dB.

**Excursion.** Sets the deviation amount of amplitude that is considered to be spurious (see Figure 3-14). A signal with amplitude higher than the set value of the Spurious Threshold and larger than the set value of Excursion is considered to be spurious. Range: 0 to 30 dB (default: 3 dB)

**Scroll Table.** Horizontally scrolls the spurious table displayed in the lower part of the screen. Up to 20 spurious signals are displayed.

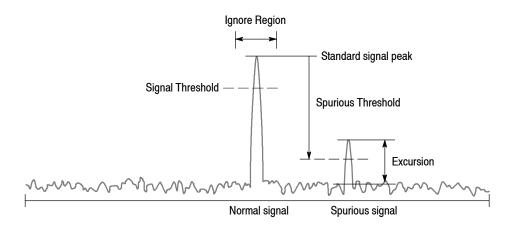



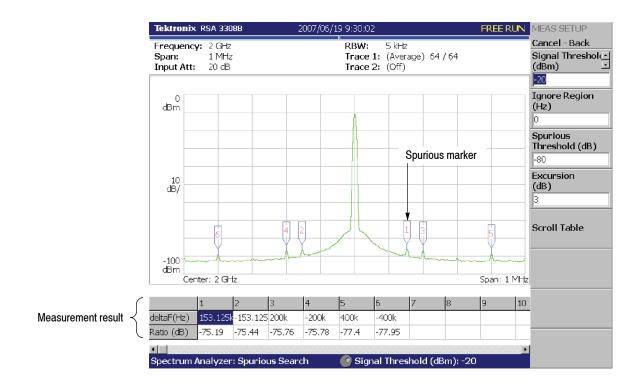


Figure 3-15 shows an example of spurious signal measurement.


Figure 3-15: Example of spurious signal measurement

Detected spurious signals are assigned markers, numbered in descending order by amplitude. The frequency difference (deltaF) and the amplitude ratio (Ratio) in relation to the standard signal are shown in the table in the lower part of the screen.

# **DPX Spectrum Display**

DPX Spectrum technology enables you to see how traces change over time, displaying transient and intermittent events that cannot be seen on a swept spectrum analyzer. A DPX Spectrum indicates how traces change in two ways. First, it uses color grading to show how often a particular signal occurs. Second, it uses persistence to keep events visible in the display, allowing newer events to be compared to older events.

The DPX Spectrum display works by using a two-dimensional array to accumulate data points for the display. As shown in the Figure 3-16, input signal data after A/D conversion go through the proprietary DPX engine to compute more than 50,000 FFTs per second. These results are written to the pixel buffer memory. Each time a trace writes to a point on the display, a counter in the array is incremented. A color is assigned to each point in the display based on the value of its counter. Thus, as acquisitions occur over time, a colored waveform, the Bitmap, develops on the display that shows how frequently a display point has been addressed.

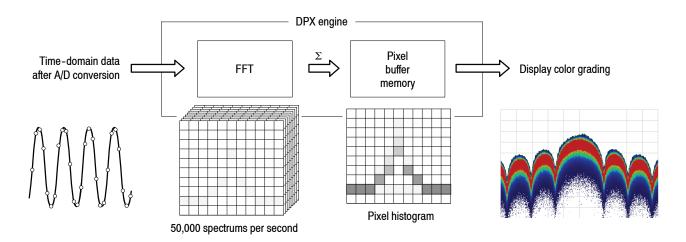

# Figure 3-16: DPX process flow

Displaying the DPX Spectrum Use the following steps to display the DPX spectrum.

- 1. Press the **DPX** key on the front panel. (Or press  $S/A \rightarrow DPX$  Spectrum.)
- 2. Press Trace/Avg  $\rightarrow$  Select Trace to select 1 (Trace 1).
- **3.** Press the **Trace 1 Type...** side key and select the type of Trace 1. Trace 1 defaults to Bitmap.
  - Bitmap. Displays the density of acquired data. The number of data points acquired at each pixel is indicated by color. Pressing the Reset Bitmap side key restarts the process.

- +Peak. Displays the maximum amplitude acquired at each frequency in each screen update.
- MaxHold. Holds the maximum amplitude at each frequency. Pressing the Reset MaxHold side key restarts the process.
- -Peak. Displays the minimum amplitude acquired at each frequency in each screen update.
- MinHold. Holds the minimum amplitude at each frequency. Pressing the Reset MinHold side key restarts the process.
- Average. Computes the average amplitude of the spectrum points within each frequency bin to create a trace for each screen update, and then average the traces with the number set in the Number Of Averages side key to display. Pressing the Reset Average side key restarts the process.

Refer to *Trace Comparison Display and Average Function* on page 3-175 for the details.

- 4. Press Select Trace to select 2 (Trace 2).
- **5.** Press the **Trace 2 Type...** side key and select the type of Trace 2. The selections are the same as Trace 1. Trace 2 defaults to +Peak.
- 6. Set the scale as needed by pressing Scale/Lines  $\rightarrow$  View Scale... Refer to *Setting Views* on page 3-187 for the details.

You can save and load a trace in a file with the extension of .*dpt* for the bitmap or .*trc* for other than bitmap using the **Save Trace 1/2** and **Load Trace 1/2** side keys in the Trace/Avg menu. Refer to *File Operations* on page 3-233 for the details.


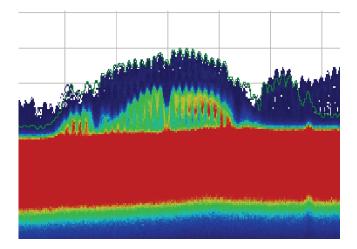


Figure 3-17: DPX spectrum display

**Persistence** You can select the persistence mode for the Bitmap trace. While the screen display is generally updated periodically to rewrite waveforms, the persistence mode enables a data point once displayed to remain visible for a decay period. This allows you to see the history of occurrence for dynamic signals while still receiving immediate feedback on signal variations over time. You can adjust the persistence to optimize display characteristics for varying signal conditions, from a live RF view of dynamic signals to the discovery of single occurrences.

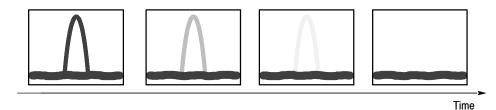

The persistence function facilitates observing single or intermittent events by fading out waveforms gradually. You can set a decay period ranging from zero to infinity.

#### Figure 3-18: Persistence display

Use the following steps to adjust the persistence.

- 1. Press the Trace/Avg key on the front panel.
- **2.** Press the side key **Trace 1** (or **2**) **Type**  $\rightarrow$  **Bitmap**.
- 3. Press the **Dot Persistence** side key to select the persistence mode.
  - Off. Disables the persistence mode. All points are reset every screen update.
  - Variable. Selects the variable persistence mode.
     Set a persistence time, that is how long a point in the display is visible before fading, using the Persistence Value side key.
  - Infinite. Selects the infinite persistence mode that prevents a point in the display from fading. Once a point in the display has been written to, it will remain visible until you press the Reset Bitmap side key.
- **4.** When selecting the variable persistence mode, press the **Persistence Value** side key to set the persistence time. Range: 1 to 1000 (no unit, default: 10)
- 5. Press the **Intensity** side key to control the visibility of events. An increased intensity level allows a single, short event to be seen. This also allows the effect of variable persistence to be seen on infrequent events. Range: 1 to 100% (default: 25%).

# **Spectrogram Display**

You can display the spectrum of an input signal concurrently with the spectrogram.

Follow these steps to display the spectrogram:

- 1. Press the S/A key on the front panel.
- 2. Press the S/A with Spectrogram side key.

The spectrum and spectrogram are displayed on the same screen as shown in Figure 3-19.


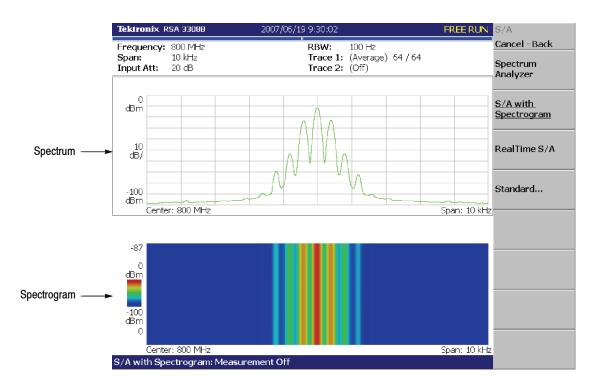


Figure 3-19: Spectrum and spectrogram concurrent display


# **Changing Display Style**

You can change display style as required.

- **1.** Press the View: **Define** key.
- 2. Press the View Orientation side key to select the view style: Wide or Tall.

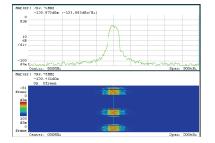


Wide display Displays spectrum and spectrogram in line-split display.

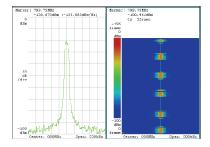

Tall display Display spectrum and spectrogram side by side.

#### Figure 3-20: View orientation

- **3.** When you display either spectrum or spectrogram only on full-screen, press the View: **Select** key to select the view. The selected view is surrounded by a light-blue frame.
- 4. Press the Show Views side key and select Single as shown in Figure 3-21.

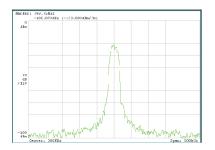

Single display Displays the selected view on one screen.

Figure 3-21: Single display

Refer to Setting Views on page 3-187 for details.

# **Real-Time Analysis**

When you select  $S/A \rightarrow Real Time S/A$ , the analyzer performs the real-time analysis, displaying the spectrogram. The screen layout is the same as that when you select  $S/A \rightarrow S/A$  with Spectrogram. Refer to Spectrogram Display on page 3-17.

**NOTE**. The TRACE/AVG menu is not available in the real-time mode.

*The Real-Time S/A mode has FFT overlap capability. Refer to* FFT and RBW *on page 3-163 for detail.* 

# Features of the Real-Time Mode

An input waveform is acquired in blocks, where a block is defined as a number of frames and a frame is defined as a number of FFT data points. The number of frames acquired at one time is called block size. In the normal spectrum analysis, the instrument acquires data for the block size determined from RBW and creates one spectrum waveform. In the real-time mode, the instrument acquires data for the block size specified with the Acquisition/Analysis menu, performs the FFT process, and creates spectrum waveforms for each frame so that you can observe spectrum variation in time seamlessly. The difference between normal and real-time modes is shown in Figure 3-22.

Frame -1

Frame 0






Frame -N

Frame -2

Time —

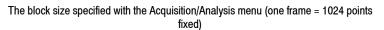

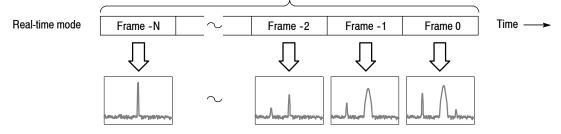




Table 3-3 shows the features of the real-time mode compared with the normal spectrum analysis.

ltem	Normal spectrum analysis	Real-time mode
Span	Up to 3 GHz (arbitrary value can be set)	Up to 36 MHz (1-2-5 sequence)
Trigger	Only the Repeat menu item available	All the Trigger menu items available
RBW/FFT	FFT points: 64 to 8192 (powers of 2) RBW: 1 Hz to 10 MHz	FFT points: 1024 fixed RBW: No RBW process FFT overlap capability
Timing	No timing parameters	Acquisition Length and Spectrum Offset can be set

Table 3-3: Features of the real-time mode

**Basic Procedure** The following steps show the basic procedure for real-time spectrum analysis with spectrogram:

- 1. Press the S/A key on the front panel.
- 2. Press the **Real Time S/A** side key. The spectrum is displayed concurrently with the spectrogram.
- **3.** Display the spectrum waveform of the measurement signal:

**NOTE**. For details on setting frequency, span, and amplitude, refer to Setting Frequency and Span on page 3-119 and Setting Amplitude on page 3-127.

For details on setting trigger, refer to Trigger on page 3-147.

- a. Set the frequency by pressing the Frequency/Channel key on the front panel.
- **b.** Set the span by pressing the **Span** key on the front panel. Table 3-4 shows the span setting range.

#### Table 3-4: Span setting range

Measurement band <sup>1</sup>	Setting range
RF	100 Hz to 10 MHz (1-2-5 sequence) and 15 MHz
Baseband	100 Hz to 20 MHz (1-2-5 sequence)
<sup>1</sup> Baseband: DC to 20 MHz: BE: 15 MHz to 3 GHz (BSA3303B) / 8 GHz (BSA3308B)	

Baseband: DC to 20 MHz; RF: 15 MHz to 3 GHz (RSA3303B) / 8 GHz (RSA3308B)

- c. Set the amplitude by pressing the Amplitude key on the front panel.
- d. Set the trigger by pressing the Trig key on the front panel.
- 4. Press the Acquisition/Analysis key on the front panel and then the Acquisition Length side key to set the time length to acquire one block.

Suppose that one block contains N frames; the acquisition length is calculated with this equation:

(One block acquisition length) =  $N \times$  (One frame acquisition length)

Where N = 1 to 16000 (standard) or 64000 (Option 02)

One frame acquisition length is determined by span internally, and shown with Spectrum Length in the setup display area (see Figure 3-23). For the details, refer to Table A-10 on page A-8 in *Appendix A: Specifications*.


Waveform data is acquired and displayed on the basis of blocks. For details about the time parameters, refer to page 3-141.

- 5. When you have acquired the measurement data, stop the data acquisition. If the analyzer is in the continuous acquisition mode, press **Run/Stop**.
- **6.** Press the **Spectrum Offset** side key and set the number of the frame to measure and display in the spectrum view (upper view) by turning the general purpose knob or using the numeric keypad.

Serial numbers are given to each frame regardless of the block size, with zero representing the latest frame. The selected frame is indicated by the marker on the spectrogram (see Figure 3-23).

Alternatively, the frame number can be set as follows:

- a. Press the View: Select key to select the spectrogram view.
- b. Press the Marker Setup key.
- c. Press the Go to page  $2 \rightarrow$  Marker X Vertical side key.
- **d.** Set the frame number by turning the general purpose knob or using the numeric keypad.

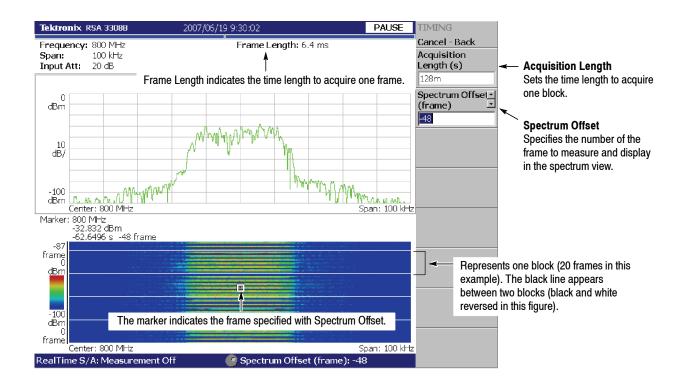

## Figure 3-23: Real-time mode

- 7. To take a spectrum measurement, press the **Measure** key on the front panel. The measurement items and procedures are the same as in the normal spectrum analysis. Refer to *Spectrum Analysis* on page 3-4.
- 8. To change the display format, press the **Define** key on the front panel. The setting procedure is the same as in the normal spectrum analysis. Refer to *Spectrogram Display* on page 3-17.

**NOTE**. In real-time mode, you cannot turn on or off the spectrogram in the View: Define menu.

#### **Zoom Function**

When you select  $S/A \rightarrow Real Time S/A$  with Zoom, the real-time analysis area can be specified on the spectrogram. Even though frequency hopping occurs, measurement can only be made by clipping the optional burst (see Figure 3-24).

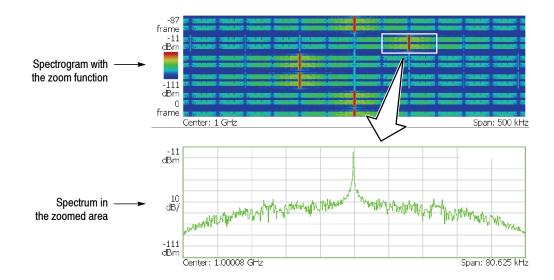


Figure 3-24: Zoom function

#### Acquisition/ Analysis Menu

In the spectrogram with the zoom function display, the zoomed area is set by the following Acquisition/Analysis menu (see Figure 3-25). The area is the analysis range.

Acquisition Length. Sets the capture time of one block.

Acquisition History. A block number that analyzes/displays data is specified. The latest block number is 0. Older blocks have larger negative numbers.


**Analysis Length.** Sets the time length of the zoomed area in the block set in the Acquisition History.

**Analysis Offset.** Sets the starting point of the zoomed area by considering the trigger output point as the reference.

Frequency Center. Sets the frequency of the center of zoomed area.

Frequency Width. Sets the frequency width of the zoomed area.

Refer to Setting Acquisition and Analysis Parameters on page 3-141.

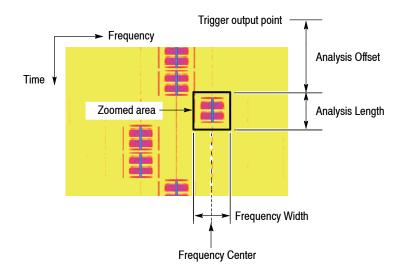

#### Figure 3-25: Zoom domain settings using the Acquisition/Analysis menu

#### Setting a Zoomed Area Using the Marker

The zoomed area can be set by using the marker and reference cursor instead of entering a numeric value in the Acquisition/Analysis menu (see Figure 3-26). The mouse can also be used in the marker operation (refer to *USB connectors* on page 2-6 for connecting the mouse). In the following procedure, [Mouse] indicates an operation where the mouse is used.

**1.** Press the View: **Select** key on the front panel and select the spectrogram with zoom function.

[Mouse] Click on the spectrogram to select it and then go to step 4.

- 2. Press the Marker Setup key on the front panel.
- 3. Press the Markers side key and select Single to display the marker.
- 4. Move the marker to the upper left point of the zoomed area using the Marker X Position and the Marker X Vertical side keys.

[Mouse] Move the mouse pointer to the upper left point of the zoomed area and left-click to put the marker on it.

5. Press the **Reference Cursor to the Marker X** side key to display the reference cursor on the marker position.

[Mouse] Right-click to display the reference cursor.

6. Move the marker to the lower right point of the zoomed area using the Marker X Position and Marker X Vertical side keys.

[Mouse] Move the mouse pointer to the lower right point of the zoomed area and left-click to put the marker on it.

7. Press the Marker → key on the front panel, and then press the Analysis Time = Marker Time side key.


The zoomed area is enclosed in a green frame.

- **8.** When you shift the zoomed area to the left or right on the frequency axis, perform these steps:
  - a. Select Marker Setup  $\rightarrow$  Marker X Position and set the new frequency center to put the marker on it.

[Mouse] Move the mouse pointer to the new frequency center and then left-click to put the marker on it.

#### **b.** Select Marker $\Rightarrow$ $\rightarrow$ Center Zoom = Marker Freq.

The zoomed area shifts to the specified frequency.

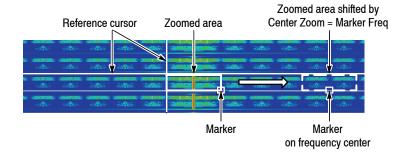


#### Figure 3-26: Zoomed area settings using the marker

**NOTE**. The zoomed spectrum is calculated based on time domain data acquired with a higher sampling rate than for obtaining a normal spectrum. So, a large zoom ratio may cause the skirt shape of the spectrum to differ from that of the spectrum without zoom. However, it does not affect carrier power measurement results.

## **Modulation Analysis (Demod Mode)**

This section describes how to perform measurements in the modulation analysis (Demod) mode.

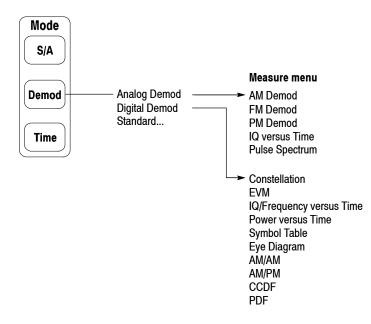

#### Figure 3-27: Demod menu structure

There are three items in the Demod menu:

- Analog Demod. Performs analog modulation analysis. Refer to page 3-37 for details.
- Digital Demod. Option 21 Only. Performs digital modulation analysis. Refer to page 3-44 for details.
- Standard... Performs modulation analysis according to a communication standard. For details, refer to the user manual included in each option as shown in Table 3–5. For RFID analysis, refer to page 3–61.

Menu item Option		Description			
W-CDMA-UL	Option 30	W-CDMA uplink analysis			
3GPP-R5-DL	Option 30	3GPP Release 5 downlink analysis			
3GPP-R5-UL	Option 30	3GPP Release 5 uplink analysis			
3GPP-R6-DL	Option 40	3GPP Release 6 downlink analysis			
3GPP-R6-UL	Option 40	3GPP Release 6 uplink analysis			
TD-SCDMA	Option 28	TD-SCDMA analysis			
GSM/EDGE	Option 24	GSM/EDGE analysis			
cdma2000-Fwd	Option 25	cdma2000 forward link analysis			
cdma2000-Rev	Option 25	cdma2000 reverse link analysis			
1xEV-DO-Fwd	Option 26	cdma2000 1xEV-DO forward link analysis			
1xEV-DO-Rev	Option 26	cdma2000 1xEV-DO reverse link analysis			
RFID	Option 21	RFID analysis			

Table 3-5: Standard menu in the Demod mode

#### **Measurement Screen Layout**

The following three views are displayed on screen by default in the Demod mode (see Figure 3-28):

- Overview: Displays all data in one acquisition block. The Timing field at the bottom of the overview indicates the contents of the main view and the subview, as well as the trigger point. For information about the trigger point, refer to *Trigger Point Indicator* on page 3-161.
- Main view: Displays the measurement results and waveform for the range specified by the green horizontal line in the overview. Some measurements show the results and waveform separately on the right and left sides of the main view.
- **Subview:** Displays the spectrum (default). The FFT processing range is indicated by the pink horizontal line shown in the overview.

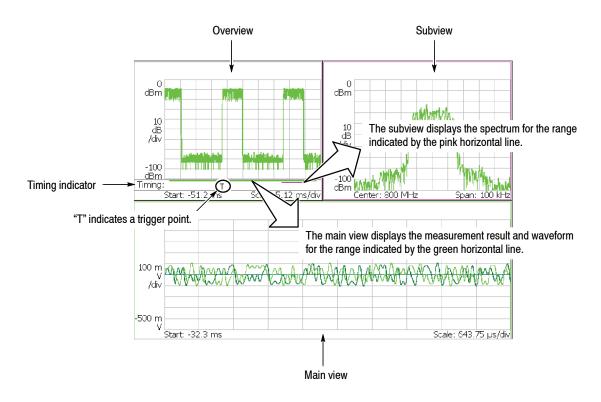


Figure 3-28: Demod mode screen

#### Setting the Analysis Range

The analysis is performed for the range specified in the overview (see Figure 3-29), and then the measurement result and waveform are displayed in the main view. Set the analysis range after data acquisition with the following procedure, using the Acquisition/Analysis menu. The range is indicated by green lines.

- 1. Press the Acquisition/Analysis key on the front panel.
- 2. Set the time length to acquire one block by pressing the Acquisition Length side key.

Suppose that one block contains N frames; the acquisition length is calculated with this equation:

(One block acquisition length) =  $N \times$  (One frame acquisition length)

One frame acquisition length is determined by span and indicated on the **Spectrum Length** side key.

- **3.** For data acquired in the continuous mode: Specify the number of the block to be analyzed by pressing the **Acquisition History** side key. 0 (zero) represents the latest block.
- **4.** Specify the time length of analysis range by pressing the **Analysis Length** side key.
- **5.** Specify the start point of the analysis range by pressing the **Analysis Offset** side key.

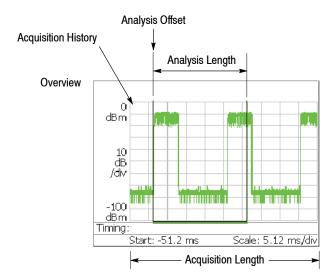


Figure 3-29: Analysis range setting in the overview

**Specifying Origin of Analysis Range Using the Marker.** You can also use the marker  $(\Box)$  in place of the **Analysis Offset** side key to specify the origin of analysis range. Do the following steps in place of step 5 on the previous page.

- 1. Press the Marker Setup key on the front panel (see Figure 3-30).
- **2.** Press the **Markers** side key to select **Single**. The marker  $(\Box)$  appears on the screen.
- **3.** Rotate the general purpose knob to move the marker to the measurement starting point.
- Press the Marker → key on the front panel and then press the Analysis Time = Marker Time side key. The green line resizes to the specified range.


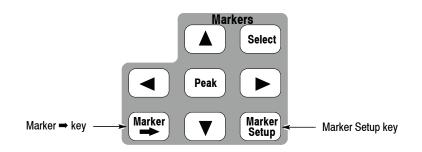


Figure 3-30: Markers keys

**Specifying Analysis Range Using the Marker and Reference Cursor.** You can use the marker and the reference cursor in place of the **Analysis Length** and the **Analysis Offset** side keys to specify the analysis range. Perform the following steps in place of steps 4 and 5 on page 3–30.

- 1. Press the View: Select key on the front panel to select the overview.
- 2. Press the Marker Setup key on the front panel.
- **3.** Press the **Markers** side key to select **Single**. The marker  $(\Box)$  appears on the screen.
- **4.** Rotate the general purpose knob to move the marker to the measurement start point.
- **5.** Press the **Reference Cursor to Marker X** side key. The reference cursor appears at the marker position (see Figure 3-31).
- **6.** Rotate the general purpose knob to move the marker to the measurement end point.
- Press the Marker → key on the front panel and then press the Analysis Time = Marker Time side key. The green line resizes to the specified range.



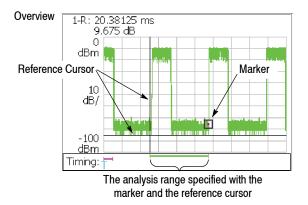

**Specifying Analysis Range Using a Mouse.** You can use a USB mouse in place of the **Analysis Length** and the **Analysis Offset** side keys to set the analysis range. Perform the following steps in place of steps 4 and 5 on page 3-30. For connecting the mouse, refer to page 2-6.

- 1. Click on the overview to select it.
- **2.** Move the mouse pointer to the measurement start point, and then left-click to put the marker on it. See Figure 3-32.
- 3. Right-click to display the reference cursor.
- **4.** Move the mouse pointer to the measurement end point, and then left-click to put the marker on it.
- Press the Marker ⇒ key on the front panel and then press the Analysis Time = Marker Time side key. The green line resizes to the specified range.

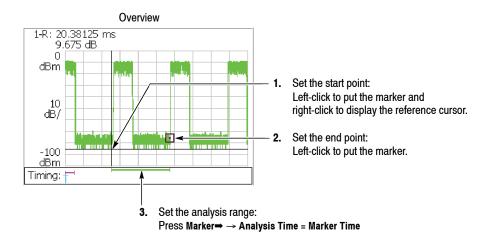


Figure 3-32: Specifying the analysis range using the mouse


#### Setting FFT Processing Range for the Subview

Set the FFT processing range (shown in Figure 3-33) for the spectrum displayed in the subview with the following procedure using the Acquisition/Analysis menu after data acquisition. The range is indicated by a pink line.

1. Press the Acquisition/Analysis key on the front panel.

The **Spectrum Length** side key shows the time for creating the FFT available in the subview. It is determined by the span internally.

2. Press the **Spectrum Offset** side key and specify the beginning of the range using the general purpose knob or numeric input keypad.

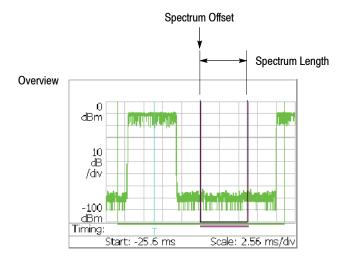

#### Figure 3-33: FFT processing range setting on the overview

## Changing the Overview and Subview

By default, the overview shows the waveform representing the signal level change over time, and the subview shows the spectrum waveform. To change the views, use the following procedure:

**NOTE**. You can change the subview only in the digital modulation analysis (*Option 21*).

- 1. Press the View: Define key on the front panel.
- 2. Press the **Overview Content...** side key and select the overview:
  - Waveform (Amplitude versus Time, default display)
  - Spectrogram

The overview displays all data in one acquisition block.

In the RFID analysis (Option 21, refer to page 3-61), you can also select Zoom (spectrogram with zoom function) for the overview.

**3.** *Digital modulation analysis (Option 21) only* Press the **Subview Content...** side key and select a view.

The subview content depends on measurement. Refer to each measurement description.


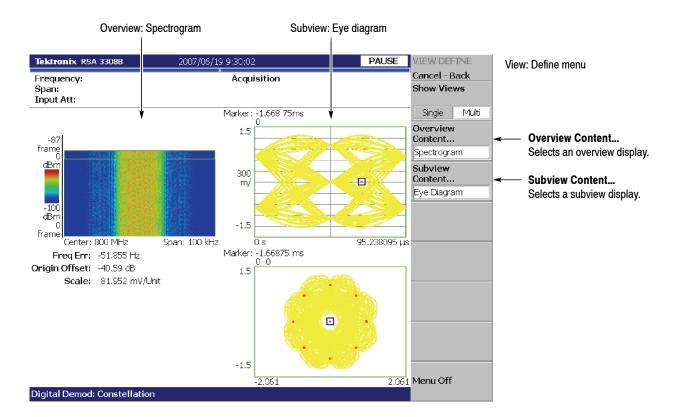


Figure 3-34: Changing the overview and subview

#### **One-View Display**

Three views are displayed on screen by default. To display one view, perform the following steps:

1. Press the View: **Define** key on the front panel.

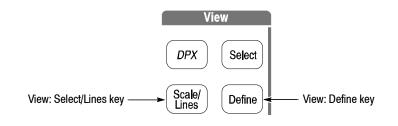

- 2. Select a view to make the single display by pressing the **Select** key in the View area. The selected view is surrounded by a white box.
- **3.** Press the **Show Views** side key to select **Single**. Figure 3-35 shows the one view display.

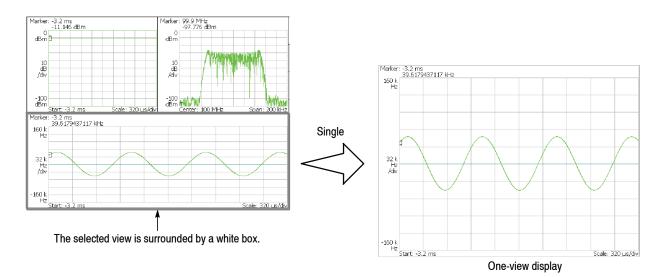


Figure 3-35: One-view display

#### **Analog Modulation Analysis**

You can select the following measurement items with the **Measure** key when you select **Analog Demod** in the Demod menu.

	AM Demod	page 3-39
--	----------	-----------

FM Demod	page 3-40
PM Demod	page 3-41
IQ versus Time	page 3-42

■ Pulse Spectrum page 3-43

**Basic Procedure** Follow this procedure to perform analog modulation analysis:

- 1. Press the Mode: **Demod** key on the front panel.
- 2. Press the Analog Demod side key.
- 3. Select a measurement item.
- 4. Display the measurement waveform:

**NOTE**. You should set appropriate frequency and span. It is important to set frequency and span as close to the measurement signal band as possible and finely adjust them. The modulated signal will not be recognized unless frequency and span are set appropriately.

For information about setting frequency, span, and amplitude, refer to Setting Frequency and Span on page 3-119 and Setting Amplitude on page 3-127.

- **a.** Set frequency by pressing the **Frequency/Channel** key on the front panel.
- **b.** Set span by pressing the **Span** key on the front panel.
- c. Set amplitude by pressing the Amplitude key on the front panel.
- **5.** Set analysis range by pressing the **Acquisition/Analysis** key on the front panel. Refer to page 3–30 for information about setting the analysis range.
- 6. Set the measurement parameters by pressing the **Meas Setup** key on the front panel. For the Meas Setup menu, refer to the following sections.

**ASK/FSK Modulation Analysis.** Basic ASK and FSK modulation analysis can be performed using Analog Modulation Analysis. This analysis includes frequency deviation and depth of modulation. Typical settings are as follows:

■ ASK signal measurement

Measurement item ... AM Demod Span ..... 500 kHz Acquisition Length ... 10.24 ms

■ FSK signal measurement

Measurement item ... FM Demod Span ...... 500 kHz Acquisition Length ... 10.24 ms Vertical Scale ...... 1.6 MHz (deviation scale: 800 kHz)

The following pages show examples for each measurement item. Refer to page 3-187 for information about setting scale and formatting each view. Refer to page 3-34 for information about changing the overview and subview.

#### **AM Demod Measurement**

Demodulates and measures an AM signal. Figure 3-36 shows an example. The main view displays the measurement results and the modulation factor versus time graph. **Measure**  $\rightarrow$  **Show Measurements** is selected in Figure 3-36.

The modulation factor is defined by the following equation:

(modulation factor) =  $(A_m - A_0) / A_0$ 

Where

A<sub>m:</sub>Amplitude of the modulation signal

A<sub>0</sub>: Amplitude of the carrier with no modulation (0% modulation)

**Meas Setup Menu.** The Meas Setup menu for the AM demod measurement contains the following control:

**Carrier Amplitude Detection.** Selects how to calculate the amplitude of the carrier with no modulation ( $A_0$  in the above equation).

- Average. *Default*. Defines A<sub>0</sub> as the average amplitude in the analysis range.
- Median. Defines A<sub>0</sub> as the median amplitude ([(maximum)+(minimum)]/2) in the analysis range.

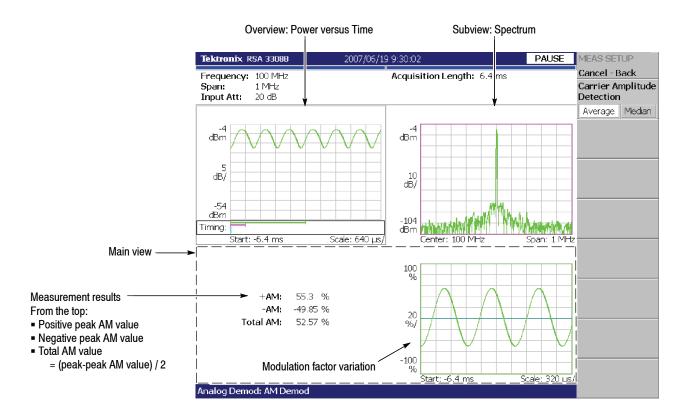


Figure 3-36: AM demod measurement

## **FM Demod Measurement** Demodulates and measures an FM signal. Figure 3-37 shows an example. The main view displays the measurement results and the frequency deviation versus time graph. **Measure** $\rightarrow$ **Show Measurements** is selected in Figure 3-37.

**Meas Setup Menu.** The Meas Setup menu for the FM demod measurement contains the following controls:

Auto Carrier. Determines whether to detect the carrier automatically.

- On. *Default*. Automatically detects the carrier for every frame.
   The error from the center frequency is shown on the Freq Error side key.
- Off. Sets the carrier frequency using Frequency Offset described below.

**Frequency Offset.** Sets the carrier frequency when Off is selected in Auto Carrier. Enter the carrier offset from the center frequency. Range: -30 to +30 MHz

**Threshold.** Sets the threshold level above which the input signal is determined to be a burst in time domain. The burst detected first is used for the analysis. Range: -100.0 to 0.0 dB


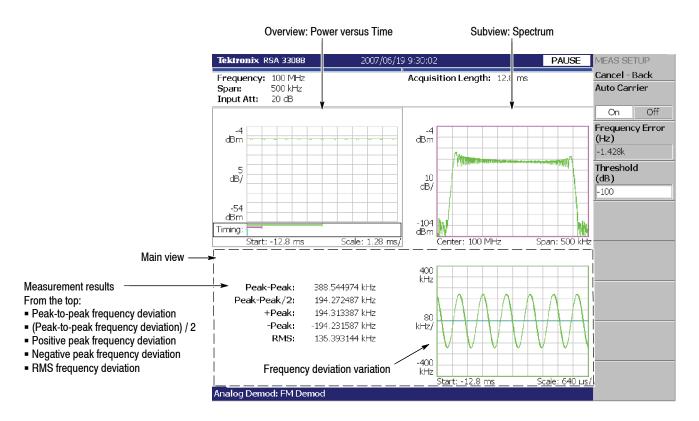



Figure 3-37: FM demod measurement

**PM Demod Measurement** Demodulates and measures a PM signal. Figure 3-38 shows an example. The main view displays the phase error versus time graph.

**Meas Setup Menu.** The Meas Setup menu for the PM demod measurement contains the following control:

**Threshold.** Sets the threshold level above which the input signal is determined to be a burst in time domain. The burst detected first is used for the measurement.

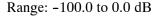

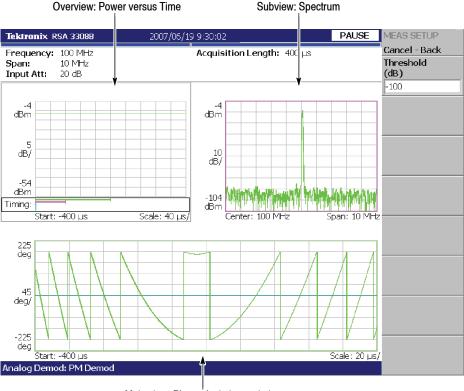

Main view: Phase deviation variation


Figure 3-38: PM demod measurement

You can select degree or radian for the angular unit by pressing System  $\rightarrow$  Instrument Setup... $\rightarrow$  Angular Units.

#### IQ versus Time Measurement

Measures I and Q level over time. Figure 3-39 shows an example. The main view shows the I/Q voltage versus time graph. I and Q are indicated in yellow and green, respectively.

The IQ level measurement has no Meas Setup menu items.

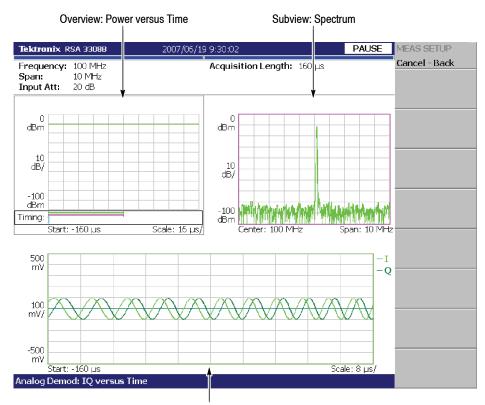


Main view: I/Q voltage versus Time (I and Q are indicated in yellow and green, respectively)

Figure 3-39: IQ versus time measurement

# Pulse Spectrum<br/>MeasurementPerforms FFT for the analysis range specified in the overview. On specifying the<br/>analysis range, refer to Setting the Analysis Range on page 3-30. Figure 3-40<br/>shows the pulse spectrum measurement. The main view shows the pulse<br/>spectrum (FFT for the analysis range).

The pulse spectrum measurement has no Meas Setup menu items.

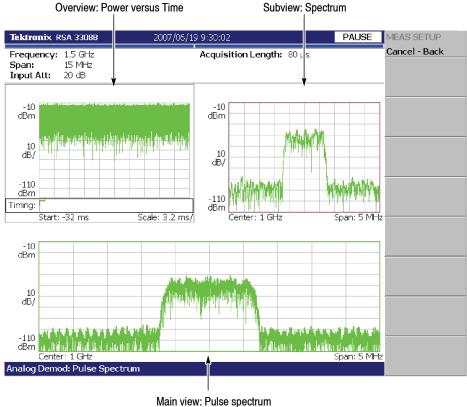

(FFT for the analysis range specified in the overview)

Figure 3-40: Pulse spectrum measurement

#### **Digital Modulation Analysis (Option 21 Only)**

You can select the following measurement items using the **Measure** key when you select **Digital Demod** in the Demod menu.

	Constellation	page 3-50
•	EVM	page 3-51
•	IQ/Frequency versus Time	page 3-52
•	Power versus Time	page 3-53
•	Symbol Table	page 3-54
•	Eye Diagram	page 3-55
•	AM/AM	page 3-56
•	AM/PM	page 3-58
•	CCDF	page 3-59
	PDF	page 3-60

Refer to Appendix C for digital demodulation symbol mapping.

**Meas Setup Menu** This subsection describes the Meas Setup menu items common to all the measurements in the digital modulation analysis. For the Meas Setup menu items specific to a measurement, refer to that measurement description.

**Parameter Presets...** Selects a communication standard. Parameters are set according to the standard you have selected (refer to Table 3–6).

Standard	Modulation	Symbol rate	Filter	α/BT	
802.15.4/OQPSK	OQPSK	1 Msps	None	0	
NADC	1/4πQPSK	24.3 ksps	Root Raised Cosine	0.35	
PDC	1/4πQPSK	21 ksps	Root Raised Cosine	0.5	
PHS	1/4πQPSK	192 ksps	Root Raised Cosine	0.5	
TETRA	1/4πQPSK	18 ksps	Root Raised Cosine	0.35	
GMS	GMSK	270.833 ksps	None	0.3	
CDPD	GMSK	19.2 ksps	None	0.5	
Bluetooth	GFSK	1 Msps	None	0.5	
P25_C4FM	P25_C4FM	4.8 ksps	Raised Cosine	0.2	

Table 3-6: Communication standard and parameters

**Modulation Type...** Selects a modulation method. Available measurement items depend on the modulation method, as shown in Table 3–7.

Measurement item	1/4πQPSK	PSK/QAM <sup>1</sup>	GMSK	GFSK	ASK	FSK	OQPSK	P25_C4FM
Constellation	~	1	~	1	~	~		~
EVM	~	1	~				~	
IQ/Frequency versus Time	~	$\checkmark$	~			~		~
Power versus Time					~			
Symbol Table	~	$\checkmark$	~		~	~		~
Eye Diagram	~	$\checkmark$	~	$\checkmark$	~	~	~	~
AM/AM		$\checkmark$						
AM/PM		$\checkmark$						
CCDF	$\checkmark$	$\checkmark$	~				$\checkmark$	
PDF	1	1	1				1	

 Table 3-7: Modulation type and available measurement item

<sup>1</sup> Includes BPSK, QPSK, 8PSK, D8PSK, 16QAM, 32QAM, 64QAM, 128QAM, and 256QAM.

#### **Restrictions of Use.**

- When the modulation type is P25 C4FM, the analysis is performed only for the span of 20 kHz and 50 kHz. For the 50 kHz span, Acquisition Length in the Acquisition/Analysis menu must be set to greater than 25.6 ms.
- When the modulation type is OQPSK, set the span so that the sampling rate is more than three times as fast as the symbol rate. The sampling rate varies according to span setting. For the detail, refer to Table A-9 in *Appendix A*.

Modulation Parameters... Sets the following modulation parameters:

**Symbol Rate.** Enters the symbol rate for demodulating digitally modulated signals. Symbol rate and bit rate are related as follows:

(Symbol rate) = (Bit rate)/(Number of bits per symbol)

For example, the number of bits per symbol is 3 for 8PSK.

**Measurement Filter...** Selects a filter for demodulating a digitally modulated signal:

- None (no filter)
- Root Raised Cosine

**Reference Filter...** Selects a filter for creating a reference signal:

- None (no filter)
- Raised Cosine
- Gaussian
- Half Sine

Refer to *Process Flow for Digitally-Modulated Signals* on page 3-48 for the measurement and reference filters.

**Filter Parameter.** Enters an  $\alpha$ /BT value for the Measurement Filter and Reference Filter. Range: 0.0001 to 1

**Auto Modulation Depth.** Valid when the modulation type (refer to page 3-45) is ASK. Determines whether to automatically detect or manually set the modulation depth used to distinguish between the two states of an ASK signal.

- **On.** *Default*. Automatically calculates the modulation depth for the analysis range and displays the value in the Modulation Depth side key.
- Off. Sets the modulation depth using the Modulation Depth side key.

**Modulation Depth.** Sets the modulation depth when Auto Modulation Depth described above is Off. Range: 0 to 100%

**Auto Frequency Deviation.** Valid when the modulation type (refer to page 3-45) is FSK or GFSK. Determines whether to automatically detect or manually set the frequency deviation used to distinguish between the two states of an FSK or GFSK signal.

- **On.** *Default.* Automatically calculates the frequency deviation for the analysis range and displays the value in the Frequency Deviation side key.
- Off. Sets the frequency deviation using the Frequency Deviation side key.

**Frequency Deviation.** Sets the frequency deviation when Auto Frequency Deviation is Off. Range: 0 to Span/2 Hz

**Decoding Format.** Valid when the modulation type (refer to page 3-45) is ASK, FSK or GFSK. Selects the method used to decode data bits from each symbol.

- NRZ
- Manchester
- Miller

Auto Carrier. Selects whether to use automatic carrier detection.

- **On.** *Default.* Automatically detects carriers for every frame and displays errors from center frequency as "Freq Error" on the screen.
- Off. Sets the carrier frequency using Frequency Error.

**Frequency Error.** Sets the carrier frequency when Off is selected in Auto Carrier. Enter the carrier offset from the center frequency.

## Shifting Q Data for<br/>OQPSKWhen the modulation type is OQPSK, you can shift Q data by half a symbol<br/>relative to I data by pressing View: Define $\rightarrow$ Q Data Half Symbol Shift.

The Q Data Half Symbol Shift side key has the following selections:


- +. Shifts Q data by half a symbol in the positive direction on the time axis.
- 0. Does not shift Q data (default).
- -. Shifts Q data by half a symbol in the negative direction on the time axis.

This function is available in the following views:

- Constellation
- EVM
- Eye diagram
- IQ versus Time

#### Process Flow for Digitally-Modulated Signals

To determine the necessary settings in digital modulation analysis, you must know the digitally-modulated signal process in the analyzer. Figure 3-41 outlines the process.

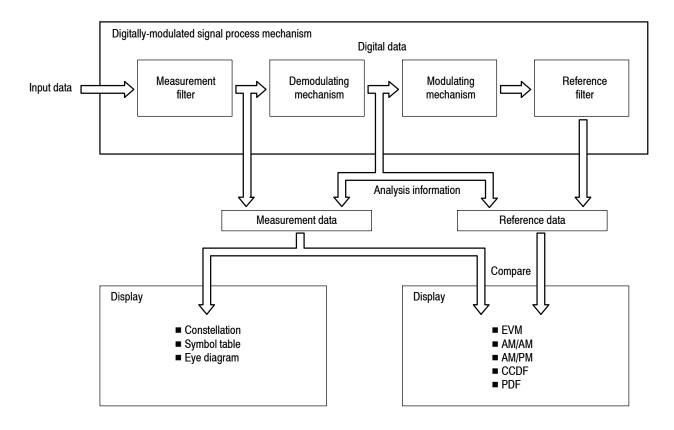

#### Figure 3-41: Process flow for digitally-modulated signals

The input signal goes through the measurement filter after being converted into a digital signal, and then is stored as measurement data and is concurrently demodulated. The demodulated signal is modulated again, goes through the reference filter, and is stored as reference data. The Constellation, Symbol table, and Eye diagram displays are based on measurement data, and the EVM, AM/AM, AM/PM, CCDF, and PDF displays are based on the comparison of the measurement data and the reference data.

**Basic Procedure** To perform digital modulation analysis, follow these steps:

- 1. Press the Mode: **Demod** key on the front panel.
- 2. Press the Digital Demod side key and then select the measurement item.
- 3. Display the measurement waveform:

**NOTE**. You should set appropriate frequency and span. It is important to set frequency and span as close to the measurement signal band as possible and finely adjust them. The modulated signal will not be recognized unless frequency and span are set appropriately.

For information about setting frequency, span, and amplitude, refer to Setting Frequency and Span on page 3-119 and Setting Amplitude on page 3-127.

- **a.** Set the frequency by pressing the **Frequency/Channel** key on the front panel.
- **b.** Set the span by pressing the **Span** key on the front panel.
- c. Set the amplitude by pressing the Amplitude key on the front panel.
- 4. Set analysis range by pressing the Acquisition/Analysis key on the front panel. Refer to page 3-30 for details.

**NOTE**. When the modulation type (refer to page 3-45) is ASK, FSK or GFSK, the measurement requires at least 16 symbols in the analysis length (refer to page 3-30 for setting the analysis length).

5. Set the measurement parameters by pressing the Meas Setup key on the front panel. Refer to page 3-45 for the Meas Setup menu description.

**If No Measurement Result or Waveform is Displayed on the Main View.** Unless the analyzer obtains valid data for the measurement, neither the result nor waveform appears on the main view. In this case, try the following steps:

- Set the center frequency to the middle of the measurement signal bandwidth.
- Set the span near the measurement signal bandwidth.
- Set the data acquisition length (Acquisition/Analysis → Acquisition Length) larger to increase the amount of data.


The following pages show examples for each measurement item. Refer to page 3-187 for information about setting scale and formatting each view. Refer to page 3-34 for information about changing the overview and subview.

#### Constellation Measurement


Performs digital demodulation processing and displays a constellation diagram. Figure 3-42 shows an example. The main view displays the measurement results and the constellation diagram. When the modulation type (Meas Setup  $\rightarrow$ Modulation Type...) is P25 C4FM, the modulation fidelity is displayed as the measurement result.

**NOTE**. In the constellation view, I and Q signals are normalized to prevent the scale from changing when signal attenuation changes.

For the Meas Setup menu, refer to page 3-45. For setting views, refer to page 3-187.

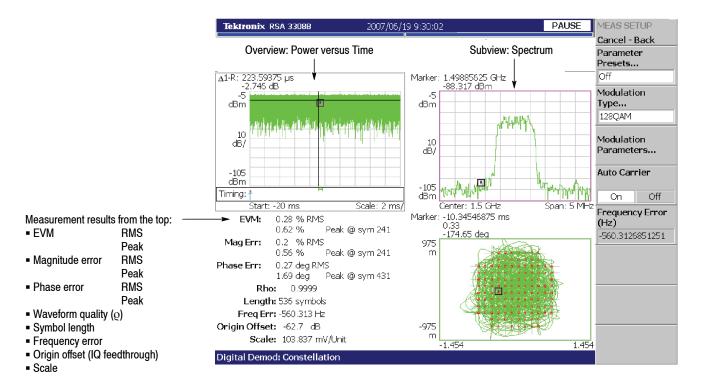


#### P25 C4FM modulation only

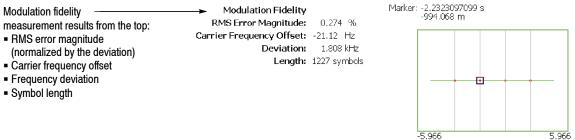

#### Figure 3-42: Constellation analysis

**EVM Measurement** Measures EVM (Error Vector Magnitude). Figure 3-43 shows an example. The main view displays the measurement results and the EVM versus time graph.

For the Meas Setup menu, refer to page 3-45. For setting views, refer to page 3-187.

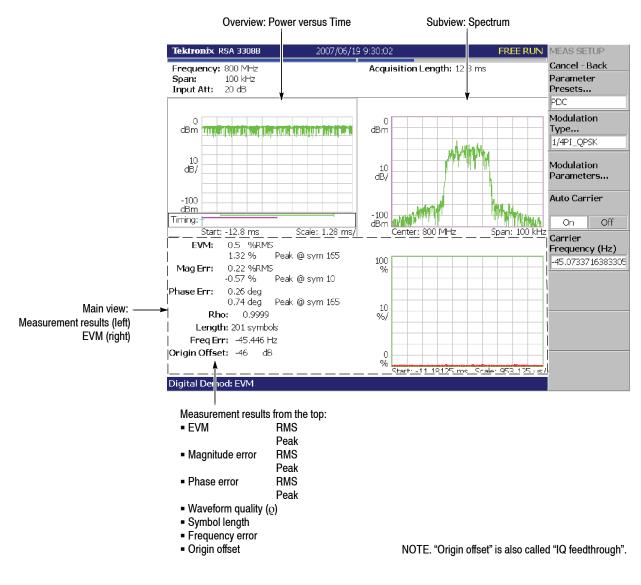



You can select degree or radian for the angular unit by pressing System  $\rightarrow$  Instrument Setup... $\rightarrow$  Angular Units.

#### IQ/Frequency versus Time Measurement

Observes I/Q signal voltage change over time. When the modulation type (Meas Setup  $\rightarrow$  Modulation Type...) is FSK, GFSK, or P25 C4FM, it displays frequency change over time. Figure 3-44 shows an example. The main view shows I/Q voltage versus time graph where I and Q are indicated in yellow and green, respectively, or shows frequency deviation versus time graph for the FSK, GFSK, or P25 C4FM modulation. For the P25 C4FM modulation, the main view also shows the modulation fidelity measurement results and the waveforms with the red dots indicating the symbol points.

For the Meas Setup menu, refer to page 3-45. For setting views, refer to page 3-187.

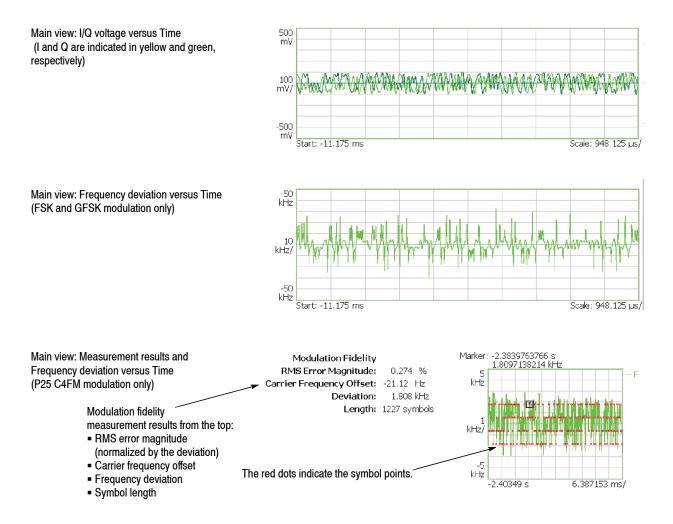


Figure 3-44: IQ/Frequency versus time measurement

#### Power versus Time Measurement

Observes signal power change over time. Figure 3-45 shows an example. The main view displays a power versus time graph.

For the Meas Setup menu, refer to page 3-45. For setting views, refer to page 3-187.




Main view: Power versus Time

Figure 3-45: Power variance measurement

## **Symbol Table Analysis** Performs digital demodulation processing and displays a symbol table in the main view. Figure 3-46 shows an example.

For the Meas Setup menu, refer to page 3-45. For setting views, refer to page 3-187.

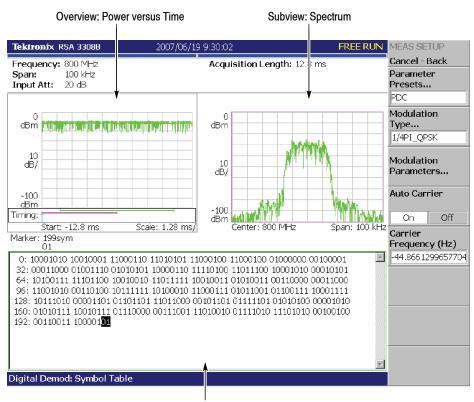


Main view: Symbol table

Figure 3-46: Symbol table analysis

### **Eye Diagram Analysis** Performs digital demodulation processing and displays an eye diagram in the main view. Figure 3-47 shows an example.

For the Meas Setup menu, refer to page 3-45. For setting views, refer to page 3-187.

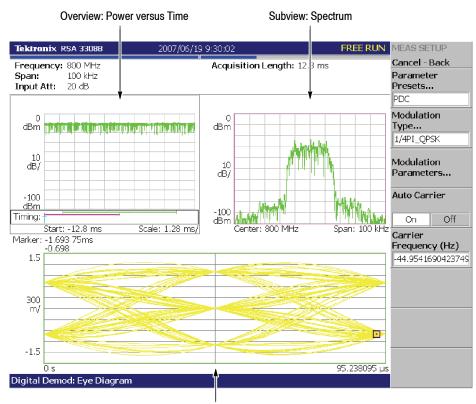

Main view: Eye diagram

Figure 3-47: Eye diagram analysis

You can select degree or radian for the angular unit by pressing System  $\rightarrow$  Instrument Setup... $\rightarrow$  Angular Units.

AM/AM Measurement Performs digital demodulation processing and displays an AM/AM distortion graph for measuring non-linearity characteristics of a DUT (device under test) such as an RF amplifier. Figure 3-49 on page 3-57 shows an example. The main view displays the XY graph of recovered reference amplitude versus measured signal amplitude and the table of AM/AM coefficients extracted from curve fitting.

**Meas Setup Menu.** The Meas Setup menu for the AM/AM measurement contains the following controls. For information about the common controls for the digital modulation analysis, refer to *Meas Setup Menu* on page 3-45.

Linear Signal Region Unit. Selects the unit to set the Linear Signal Region.

- **dB.** *Default.* Specifies the Linear Signal Region with a value relative to the maximum power measured in the analysis range.
- **dBm.** Specifies the Linear Signal Region with an absolute power.

**Linear Signal Region.** Sets a region that should have an ideal AM/AM characteristic. Generally, signal distortion increases with amplitude. Linear Signal Region sets the region where the trend is assumed to be linear (see Figure 3-48). The analysis is performed with this assumption. The region is indicated by blue oblique lines in the main view.

Range: -100 to 50 dB (the default is -10 dB)

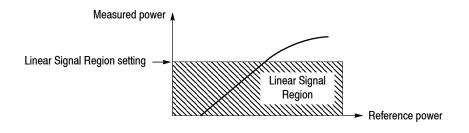


Figure 3-48: Setting Linear Signal Region

**Maximum Coefficient.** Sets the maximum order of the best-fit curve polynomial. The polynomial is expressed as follows (n: Max 15):

$$f(x) = a_0 x^0 + a_1 x^1 + a_2 x^2 + \dots + a_n x^n$$

The values of  $a_n$  are shown in the main view (see Figure 3-49 on page 3-57).

**Display Reference Line.** Determines whether the reference line is visible. The reference line represents an ideal AM/AM trend line.

Display Best-Fit Line. Determines whether the best-fit line is visible.

**Linear Signal Region Mask.** Determines whether the linear signal region is visible.

Scroll Table. Scrolls the coefficient table when it has many rows.

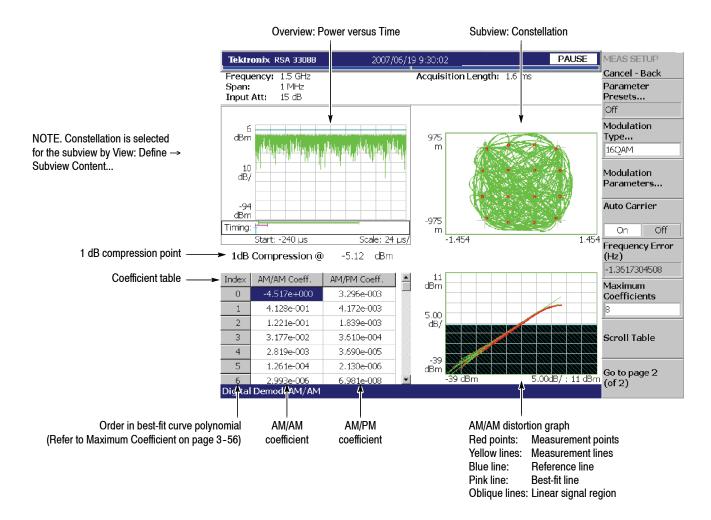


Figure 3-49: AM/AM measurement

**AM/PM Measurement** Performs digital demodulation processing and displays an AM/PM distortion graph for measuring non-linearity characteristics of a DUT (device under test) such as an RF amplifier. Figure 3–50 shows an example. The main view displays the XY graph of recovered reference amplitude versus measured phase error and the table of AM/PM coefficients extracted from curve fitting.

**Meas Setup Menu.** The Meas Setup menu for the AM/PM measurement has the same controls as in the AM/AM measurement, except for the Linear Signal Region Mask. Refer to page 3-56.

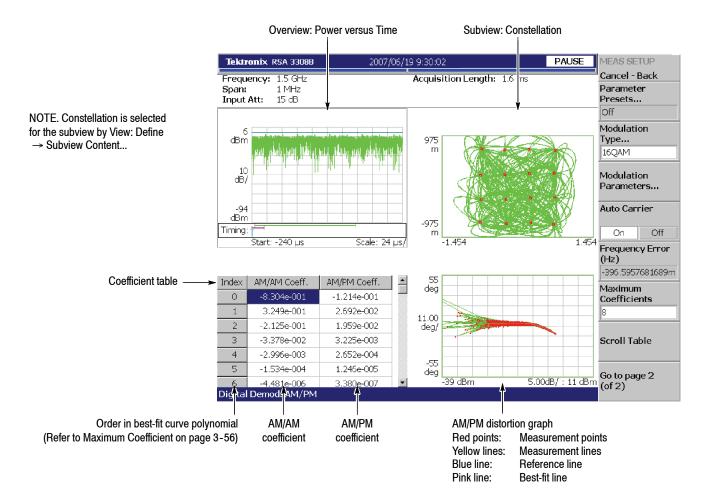


Figure 3-50: AM/PM measurement

**CCDF Measurement** Performs digital demodulation processing and displays a CCDF graph. Figure 3-51 shows an example. The main view displays the XY graph of power level versus the probability that the instantaneous power of the input is above that power level. For the CCDF measurement, refer to page 3-88 in the Time mode.

**Meas Setup Menu.** The Meas Setup menu for the CCDF measurement contains the following controls. (For information about the common controls for the digital modulation analysis, refer to page 3–45.)

**Linear Signal Region Unit.** *Display only.* You can set it in the AM/AM or AM/PM measurement (refer to page 3-56).

**Linear Signal Region.** *Display only.* You can set it in the AM/AM or AM/PM measurement (refer to page 3-56).

**Horizontal Division.** Sets the horizontal interval between display points. Range: 0.01 to 1 dB (the default is 0.1 dB)

Display Gaussian Line. Determines whether to display the Gaussian line.

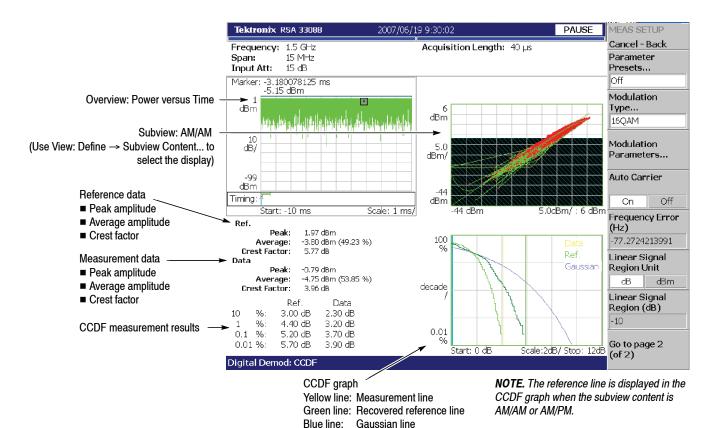


Figure 3-51: CCDF measurement

**PDF Measurement** Performs digital demodulation processing and displays a PDF (Probability Distribution Function) graph. Figure 3-52 shows an example. The main view displays the XY graph of power level versus occurrence probability. The horizontal axis represents power level relative to the mean value at the center.

**Meas Setup Menu.** The Meas Setup menu for the PDF measurement contains the following controls. On the common controls for the digital modulation analysis, refer to *Meas Setup Menu* page 3-45.

**Linear Signal Region Unit.** *Display only.* You can set it in the AM/AM or AM/PM measurement (refer to page 3-56).

**Linear Signal Region.** *Display only.* You can set it in the AM/AM or AM/PM measurement (refer to page 3-56).

**Horizontal Division.** Sets the horizontal interval between display points. Range: 0.01 to 1 dB (the default is 0.1 dB)

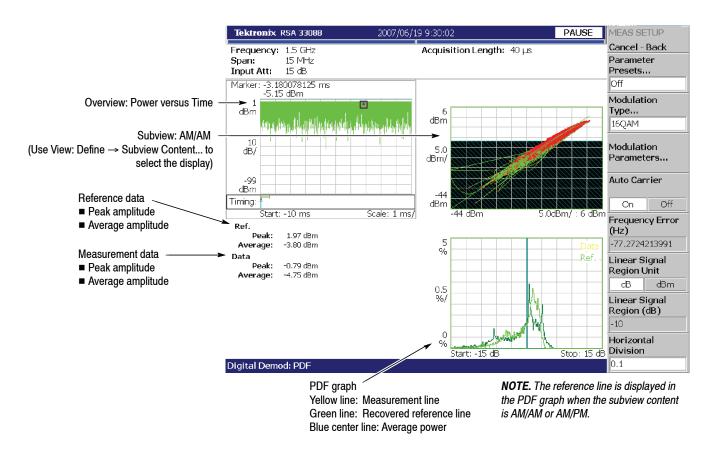


Figure 3-52: PDF measurement

## **RFID Analysis (Option 21 Only)**

This section describes the method of RFID (Radio Frequency Identification) analysis. RFID is an automatic individual identification technology that uses a radio. The IC microprocessors in which data is recorded and the antenna module are considered to be the "Tag". A tag can be incorporated in various objects or living beings. The tag information is read by the radio wave from the "Interrogator" (also called Reader/Writer) and if necessary, the data is written in the tag. This instrument performs an analysis based on the ISO/IEC and EPCglobal standards related to RFID.

#### **Compatible Standards** RFID analysis is based on the standards listed in Table 3-8.

Standard	Description
ISO/IEC 18000-4	Information technology - Radio frequency identification for item management - Part 4: Parameters for air interface communications at 2.45 GHz Mode-1 <sup>1</sup>
ISO/IEC 18000-6	Information technology - Radio frequency identification for item management - Part 6: Parameters for air interface communications at 860 MHz to 960 MHz Type-A, Type-B
ISO/IEC 18000-7	Information technology - Radio frequency identification for item management - Part 7: Parameters for active air interface communications at 433MHz
ISO/IEC 18092	Information technology - Telecommunications and information exchange between systems - Near Field Communication - Interface and Protocol (NFCIP-1)
ISO/IEC 14443-2	Identification cards - Contactless integrated circuit(s) cards - Proximity cards - Part 2:Radio frequency power and signal interface Type-A, Type-B
ISO/IEC 15693-2	Identification cards - Contactless integrated circuit(s) cards - Vicinity cards - Part 2: Air interface and initialization
EPCglobal	EPC Radio-Frequency Identity Protocols Class-1 Generation-2 UHF RFID Protocol for Communications at 860 MHz-960 MHz Version 1.0.9 (ISO/IEC 18000-6 Type-C)
C0G1 (EPCglobal Gen1 Class0)	Auto-ID Center Draft protocol specification for a 900 MHz Class 0 Radio Frequency Identification Tag
C1G1 (EPCglobal Gen1 Class1)	Auto-ID Center TECHNICAL REPORT 860MHz-930MHz Class 1 Radio Frequency Identification Tag Radio Frequency & Logical Communication Interface Specification Candidate Recommenda- tion, Version 1.0.1

#### Table 3-8: Compatible RFID standards

<sup>1</sup> This instrument does not support ISO/IEC 18000-4 Mode-2.

#### **Measurement Items**

Table 3-9 shows the measurement items in the RFID analysis. (Select the items using the **Measure** key.)

Table 3-9: Measurement items for RFID analysis

Measure menu	Measurement items	
Carrier	<ul> <li>Carrier frequency</li> <li>OBW (Occupied Bandwidth)</li> <li>EBW (Emission Bandwidth)</li> <li>Maximum EIRP (Effective Isotropically Radiated Power)</li> </ul>	
Spurious	■ Spurious	
ACPR	ACPR (Adjacent Channel Leak Power Ratio)	
Power On/Down	<ul> <li>Transmission Power on/down rise/fall time</li> <li>Settling time</li> <li>Overshoot</li> <li>Undershoot</li> </ul>	
RF Envelope / FSK Pulse	<ul> <li>On/Off width</li> <li>Duty cycle</li> <li>On/Off ripple</li> <li>Rise/Fall time</li> </ul>	
Constellation	Modulation depth	
Eye Diagram	<ul> <li>Modulation index</li> <li>Frequency error</li> </ul>	
Symbol Table	■ Bit rate or Tari	

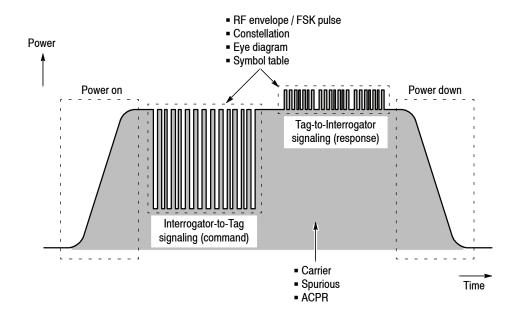


Figure 3-53: RFID waveform and measurement items

Figure 3-53 illustrates the measurement items with a typical RFID waveform. The analyses are performed on the carrier in the Carrier, Spurious, and ACPR measurements, on the transmission power on/down in the Power On/Down measurement, and on the Interrogator-to-Tag or Tag-to-Interrogator signaling in the RF Envelope, Constellation, Eye Diagram, and Symbol Table measurements.

**Basic Procedure** Use the following procedure to perform the RFID analysis.

**NOTE**. If you use a signal input (Input  $\rightarrow$  Signal Input Port...) other than RF, the measurements are not guaranteed.

For Option 02 (256 MB memory), the acquisition length (Acquisition/Analysis  $\rightarrow$  Acquisition Length) and the analysis length (Acquisition/Analysis  $\rightarrow$  Analysis Length) are both up to 1.28 s at the maximum real-time span (15 MHz). The length of a burst is up to 40 ms. A burst more than 40 ms is analyzed every 40 ms as one burst each. The time length varies according to span settings.

- 1. Press the **Demod** key on the front panel.
- 2. Select Standard...  $\rightarrow$  RFID using the side keys.
- **3.** Select a measurement item by pressing the side key. For example, press the **Carrier** side key to perform the carrier measurement.
- 4. Adjust the amplitude and frequency to display the measurement waveform.

Refer to page 3–119 for setting frequency and span. Refer to page 3–127 for setting amplitude.

5. Set the analysis range in the overview. (This is not necessary in the spurious and ACPR measurements.)


Refer to page 3-30 for details of the analysis range settings.

If you use the zoom function, do these steps:

- Press the View: **Define** key.
- Press the **Overview Content...** side key and select **Zoom**.

Refer to *Zoom Function* on page 3-23 for setting the analysis area on the spectrogram with the zoom function.

<ul> <li>6. Press the Meas Setup key and set the measurement parameters. Refer to Meas Setup Menu in each measurement description.</li> <li>7. Press the Analyze side key and carry out the analysis. (It is not necessary in the spurious and ACPR measurements.)</li> <li>Carrier Measurement Measures the carrier frequency, OBW (Occupied Bandwidth), EBW (Emission Bandwidth), and the maximum EIRP (Effective Isotropically Radiated Power).</li> <li>Meas Setup Menu. The following items are provided in the Meas Setup menu of the carrier measurement.</li> <li>Analyze. Perform measurements for acquired data in the analysis range.</li> <li>NOTE. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.</li> <li>Counter Resolution. Sets the resolution of the carrier frequency measurement. Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)</li> <li>Refer to page 3-10 for the carrier frequency measurement in the S/A mode.</li> <li>Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (MB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the CBW measurement. Range: 80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset or the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> <li>Figure 3-54 shows an example of the carrier measurement.</li> </ul>				
<ul> <li>7. Press the Analyze side key and carry out the analysis. (It is not necessary in the spurious and ACPR measurements.)</li> <li>Carrier Measurement Measures the carrier frequency, OBW (Occupied Bandwidth), EBW (Emission Bandwidth), and the maximum EIRP (Effective Isotropically Radiated Power).</li> <li>Meas Setup Menu. The following items are provided in the Meas Setup menu of the carrier measurement.</li> <li>Analyze. Perform measurements for acquired data in the analysis range.</li> <li>NOTE. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.</li> <li>Counter Resolution. Sets the resolution of the carrier frequency measurement. Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)</li> <li>Refer to page 3-10 for the carrier frequency measurement in the S/A mode.</li> <li>Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 Hz).</li> </ul>		6. Press the Meas Setup key and set the measurement parameters.		
(It is not necessary in the spurious and ACPR measurements.)         Carrier Measurement       Measures the carrier frequency, OBW (Occupied Bandwidth), EBW (Emission Bandwidth), and the maximum EIRP (Effective Isotropically Radiated Power).         Meas Setup Menu. The following items are provided in the Meas Setup menu of the carrier measurement.       Analyze. Perform measurements for acquired data in the analysis range.         NOTE. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.       NOTE. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.         Counter Resolution. Sets the resolution of the carrier frequency measurement. Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)         Refer to page 3-10 for the carrier frequency measurement in the S/A mode.         Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.         Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).         Power Ratio for OBW (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).         Power Ratio for CBW (dB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 1000 dB (default: -20 dB).         Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Ra		Refer to Meas Setup Menu in each measurement description.		
<ul> <li>Bandwidth), and the maximum EIRP (Effective Isotropically Radiated Power).</li> <li>Meas Setup Menu. The following items are provided in the Meas Setup menu of the carrier measurement.</li> <li>Analyze. Perform measurements for acquired data in the analysis range.</li> <li>NOTE. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.</li> <li>Counter Resolution. Sets the resolution of the carrier frequency measurement. Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)</li> <li>Refer to page 3-10 for the carrier frequency measurement in the S/A mode.</li> <li>Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the CBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>				
<ul> <li>the carrier measurement.</li> <li>Analyze. Perform measurements for acquired data in the analysis range.</li> <li><b>NOTE</b>. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.</li> <li><b>Counter Resolution.</b> Sets the resolution of the carrier frequency measurement. Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)</li> <li>Refer to page 3-10 for the carrier frequency measurement in the S/A mode.</li> <li><b>Power Ratio Unit for OBW or EBW.</b> Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to <i>OBW Measurement</i> (page 3-8) and <i>EBW Measurement</i> (page 3-10) in the spectrum analysis (S/A mode).</li> <li><b>Power Ratio for OBW</b> (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li><b>Power Ratio for EBW (dB).</b> Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li><b>Antenna Gain.</b> Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li><b>Channel Bandwidth.</b> Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>	Carrier Measurement			
<ul> <li>NOTE. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.</li> <li>Counter Resolution. Sets the resolution of the carrier frequency measurement. Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)</li> <li>Refer to page 3-10 for the carrier frequency measurement in the S/A mode.</li> <li>Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit.</li> <li>Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit.</li> <li>Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>				
<ul> <li>side key to perform the measurement for the modified settings.</li> <li>Counter Resolution. Sets the resolution of the carrier frequency measurement. Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)</li> <li>Refer to page 3-10 for the carrier frequency measurement in the S/A mode.</li> <li>Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>		Analyze. Perform measurements for acquired data in the analysis range.		
<ul> <li>Setting range: 1 mHz to 1 MHz (switching by 10 times, default: 1 Hz)</li> <li>Refer to page 3-10 for the carrier frequency measurement in the S/A mode.</li> <li>Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit.</li> <li>Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit.</li> <li>Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>		side key to perform the measurement for the modified settings. Counter Resolution. Sets the resolution of the carrier frequency measurement.		
<ul> <li>Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB.</li> <li>Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>				
<ul> <li>measurement: % or dB.</li> <li>Refer to <i>OBW Measurement</i> (page 3-8) and <i>EBW Measurement</i> (page 3-10) in the spectrum analysis (S/A mode).</li> <li><b>Power Ratio for OBW</b> (%). <i>Valid when % is selected for the unit</i>. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li><b>Power Ratio for EBW (dB).</b> <i>Valid when dB is selected for the unit</i>. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li><b>Antenna Gain.</b> Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li><b>Channel Bandwidth.</b> Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>		Power Ratio Unit for OBW or EBW. Selects the unit for the OBW or EBW measurement: % or dB. Refer to OBW Measurement (page 3-8) and EBW Measurement (page 3-10) in		
<ul> <li>the spectrum analysis (S/A mode).</li> <li>Power Ratio for OBW (%). Valid when % is selected for the unit. Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>				
<ul> <li>Sets the power ratio of the carrier and span region of the OBW measurement. Range: 80 to 100% (default: 99%).</li> <li>Power Ratio for EBW (dB). Valid when dB is selected for the unit. Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>				
<ul> <li>Sets the power ratio of the carrier and span region of the EBW measurement. Range: -80 to 100 dB (default: -20 dB).</li> <li>Antenna Gain. Sets the amplitude offset for the maximum EIRP measurement. The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB).</li> <li>Channel Bandwidth. Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).</li> </ul>		Sets the power ratio of the carrier and span region of the OBW measurement.		
The maximum EIRP value to which the amplitude offset value is added is displayed as the measurement result. Range: -100 to 100 dB (default: 0 dB). <b>Channel Bandwidth.</b> Sets the channel bandwidth for the maximum EIRP measurement. Range: 0 to Span (default: 0 Hz).		Sets the power ratio of the carrier and span region of the EBW measurement.		
measurement. Range: 0 to Span (default: 0 Hz).		The maximum EIRP value to which the amplitude offset value is added is		
Figure 3-54 shows an example of the carrier measurement.				
		Figure 3-54 shows an example of the carrier measurement.		

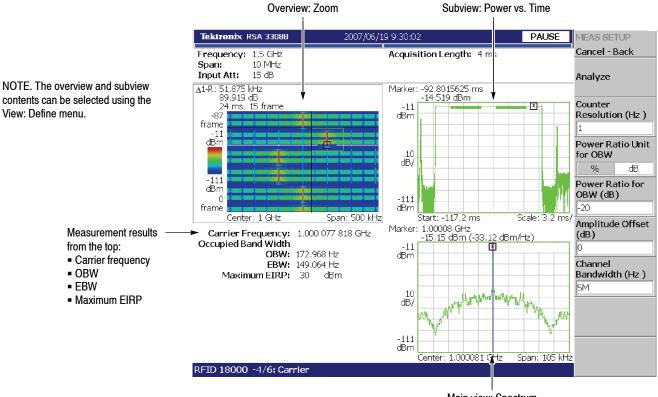

Main view: Spectrum



**Changing the View.** The view contents can be changed in the next item by pressing the View: **Define** key.

Show Views. Select whether to display a single view or multiple views.

- Single. Only displays the view selected using the View: Select key.
- Multi. Displays the overview, sub view and main view.

Overview Content... Selects the overview information.

- Spectrogram
- Waveform (Power versus Time waveform)
- Zoom (Spectrogram with the zoom function)

Selects Waveform or Zoom while selecting the analysis range.

Refer to *Zoom Function* on page 3-23 for setting the analysis area on the spectrogram with the zoom function.

Subview Content... Selects the subview information.

- Spectrum
- Power versus Time
- Frequency versus Time
- Zoomed Spectrum (Spectrum of a zoomed area in the Zoom display)

Setting the View Scale. Refer to page 3-187 for information on setting the scale for each view. Note that when the subview is Power versus Time or Frequency versus Time, Default Scale is used instead of Full Scale in the View Scale menu. The time-domain waveform is processed in units of 256K data points. When the number of points equivalent to the analysis length (Acquisition/Analysis  $\rightarrow$  Analysis Length) is 256K or less, pressing the **Default Scale** side key sets the horizontal scale so that the whole waveform is displayed in the analysis range. When it is greater than 256K, the horizontal full scale is set to the time equivalent of 256K points, displaying part of the waveform in the analysis range.

# **Spurious Measurement** Measures the spurious level of the RF signal. It is the same as the spurious measurement in the S/A mode (refer to page 3-12).

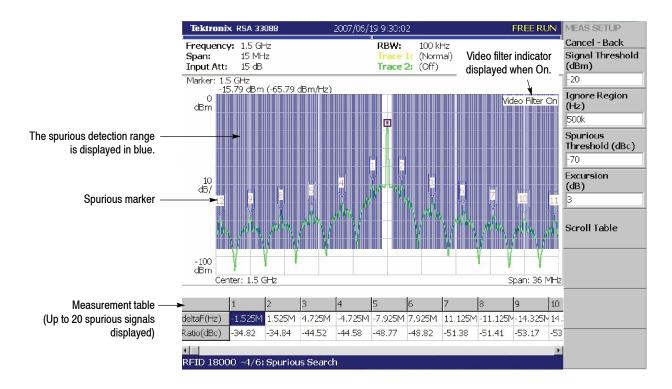


Figure 3-55: Spurious measurement

**Meas Setup Menu.** The Meas Setup menu of the spurious measurement is same as that of the S/A mode. Refer to page 3-12.

Figure 3-55 (on page 3-66) shows an example of the spurious measurement. For the detected spurious signals, the spurious marker is displayed on the waveform by assigning the numbers from 1 in the maximum order of amplitude. Furthermore, the frequency difference (deltaF) and amplitude ratio (Ratio) with the carrier is displayed in the table, in the lower part of the screen.

**Using Video Filter.** The spurious measurement in the RFID analysis has a video filter function for compatibility with measurement data from conventional swept spectrum analyzers. Normally swept spectrum analyzers smooth the detected signal with a video filter having a bandwidth equal to the resolution bandwidth. However, you can set the bandwidth and sweep time on this analyzer. When you use the video filter, press the **RBW/FFT** key on the front panel to set the following items:

**Video Filter.** Determines whether to use the video filter. When you use the filter, select On. Then the following items appear:

**VBW.** Sets the frequency bandwidth of the video filter. Range: 1 Hz to 1 GHz. The setting value may be limited by the sweep time setting.

**Sweep Time for VBW.** Sets the sweep time to scan a set span. Range: 1 µs to 100 s. **ACPR Measurement** The ACPR (Adjacent Channel Leak Power Ratio) measurement is similar to that of the S/A mode (refer to page 3–6). However, in the ACPR measurement of the RFID analysis, the adjacent channel is measured for up to 25 channels from up and down.

**Meas Setup Menu.** The Meas Setup menu of the ACPR measurement is the same as that of the S/A mode except that Scroll Table is added. Refer to page 3-6.

**Scroll Table.** Selects an adjacent channel (1 to 25 in the upper/lower side) from the measurement table displayed in the bottom of the screen.

Figure 3-56 shows an example of the ACPR measurement.

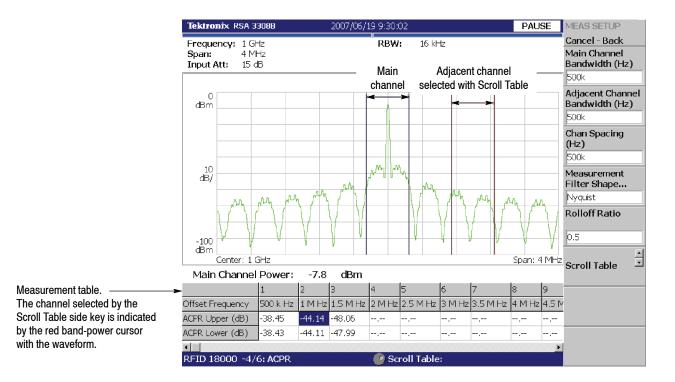


Figure 3-56: ACPR measurements

# Power On/Down<br/>MeasurementMeasures the rise/fall time, settling time, overshoot, and undershoot of the<br/>transmission power pulse.

**NOTE**. In the analysis range, there must be at least one pulse rise or fall. If there is not, the analysis is not performed.

**Meas Setup Menu.** The following menu items are common to the power on/down, RF envelope, constellation, eye diagram, and symbol table measurements.

Analyze. Perform measurements for acquired data in the analysis range.

**NOTE**. When you change settings in the Meas Setup menu, press the **Analyze** side key to perform the measurement for the modified settings.

**Standard Type...** Sets the parameters according to the standard as shown in Table 3-10.

#### Table 3-10: Standard settings

Standard type	Link	Modulation type	Decoding format	Bit rate / Tari <sup>1</sup>
18000-4-1	Interrogator to Tag	ASK	Manchester / NRZ	30 kbps
	Tag to Interrogator		FM0 / NRZ	
18000-6-A	Interrogator to Tag	ASK	PIE (Type A)	20 µs (Tari)
			NRZ	33 kbps
	Tag to Interrogator		FM0 / NRZ	33 kbps
18000-6-B	Interrogator to Tag	ASK	Manchester / NRZ	10 kbps
	Tag to Interrogator		FM0 / NRZ	
18000-6-C	Interrogator to Tag	DSB-ASK	PIE (Type C)	6.25 μs (Tari)
	SSB-ASK PR-ASK		NRZ	26.7 kbps
	Tag to Interrogator	ASK	FM0 Miller (M_2 / M_4 / M_8) NRZ	26.7 kbps
18000-7	-	FSK	Manchester / NRZ	27.7 kbps
14443-2-A	Interrogator to Tag	ASK	Modified Miller / NRZ	105.9375 kbps
	Tag to Interrogator	Subcarrier OOK	Manchester / NRZ	105.9375 kbps
		Subcarrier BPSK	NRZ-L (4 periods)	211.875 kbps
			NRZ-L (2 periods)	423.75 kbps

Standard type	Link	Modulation type	Decoding format	Bit rate / Tari <sup>1</sup>	
14443-2-B	Interrogator to Tag	ASK	NRZ	105.9375 kbps	
	Tag to Interrogator	Subcarrier BPSK	NRZ-L (8 periods)	105.9375 kbps	
			NRZ-L (4 periods)	211.875 kbps	
			NRZ-L (2 periods)	423.75 kbps	
18092(424k)	-	ASK	Manchester / NRZ	212 kbps	
15693-2	Interrogator to Tag	ASK	1 out of 4 / 1 out of 256	Not settable	
	Tag to Interrogator		SSC High / SSC Low DSC High / DSC Low	(fixed internally)	
C0G1	Interrogator to Tag	ASK	PWM	40 kbps	
Та	Tag to Interrogator	FSK	NRZ		
C1G1	Interrogator to Tag	ASK	PWM	70.18 kbps	
	Tag to Interrogator		Bit Cell	140.35 kbps	
Manual	-	ASK DSB-ASK SSB-ASK	Manchester / FM0 / NRZ Miller / Modified Miller Miller (M_2 / M_4 / M_8)	40 kbps	
		PR-ASK OOK	PIE (Type A)	20 µs (Tari)	
	UUK		PIE (Type C)	6.25 μs (Tari)	
		FSK	Manchester / FM0 / NRZ Miller / Modified Miller	9.6 kbps	

#### Table 3-10: Standard settings (Cont.)

[Abbrev.] DSB: Double sideband; SSB: Single sideband; PR: Phase reversal; OOK: On/Off key; PIE: Pulse interval encoding; PWM: Pulse Width Modulation SSC: Single Subcarrier; DSC: Double Subcarrier

<sup>1</sup> Default values are shown. It can be set from 1 bps to 51.2 Mbps (bit rate) or from 1 ns to 1 s (Tari).


**NOTE**. The standard types of C0G1 and C1G1 disable the constellation measurement.

Tari is the reference time interval for a data-0 in Interrogator-to-Tag signaling, defined in ISO/IEC 18000-6 Part A specification (an abbreviation for Type A Reference Interval).

The period of the NRZ-L decoding format (2, 4, and 8 periods) indicate the number of subcarrier cycles per symbol.

**NOTE**. This instrument does not support one-period NRZ-L.

The M values (M\_2, 4, and 8) of the Miller decoding format indicate the number of subcarrier cycles per symbol. Refer to the ISO/IEC and EPCglobal documents listed in *Compatible Standards* on page 3-61 for details on the encoding/decoding formats except the modified Miller code illustrated in Figure 3-57.

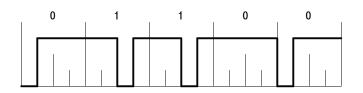

#### Figure 3-57: Modified Miller code

**Link...** Selects whether to measure Interrogator-to-Tag signaling or Tag-to-Interrogator signaling.

- Interrogator. Detects the interrogator preamble from a measurement signal and decodes the signal with the interrogator decoding format.
- **Tag.** Detects the tag preamble from a measurement signal and decodes the signal with the tag decoding format.

Modulation Type. Selects the modulation type. Refer to Table 3-10.

Decoding Format. Selects the decoding format. Refer to Table 3-10.

**Auto Bit Rate** or **Auto Tari.** Selects whether to set the bit rate or Tari automatically. Tari is the parameter for the PIE Type A and C decoding formats.

- **On.** Detects the bit rate or Tari automatically based on the Bit Rate or Tari set value.
- **Off.** *Default.* Sets the bit rate or Tari using the Bit Rate or Tari side key.

**NOTE.** The Auto Bit Rate or Auto Tari can be set to On or Off for the four measurements of RF Envelope, Constellation, Eye Diagram, and Symbol Table, but it is fixed to Off for the Power On/Down measurement. When you run the Power On/Down measurement, these four measurements are also performed at the same time with the Auto Bit Rate or Auto Tari still set to Off. If you want to turn it On to perform the measurements, select **On** for **Auto Bit Rate** or **Auto Tari** and then press the **Analyze** side key in each measurement.

Auto Bit Rate, Auto Tari, Bit Rate, and Tari are invalid for the 15693-2 standard.

Bit Rate. Sets the bit rate. Range: 1 bps to 51.2 Mbps. Refer to Table 3-10.

Tari. Sets the Tari. Range: 1 ns to 1 s. Refer to Table 3-10.

**NOTE**. You must set the bit rate or Tari manually to perform the analysis properly. Refer to Setting the Bit Rate / Tari on page 3-74.

**Settling Error Width +-.** Sets an error range for determining the settling time. See Figure 3-58. Range: 1 to 100% (default: 5%)

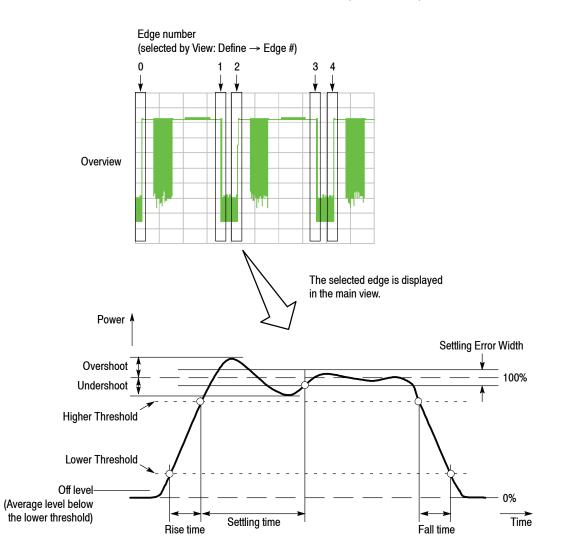


Figure 3-58: Transmission power on/down measurement parameters

**Lower Threshold.** Sets the lower threshold value for measuring the rise/fall time of the pulse. See Figure 3–58. Range: 1% to (Middle Threshold). Default value: 10%, or 5% only for the 14443-2-A, 18092(424k), and 15693-2 standards.

**Middle Threshold.** Sets the middle threshold value for measuring the width of the pulse. See Figure 3-58. Range: (Lower Threshold) to (Higher Threshold). Default value: 50%, or 60% only for the 14443-2-A, 18092(424k), and 15693-2 standards.

**Higher Threshold.** Sets the higher threshold value for measuring the rise/fall time of the pulse. See Figure 3-58. Range: (Middle Threshold) to 99%. Default value: 90%, or 95% only for the 15693-2 standard.

**Interpolation Points.** When the sample rate is lower than the bit rate, because the data cannot be analyzed correctly, the shortage of the sample rate is calculated and interpolated by the spline method.

Range: 0 to 7 (default: 0). Zero means no interpolation.

(Effective sample rate) = (Actual sample rate) × [(Interpolation Points setting value) + 1]

The actual sample rate and effective sample rate are displayed in the main view of the power on/down and RF envelope measurements (see Figure 3–59 on page 3–75 and Figure 3–60 on page 3–77).

Advanced Parameter... For the standard types of 14443-2-A and 14443-2-B only. Sets the following three parameters:

Channel Filter... Selects the filter for demodulating an input signal.

- None (no filter, default for Interrogator)
- Raised Cosine (default for Tag)

**Preamble.** Determines whether to search for the preamble.

- **On.** Searches for the preamble while analyzing data. The preamble is displayed in yellow in the symbol table.
- Off. Analyzes data without searching for the preamble.

Side Band. Selects the sideband to analyze.

- **Upper.** *Default.* Analyzes the upper sideband.
- Lower. Analyzes the lower sideband.

Setting the Bit Rate / Tari. You must set the bit rate or Tari near the true value (within  $\pm 20\%$  of the true value approximately depending on waveform) to perform the analysis properly in the power on/down, RF envelope, constellation, eye diagram, and symbol table measurements.

**NOTE**. If the bit rate or Tari setting value is out of range, the analyzer cannot distinguish between the rising edge of the power-on CW and that of the RF envelope pulse.

In the measurement setup, enter the bit rate or Tari and perform the analysis using the following steps:

- 1. Press the Meas Setup key on the front panel.
- 2. Press the Standard Type... side key and select the standard.
- 3. Select Link, Modulation Type, and Decoding Format as appropriate.
- 4. Select Off in the Auto Bit Rate or Auto Tari side key.
- 5. Press the **Bit Rate** side key and set the bit rate, or press the **Tari** side key and set the Tari. Enter a value near the true value (within about  $\pm 20\%$  of the true value).
- 6. Press the Auto Bit Rate or Auto Tari side key to select On.


**NOTE**. When you set Auto Bit Rate or Auto Tari to On, the analyzer estimates the bit rate or Tari from your input value and performs measurements based on the estimate. If you want to use your input value directly in measurements without using an estimate, leave the Auto Bit Rate or Auto Tari setting Off.

7. Press the **Analyze** side key to perform the analysis. The analysis is performed based on the bit rate or Tari set value.

After the analysis has been completed, the estimated bit rate or Tari is displayed in the **Bit Rate** side key.

In the constellation, eye diagram, and symbol table measurements, the calculated bit rate or Tari is also displayed in the main view (see Figure 3-62 on page 3-80).

Figure 3-59 shows an example of the power on/down measurements. You can save the measurement table to a file (\*.csv) by **Save**  $\rightarrow$  **Save Table**. For file operations, refer to page 3-233.

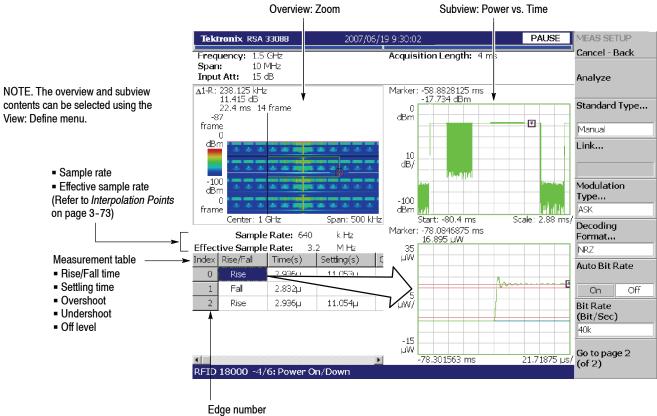

(selected by View: Define  $\rightarrow$  Edge #)

#### Figure 3-59: Power on/down measurements

**Changing the View.** The contents of the view can be changed in the following item by pressing the View: **Define** key.

The following items are the same as the RF carrier measurement. Refer to page 3-64:

Show Views Overview Content... Subview Content... **Edge #.** Selects an edge number (index number of the measurement table) for displaying a waveform on the main screen. The edge number is assigned to the rising/falling edge of each pulse. (Refer to page 3-72, Figure 3-58.)

**Guidelines.** Selects whether to display the guideline (red) along with the waveform in the main view.

- **On.** *Default*. Displays the guideline.
- **Off.** No guideline is displayed.

Scroll Table. Scrolls the measurement table in the main view horizontally.

**Changing the Vertical Unit.** The unit is W (watts) by default for the vertical axis of the graph displayed in the main view. You can also select dBm and V (volts) using View: Scale/Lines  $\rightarrow$  Scale  $\rightarrow$  Unit.

**NOTE**. The measurement table in the main view always displays the results for waveform data with the default unit of W. Selecting the unit does not affect the table contents.

#### RF Envelope / FSK Pulse Measurement

Measures time attributes for each envelope on the Interrogator-to-Tag or Tag-to-Interrogator signaling and displays the measurement table with the amplitude versus time waveform.

The RF Envelope measurement switches to the FSK Pulse measurement with one of the following settings in the Meas Setup menu:

- Standard Type is set to 18000-7
- Standard Type is set to Manual and Decoding Format FSK

The vertical axis changes from amplitude in the RF Envelope measurement to frequency in the FSK Pulse.

**Meas Setup Menu.** The Meas Setup menu of the RF envelope / FSK pulse measurement is the same as the power on/down measurement. Refer to page 3–69.

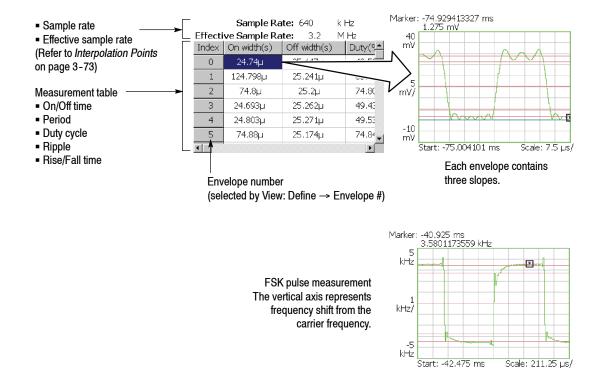


Figure 3-60: RF envelope / FSK pulse measurement (main view)

Figure 3-60 shows an example of the RF envelope / FSK pulse measurement main view. The following items are provided in the measurement table (see also Figure 3-61).

ltem	Description	
Index	Envelope # (Select by View: Define $\rightarrow$ Envelope)	
On-width / Off-width	Pulse width for pulse on or off	
Period	Period ( = (On-width) + (Off-width) )	
Duty	Duty cycle ( = (On-width) / (Period) )	
On-ripple / Off-ripple	Ripple of on-time or off-time	
Slope 1, 2, and 3	Rise/Fall time of slope 1, 2, and 3 (Slope 3 is Slope 1 of the envelope on the right)	
T1, T2, T3, and T4 <sup>1</sup>	Pulse characteristic values defined in the 14443-2-A, 18092(424k), and 15693-2 standards	
<sup>1</sup> For the definitions of T1 T2 T3 and T4 refer to the ISO/IEC 18092 14443-2 and		

Table 3-11: Details of the RF envelope / FSK pulse measurement table

<sup>1</sup> For the definitions of T1, T2, T3, and T4, refer to the ISO/IEC 18092, 14443-2, and 15693-2 standard specifications shown in Table 3-8 on page 3-61.

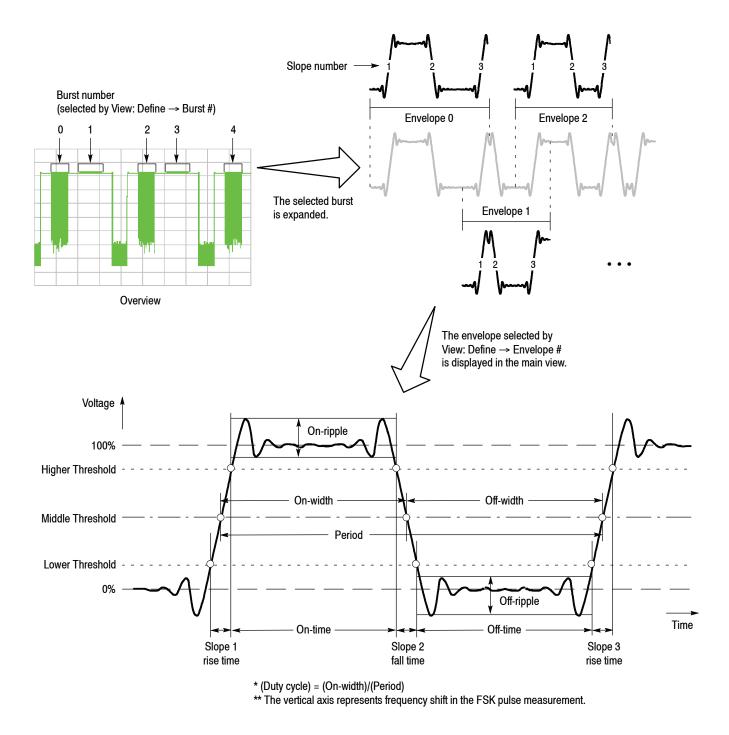


Figure 3-61: RF envelope / FSK pulse measurement parameters

You can save the measurement table to a file (\*.csv) by **Save**  $\rightarrow$  **Save Table**. For file operations, refer to page 3-233. **Changing the View.** The contents of the view can be changed in the next item by pressing the View: **Define** key.

The next items are the same as the carrier measurement. Refer to page 3-65.

Show Views Overview Content...

Subview Content... Selects the information displayed in the subview.

- Spectrum
- Power versus Time
- Frequency versus Time
- Zoomed Spectrum (Spectrum for a zoomed area in the Zoom display)
- RF Envelope (Waveform same as the RF envelope in the main view)
- Constellation (Waveform same as the constellation in the main view)
- Eye Diagram (Waveform same as the eye diagram in the main view)
- Symbol Table (Same as the symbol table in the main view)

Refer to page 3-23 for the zoom function. Refer to page 3-66 for setting the view scale.

**Burst #.** Selects the burst number to be measured. (Refer to page 3-78, Figure 3-61.) The setting range is limited as follows (default: 0):

(Burst # setting) + (Edge # setting)  $\leq$  20000 (for Edge #, refer to page 3-76)

**Envelope #.** Selects the envelope number (the index number in the RF envelope measurement table) that displays the waveform in the main view. Range: 0 to 65536. (Refer to page 3-78, Figure 3-61.)

**Guidelines.** Selects whether to display the guidelines (red color) along with the waveform in the main view.

- **On.** *Default*. Displays the guidelines.
- Off. No guideline is displayed.

Scroll Table. Scrolls the measurement table in the main view horizontally.

**Display Area.** *Valid for the constellation and the eye diagram views.* Sets the percentage of display area (sample points) from the beginning in the selected burst.

Range: 0.1 to 100% (default: 100%). 100% represents the whole burst.

**Changing the Vertical Unit.** The unit for the vertical axis of the graph displayed in the main view is V (volts) by default. You can also select dBm and W (watts) using View: Scale/Lines  $\rightarrow$  Scale  $\rightarrow$  Unit.

**NOTE**. The measurement table in the main view always displays the results for waveform data with the default unit of W. Selecting the unit does not affect the table contents.

Constellation, Eye Diagram, and Symbol Table Measures modulation attributes on the Interrogator-to-Tag or Tag-to-Interrogator signaling and displays the constellation, eye diagram, or symbol table (see Figure 3-62).

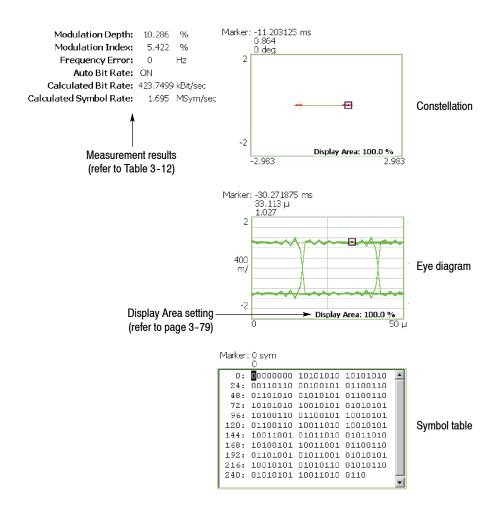


Figure 3-62: Constellation, eye diagram, and symbol table (main view)

The constellation, eye diagram, and symbol table are based on the digital modulation analysis function. Refer to *Digital Modulation Analysis* on page 3-44.

**Meas Setup Menu.** The Meas Setup menu for the constellation, eye diagram, and symbol table is the same as the power on/down measurement. Refer to page 3-69.

The measurement result display items depend on standards (Meas Setup  $\rightarrow$  Standard Type) as shown in Table 3-12.

#### Table 3-12: Measurement result display items

18000-4-1, 18000-6-A/B/C, 14443-2-A/B, 18092(424k), C0G1, and C1G1 standards and Manual	18000-7 standard	15693-2 standard	Manual (FSK modulation only)
<ul> <li>Modulation depth</li> <li>Modulation index</li> <li>Frequency error</li> <li>Auto bit rate (on/off) <sup>1</sup></li> <li>Calculated bit rate <sup>2</sup></li> <li>Calculated symbol rate <sup>2</sup></li> <li>Subcarrier jitter (rms) <sup>3</sup></li> <li>Frequency offset <sup>3</sup></li> <li>Calculated Tari data-0 <sup>4</sup> (duration of a data-0)</li> <li>Calculated Tari data-1 <sup>4</sup> (duration of a data-1)</li> </ul>	<ul> <li>Frequency deviation</li> <li>Preamble detection (Interrogator / Tag)</li> <li>Frequency error</li> <li>Auto bit rate (on/off) <sup>1</sup></li> <li>Calculated bit rate</li> <li>Calculated symbol rate</li> </ul>	<ul> <li>Modulation depth</li> <li>Modulation index</li> <li>Frequency error</li> <li>Calculated bit rate</li> </ul>	<ul> <li>Frequency deviation</li> <li>Frequency error</li> <li>Auto bit rate (on/off) <sup>1</sup></li> <li>Calculated bit rate</li> <li>Calculated symbol rate</li> </ul>

<sup>1</sup> Set in the Meas Setup menu.

<sup>2</sup> For other than PIE decoding.

<sup>3</sup> For Subcarrier BPSK only.

<sup>4</sup> For PIE decoding only.

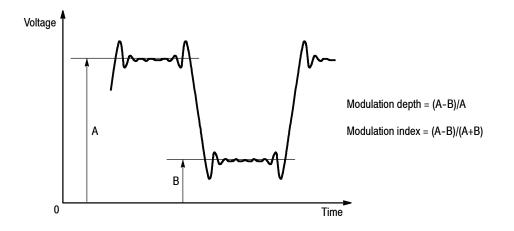


Figure 3-63 shows the definition of the modulation depth and index.

#### Figure 3-63: Definition of the modulation depth and index

Table 3-13 shows the definition of symbol values displayed in the symbol table.

Symbol value	Definition
0	0
1	1
Х	Don't care
Р	Preamble
S	Frame Sync or SOF <sup>1</sup>
Ν	Null
I	Interrogator
Т	Tag
E	EOF <sup>1</sup>


Table 3-13: Symbol value definition

<sup>1</sup> 15693-2 standard only.

**Changing the View.** The view controls for the constellation, eye diagram, and symbol table are the same as the RF envelope measurement except for the Scroll Table side key. Refer to page 3-79.

# Time Analysis (Time Mode)

This section describes how to perform measurements in the time analysis (Time mode).

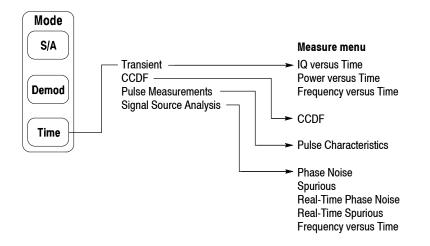

#### Figure 3-64: Time menu structure

There are three items in the Time menu:

- Transient. Performs time variation measurement. Refer to page 3-84 for details.
- CCDF. Performs CCDF measurement. Refer to page 3-88 for details.
- Pulse Measurements. Performs pulse measurements. Refer to page 3-93 for details.
- Signal Source Analysis. Option 21 Only. Performs signal source analysis. Refer to page 3-103 for details.

#### Measurement Screen Layout

The measurement screen layout in the time analysis (Time mode) is the same as in the modulation analysis (Demod mode). Refer to *Measurement Screen Layout* on page 3-29 for the view functions, setting an analysis range, and changing the view content.

### **Time Variation Measurement**

You can select the following time variation measurement items using the **Measure** key when you select **Transient** in the Time menu.

- IQ versus Time
- Power versus Time
- Frequency versus Time

The time variation measurements have no Meas Setup menu items.

**Basic Procedure** Follow these steps to perform the time variation measurement:

- 1. Press the Mode: Time key on the front panel.
- 2. Press the **Transient** side key and select a measurement item.
- 3. Display the measurement waveform:

**NOTE**. You should set appropriate frequency and span. It is important to set frequency and span as close to the measurement signal band as possible and finely adjust them. The modulated signal will not be recognized unless frequency and span are set appropriately.

For details on setting frequency, span, and amplitude, refer to Setting Frequency and Span on page 3-119 and Setting Amplitude on page 3-127.

- **a.** Set the frequency by pressing the **Frequency/Channel** key on the front panel.
- **b.** Set the span by pressing the **Span** key on the front panel.
- c. Set the amplitude by pressing the Amplitude key on the front panel.
- **4.** Set the analysis range by pressing the **Acquisition**/**Analysis** key on the front panel.

Refer to page 3-30 for details about setting analysis range.

The following pages show the examples for each measurement item. Refer to page 3-187 for information about setting scale and formatting each view. Refer to page 3-34 for information about changing the overview.

## **IQ versus Time Measurement** Observes I and Q signal level over time. Figure 3-65 shows an example. The main view displays I/Q voltage versus time. I and Q are indicated in yellow and green, respectively.




yellow and green, respectively)

Figure 3-65: IQ versus Time measurement

#### Power versus Time Measurement

Observes signal power change over time. Figure 3-66 shows an example. The main view displays power level versus time.

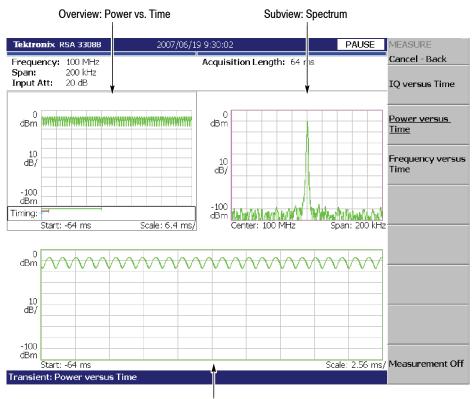


Main view: Power vs. Time

Figure 3-66: Power variance measurement

#### Frequency versus Time Measurement

Observes signal frequency change over time. Figure 3–67 shows an example. The main view displays the frequency deviation from the center frequency versus time.

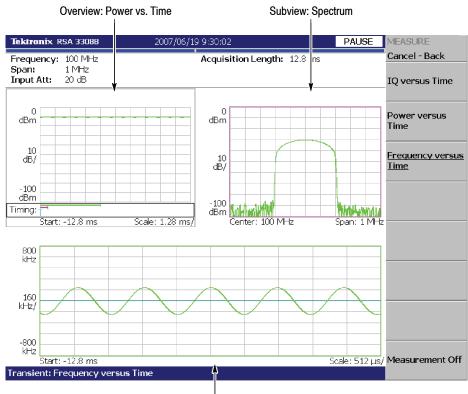

Main view: Frequency deviation vs. Time

Figure 3-67: Frequency measurement

## **CCDF Measurement**

CCDF (Complementary Cumulative Distribution Function) represents the probability that the peak power above average power of a measured signal exceeds a threshold. The analyzer displays the ratio of peak power to average power along the horizontal axis, and the probability that the ratio is exceeded along the vertical axis.

This CCDF analysis function and the real-time analysis function allow you to measure the time-varying crest factor in the time series for code-multiplexing signals such as CDMA/W-CDMA signals, and multi-carrier signals such as OFDM signals. This function is useful in designing amplifiers for CDMA/W-CDMA and OFDM.

# **CCDF Calculation Process** The CCDF analysis obtains the distribution of the amplitude of observed signals and makes a graph of accumulation from the threshold. CCDF is calculated using this formula:

$$SP(X) = \int_{X}^{Max} P(Y) \, dY$$
  
where  
$$CCDF(X) = SP(X + Average)$$
$$P = Probability density of amplitudeMax = Maximum of amplitudeAverage: = Average of amplitude$$

The analyzer processes input signals internally using the following procedure (see Figure 3-68):

- 1. Measure the amplitude of the input signal over time.
- **2.** Determine the amplitude distribution.
- **3.** Obtain CCDF using the above formula.

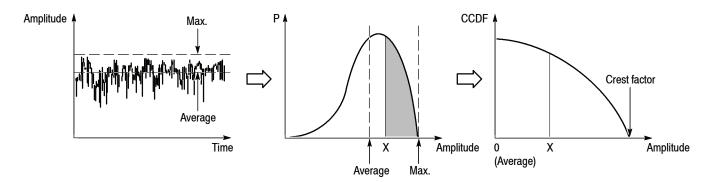


Figure 3-68: CCDF calculation process

**Basic Procedure** 

Follow these steps to perform the CCDF measurement:

- 1. Press the Time key on the front panel.
- 2. Press the CCDF side key.
- **3.** Display the measured waveform.

**NOTE**. Make sure to set the proper frequency and span. It is important to set the frequency and the span as close to the measured signal bandwidth as possible and adjust precisely. If not, modulated signals are not recognized correctly.

For details about setting frequency, span, and amplitude, refer to Setting Frequency and Span on page 3-119 and Setting Amplitude on page 3-127.

- **a.** Press the **Frequency/Channel** key on the front panel to set the frequency.
- **b.** Press the **Span** key on the front panel to set the span.
- c. Press the Amplitude key on the front panel to set the amplitude.
- 4. Set the analysis range by pressing the Acquisition/Analysis key on the front panel. Refer to page 3-30 for details about setting the analysis range.

**NOTE**. The CCDF measurement does not have Analysis Length and Analysis Offset controls in the Acquisition/Analysis menu.

5. Set the measurement parameters described just below by pressing the **Meas Setup** key on the front panel.

**Single Acquisition.** In single acquisition mode, you can set the total number of accumulated data points with **Trace/Avg**  $\rightarrow$  **Maximum Points**, ranging from 2048 to 10<sup>15</sup>.

**Meas Setup Menu** The Meas Setup menu for the CCDF measurement contains the following controls:

Display Gaussian Line. Determines whether to display Gaussian line on screen.

**Display Reference Line.** Determines whether to display the most recently stored reference line on screen.

**Store Reference Line.** Stores the current CCDF curve as a new reference line. Automatically enables the reference line display.

**Reset Measurement.** Performs the CCDF calculation again from the beginning. The calculation is performed accumulatively until this side key is pressed.

**CCDF Auto-Scaling.** Selects whether to fix the scale of the horizontal axis (amplitude) of CCDF graph display.

- **Off.** *Default.* The scale of the horizontal axis is set to a fixed value using CCDF Scale.
- **On.** Displays the graph with the signal peak value as the maximum value (right edge) of the horizontal axis.

**CCDF Scale.** Sets the full-scale of the horizontal axis of CCDF graph display when CCDF Auto Scaling is Off. Range: 1 to 100 dB.

**Threshold.** Sets the threshold which defines the samples to be included in the CCDF calculation. Range: -250 to 130 dBm.

#### **Measurement Display**

Figure 3-69 shows an example of the CCDF measurement.

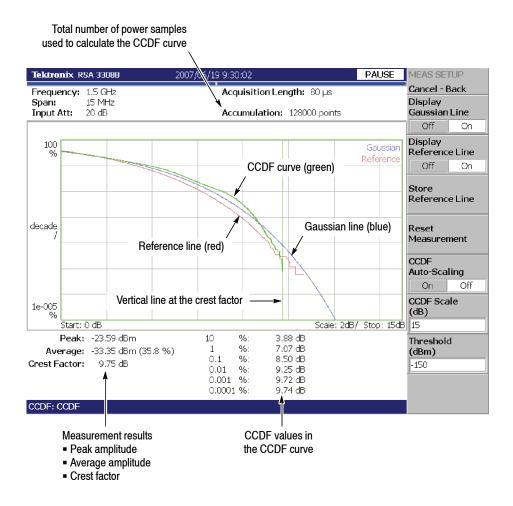


Figure 3-69: CCDF single view

**Changing the View** You can change the contents of the views by pressing the View: **Define** key with the following controls:

Show Views. Selects the view style:

- Single. Displays only the view selected by the View: Select key.
- **Multi.** Displays the overview, the subview, and the main view.

Overview Content... Selects the content of the overview:

- Spectrogram
- Waveform (power versus time)

Refer to page 3-187 for information about setting the view scale and format. Figure 3-70 shows an example of the CCDF multi-view.

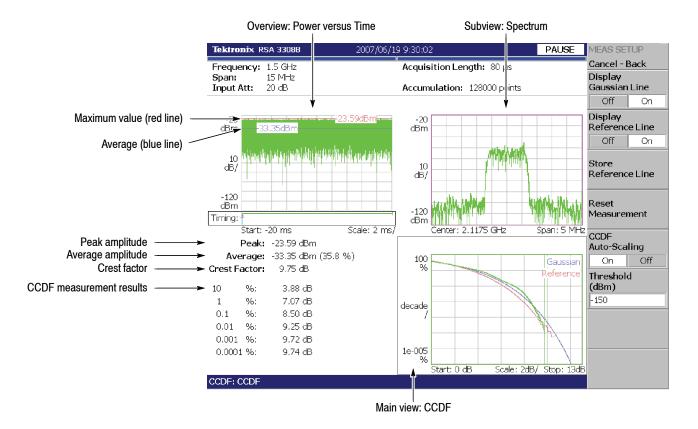


Figure 3-70: CCDF multi-view

## **Pulse Measurements**

This section describes how to perform the pulse measurements. Select **Time**  $\rightarrow$  **Pulse Measurements**  $\rightarrow$  **Pulse Characteristics** to measure various pulse characteristics.

# **Measurement Items** The following lists shows the measurement items (Measure menu) for pulse characteristics and their definitions (see Figure 3-71).

- Pulse Width. Measures the time from rising edge to falling edge at the level of -3 dB (50%) of the mean power for the pulse-on.
- **Peak Power.** Measures the maximum power during the pulse-on.
- **On/Off Ratio.** Measures the ratio of the mean power during the pulse-on to that during the pulse-off time.
- **Pulse Ripple.** Measures the difference between the maximum and the minimum power during the pulse-on.
- **Repetition Interval.** Measures the time from a pulse rising edge to the next pulse rising edge.
- **Duty Cycle.** Measures the ratio of the pulse width to the repetition interval.
- Pulse-Pulse Phase. Measures the phase difference between the first pulse and the selected pulse in the analysis window.
- **Channel Power.** Measures the channel power for the pulse-on spectrum.
- **OBW.** Measures the OBW (Occupied Bandwidth) for the pulse-on spectrum.
- **EBW.** Measures the EBW (Emission Bandwidth) for the pulse-on spectrum.
- **Frequency Deviation.** Measures the signal frequency change from the center frequency over time during the pulse-on.

For information about Channel Power, OBW and EBW measurements, refer to the S/A (Spectrum Analysis) mode.

Channel Power	page 3-5
OBW	page 3-8
EBW	page 3-10

The frequency deviation measurement is based on the Frequency versus Time measurement in the Transient mode. Refer to page 3-87.


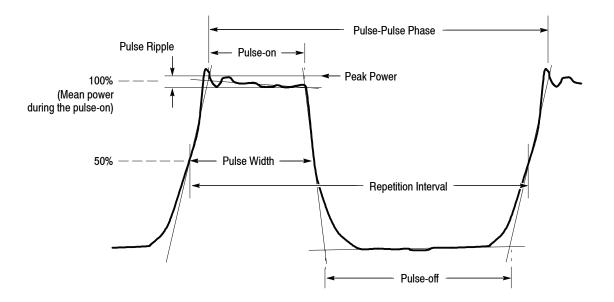


Figure 3-71: Definition of pulse characteristics

#### **Definition of Pulse-On/Off Time** The pulse-on and -off times are fundamental parameters for pulse measurements. These parameters are defined in Figure 3-72. The threshold is a level to detect a pulse, relative to the maximum peak in the acquired data. You can set it by using

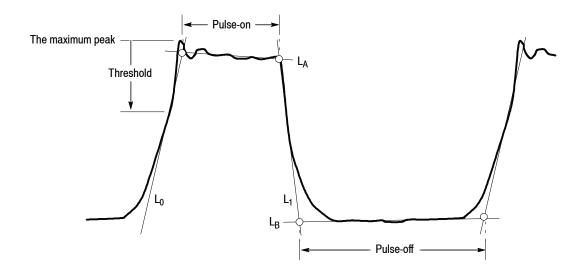

Detection Threshold in the Meas Setup menu (refer to page 3-100).

Figure 3-72: Definition of pulse-on/off

The L<sub>0</sub>, L<sub>1</sub>, L<sub>A</sub>, and L<sub>B</sub> lines represent the following:

- L<sub>0</sub>: Tangent line through the threshold point on the rising edge
- L<sub>1</sub>: Tangent line through the threshold point on the falling edge
- L<sub>A</sub>: Regression line calculated on the pulse top
- L<sub>B</sub>: Regression line calculated on the pulse bottom

These lines are displayed in red in the subview. The pulse-on and -off times are determined by the nodes of these lines as shown in Figure 3-72.

## Basic Measurement Procedure

Use this procedure to perform the pulse measurements:

- 1. Press the Time key on the front panel.
- 2. Press the side key Pulse Measurement  $\rightarrow$  Pulse Characteristics.
- **3.** Display the measurement waveform:
  - **a.** Press the **Frequency/Channel** key on the front panel to set the frequency.
  - **b.** Press the **Span** key on the front panel to set the span.
  - c. Press the Amplitude key on the front panel to set the amplitude.

For information about setting frequency, span, and amplitude, refer to *Setting Frequency and Span* on page 3-119 and *Setting Amplitude* on page 3-127.

- **4.** For the channel power, OBW, and EBW measurements only: You can select an FFT window for the measurement. Refer to *RBW/FFT Menu* on page 3-97.
- Press the Acquisition/Analysis key on the front panel to set the analysis range. Refer to page 3-30 for details on setting analysis range.
- 6. Press the View: **Define** key on the front panel to select the measurement item(s) for display. Refer to page 3-97 for the View Define menu.

By default, only the pulse-width measurement result is displayed on screen.

7. Press the **Meas Setup** key on the front panel to set the measurement parameters. Refer to page 3-100 for the Meas Setup menu.

- **8.** After you capture a measurement waveform, stop the data acquisition. If you are acquiring data in the continuous mode, press the **Run/Stop** key.
- **9.** Press the **Meas Setup** key and then **Analyze** side key to start the analysis of the acquired data.
- **10.** Press the View: **Define** key to select a pulse to observe from the measurement results.

## Hints for Taking Measurements.

- It takes longer to detect pulses when the analysis length in the Acquisition/ Analysis menu is set to a larger value.
- The results are not displayed in the Channel Power, OBW, or EBW measurement.

The number of data samples in Pulse-On must be 20 to 16,384 for FFT calculation. Change the span to display the results.

The results are not displayed in the Pulse Repetition Interval or the Duty Cycle measurement.

The observed pulse and the next one must be contiguous. If either pulse is not recognized due to noise and so on, the measurement results are not shown.

### Error messages (displayed on the bottom of the screen)

Message: Filter bandwidth is too wide or Channel bandwidth is too wide. Indicates that you set Filter Bandwidth or Channel Bandwidth inappropriately in the Meas Setup menu (refer to page 3-100).

In this case, change the setting of Filter Bandwidth or Channel Bandwidth until the message disappears.

Message: Too long pulse.
 Indicates that the number of data samples of a pulse width exceeds about 260,000.

In this case, change the span.

 Message: Too long repetition interval. Indicates that the number of data samples of a pulse repetition interval exceeds about 260,000.

In this case, change the span.

**RBW/FFT Menu** The RBW/FFT menu for the pulse measurements contains the following controls:

**FFT Window...** Selects an FFT window for the channel power, OBW, and EBW measurements from these two:

- Nyquist (default)
- Blackman-Harris 4B

**Rolloff Ratio...** Sets the roll-off ratio for the Nyquist window. Range: 0.0001 to 1 (default: 0.2)

For details on FFT and RBW, refer to page 3-163.

**View Define Menu** The View Define menu for the pulse measurements contains the following controls:

Show Views. Selects a single view or multi-view display.

- Single. Displays only the view selected by the View: Select key on the front panel.
- Multi. Default. Displays the overview, subview, and main view.

**Overview Content...** Selects a view displayed in the overview.

- Waveform (Power versus Time)
- Spectrogram

Refer to page 3-187 for information about setting views.

**Select Measurement...** Selects a measurement item to display in the subview from the items selected in Displayed Measurement.

**Select Pulse.** Selects a pulse to measure when you select A Single Pulse in View Results For... described below. 0 (zero) represents the most current pulse while larger negative numbers represent subsequent pulses.

**NOTE**. The main view shows the measurement results for up to 1000 pulses (the pulse number of -999 to 0).

**View Results For...** Selects how to display the measurement results in the subview (see Figure 3-73).

- A Single Pulse. Displays the measurement results for one pulse specified with Select Pulse.
- All Pulses. Displays the measurement results for all pulses in the analysis range. The horizontal axis represents the pulse number and the vertical axis represents the measurement results.

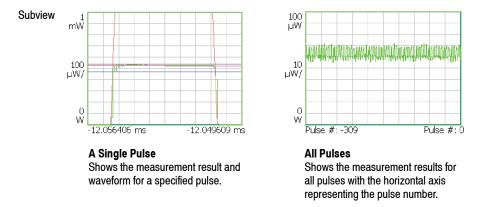


Figure 3-73: Setting for View Results For...

**Displayed Measurement...** Determines whether to show the measurement result in the main view for each measurement item. Select On to show the result. For the measurement items, refer to page 3–93.

**NOTE**. You must select at least one measurement item in Displayed Measurement. If you turn all items off, an error message "One or more measurement items must be selected" will appear on the bottom of the screen.

Display Time Range... Selects how to display a pulse in the subview.

- Adaptive. *Default*. Adjusts the horizontal scale for each pulse to fit the pulse width to the subview.
- Max. Adjusts the horizontal scale to fit the maximum pulse width in the analysis range to the subview.

Guidelines... Determines whether to display the guidelines in the subview.

- **On.** *Default*. Displays the guidelines in the subview.
- Off. Hides the guidelines in the subview.

For the guidelines, refer to Definition of Pulse-On/Off Time on page 3-94.

**Menu Off.** Turns off the side menu to expand the measurement results and waveform view area on screen. To recall the side menu, press the key on the front panel.

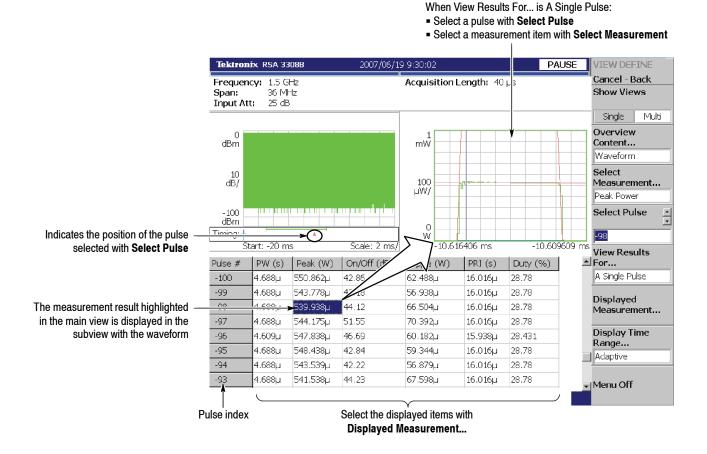


Figure 3-74: Setting the View Define menu

You can save the measurement table to a file (\*.csv) by **Save**  $\rightarrow$  **Save Table**. For file operations, refer to page 3-233.


Meas Setup Menu The Meas Setup menu for the pulse measurements contains the following controls: Analyze. Starts analysis on acquired data in the analysis range. Detection Threshold. Sets the threshold level to detect pulses in acquired data. This level is shown with a green arrow in the subview (see Figure 3-75). Range: -100 to 0 dBc (default: -3 dBc) Channel Bandwidth. For the Channel Power measurement only. Sets the measurement frequency range. Range: (Bin bandwidth) × 8 to set span. Refer to Trace Compression on page 3-183 for the bin bandwidth. **OBW Power Ratio.** For the OBW measurement only. Specifies the power ratio of the carrier and span regions for calculating OBW. Range: 80 to 99.9% (default: 99% as defined in T-53 or IS-95 standard). EBW Meas. Level. For the EBW measurement only. Specifies how far down from the peak level the bandwidth is measured. Range: -100 to -1 dB (default: -30 dB) **Counter Resolution.** For the Frequency Deviation measurement only. Sets the frequency counter resolution. It is the same as in the Carrier Frequency measurement in the S/A mode (refer to page 3-10). Range: 1 mHz to 1 MHz in a tenfold sequence (default: 1 Hz) P-P Phase Time Offset. For the Pulse-Pulse Phase measurement only. Sets the time offset for the measurement point. The default value is 0 (zero) second; the measurement point is at the beginning of the pulse-on time. **Measurement Filter.** Selects the measurement filter applied to time domain I/Q data for each pulse. None (no filter) Gaussian Filter Bandwidth. Sets the bandwidth for the Gaussian measurement filter. Range: Span/10 to full span (default: 3.6 MHz) **Filter Parameter.** Sets the  $\alpha/BT$  value for the Gaussian measurement filter. Range: 0.0001 to 1.0 (default: 0.35)


**Frequency Recovery.** For the Pulse-Pulse Phase and the Frequency Deviation measurements only. Selects the frequency correction method.

- 1st. Sets the correction value automatically, using the first pulse in the analysis range. The value is displayed in the Frequency Offset side key.
- User. Sets the correction value using the Frequency Offset side key described below.
- **Off.** *Default*. Disables the frequency correction.

**Frequency Offset.** Sets the frequency correction value when you select User in Frequency Recovery described above.

**Measurement Screen** The measurement waveform is displayed in the subview with the yellow lines (bold lines in Figure 3-75 on page 3-102) indicating the measurement points. Also, the pulse detection threshold (**Meas Setup**  $\rightarrow$  **Detection Threshold**) is shown with the green arrow. These guidelines can be turned on or off by View: **Define**  $\rightarrow$  **Guidelines**.

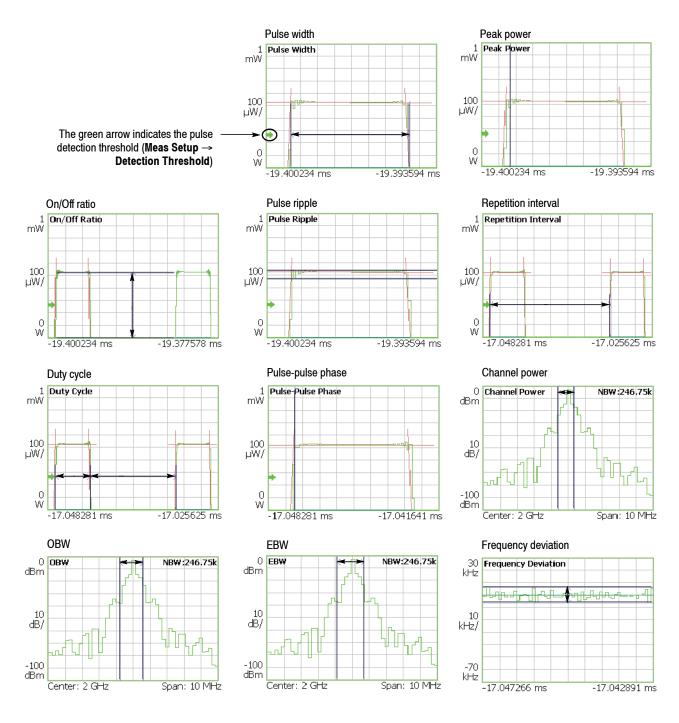

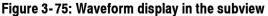

# Signal Source Analysis (Option 21 Only)

The signal source analysis function allows you to measure phase noise, jitter, and settling time in a signal source such as a PLL (Phase-Locked Loop) system.

Measurement Items	The following items are measured in the signal source analysis.
	(Select the items using the Measure key.)

Table 3-14: Signal source measurement items

Measure menu	Measurement items
Phase Noise Real-Time Phase Noise	<ul> <li>Carrier frequency</li> <li>Channel power</li> <li>Phase noise</li> <li>Random jitter</li> </ul>
	<ul> <li>Periodic jitter</li> </ul>
Spurious	Spurious
Real-Time Spurious	
Frequency versus Time	Frequency settling time

For the "Real-Time" concept, refer to: Real-Time Analysis page 1-3

2	10
Features of the Real-Time Mode	page 3-19

## **Basic Procedure** Use the following procedure to perform the signal source analysis.

**NOTE**. If you use a signal input (Input  $\rightarrow$  Signal Input Port...) other than RF, the measurements are not guaranteed.

- 1. Press the **Time** key on the front panel.
- 2. Select Signal Source Analysis using the side keys.
- **3.** Select a measurement item by pressing the side key. For example, press the **Phase Noise** side key to perform the phase noise measurement.
- 4. Adjust the amplitude and frequency to display the measurement waveform.

Refer to page 3-119 for setting frequency and span. Refer to page 3-127 for setting amplitude.

5. Set the analysis range in the overview. (This is not necessary in the Phase Noise and the Spurious measurements.)

Refer to page 3-30 for details of the analysis range settings.

If you use the zoom function, do these steps:

- Press the View: **Define** key.
- Press the **Overview Content...** side key and select **Zoom**.

Refer to *Zoom Function* on page 3-23 for setting the analysis area on the spectrogram with the zoom function.

6. Press the Meas Setup key and set the measurement parameters.

Refer to Meas Setup Menu in each measurement description.

7. Press the **Analyze** side key and carry out the analysis. (This is for the Real-Time Phase Noise and Real-Time Spurious measurements only.)

**Selecting the Phase Unit.** The phase unit is set to degrees by default. It is possible to select degrees or radians as the unit using the following procedure:

- 1. Press the System key on the front panel.
- 2. Press the Instrument Setup... side key.
- 3. Press the Angular Units side key to select degrees or radians.

# **Phase Noise Measurement** C/N (Carrier to Noise ratio) is measured at the offset frequency from the carrier. Jitter is also calculated from the phase noise data.

**Meas Setup Menu.** The following items are provided in the Meas Setup menu for Phase Noise measurement (see Figure 3-76).

**Carrier Threshold Level.** Sets the threshold for carrier detection. A signal with amplitude above the threshold is detected as a carrier.

Setting range: -100 to +30 dBm (default: -20 dBm)

Carrier Bandwidth. Sets the carrier bandwidth.

Setting range: Span/100 to Span/2 (default: Span/100)

C/N Sideband. Selects the sideband for measuring phase noise.

- Upper. *Default*. Measures the upper sideband.
- Lower. Measures the lower sideband.

**Minimum Offset Frequency...** Sets the minimum frequency in the phase noise measurement range as the offset from carrier frequency.

Setting value:

10 Hz (default), 100 Hz, 1 kHz, 10 kHz, 100 kHz, 1 MHz, or 10 MHz.

**Maximum Offset Frequency...** Sets the minimum frequency in the phase noise measurement range as the offset from carrier frequency.

Setting value:

100 Hz, 1 kHz, 10 kHz, 100 kHz, 1 MHz, 10 MHz, or 100 MHz (default)

**Rj Start Offset Frequency.** Sets the random jitter measurement start frequency as the offset from carrier frequency.

Setting range: 10 Hz to Rj Stop Offset Frequency set value (default: 10 Hz)

**Rj Stop Offset Frequency.** Sets the random jitter measurement stop frequency as the offset from carrier frequency.

Setting range: Rj Start Offset Frequency set value to 100 MHz (default: 100 MHz)

**NOTE**. The integrated phase noise and random jitter will be calculated if the frequency range from Rj Start Offset Frequency to Rj Stop Offset Frequency is within the range from Minimum Offset Frequency to Maximum Offset Frequency.

**Max Pj Threshold.** Sets the threshold level to determine periodic jitter by using the C/N value averaged by the Median filter (the yellow trace in the lower view) as the reference. A signal whose C/N value [dBc/Hz] is larger than the threshold is regarded as a periodic jitter. The maximum periodic jitter in the analysis range is shown in the measurement results (see Figure 3-77).

Setting range: 1 to 50 dB (default: 10 dB)

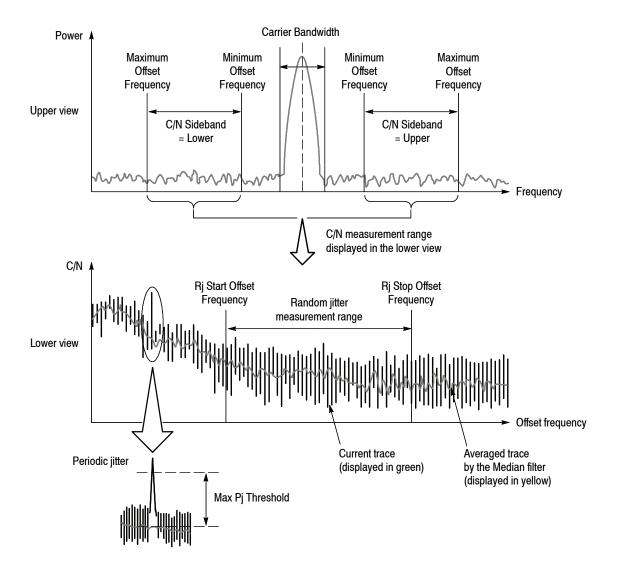


Figure 3-76: Phase noise measurement parameters

Figure 3-77 shows an example of the phase noise measurement. The *Integrated Phase Noise* in the measurement results indicates a phase noise value obtained by integrating C/N in the random jitter measurement range illustrated in Figure 3-76. The *Random Jitter* in the measurement results indicates a jitter value equivalent to the integrated phase noise.

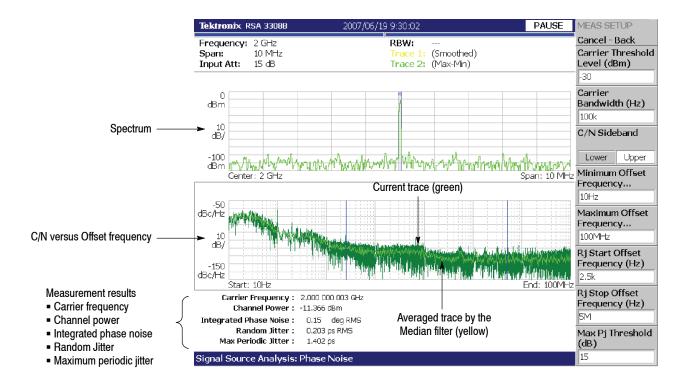


Figure 3-77: Phase noise measurement

**Selecting Trace 2.** In the C/N versus Offset frequency graph, Trace 1 (yellow) represents a trace averaged by the Median filter and Trace 2 (green) represents a Max-Min waveform by default. You can change Trace 2 using the Trace/Avg menu.

To change Trace 2, press the Trace/Avg key and set the following items:

Trace 2... Selects Trace 2.

- Max-Min. Default. Shows the maximum and minimum C/N values at alternate frequency points; its advantage is its resemblance to an analog display.
- **Reference.** When loaded as Trace 2, displays the waveform as the reference, which has been saved using **Save Trace 1**.
- Off. Displays no waveform.

**Load Trace 2.** Loads the waveform that was saved using **Save Trace 1** for Trace 2 as the reference.

Save Trace 1. Saves Trace 1 waveform data as a reference waveform.

Refer to page 3-233 for file operations.

## Hints for Taking Measurements.

- Widening the span expands carrier detection range, but loses accuracy. Narrowing the span increases measurement time.
- The phase noise measurement is performed within one of the four frequency bands listed in Table 3-15. A measurement outside the band is not performed.

Measurement band	Frequency range	Note
Baseband	0 Hz to 20 MHz	Set a value of [(center frequency) $\pm$ (span/2)] within the frequency range.
RF1	15 MHz to 3.5 GHz	Set a center frequency within the
RF2	3.5 GHz to 6.5 GHz	frequency range.

### Table 3-15: Phase noise measurement frequency band

5 GHz to 8 GHz

■ Tables 3-16 and 3-17 show the bin width and the number of waveform points for C/N versus Offset frequency measurement.

## Table 3-16: Bin width for each decade

RF3

Decade	Bin width
10 Hz to 100 Hz	0.24414 Hz
100 Hz to 1 kHz	2.4414 Hz
1 kHz to 10 kHz	2.4414 Hz
10 kHz to 100 kHz	19.531 Hz
100 kHz to 1 MHz	195.31 Hz
1 MHz to 10 MHz	1.5625 kHz
10 MHz to 100 MHz	12.5 kHz

(RF2 and 3 are RSA3308B only)

Point placement

		•	•	
	Averaged waveform	100 per decade	Linear on the logarithmic frequency axis.	
	Max-Min waveform	368 per decade	Linear on the linear frequency axis.	
	Error messages			
	<ul> <li>Message: No Carrier.</li> </ul>			
	In this case, lower Carrier Threshold Level in the Meas Setup menu.			
	■ Message: Out o	of Span.		
		crease Carrier Bandwid	th in the Meas Setup menu, or match equency.	
Spurious Measurement	pages 3-12). Howe	ever, for the spurious me	me way as in the S/A mode (refer to easurements in the signal source symmetrical spurious signals only.	
	•	The next items are ident ode. Refer to pages 3-12	ical to those in the spurious measure- 2.	
	<b>Carrier Threshold Level</b> ("Signal Threshold" in the S/A mode) <b>Ignore Region</b> <b>Spurious Threshold</b> <b>Excursion</b> <b>Scroll Table</b>			
	The following item	s are added for the sign	al source analysis.	
	Symmetrical Filte spurious signals on		the filter for extracting symmetrical	
	■ <b>On.</b> <i>Default</i> . D	visplays symmetrical spu	urious signals only.	
	<ul> <li>Off. Displays all spurious signals.</li> </ul>			
<b>Carrier Tracking</b> . Selects whether carrier tracking is executed or not. Carrier tracking ensures that the carrier frequency is always positioned c even when the signal drifts (it does not affect the waveform display).			quency is always positioned centrally,	
	•	arrier tracking is execut frequency always posit	ed. The measurement is conducted ioned centrally.	
■ Off. Carrier tracking is not executed. The measurement is conducted central frequency that was set using Frequency/Channel → Center				

Table 3-17: The number of waveform points

The number of points

Waveform

Figure 3-78 shows an example of the spurious measurement. Detected spurious signals are numbered sequentially from No. 1 in descending order of amplitude and spurious markers are displayed on the waveform. A maximum of 20 spurious signals are displayed. In addition, frequency difference (deltaF) and amplitude ratio (Ratio) in relation to the carrier are displayed in the measurement table at the bottom of the screen. You can save the table to a file (\*.csv) by **Save**  $\rightarrow$  **Save Table**. For file operations, refer to page 3-233.

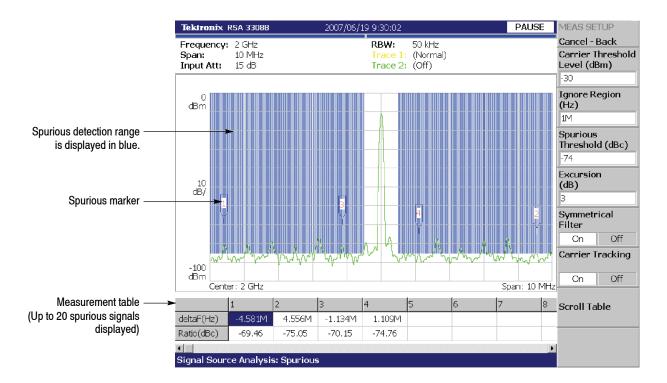


Figure 3-78: Spurious measurement

# Real-Time Phase Noise Measurement

Measures the phase noise in real time. The measurement results can be displayed in three dimensions like a spectrogram, called "Noisogram": the color axis shows the C/N (dBc/Hz), the horizontal axis shows the offset frequency (Hz) and the vertical axis shows the time (frame number). Jitter equivalent to the RMS value of the phase noise and the settling time of the phase noise and jitter are also calculated from the phase noise data.

**Meas Setup Menu.** The following items are provided in the Meas Setup menu for the Real-Time Phase Noise measurement.

Analyze. Performs measurement for acquisition data in the analysis range.

**NOTE**. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.

The next items are the same as for the phase noise measurement. Refer to page 3-105.

Carrier Threshold Level Carrier Bandwidth C/N Sideband Rj Start Offset Frequency Rj Stop Offset Frequency Max Pj Threshold

The following items are added for the real-time measurement (see Figure 3-79). They depend on the subview content (refer to *Selecting the Subview Content* on page 3-113).

**Rj Settling Threshold.** Enabled when Random jitter versus Time is displayed in the subview. Sets the threshold value for obtaining the random jitter settling time.

Setting range: 0 to 1 s (default: 0 s).

**C/N Offset Frequency.** Enabled when C/N versus Time is displayed in the subview. Sets the value of the frequency displaying the C/N versus Time in the subview as the offset from carrier frequency.

Setting range: The upper limit is Span/2. The lower limit is determined by the span and FFT Points for C/N (refer to *Setting the Number of FFT Points* on page 3-112) settings.

**C/N Settling Threshold.** Enabled when C/N versus Time is displayed in the subview. Sets the threshold value for obtaining the phase noise settling time.

Setting range: -200 to 0 dBc/Hz (default: 0 dBc/Hz).

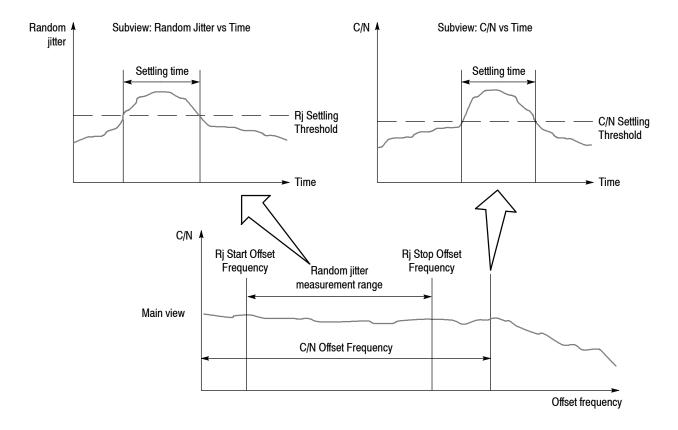


Figure 3-79: Real-time phase noise measurement parameters

Setting the Number of FFT Points. The real-time phase noise measurement uses 1024 FFT sample points by default for the upper and lower side bands to obtain a C/N versus offset frequency graph. You can change the number of FFT sample points in **RBW/FFT**  $\rightarrow$  FFT Points for C/N.

**FFT Points for C/N**. Sets the number of FFT sample points for the C/N versus offset frequency measurement.

Setting range: 64 to 65536 (2<sup>n</sup>, default: 1024).

The more points, the higher the resolution. The fewer points, the faster the measurement.

Figure 3-80 shows an example of the real-time phase noise measurement. The subview displays the noisogram (Color axis: C/N in dBc/Hz, Horizontal axis: Frequency in Hz, Vertical axis: Time in frame number) selected in View: **Define**  $\rightarrow$  **Subview Content...** The main view is the same as in the phase noise measurement (see Figure 3-77 on page 3-107), except that this main view is real time.

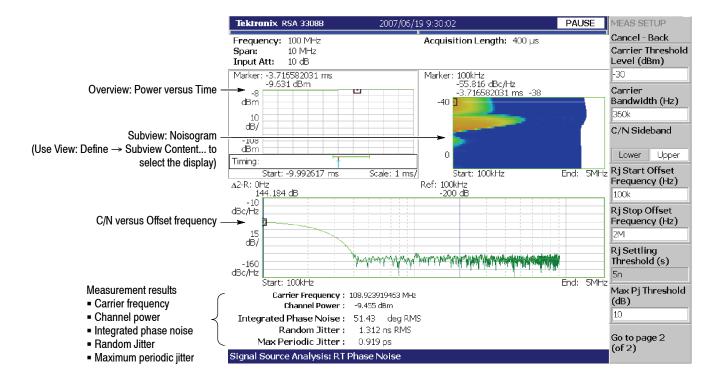



Figure 3-80: Real-time phase noise measurement

Selecting the Subview Content. The real-time phase noise measurement subview has displays that are specific to the signal source analysis. The subview is selected by pressing View: Define  $\rightarrow$  Subview Content...

**Subview Content...** Selects the information that is displayed in the subview. Refer to *Setting Views* on page 3–187 for scaling a view.

Spectrum

- Noisogram Color axis: C/N, Horizontal axis: Frequency, Vertical axis: Time (Frame No.) It is like a spectrogram but the color axis represents C/N in dBc/Hz.
- Random Jitter vs Time See Figure 3–81. The random jitter settling time is also indicated.
- Integrated Phase Noise vs Time
- C/N vs Time See Figure 3-81. The C/N settling time is also indicated.

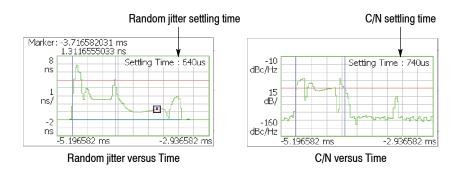

## Figure 3-81: Settling time displayed in the subview

## Hints for Taking Measurements.

If the subview displays nothing, check the following:

Subview: Random Jitter vs Time or Integrated Phase Noise vs Time

Make sure that the jitter measurement range set by Rj Start/Stop Offset Frequency in the Meas Setup menu is within the display range of the C/N versus Offset frequency in the main view.

■ Subview: C/N vs Time

Make sure that C/N Offset Frequency in the Meas Setup menu is within the display range of the C/N versus Offset frequency in the main view.

## Error messages

■ Message: No Carrier.

In this case, lower Carrier Threshold Level in the Meas Setup menu.

■ Message: Out of Span.

In this case, decrease Carrier Bandwidth in the Meas Setup menu, or match the carrier frequency with the center frequency.

# Real-Time Spurious Measurement

Spurious measurements are made in real time. The main view shows the spurious display the same as in the spurious measurement described on page 3-109, except that this spurious display is real time. You can observe the noisogram or C/N versus offset frequency in the subview.

**NOTE**. In the real-time spurious measurement, the number of FFT points is fixed to 1024.

**Meas Setup Menu.** The following items are provided in the Meas Setup menu for the Real-Time Spurious measurement.

Analyze. Performs measurement for acquisition data in the analysis range.

**NOTE**. When you change settings in the Meas Setup menu, press the Analyze side key to perform the measurement for the modified settings.

The following menu items are identical to those in the spurious measurement. Refer to page 3-109.

Carrier Threshold Level Ignore Region Spurious Threshold Excursion Symmetrical Filter Carrier Tracking Scroll Table

The following items are added for the real-time measurement. They depend on the subview content (refer to *Selecting the Subview Content* on page 3-116).

**Carrier Bandwidth**. This value is used to display C/N versus Offset frequency in the subview. Sets the frequency bandwidth for calculating channel power.

Setting range: Span/100 to Span/2 (default: Span/100).

C/N Sideband. This value is used to display C/N versus Offset frequency in the subview. Selects the sideband for measuring phase noise.

- Upper. *Default*. Measures upper sideband (half span).
- Lower. Measures lower sideband (half span).

Figure 3-82 shows an example of the real-time spurious measurement. The main view is the same as the spurious measurement (refer to page 3-109), except that this main view is real time.

## Time Analysis (Time Mode)

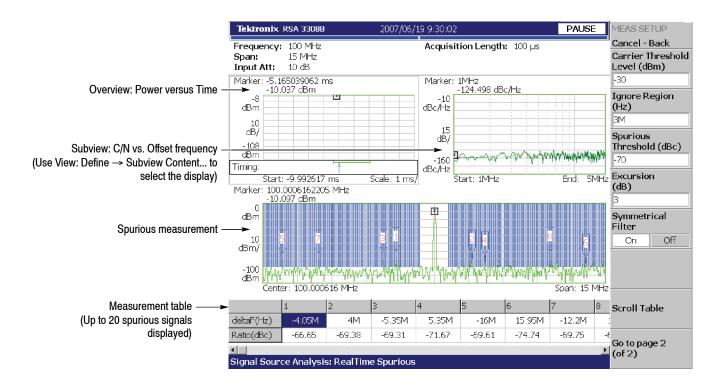


Figure 3-82: Real-time spurious measurement

The main view displays the spurious measurement results when the carrier level exceeds Carrier Threshold Level in the Meas Setup menu. The subview displays the measurement results in C/N versus Offset frequency and Noisogram when the carrier channel power exceeds Carrier Threshold Level.

Selecting the Subview Content. The real-time spurious measurement subview includes displays specific to the signal source analysis. The subview is selected by pressing View: Define  $\rightarrow$  Subview Content...

Subview Content... Selects the information displayed in the subview.

- Spectrum
- Noisogram

Color axis: C/N, Horizontal axis: Frequency, Vertical axis: Time (Frame No.) This is like a spectrogram except that the color axis represents C/N in dBc/Hz.

 C/N vs Offset Freq This is the same as the main view in the real-time phase noise measurement (see Figure 3-80 on page 3-113).

# Frequency versus Time Measurement

Measures changes in frequency. A running average of the frequency versus time waveform is displayed in the main view and the frequency settling time is obtained. However, the measurement signal must meet the following conditions:

- The frequency shift does not exceed the threshold at the start and end of the analysis range.
- Frequency hopping occurs only once within the analysis range.

**NOTE**. The frequency versus time measurement requires at least 2048 samples. If the analysis length is set to less than 2048, an error message "Analysis Length is too short" is displayed. Also, 1024 samples stable in frequency are required at both the start and end of the analysis range to calculate the average frequency.

**Meas Setup Menu.** The following items are provided in the Meas Setup menu for the frequency versus time measurement.

**Freq Settling Threshold**. Sets the threshold for judging the frequency settling time (see Figure 3-83). When frequency deviation exceeds this value, hopping is considered to have occurred.

Setting range: 10 Hz to Span (default: 10 Hz).

**Smoothing Factor**. Sets the number of points from which the running average is calculated.

Setting range: 1 to (analysis range)/2 or 9999 points (default: 1).

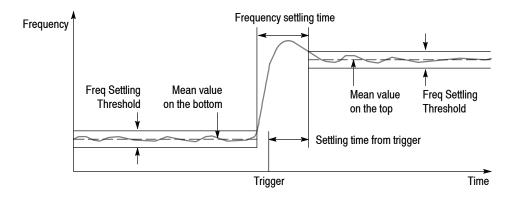


Figure 3-83: Frequency versus Time measurement parameters

**NOTE**. In the frequency versus time measurement, the analysis range is approximately 500 frames (512,000 points) maximum. For example, it is about 40 ms for the span of 10 MHz.

Figure 3-84 shows an example of the frequency versus time measurement. The subview displays only spectrum. The frequency settling time is indicated on the bottom of the screen. The value from trigger is also indicated when the trigger occurs during or before the frequency settling time in the analysis range.


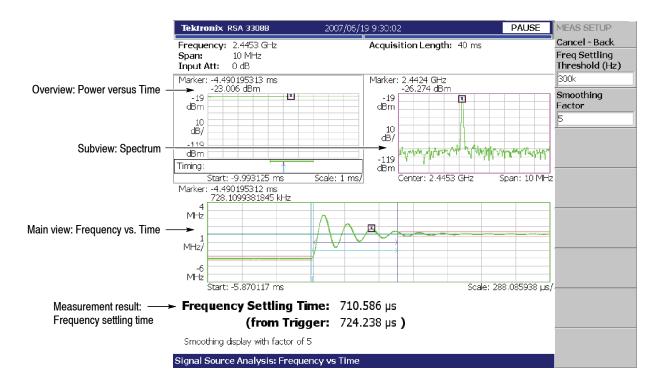


Figure 3-84: Frequency versus Time measurement

# **Setting Frequency and Span**

This section describes frequency and span, which are fundamental settings for observing the spectrum. These items are set using the general purpose knob and the numeric input keypad. You can also use the marker peak search function to set frequency.

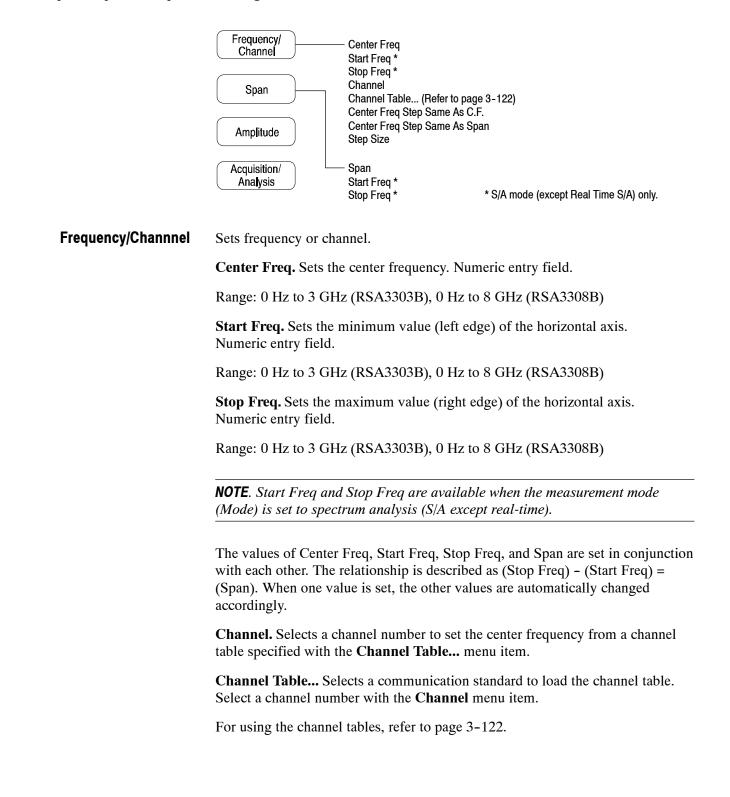

# **Frequency and Span Setting Menu**

**Center Freq Step Same As C.F.** Sets the center frequency step size equal to the center frequency.

**Center Freq Step Same As Span.** Sets the center frequency step size equal to the span.

Step Size. Sets the step size (amount per press by which the up or down key  $(\blacktriangle \nabla)$  changes the setting value) for setting frequency.

**Span** Controls span. Span and frequency settings are shown in Figure 3–85.

**Span.** Sets span. The available range depends on the measurement frequency band and mode as shown in Table 3-18 on page 3-124.

Turn the general purpose knob to set span in the specified sequence. In the S/A mode (except real-time), you can set an arbitrary span within the limit using the numeric keypad.

Start Freq. Same as *Start Freq* on page 3-120.

Stop Freq. Same as *Stop Freq* on page 3-120.

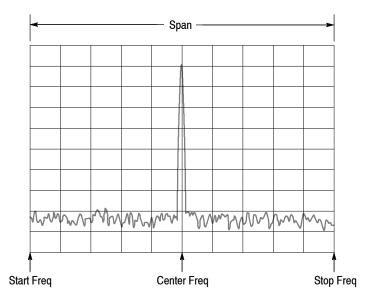


Figure 3-85: Setting frequency and span

# **Using the Channel Table**

The channel table contains channel numbers and the corresponding frequencies for a communication system. When you use W-CDMA standard signals, for example, you can set a center frequency by selecting a channel number from the W-CDMA table.

- 1. Press the Frequency/Channel key on the front panel.
- 2. Press the Channel Table... side key and select one of these items:

None selects no channel table.

You can select one of the following communication standards:

CDMA2000 EU PAMR400-FL	CDMA2000 NMT450 25k-FL
CDMA2000 EU PAMR400-RL	CDMA2000 NMT450 25k-RL
CDMA2000 EU PAMR800-FL	CDMA2000 SMR800-FL
CDMA2000 EU PAMR800-RL	CDMA2000 SMR800-RL
CDMA2000 GSM BAND 1-FL	CDMA2000 TACS BAND-FL
CDMA2000 GSM BAND 1-RL	CDMA2000 TACS BAND-RL
CDMA2000 GSM BAND 2-FL	DCS1800-DL
CDMA2000 GSM BAND 2-RL	DCS1800-UL
CDMA2000 IMT2000-FL	GSM850-DL
CDMA2000 IMT2000-RL	GSM850-UL
CDMA2000 JTACS BAND-FL	GSM900-DL
CDMA2000 JTACS BAND-RL	GSM900-UL
CDMA2000 KOREA PCS-FL	IEEE802.11a
CDMA2000 KOREA PCS-RL	IEEE802.11b/g
CDMA2000 N.A. 700MHz Cellular-FL	NMT450-DL
CDMA2000 N.A. 700MHz Cellular-RL	NMT450-UL
CDMA2000 N.A. Cellular-FL	PCS1900-DL
CDMA2000 N.A. Cellular-RL	PCS1900-UL
CDMA2000 N.A. PCS-FL	TD-SCDMA
CDMA2000 N.A. PCS-RL	W-CDMA-DL
CDMA2000 NMT450 20k-FL	W-CDMA-UL
CDMA2000 NMT450 20k-RL	

\* FL: Forward link; RL: Reverse link; UL: Uplink; DL: Downlink

3. Press the Channel side key and select a channel number.

For example, when selecting channel 10551 in the W-CDMA downlink table, the center frequency is automatically set to 2.1102 GHz.

# **Using the Marker and Peak Search**


You can use the search function to position the marker on the spectrum peak, and then set the center frequency to the frequency at the marker, as shown in Figure 3-86.

**In S/A Mode.** You can set the peak spectrum to the center frequency using the marker search functions when the measurement mode is S/A (spectrum analysis).

- **1.** Display the spectrum on the screen.
- 2. Press the **Peak** key on the front panel. The maximum peak spectrum is detected and the marker is moved to that point.

Use the arrow keys  $(\blacktriangleleft \blacktriangleright \blacktriangle)$  to move the marker to another peak.

3. Press the Marker → key and then the Center Freq = Marker Freq side key. The center frequency is set to the frequency at the marker position.

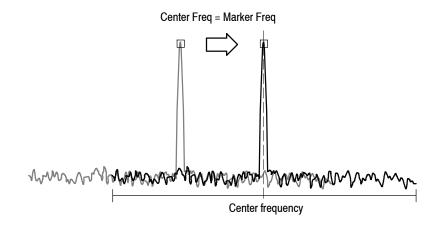

# Figure 3-86: Setting the center frequency using Marker →

The frequency set in this procedure may not take effect with all span settings. Refer to *Frequency Setting Range* on page 3-124.

In Demod and Time Modes. The Marker  $\Rightarrow$  key is used to set the analysis range when the measurement mode is Demod or Time (modulation and time analyses). Refer to *Setting the Analysis Range* on page 3-30 for details.

In the overview of spectrogram or the subview of spectrum, you can use the **Center Freq = Marker Freq** side key as in the S/A mode.

# **Frequency Setting Range**

Frequency bands are defined as shown in Table 3–18, based on the analyzer hardware architecture. The frequency band is switched automatically by a frequency setting. Span setting range depends on the frequency band and the measurement mode (S/A, Demod, or Time).

Measurement mode	Frequency band		Frequency range	Span setting range
S/A (except real-time)	Baseba	and	DC to 20 MHz	50 Hz to 20 MHz (1-1.2-1.5-2-2.5-3-4-5-6-8 sequence)
	RF	(RSA3303B)	15 MHz to 3 GHz	
	RF1	(RSA3308B)	15 MHz to 3.5 GHz	
	RF2	(RSA3308B)	3.5 to 6.5 GHz	50 Hz to 3 GHz (1-1.2-1.5-2-2.5-3-4-5-6-8 sequence)
	RF3	(RSA3308B)	5 to 8 GHz	
Real Time S/A	Baseba	and	DC to 20 MHz	100 Hz to 20 MHz (1-2-5 sequence)
Demod, Time	RF	(RSA3303B)	15 MHz to 3 GHz	
	RF1	(RSA3308B)	15 MHz to 3.5 GHz	
	RF2	(RSA3308B)	3.5 to 6.5 GHz	- 100 Hz to 10 MHz (1-2-5 sequence) and 15 MHz
	RF3	(RSA3308B)	5 to 8 GHz	

Table 3-18: Frequency and span setting range

In Demod and Time modes, the frequency and span settings must meet the following conditions:


(Center frequency) + (Span)/2

- $\leq$  Upper limit of the frequency setting range (RF mode)
- $\leq$  40 MHz (Baseband)

(Center frequency) - (Span)/2

- $\geq$  Lower limit of the frequency setting range (RF mode)
- $\geq$  0 Hz (Baseband)

When you enter an out-of-range value, the value is limited based on these conditions. In spectrum analyzer mode, the out-of-range values are allowed, but part of the trace may not be displayed because the analyzer cannot acquire that part of the waveform (see Figure 3–87).

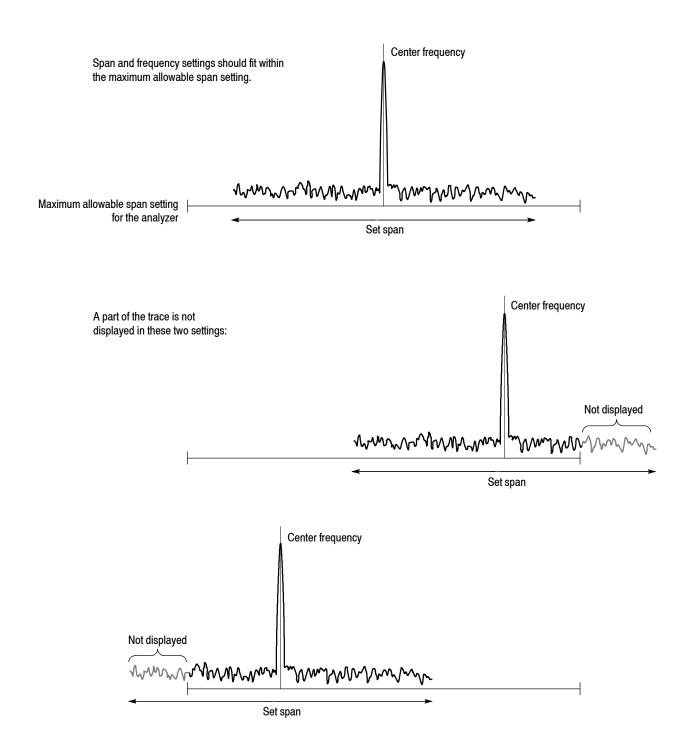



# **Vector Span**

The input signals are scanned in frame units (one frame = 1024 points). There are two frames: a physical frame to store scanned data and a logical frame to store display data, as illustrated in Figure 3–88.

In the baseband, one logical frame is acquired at one scan regardless of the span setting. You can acquire one logical frame with one scan of a span setting below 15 MHz in the RF band. In larger spans, one logical frame is built up by acquiring data with multiple physical frames. For example, when the span is 30 MHz, one logical frame is composed from two (30 MHz/15 MHz) scans.

When the span is below 15 MHz in the baseband and RF band, one physical frame corresponds to one logical frame; this is referred to as Vector mode, and the span in vector mode is called the Vector span. In the other cases, one logical frame is composed of multiple physical frames, and this is called Scalar mode.

#### Span $\leq$ 15 MHz: Vector mode

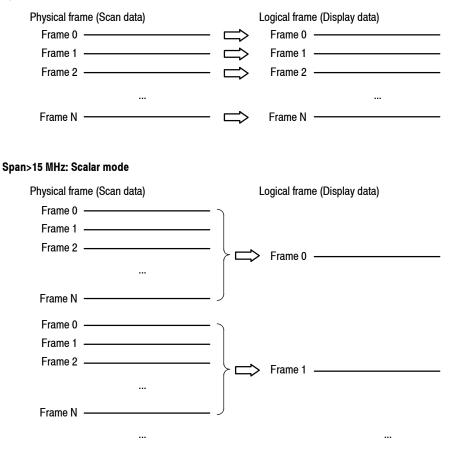


Figure 3-88: Vector mode and scalar mode

# **Setting Amplitude**

This section describes the fundamental amplitude settings for observing a spectrum. You can add correction to the waveform display by considering frequency characteristics of external devices such as antennas or preamplifiers.

# **Amplitude Menu**

Press the **Amplitude** key to set amplitude. Figure 3-89 shows the Amplitude menu structure.

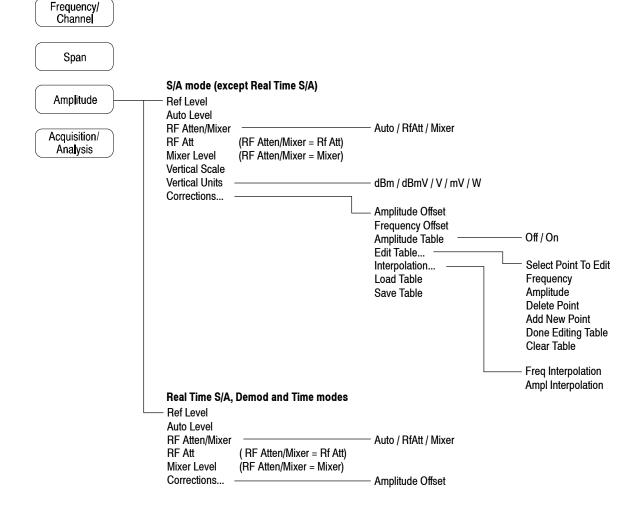


Figure 3-89: Amplitude menu structure

**Amplitude** Sets the scale of the amplitude of the displayed waveform.

**Ref Level.** Sets the maximum (top) edge of the vertical axis. The available range depends on the measurement frequency band as shown in Table 3-19.

Table 3-19: Reference level setting range

Frequency band		Setting range
Basebar	nd	-30 to +20 dBm (2 dB step)
RF RF1	(RSA3303B) (RSA3308B)	-51 to +30 dBm (1 dB step)
RF2, 3	(RSA3308B)	-50 to +30 dBm (1 dB step)
IQ (Option 03 only)		-10 to +20 dBm (10 dB step)

**Auto Level.** Automatically adjusts reference level for the best system performance based on power measurement within the set span.

**NOTE**. An input signal within 10 MHz of the center frequency may cause Auto Level to choose the wrong reference level, even if the span setting does not make the signal appear on sereen. To prevent nonlinear operation, the reference level may need to be adjusted manually.

Auto Level will force the analyzer to acquire new data, which will overwrite the existing data. To scale the existing displayed waveform, do not use Auto Level. Instead, use the View: Scale menu.

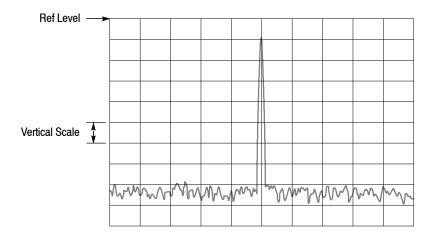


Figure 3-90: Setting the amplitude

**RF Atten/Mixer.** Input signals are attenuated through the attenuator and converted to IF signals through the mixer in the down-converter (refer to *Architecture* on page 1–7). The attenuation level and the mixer level are normally set automatically. Select **RF Att** or **Mixer** to manually set either parameter, if necessary. RF Atten/Mixer is unavailable in the baseband (DC to 20 MHz).

- Auto. Sets the mixer level and the RF attenuation level automatically.
- **RF** Att. Sets the RF attenuation level with **RF** Att described below.
- Mixer. Sets the mixer level using the Mixer Level setting.

**NOTE**. *RF* Atten/Mixer is set to Auto by default. In this setting, the mixer level is fixed to -25 dBm.

**RF** Att. Changes the RF attenuation level when **RF** Att is selected in **RF** Atten/Mixer above.

#### Table 3-20: RF attenuation level setting

Frequency band		Attenuation level (dBm)
RF RF1	(RSA3303B) (RSA3308B)	0 to 50 (2 dB step)
RF2, 3	(RSA3308B)	0, 10, 20, 30, 40, or 50

Reducing the attenuation may increase the signal level compared to the noise floor, because the noise is generated after the RF attenuator.

**Mixer Level.** Selects the input level of the primary mixer when **Mixer** is selected in **RF Atten/Mixer**.

#### Table 3-21: Mixer level setting

Frequency band		Mixer level (dBm) <sup>1</sup>
RF RF1	(RSA3303B) (RSA3308B)	-5, -10, -15, -20, or -25
RF2, 3	(RSA3308B)	-5, -15, or -25

- 25 dBm fixed when RF Atten/Mixer is set to Auto.

Select the level according to measurement type. The default value is -25 dBm. Use the default value in most cases. This level can be increased up to -5 dBm when a high dynamic range is required for measurements such as the ACPR (Adjacent Channel Leakage Power Ratio) measurement.

1

**NOTE**. As the mixer level increases, the distortion also increases.

**Vertical Scale.** *S*/*A mode (except Real Time) only.* Sets the vertical scale (per division). The setting range depends on Vertical Units described below, as shown in Table 3–22.

 Table 3-22: Vertical scale setting range

Vertical units	Scale setting <sup>1</sup>
dBm, dBμV	1 to 10
V	223.6 n to 22 m
μV	223.6 m to 22 k
W	100 p to 100 μ

<sup>1</sup> 1-2-5 sequence with the general purpose knob. Arbitrary value with the numeric keypad.

Vertical Units. Selects the unit of amplitude scale: dBm, dBµV, V, µV, or W.

**Corrections...** Sets amplitude correction. Refer to page 3-132 for details of amplitude correction.

**Basic Setting Procedure** To set amplitude, follow these steps: 1. Press the Amplitude key on the front panel. 2. Set the reference level using the **Ref Level** side key. 3. To best display the waveform automatically, press the Auto Level side key. 4. To set the attenuation level or the mixer level manually, use the RF Atten/Mixer side key to select RF Att or Mixer. When you select **RF** Att: Select the attenuation level with the **RF** Att side key. When you select **Mixer**: Select the first mixer input level using the **Mixer Level** side key. 5. S/A mode (except Real Time) only. Use the Vertical Scale side key to set the vertical scale (per division). Use the Vertical Units side key to select the unit. 6. To apply amplitude correction, press the **Corrections...** side key to set the parameters. Refer to Amplitude Correction on page 3-132 for details.

### **Over-Voltage Input**

Set the reference level (Ref Level) according to the input signal level. The default setting is 0 dBm. If the signal level gets too high or the reference level is set too low, over-voltage input may occur. If an over-voltage input occurs, the status indicator "Overrange – increase RefLev or Atten" is displayed in the red box (see Figure 3-91).



**CAUTION.** If a signal exceeding +30 dBm (1 W) is applied, it can damage the analyzer. Be sure to limit input signals to +30 dBm or below.

When the input signal level is too high, overrange is indicated in the red box.

	*					
Tektronix R	Overrange - increase RefLev o	r Atten 1 AM	INPUT: CAL	FREE RUN	AMPLITUDE	
Frequency:	50 MHz	RBW: 50 kHz			Cancel - Back	
Span:		Trace 1: (Normal	)		Ref Level	-
Input Att:	0 dB	Trace 2: (Off)			(dBm)	-
					-30	

### Figure 3-91: Overrange indicator

**NOTE**. If "Overrange – increase RefLev or Atten" is displayed, it indicates that the A/D converter in the subsequent part of the downconverter inside this instrument is overloaded. In this case, data display is distorted and the measurement is not accurate.

If a signal larger than the reference level set value by 20 dB or more is applied continuously, however, the limiter of the IF amplifier in the downconverter is automatically activated to prevent large level signals from passing through the A/D converter. Even if a signal exceeding the reference level is applied, therefore, "Overrange – increase RefLev or Atten" may not be displayed. Auto Leveling the signal can usually eliminate all overload conditions.

Take sufficient care about the input signal level.

The overflow indicator is updated at every physical frame acquisition. When a high level signal is applied under a setting that uses multiple physical frames for one scan, "Overrange – increase RefLev or Atten" may be displayed for a moment and then turned off immediately. When a single high-level signal is applied under a setting using one physical frame in one scan, the same phenomenon may be observed.

## **Amplitude Correction**

If an external device such as an antenna or preamplifier is connected to the analyzer, you can amplitude correct the waveform by considering the amplitude characteristics of the external device.

**NOTE**. The amplitude correction function is fully controlled only in the S/A mode (except real-time mode). In the other modes (Real Time S/A, Demod, and Time), amplitude offset control is available. Refer to page 3-139 for setting the offset.

Figure 3-92 shows the concept of amplitude correction. In this example, a signal with -80 dBm is sent to a preamplifier with a gain of +20 dB around 1 GHz. In the normal display, without any amplitude correction, the signal peak becomes -80+20 = -60 dBm. If amplitude correction is performed, the peak becomes -60-20 = -80 dBm to obtain the original peak value of the input signal.

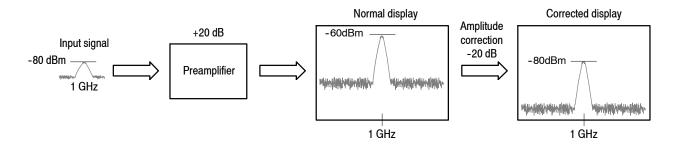


Figure 3-92: The amplitude correction concept

### Amplitude Correction Menu

Press **Amplitude**  $\rightarrow$  **Corrections...** to set the correction parameters. The correction menu has the following controls:

**Amplitude Offset.** Sets the amplitude offset. The amplitude of an entire waveform decreases by the offset value.

**NOTE**. The following Corrections menu items are available in the S/A (spectrum analysis) mode except Real Time S/A.

**Frequency Offset.** Sets the frequency offset. The effective correction range of the amplitude correction table shifts by the offset value.

**Amplitude Table.** Enables or disables amplitude correction. Select On to enable the correction.

**Edit Table...** Creates a correction table. Input pairs of frequency and amplitude correction value.

- Select Point To Edit. Selects a row to be edited.
- **Frequency.** Enters the frequency of a correction point.
- Amplitude. Enters an amplitude correction value for the specified frequency.
- **Delete Point.** Deletes the selected row.
- Add New Point. Adds a row with the values in the previous row copied as initial values.
- **Done Editing Table.** Confirms the input and adds a new row.
- Clear Table. Deletes the correction data from the memory.

**Interpolation...** Selects the horizontal and the vertical scale for interpolating correction data.

- Freq Interpolation. Selects the horizontal scale for interpolating correction data: Linear or Logarithmic.
- Ampl Interpolation. Selects the vertical scale for interpolating correction data: Linear or dB.

Load Table. Reads an amplitude correction table from a file.

Save Table. Writes a created correction table to a file.

Amplitude Correction FileBefore performing amplitude correction, you must describe the frequency<br/>characteristics of the external device in an amplitude correction file (\*.cor).

You can use word processing software to create the file on a computer and then use a disk or network to load the file into the analyzer. Alternatively, you can create correction data on the analyzer screen.

The rules for creating an amplitude correction file are given on page 3-135 and the procedures are given on page 3-136.

**File Format.** A correction file contains amplitude correction values at each frequency, in the following format:

<Frequency 1> = <Amplitude correction value 1> <Frequency 2> = <Amplitude correction value 2> <Frequency 3> = <Amplitude correction value 3>

*Example:* The following amplitude correction file has correction data for three points:

Correction data:	Correction file description:
10 dB at 10 MHz	10M=10
5 dB at 100 MHz	100M=5
0 dB at 1 GHz	1G=0

In this example, only the data from 10MHz to 1GHz is corrected (see Figure 3-93). The correction value in the display range is obtained using linear interpolation between input points. A waveform obtained by subtracting the correction value from the input waveform is displayed.

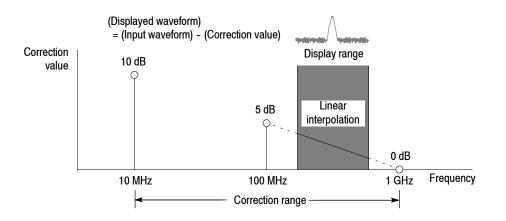


Figure 3-93: Amplitude correction example

Linear or Logarithmic (Log/dB) interpolation can be selected using the following menus:

### ■ Amplitude → Corrections...→ Interpolation...→ Freq Interpolation

**Lin.** Linear interpolation of correction values is performed on the frequency axis using a linear scale.

**Log.** Linear interpolation of correction values is performed on the frequency axis using a logarithmic scale.

### ■ Amplitude → Corrections...→ Interpolation...→ Ampl Interpolation

**Lin.** Linear interpolation of correction values is performed on the amplitude axis using a linear scale.

**dB.** Linear interpolation of correction values is performed on the amplitude axis using a logarithmic scale.

#### **Rules for Creating an Amplitude Correction File.**

- Create the file as a text file and save it with the extension ".cor".
- The maximum number of input lines is 3000.
- The order of correction data input is unimportant because the file is sorted when it loads. However, it may be easier to follow if you enter data in the order of ascending frequency.
- Numeric values are described without the unit of frequency or amplitude (Hz, dB, W, etc.). For example, 5 MHz frequency is expressed as 5M.
- Frequency can be expressed as a floating point number or with the SI unit (k, M, or G). For example, the following lines show three different ways to express the same value:

1000, 1E+3, 1k 1230000, 1.23E+6, 1.23M 1000000000, 1.0E+9, 1.0G

- Amplitude is expressed as a decimal number or integer (such as 1.23 or 10).
- Do not use a space in a numeric value; however, you can space before and after "=".

*Correct:* 10M = 10 (A space is inserted before and after "=".) *Incorrect:* 10 M=10 (A space is inserted between "10" and "M".) **Creating an Amplitude Correction File on a Computer.** Use word processing software on a computer to create a text file with the extension ".cor".

Be sure to use the correct file format (refer to *File Format* on page 3-134) and follow the *Rules for Creating an Amplitude Correction File* on page 3-135.

For information about performing amplitude correction, refer to page 3-140.

**Creating Correction Data on the Analyzer Screen.** The procedure for entering new correction data or modifying existing data on the screen is as follows. Refer also to *Rules for Creating an Amplitude Correction File* on page 3–135.

- 1. Press the Amplitude key on the front panel.
- 2. Press the Corrections... side key.
- 3. To edit an existing file: Refer to page 3-233 for file handling.
  - a. Press the Load Table side key to load the file.
  - **b.** Press the **Edit Table...** side key.
- 4. To enter new data: See Figure 3-94.
  - **a.** Press the **Frequency** side key and enter the frequency of a correction point.
  - **b.** Press the **Amplitude** side key and enter the amplitude correction value at the correction point.
  - c. Press the Add New Point side key.

A new row is added, containing the same frequency and amplitude correction values as the previous row. Modify the values appropriately.

**NOTE**. If you set the same frequency as the previous row, the amplitude correction value in the previous row will be overwritten the next time you press the Add New Point side key.

It is not necessary to enter correction data in order by frequency. The rows are automatically reordered in ascending order of frequency.

- **d.** Repeat step c to enter the frequency and amplitude correction values for all points.
- e. Press the Done Editing Table side key.

The input is accepted and a new row is added.

- 5. To add correction data:
  - **a.** Press the **Select Point To Edit** side key and turn the general purpose knob to move the cursor to the last row (empty) in the table.
  - **b.** Perform step 4 to enter the frequency and amplitude correction values.

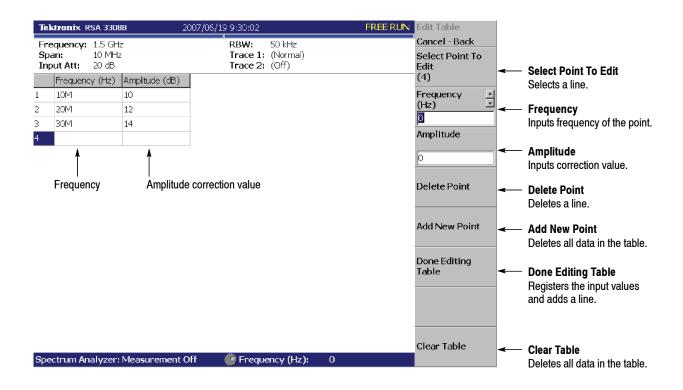


Figure 3-94: Amplitude correction data input

- 6. *To modify input data:* Press the **Select Point To Edit** side key and turn the general purpose knob to move the cursor to the row to be modified. Use the following side keys as needed:
  - To modify a frequency, use the **Frequency** side key.
  - To modify an amplitude, use the **Amplitude** side key.
  - To delete a row, press the **Delete Point** side key.
  - To delete all data in the table, press the **Clear Table** side key.
- 7. Repeat steps 5 through 6 as needed.

- 8. When you have finished the data input, save the file as follows:
  - a. Press the Amplitude key on the front panel.
  - **b.** Press the **Corrections...** side key.
  - **c.** Press the **Save Table** side key to specify the save file. Refer to page 3-233 for file handling.

For information about performing amplitude correction, refer to page 3-140.

### Setting the Offset

The amplitude correction function includes amplitude and frequency offset.

**NOTE**. The amplitude offset is available in all the measurement modes of S/A, Demod, and Time.

**Amplitude Offset.** Subtracts the specified offset from the original amplitude (see Figure 3-95). The offset value is set using the menu item **Amplitude**  $\rightarrow$  **Corrections...** $\rightarrow$  **Amplitude Offset**.

The operation varies slightly between the measurement modes as follows:

- In the S/A mode (except real-time mode): The amplitude offset is effective whenever the amplitude correction is turned on (Amplitude → Corrections...→ Amplitude Table → On) even if the table is not set.
- In the other modes (Real Time S/A, Demod, and Time): The amplitude offset is always effective. The default value is zero. If you set the offset to a non-zero value, a waveform shifts vertically as shown in Figure 3-95.


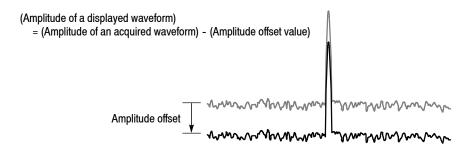


Figure 3-95: Amplitude offset

**Frequency Offset.** Shifts the correction range by the specified offset with the same correction table (see Figure 3-96). The offset value is set using the **Amplitude**  $\rightarrow$  **Corrections...** $\rightarrow$  **Frequency Offset** menu item.



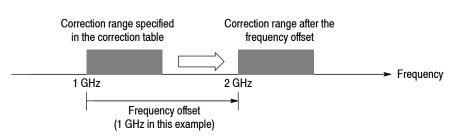


### Performing Amplitude Correction

Load an amplitude correction file and acquire input signals:

- 1. Press the Amplitude key on the front panel.
- 2. Press the Corrections... side key.
- **3.** Do either of these two steps:
  - Press the Edit Table... side key to create an amplitude correction file.
     Refer to page 3-134 for creating the correction file.
  - Press the Load Table side key to specify an amplitude correction file. Refer to page 3-233 for file operation.
- 4. Press the Amplitude Offset side key to set the amplitude offset, if necessary.
- 5. Press the Frequency Offset side key to set the frequency offset, if necessary.
- 6. Press the Interpolation... side key to select scaling for the interpolation:
  - **a.** Press the **Freq Interpolation** side key to select the scale for frequency interpolation: Lin (linear) or Log (logarithmic).
  - **b.** Press the **Ampl Interpolation** side key to select the scale for amplitude interpolation: Lin (linear) or dB (logarithmic).
- 7. Press the **Amplitude Table** side key to select **On**. The amplitude correction is now applied to the waveform.

The amplitude correction is applied to the acquired data and the corrected waveform is shown. When the amplitude correction is turned on, "Correction" is displayed in the setup display area at the upper right of the screen, as shown in Figure 3-97.

Tektronix RSA 33088	2007/06/19 9:30:02	FREE RUN	Corrections
Frequency: 800 MHz Span: 10 MHz Input Att: 20 dB	RBW: 50 kHz Trace 1: (Normal) Trace 2: (Off)	Confection	Sancel - Back Amplitude Offset (dB)
0 dBm	Indicates the amplitude correction is or	/	0 Frequency Offset (Hz)

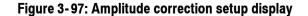

**Erasing Correction Data** The analyzer auto

The analyzer automatically saves the correction data in use when it is powered down. The data will be erased when you press these keys:

- Amplitude → Corrections...→ Edit Table...→ Clear Table
- Preset

# **Setting Acquisition and Analysis Parameters**

You can adjust the length of various time parameters and the relationship between them in the Real Time S/A, Demod, and Time modes using the Acquisition/Analysis menu.

**NOTE**. The Acquisition/Analysis menu is not available in the S/A mode (except Real Time S/A).

### **Acquisition/Analysis Menu**

Figure 3-98 shows the Acquisition/Analysis menu structure.

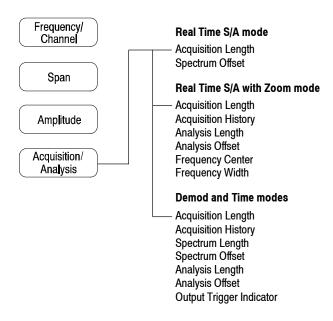


Figure 3-98: Acquisition/Analysis menu structure

### Parameters in Demod and Time Modes

The Acquisition/Analysis menu in Demod (modulation analysis) and Time (time analysis) modes contains the following items (see Figure 3-99):

**Acquisition Length.** Sets the time length to acquire one block (= M frames). The acquisition length is calculated using this equation:

(One block acquisition length) =  $M \times$  (One frame acquisition length)

where M is the number of frames in a block.

One frame acquisition length is determined by span internally. For details, refer to *Appendix A: Specifications*.

Acquisition History. Specifies the number of the block to display and analyze. The latest block is number zero. Older blocks have larger negative numbers. In most cases the analyzer retains previous acquisitions. You can view previous acquisitions by selecting the acquisition history by block number.

**Spectrum Length.** Displays the time length for FFT processing of the spectrum displayed in the subview. It is equivalent to one frame acquisition length determined by Span, RBW, and FFT points setting.


**Spectrum Offset.** Sets the beginning of Spectrum Length with respect to the trigger output point.

**Analysis Length.** Sets the time length of the analysis range in the block specified by Acquisition History.

Analysis Offset. Sets the beginning of Analysis Length with respect to the trigger output point.

**Output Trigger Indicator.** Determines whether to display the output trigger indicator ("O") on the overview. Refer to *Trigger Point Indicator* on page 3-161.

For procedures to set these parameters, refer to *Setting the Analysis Range* on page 3-30 and *Setting FFT Processing Range for the Subview* on page 3-34.

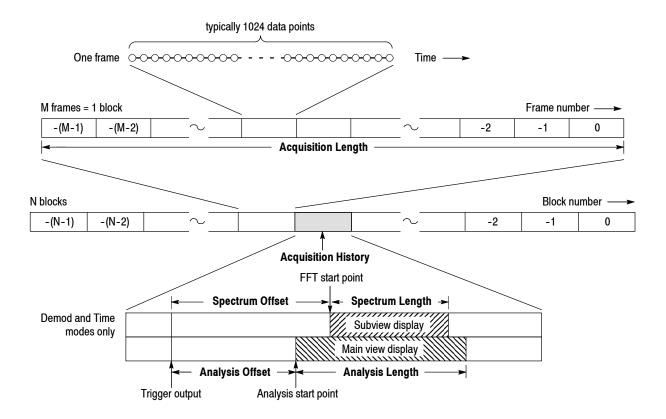

### Figure 3-99: Acquisition and analysis parameters

# Parameters in Real Time The S/A Mode anal

The Acquisition/Analysis menu in the Real Time S/A (real-time spectrum analysis) mode contains the following items:

Acquisition Length. Same as in the Demod and Time modes.

**Spectrum Offset.** Specifies the number of the frame in the spectrogram to display the spectrum. The latest frame is number zero. Older frames have larger negative numbers.

### Parameters in Real Time S/A with Zoom Mode

The Acquisition/Analysis menu in the Real Time S/A with Zoom mode contains the following items (see Figure 3-100):

Acquisition Length. Same as in the Demod and Time modes.

Acquisition History. Same as in the Demod and Time modes.

Analysis Length. Same as in the Demod and Time modes.

Analysis Offset. Same as in the Demod and Time modes.

Frequency Center. Sets the frequency of the center of analysis area.

Frequency Width. Sets the frequency width of the analysis area.


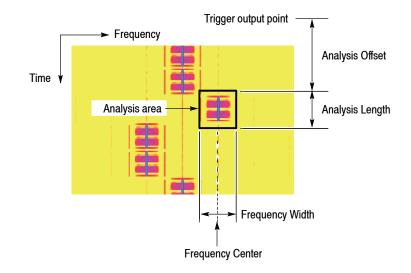


Figure 3-100: Parameters in the zoom mode


### **Seamless Acquisition**

The frame data is acquired every designated time. The interval between one acquisition and another is called a "frame cycle," as shown in Figure 3-101.

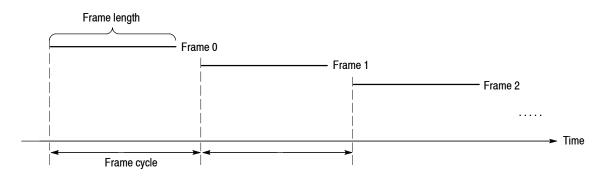

### Figure 3-101: Frame cycle

If the frame cycle is longer than the frame length, time gaps appear between frames. The shorter the frame cycle, the more precisely you can observe the spectrum waveform variation in time. You can acquire frames without gaps when the span setting is below 15 MHz. Acquiring frame data without a gap is called Seamless Acquisition. See Figure 3-102.

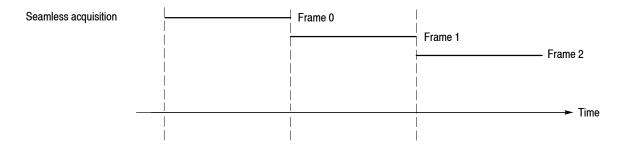

### Figure 3-102: Seamless acquisition

When the span is set greater than 15 MHz, one display frame is reproduced from multiple acquired frames and the frame cycle becomes meaningless. Refer to *Vector Span* on page 3-126 for details.

# Trigger

Triggers determine when the analyzer stops acquiring an input signal and displays measurement results. To properly acquire data that you want to measure, you need to set up the trigger conditions. This section describes the following trigger parameters:

- Mode: Selects free run or triggered acquisition.
- **Repeat:** Determines whether to acquire data continuously or singly.
- **Source:** Selects a trigger signal source.
- Level: Sets a trigger level.
- Slope: Selects the rising or falling edge of a trigger signal.
- **Position:** Specifies a trigger position.

You must create a trigger mask when triggering in the frequency domain with Option 02. Refer to page 3-155 for creating a trigger mask.

A "T" mark (with "O" optionally) indicating a trigger point is displayed on the overview in the Demod and Time modes. Refer to page 3-161 for information about displaying the trigger point.

# Trigger Menu

Use the Trigger menu to set trigger conditions. Press the **Trig** key on the front panel and set the parameters using the following menu items.

**NOTE**. Except for the Repeat menu item, the Trigger menu is available when the measurement mode is set to Real Time S/A, Demod, or Time.

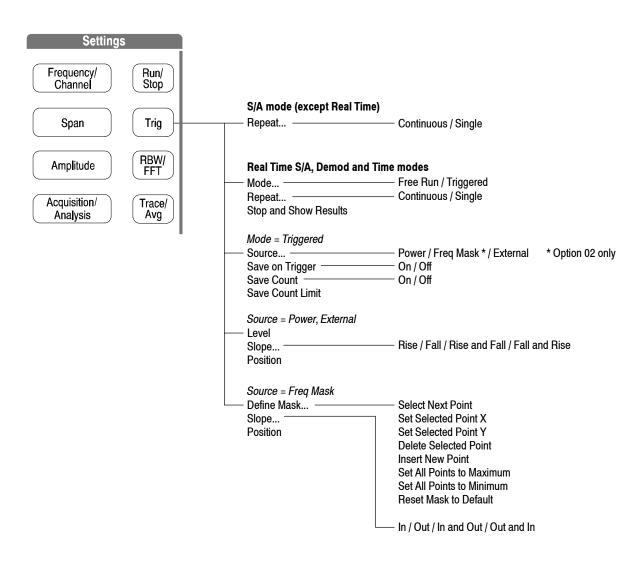


Figure 3-103: Trigger menu structure

**Mode...** Selects the trigger mode. The trigger mode determines whether to acquire data with or without triggering.

- Free Run. Acquires and displays the waveform without triggering. Press the Run/Stop key to start data acquisition. To stop the acquisition, press the Run/Stop key again.
- Triggered. Set the trigger conditions (level, slope, and position) before starting data acquisition by pressing the Run/Stop key. When the trigger occurs, data is acquired and displayed. To stop data acquisition when the trigger does not occur, press the Run/Stop key again.

Repeat... Selects whether to acquire data continuously or singly.

- **Continuous.** Repeatedly acquires and displays the waveform. Acquisition data will be overwritten while waiting for a new trigger event. Use the Single setting if you will need to examine or reanalyze the measurement results.
- Single. Acquires and displays one waveform. After the first waveform display, you must press Run/Stop to acquire and display each waveform.

Table 3-23 shows the acquisition method with the trigger and repeat modes. The concept is illustrated in Figure 3-107 on page 3-153.

Trigger mode (Trig → Mode)	Repeat mode (Trig → Repeat)	Description
Free run	Continuous	Press Run/Stop to acquire data repeatedly. Press this key again to stop acquisition.
	Single	Press Run/Stop to acquire one waveform.
Triggered	Continuous	Press Run/Stop to acquire data every trigger event.
	Single	Press Run/Stop to acquire one waveform every trigger event.

#### Table 3-23: How to acquire a waveform

**Stop and Show Results.** Stops data acquisition and shows measurement results. Instead of the **Run/Stop** key, you can stop data acquisition by pressing the **Stop and Show Results** side key when the measurement mode is Real Time S/A, Demod, or Time. However, there is the following difference in function:

- Pressing the Run/Stop key to stop data acquisition: The block data acquired at the moment of pressing the key is discarded and the measurement results are shown for the previous block of data.
- Pressing Trig → Stop and Show Results to stop data acquisition: Shows measurement results even for the block of data acquired at the moment of pressing the side key, which might not fill a whole block.

In either case, press Run/Stop again to restart acquisition.

Source. Selects a trigger source when the trigger mode is Triggered.

- **Power (Span BW).** *Default.* Uses time domain I/Q data of the input signal as the trigger source. You can set the trigger level, slope, and position.
- Freq Mask. *Option 02 only*. Uses a trigger mask as the trigger source. Refer to page 3-155 for information about creating the trigger mask.
- External. Uses the external signal input from the TRIG IN connector on the rear panel as the trigger source. You can set the trigger level, slope, and position. For the specification of the external trigger, refer to Table A-13 on page A-10 in *Appendix A*.

**Position.** Available when the trigger mode is set to Triggered. Specifies the position of the trigger in a block as a percentage of all the frames in the block. For example, if you set the trigger position to 50%, the trigger occurs at the center of the block, as shown in Figure 3-104.

Range: 0 to 100% in 1% step.

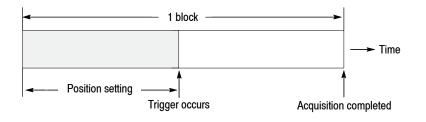


Figure 3-104: Trigger position

**Level.** Sets the trigger level. Available when the trigger mode is set to Triggered and the source is Power or External. Table 3-24 shows the setting range.

Trigger sourceTrigger levelPower (Span BW)-40 to 0 dBfs in 1 dBfs steps (in time domain)External-1.5 to +1.5 V in 0.1 V steps

Table 3-24: Trigger level setting range

*Hint for the Power trigger.* An input signal is FFT-processed every frame (= 1024 points). However, the displayed data is less than 1024 points (refer to *Relation between Frame, Bin, and Pixel* on page 3-183). For example, 801 bins are displayed when the span is 10 MHz. The Power trigger always observes one whole frame (= 1024 points) to determine trigger generation. Note that, as shown in Figure 3-105, any signal exceeding the trigger level outside the range of 10 MHz of span activates the Power trigger even if no signal exceeds the trigger level within the range of 10 MHz of span.

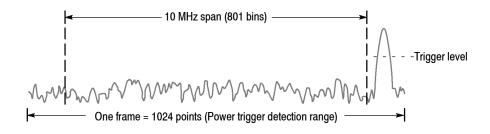


Figure 3-105: Power trigger detection range

**Define Mask...** Option 02 only. Creates a trigger mask when the trigger mode is set to Triggered and the source is set to Freq Mask. The trigger mask is created in the spectrum view. Refer to page 3–155 for details on creating a trigger mask.

**Slope.** Selects the trigger slope when the trigger mode is Triggered and the trigger source is External.

- **Rise.** The trigger occurs at the rising edge of the trigger signal.
- **Fall.** The trigger occurs at the falling edge of the trigger signal.
- Rise and Fall. The trigger occurs at the rising edge of the trigger signal to acquire the first block, and at the falling edge to acquire the next block. The rising and falling edges alternate at each block acquisition.
- Fall and Rise. The trigger occurs at the falling edge of the trigger signal to acquire the first block, and at the rising edge to acquire the next block. The rising and falling edges alternate at each block acquisition.

*When using a trigger mask (Option 02 only).* When the trigger source is set to Freq Mask for using a trigger mask, the following selection items are available:

- In. The trigger occurs when an input signal exits the blue area of a trigger mask and enters the black area.
- **Out.** The trigger occurs when an input signal exits the black area and enters the blue area of a trigger mask.
- In and Out. The analyzer triggers with In to acquire the first block, and Out to acquire the second. In and Out alternate at each block acquisition. Use the acquisition history control to view alternate In/Out blocks.
- Out and In. The analyzer triggers with Out to acquire the first block, and In to acquire the second. In and Out alternate at each block acquisition. Use the acquisition history control to view alternate In/Out blocks.

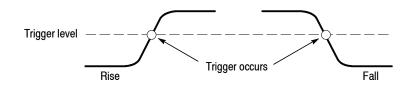


Figure 3-106: Trigger level and slope

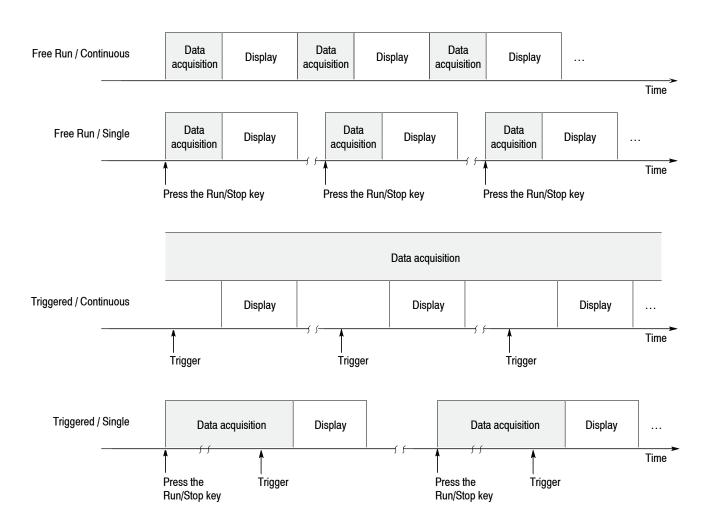


Figure 3-107: Acquiring and displaying data by Trigger and Repeat modes

Data Saving on Trigger	A Save-on-Trigger function is provided to save one block of input data to the IQT file each time a trigger occurs. This function is enabled when the trigger mode is set to Triggered. For details of file operations, refer to page 3-233.
	The file name and location are as follows:

■ File name: yyyymmdd-hhmmss-*ticks*.IQT

Where

yyyy is the year, mm the month, and dd the date. hh is the hour, mm the minutes, and ss the seconds. *ticks* is the time stamp of the final frame. (Refer to *unsigned long ticks* on page 3–251.)

Example: 20050721-140120-0000154761.IQT

 Directory: C:\Documents and Settings\<username>\My Documents\Save-on-Trigger\

Use the following Trigger menu items to control the Save-on-Trigger function:

Save On Trigger. Enables and disables the Save-on-Trigger function.

- On. Enables the Save-on-Trigger function
- **Off.** *Default.* Disables the Save-on-Trigger function

**Save Count**. Selects whether to set a limit on the number of times that data is saved.

- **On.** *Default*. When the number of data save operations reaches the Save Count Limit indicated below, data saving is halted.
- Off. No limit on data save operations is set. In this case, data saving is halted using the **Run/Stop** key on the front panel or the GPIB command.

**NOTE**. When the internal hard disk becomes full, data saving is halted and the "Media full" error message appears.

For deleting a file, refer to page 3-244.

**Save Count Limit**. Sets a limit on the number of times that data is saved. This setting is enabled when the above-mentioned Save Count is On.

Setting range: 1 to 16383 (default: 100)

## Creating a Trigger Mask (Option 02 Only)

**NOTE**. The trigger mask function is available for Option 02 in the measurement mode (Mode) of the real-time spectrum analysis (Real Time S/A), the modulation analysis (Demod), and the time analysis (Time).

Trigger masks are limited to -60 dBfs.

The trigger mask (shown in Figure 3-108) is an area created on the graticule of the spectrum view (subview in the Demod and Time modes). The trigger occurs when the input signal exits or enters the area.

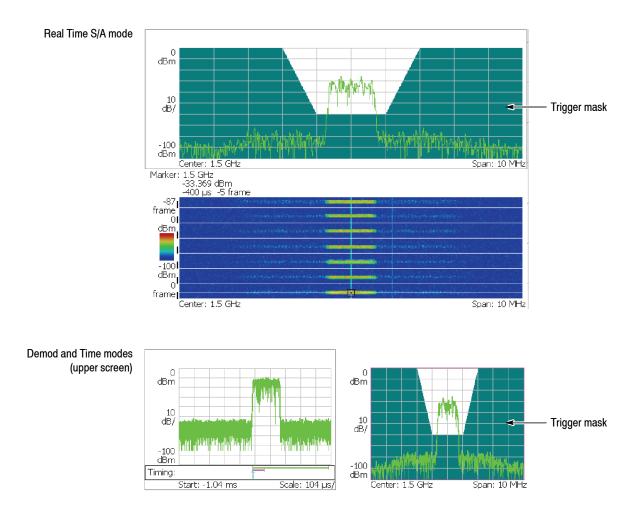


Figure 3-108: Trigger mask

**Mask Creation Menu** Use the **Trig**  $\rightarrow$  **Define Mask** menu and the marker ( $\bigcirc$ ) to create the mask.

**Select Next Point.** Selects the marker to control. The active marker is displayed in red.

Set Selected Point X. Sets the horizontal position of the active marker.

Set Selected Point Y. Sets the vertical position of the active marker.

Delete Selected Point. Deletes the point at the active marker position.

**Insert New Point.** Inserts the new point halfway between the active marker and the next marker on the right side.

**Set All Points to Maximum.** Fills the area below the maximum line (the reference level). See Figure 3-109.

**Set All Points to Minimum.** Fills the area below the minimum line (the level lower than the reference level by 60 dB). See Figure 3–109.

**Reset Mask to Default.** Displays the default mask. See Figure 3-109.

Set All Points to Maximum

Set All Points to Minimum

		1			
	- 6	50 dE	3 —		
		V			

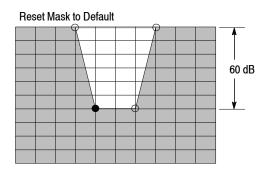


Figure 3-109: Filling operation for creating a mask


### Conditions for Creating a Mask

The required conditions for creating a trigger mask are as follows:

- Measurement mode (Mode): Real Time S/A, Demod, or Time
- Trigger mode (Trig  $\rightarrow$  Mode): Triggered
- Trigger source (Trig  $\rightarrow$  Source): Freq Mask

### Example of Mask Creation

The procedure shows you how to create the example mask shown in Figure 3-110.

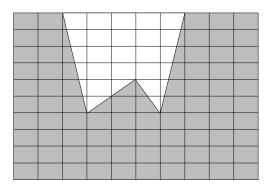


### Figure 3-110: Example mask

- **1.** Before you start creating a trigger mask, make sure that the requirements listed in *Conditions for Creating a Mask* above are satisfied.
- 2. Stop data acquisition using the Run/Stop key.
- **3.** The trigger mask is created in the spectrum view. To display only the spectrum view on the screen, follow these steps:
  - a. Press the View: Select key on the front panel to select the spectrum view.
  - **b.** Press the View: **Define** key on the front panel.
  - c. Press the Show Views side key and select Single.
- 4. Press the **Trig** key on the front panel.
- 5. Press the **Define Mask...** side key.

The default setting is for the whole area to fill with blue on the graticule.

6. Press the **Reset Mask to Default** side key.

The default mask appears (see Figure 3-111).

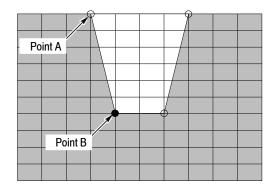

#### Figure 3-111: Default mask

**NOTE**. In the following steps, use the general purpose knob or the numeric input keypad to move the points.

- 7. Change the position of Point A:
  - a. Press the Select Next Point side key until Point A is selected.
  - **b.** Press the **Set Selected Point X** side key and set the horizontal position of Point A to two divisions from the left edge (see Figure 3-112).


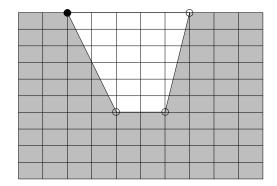


Figure 3-112: Changing the position of Point A

- **8.** Change the position of Point B:
  - a. Press the Select Next Point side key to select Point B.
  - **b.** Press the **Set Selected Point X** side key and set the horizontal position of Point B to three divisions from the left edge (see Figure 3-113).

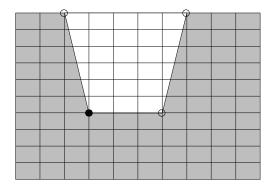

### Figure 3-113: Changing the position of Point B

- 9. Add Point C:
  - **a.** Press the **Insert New Point** side key with Point B active. The new point appears halfway between Point B and the next point in the right side.
  - **b.** Press the **Set Selected Point X** side key and set the horizontal position of Point C to five divisions from the left edge.
  - **c.** Press the **Set Selected Point Y** side key and set the vertical position of Point C to four divisions from the top edge (see Figure 3-114).

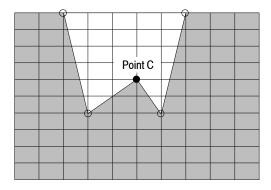


Figure 3-114: Adding Point C

The trigger mask is stored in the internal memory.

10. If you started with a multi-view display, return using the following steps:


- **a.** Press the View: **Define** key on the front panel.
- b. Press the Show Views side key and select Multi.
- 11. Set the trigger:
  - a. Press the Trig key on the front panel.
  - b. Set Slope and Position appropriately.
- 12. Start data acquisition using the **Run/Stop** key.

The analyzer stops data acquisition when the trigger event occurs.

The trigger mask is stored in the internal memory for each measurement mode of S/A: Real Time S/A, Demod: Analog Demod, Digital Demod, Time: Transient, and CCDF. When you press the **Preset** key on the front panel, the instrument settings are reset to the default only for the currently active mode.

### **Trigger Point Indicator**

In the Demod (modulation analysis) and the Time (time analysis) modes, "T" is displayed in the Timing field on the overview. "T" indicates the trigger occurrence point in the Triggered mode, and is used as a reference point for setting the acquisition and analysis parameters in the Free Run mode. However, "T" indicates the trigger output point when the trigger source is Power and the slope is *Rise and Fall* or *Fall and Rise*, or when the trigger source is Freq Mask.

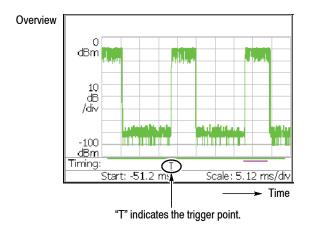



**Indicating Trigger Output** You can use the trigger output to synchronize the analyzer with other instruments. In the Demod (modulation analysis) and Time (time analysis) modes, "O" indicating the trigger output time can be also displayed in the Timing field on the overview. The default status for "O" is off. Follow these steps to display the indicator:

- 1. Press the Acquisition/Analysis key on the front panel.
- 2. Press the Output Trigger Indicator side key to select On.

"O" appears on the overview.

The trigger output timing is determined by the analyzer hardware so that you cannot change it. For the external trigger, the output timing coincides with the trigger occurrence. For other triggers, there is no relationship between the timing of the trigger occurrence and that of the trigger output.

To connect the trigger output to other instruments, use the TRIG OUT connector on the rear panel (see *Rear Panel* on page 2–4). The output specification is: H level >2.0 V, L level <0.4V, and output current <1 mA.

# Synchronous Operation with External Instruments

Use the external trigger input and output connectors on the rear panel (see Figure 2-2 on page 2-4) to synchronize the operation of the analyzer with an external instrument. By synchronizing two or more analyzers, you can acquire multiple signals in different span keeping the time correlation.

Figure 3-116 shows an example of two RSA3303B/RSA3308B analyzers performing IMT2000 signal analysis. One analyzer is used as a master to acquire an uplink (or forward link) signal and the other is used as a slave to acquire a downlink (or reverse link) signal. The master trigger output is connected to the slave trigger input. If you trigger the master with a preamble in the first PRACH signal sent from the mobile terminal, then the slave starts data acquisition. You can seamlessly acquire the uplink and downlink communication data in the time-correlated manner for 10 seconds (with 5 MHz span) to analyze them.

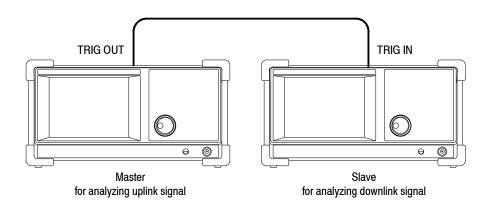


Figure 3-116: Synchronous operation of two analyzers

# **FFT and RBW**

An input signal is converted to the frequency domain data by FFT (Fast Fourier Transform) processing. RBW (Resolution Bandwidth) processing is then applied, for compatibility with measurement data from conventional swept spectrum analyzers (see Figure 3-117). Moreover, the Real Time S/A mode has FFT overlap capability.

**NOTE**. You can set the FFT and RBW parameters in the S/A (spectrum analysis) mode. In the Demod (modulation analysis) and Time (time analysis) modes, the number of FFT points is always 1024 and the window is always Blackman-Harris 4B type. However, for the pulse measurements in the Time mode, you can select Nyquist or Blackman-Harris 4B (refer to RBW/FFT Menu on page 3-97). The Demod and Time modes have no RBW process.

The adjustable parameters for FFT are:

- FFT points
- FFT window
- FFT start point (Real Time S/A only)

The adjustable parameters for RBW are:

- Resolution bandwidth (RBW)
- Filter shape
- Roll-off ratio (for the Nyquist or Root Nyquist filter)

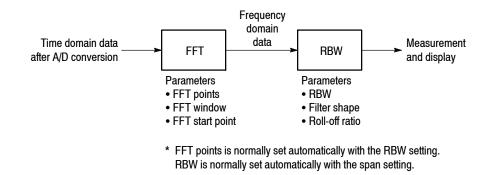


Figure 3-117: FFT and RBW process

# **RBW/FFT Menu**

Figure 3-118 shows the RBW/FFT menu structure.

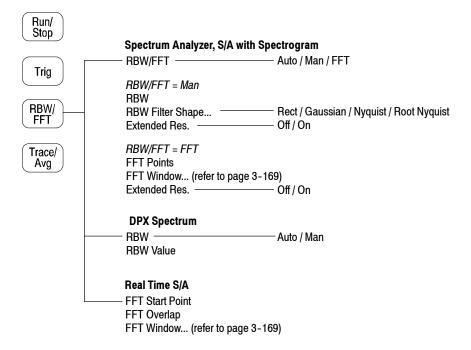


Figure 3-118: RBW/FFT menu tree

### Spectrum Analyzer and S/A with Spectrogram Modes

The RBW/FFT menu has the following controls in the Spectrum Analyzer and the S/A with Spctrogram modes:

**RBW/FFT.** Selects whether to set RBW and FFT parameters automatically or manually.

- Auto. Sets RBW automatically with the span setting. The filter shape is set to Gaussian.
- Man. Sets RBW and selects the filter manually with the RBW and RBW Filter Shape... side keys.
- FFT. Sets the FFT points and window manually with the FFT Points and FFT Window... side keys. The result of FFT processing is displayed as is, without RBW processing.

### When RBW/FFT is set to Man.

**RBW.** Sets RBW when Man is selected in RBW/FFT. Range: 2 kHz to 2 MHz (the default is 80 kHz)

**RBW Filter Shape...** Selects the filter from the following four types when RBW/FFT is set to Man:

- Rect (Rectangular)
- Gaussian (default)
- Nyquist
- Root Nyquist

**Rolloff Ratio.** Enters roll-off ratio when RBW Filter Shape is set to Nyquist or Root Nyquist. Range: 0.0001 to 1 (the default is 0.5)

**Extended Res.** The FFT point number is normally limited internally. You can eliminate the limit using the On setting. Refer to *Limit on the FFT Points* on page 3-168 for detail.

**NOTE**. It is recommended to keep Extended Res. Off as its default condition.

### When RBW/FFT is set to FFT.

**FFT Points.** Selects FFT sample point numbers per frame when RBW/FFT is set to FFT. Range: 64 to 65536 in increments of  $2^n$ . Higher numbers provide higher resolution, and lower numbers provide faster measurements.

**FFT Window...** Selects the FFT window (window function) when RBW/FFT is set to FFT. For the window type, refer to Table 3-26 on page 3-171. The default is Blackman-Harris 4B.

Extended Res. Same as when RBW/FFT is set to Man.

When RBW/FFT is set to FFT, the waveform without RBW processing is displayed on screen (see Figure 3-119).

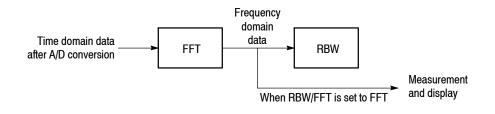


Figure 3-119: Process flow when RBW/FFT = FFT

**DPX Spectrum Mode** In the DPX spectrum mode, FFT is performed by the dedicated processor and the parameters cannot be controlled. The RBW has the following controls:

**RBW.** Selects whether to set RBW automatically or manually.

- Auto. Sets RBW automatically by the span setting using the Gaussian filter.
- Manual. Sets RBW manually using the **RBW Value** side key.

**RBW Value.** Sets the RBW when you select Manual in the **RBW** side key. Range: 1 to 10% of the span setting, rounded to the nearest settable value.

#### Overlapped FFT (Real Time S/A Mode)

The Real Time S/A mode has FFT overlap capability that computes FFT while overlapping 1024-point FFT frames by a specified number of samples. The RBW/FFT menu contains the following controls:

**NOTE**. In the Real Time S/A mode, the FFT points is always 1024 and FFTprocessed data does not go through RBW process (see Figure 3-119 on page 3-166).

**FFT Start Point.** Sets the start point of the 1024-point FFT frame by the number of samples from the previous frame (see Figure 3-120). Range: 1 to 1024 samples in increments of  $2^{n}$ .

**FFT Overlap.** Displays the amount of overlap between frames used to calculate consecutive FFTs (see Figure 3-120). Not settable. The sum of FFT Overlap and FFT Start Point is always 1024.

**FFT Window...** Selects the FFT window (window function). For the window type, refer to Table 3-26 on page 3-171. The default is Blackman-Harris 4B.

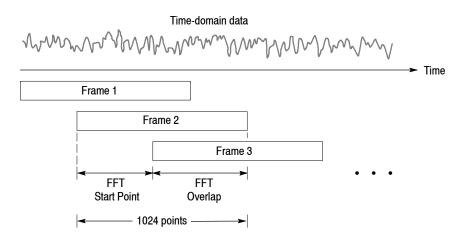


Figure 3-120: FFT overlap

FFT Points				
	The number of FFT points is normally set automatically by the RBW setting. The value is basically 1024. It can be set in powers of 2 from 64 to 8192. This value is the number of points in one physical frame in the time and frequency domains. If the number of points is reduced, the frame period is shortened and more spectrum variations can be observed in the spectrogram view. If the number of the points is increased, S/N ratio and frequency resolution are improved.			
	<b>NOTE</b> . Setting the number of FFT points greater than the limit (8192) causes the noise floor to go down and occasionally causes spurious signals to appear (see Figure 3-121). It cannot be determined whether the spurious signal comes from the input signal or the instrument.			
Limit on the FFT Points	The number of FFT points is normally limited to 8192 internally to avoid displaying internally-generated spurious signals. However, you can remove this restriction and set the number of FFT points up to 65536 using the following steps:			
	<b>1.</b> Display the measurement signal on the screen.			
	2. Press the <b>RBW/FFT</b> key on the front panel.			
	3. Press the <b>RBW/FFT</b> side key to select <b>FFT</b> .			
	4. Press the Extended Res. side key to select On.			
	<b>5.</b> Press the <b>FFT Points</b> side key and select the value with the general purpose knob. The range extends up to 65536 in powers of 2.			
	The waveform after FFT is displayed with the selected number of points.			
	Month of the second			




# **FFT Window**


	The phase of the waveform processed by FFT analysis is assumed to start at zero and end at zero, so the waveform data is an exact multiple of one cycle. If the phase at the start and end of the waveform are equal, there is no unnatural discontinuity in the signal waveform, and the frequency and the amplitude can be accurately calculated.	
	If the waveform data is not an exact multiple of one cycle, the amplitudes at the start and end of the waveform will be different. There will be a discontinuity of the waveform between the beginning and end, and a high-frequency transient phenomenon occurs. When this happens, inaccurate frequency information is recorded in the frequency domain.	
	If a window function is applied to the waveform, the start and end amplitudes will be closer, and the discontinuity is reduced. The frequency components that are calculated from the signal by FFT also become more accurate. Choose the FFT window according to your objective: accurate measurement of the frequency or accurate measurement of the amplitude of the frequency components.	
Window Characteristics	In FFT windows, the frequency resolution is inversely proportional to the amplitude accuracy. Select a proper window according to the measurement items and characteristics of the signal source. Table 3-25 shows characteristics and usages of major windows.	

FFT window	Characteristics	Usage
Rectangular	<ul> <li>Suitable for frequency measurement but not suitable for amplitude measurement.</li> </ul>	<ul> <li>Transient phenomenon or burst; signal levels before and after the event are roughly same.</li> </ul>
	Same result as measurement without a window.	<ul> <li>Sinusoidal wave with little amplitude variation and stable frequency.</li> </ul>
		<ul> <li>Wide bandwidth irregular noise; spectrum changes slowly.</li> </ul>
Hamming, Hanning	Suitable for frequency measurement.	Sinusoidal wave.
	<ul> <li>Inferior in accuracy of amplitude to the rectangular window.</li> <li>Frequency resolution of Hamming is slightly superior to that of Hanning.</li> </ul>	<ul> <li>Repeating narrow bandwidth irregular noise.</li> <li>Transient phenomenon or burst; signal levels before and after the event are remarkably different.</li> </ul>
Blackman-Harris	<ul> <li>Suitable for amplitude measurement but not suitable for frequency measurement.</li> </ul>	<ul> <li>Signal that has one-frequency ruling when a high order harmonic is detected.</li> </ul>

Table 3-25:	Characteristics ar	nd usage o	f FFT windows
		na uouge o	

Figure 3-122 outlines how frequency domain data is generated from time domain data. The FFT window serves as a bandpass filter between time and frequency domain data. The FFT frequency resolution and amplitude accuracy of each frequency component depends on the window shape.

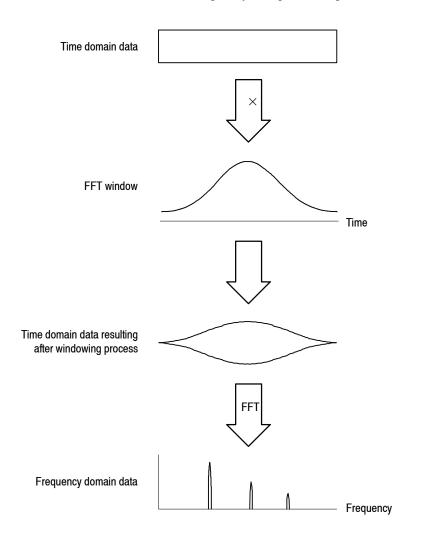

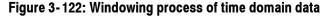

Generally, window frequency resolution is inversely proportional to the accuracy of measuring amplitude levels. For ordinal measurements, select the window capable of separating the desired frequency component. Such a window maximizes the accuracy of measuring amplitude levels and minimizes leakage error while separating each frequency component.

To select the proper window, first select the rectangular window, then sequentially switch to a window with less frequency resolution, such as Hamming, Hanning, or Blackman-Harris. Use the last window that still passes the frequency component to be separated. Suitable frequency resolution and amplitude accuracy are obtained by using the window immediately before the one from which the frequency component cannot be separated.
Consider the following characteristics when selecting a window for your purpose:
<ul> <li>Frequency resolution is improved by reducing the width of the main lobe window.</li> </ul>
Accuracy of the amplitude levels of frequency components is improved by reducing side lobe level relative to the main lobe.

**Window Type** The analyzer supports a total of 15 windows, including the major windows mentioned previously (refer to Table 3-26).

Table 3-26:	FFT window	and band	pass filter
-------------	------------	----------	-------------

Window	Bandpass filter	-3 dB bandwidth	Maximum side lobe	Equivalent noise bandwidth
Blackman-Harris 3 sample A type	0 dB -20 -40 -60 -80 -100	1.53	-62 dB	1.61075
Blackman-Harris 3 sample B type	0 dB -20 -40 -60 -80 -100	1.622	-71 dB	1.708538
Blackman-Harris 4 sample A type	0 dB -20 -40 -60 -80 -100 -120		-76 dB	1.793948
Blackman-Harris 4 sample B type (Default)	0 dB -20 -40 -60 -80 -100 -120 -100 -120		-92 dB	2.004353

Table 3-26: FFT window and bandpass filter (Cont.)

Window	Bandpass filter	-3 dB bandwidth	Maximum side lobe	Equivalent noise bandwidth
Blackman	0 dB -20 -40 -60 -80 -100	1.642	-58 dB	1.726757
Hamming	0 dB -20 -40 -60 -80	1.302	-43 dB	1.362826
Hanning	0 dB -20 -40 -60 -80 -100	1.438	-32 dB	1.5
Parzen	0 dB -20 -40 -60 -80 -100	1.27	-27 dB	1.330747
Rosenfield	0 dB -20 -40 -60 -80 -100 -120 -40 -60 -60 -100 -120 -120 -120 -120 -120 -120 -12		-48 dB	1.90989
Welch	0 dB -20 -40 -60 -80		-21 dB	1.197677
Sine Lobe	0 dB -20 -40 -40 -80	1.186	-23 dB	1.233702

Window	Bandpass filter	-3 dB bandwidth	Maximum side lobe	Equivalent noise bandwidth
Sine Cubed	0 dB -20 -40 -60 -80 -100	1.654	-39 dB	1.734891
Sine to the 4 <sup>th</sup>	0 dB -20 -40 -60 -80 -100 -120	1.85	-47 dB	1.944444
Flat Top	0 dB -20 -40 -60 -80 -100	3.182	-51 dB	3.196927
Rectangular	0 dB -20 -20 -40 -60	0.886	-13 dB	1

Table 3-26: FFT window and bandpass filter (Cont.)

FFT and RBW

# **Trace Comparison Display and Average Function**

In the spectrum analysis (S/A mode), you can display two traces concurrently on screen. Trace 1 is displayed in yellow and Trace 2 in green (see Figure 3-123). The trace can be averaged for noise reduction. You can save waveform data in files, and load them as Trace 1 or 2.

**NOTE**. Trace comparison display and average function are available in the spectrum analysis (S/A) mode except Real Time S/A.

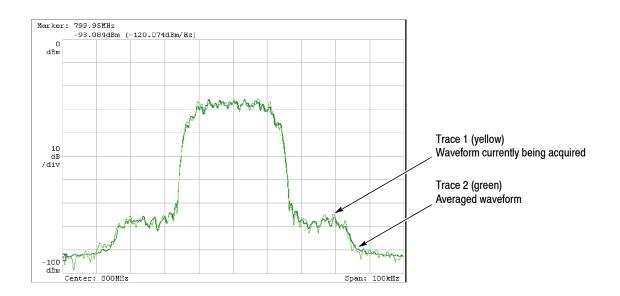


Figure 3-123: Comparison display of Trace 1 and 2

The main topics in this section are:

- Trace/Avg menu
- Displaying Trace 1 and 2
- Averaging waveform
- Saving/Loading waveform data
- Trace compression display

## Trace/Avg Menu

Controls the trace display and averaging. The menu items vary according to the the measurement mode: S/A, Demod, or Time.

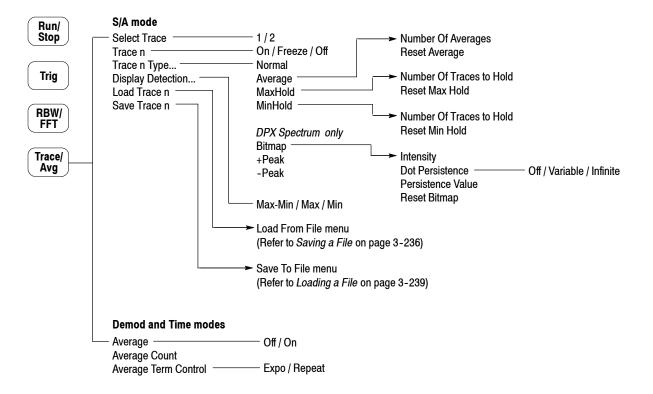


Figure 3-124: Trace/Avg menu structure

S/A Mode The Trace/Avg menu has the following controls in the S/A mode.

**NOTE**. The Trace/Avg menu is not available in the Real Time S/A mode.

Select Trace. Selects the trace to control: trace 1 or 2.

Trace 1/2. Controls the selected trace.

- **On.** Turns on the selected trace.
- **Freeze.** Freezes display of the selected trace.
- Off. Turns off the selected trace.

**Trace 1/2 Type...** Selects the type of processing for the selected trace. This menu item appears when **Trace 1/2** is set to **On**.

Table 3-27: Trace type

Trace type	Description	
Normal <sup>1</sup>	Displays normal spectrum without averaging.	
Average	<i>DPX spectrum.</i> Computes the average amplitude of the spectrum points within each frequency bin to create a trace for each screen update, and the average the traces with the number set in the <b>Number Of Averages</b> side key to display.	
	Other than DPX spectrum. Averages the traces with the number set in the <b>Number Of Averages</b> side key.	
Max Hold	Holds the maximum level at each frequency. Pressing the <b>Reset MaxHold</b> side key restarts the process.	
Min Hold	Holds the minimum level at each frequency. Pressing the <b>Reset MinHold</b> side key restarts the process.	
Bitmap <sup>2</sup>	Displays the density of acquired data. The number of data points acquired at each pixel is indicated by color.	
+Peak <sup>2</sup>	Displays the maximum values acquired in each screen update.	
-Peak <sup>2</sup>	Displays the minimum values acquired in each screen update.	

<sup>1</sup> Available in the S/A mode except DPX spectrum.

<sup>2</sup> DPX spectrum only.

**Number Of Averages.** Specifies how many traces are accumulated to create the averaged value. Setting range: 1 to 100,000 (default: 20).

Averaging is controlled in two ways as shown in Table 3-28.

Table 3-28: Averaging method

Data acquisition	Average type	Number of Averages
Free run (Continuous mode only)	Exponential RMS	Continues the average with an exponential weighting applied to old values using Number of Averages as the weighting factor.
Triggered and Single mode	RMS	Averages traces accumulating up to Number of Averages, then stops the acquisition until the next trigger event occurs.

**Intensity.** *DPX spectrum only.* Sets the display intensity for the Bitmap trace. An increased intensity level allows a single, short event to be seen. This also allows the effect of variable persistence to be seen on infrequent events. Range: 1 to 100% (default: 25%)

**Dot Persistence.** *DPX spectrum only.* 

Selects the persistence mode for the Bitmap trace.

- Off. Turns off the persistence mode. All points are reset every screen update.
- Variable. Selects the variable persistence mode.
   Set a persistence time, that is how long a point in the display is visible before fading, using the **Persistence Value** side key.
- **Infinite.** Selects the infinite persistence mode that prevents a point in the display from fading. Once a point in the display has been written to, it will remain visible indefinitely.

**Persistence Value.** *DPX spectrum only.* Sets the persistence time for the variable persistence mode. Range: 1 to 1000 (no unit, default: 10).

**Reset Average.** *Available when Trace 1 or 2 Type is set to Average.* Causes trace accumulation to start over for the average process.

**Reset MaxHold.** *Available when Trace 1 or 2 Type is set to MaxHold.* Causes trace accumulation to start over for the Max Hold process.

**Reset MinHold.** *Available when Trace 1 or 2 Type is set to MinHold.* Causes trace accumulation to start over for the Min Hold process.

**Reset Bitmap.** *Available when Trace 1 or 2 Type is set to Bitmap.* Causes trace accumulation to start over for the Bitmap display.

**Display Detection...** The bin data is compressed in the display because the number of pixels in the horizontal direction of the screen is generally smaller than the number of bins in the FFT. Display Detection selects the method to be used for decimating traces to fit the available space. It is related only to the decimation of displayed data.

- **Max-Min.** Draws a line connecting the maximum and the minimum values of data corresponding to each pixel.
- Max. Displays the maximum value of data corresponding to each pixel.
- Min. Displays the minimum value of data corresponding to each pixel.

The analyzer performs RMS detection using the linear A/D converter regardless of the Display Detection setting.

Refer to page 3-183 for details on compression of waveform display data.

Load Trace. Loads trace data from a file.

Save Trace. Saves trace data to a file.

Refer to *File Operations* on page 3-233 for details on saving/loading data to/from a file.

**Demod and Time Modes** The Trace/Avg menu has the following controls in the Demod and Time modes.

Average. Determines whether to perform averaging.

- **On.** Turn averaging on.
- **Off.** Turn averaging off.

**NOTE**. For Demod and Time modes, data is always acquired without averaging.

**Average Count.** Specifies the number of measurements to combine. Setting range: 1 to 10000 (default: 20).

**Average Term Control.** Specifies the action when more than Average Count measurement results are generated.

- **Expo.** Continues the average with an exponential weighting applied to old values.
- **Repeat.** Clears average data and counter, and restarts the average process.

#### **Displaying Trace 1 and 2**

Only Trace 1 (yellow) is displayed by default. Use the following steps to select Trace 1 and 2:

- 1. Press the Trace/Avg key on the front panel.
- 2. Press the Select Trace side key to select the trace (1 or 2) to be controlled. For example, to control Trace 2, select 2.
- **3.** Press the **Trace 1** (or **2**) side key to select the display method: On, Freeze, or OFF.
- 4. Press the **Trace 1** (or **2**) **Type...** side key to select the trace type: Normal, Average, MaxHold or MinHold (refer to *Trace 1/2 Type* on page 3-177).
- 5. Repeat steps 2 through 4 for Trace 1 and 2.

#### Averaging the Waveform

The averaging technique is generally used to average and reduce the noise of the waveform. The average function includes the peak hold which maintains the maximum and minimum values as well as the averaging process.

**Average Type** There are four averaging types as follows, where these variables are used:

 $X(p)_n$ : Display data for the n<sup>th</sup> frame

 $x(p)_n$ : Active data for the n<sup>th</sup> frame

*P*: Frame point

*N*: Value of "Number of Averages"

**RMS.** Root-mean-square. Used when acquiring data in single mode. N frames average, then the acquisition stops.

$X(p)_n = x(p)_n$	for : n = 1
$X(p)_n = \frac{(n-1) \times X(p)_{n-1} + x(p)_n}{n}$	for : $2 \leq n \leq N$
$X(p)_n = x(p)_N$	for : n > N

**Exponential RMS.** Exponential root-mean-square. Used when acquiring data in continuous mode. Continuous averaging weights older sweeps so that they have a progressively smaller effect on the average.

$$X(p)_n = x(p)_n \qquad \text{for : } n = 1$$

$$X(p)_n = \frac{(n-1) \times X(p)_{n-1} + x(p)_n}{n} \qquad \text{for : } 2 \le n \le N$$

$$X(p)_n = \frac{(N-1) \times X(p)_{n-1} + x(p)_n}{N} \qquad \text{for : } n > N$$

Max Hold. Displays only the maximum value at each data point.

$$\begin{split} X(p)_n &= x(p)_n & \text{for : } n = 1 \\ X(p)_n &= \max(X(p)_{n-1}, x(p)_n) & \text{for : } n \ge 2 \end{split}$$

Min Hold. Displays only the minimum value at each data point.

$$\begin{split} X(p)_n &= x(p)_n & \text{for : } n = 1 \\ X(p)_n &= \min(X(p)_{n-1}, x(p)_n) & \text{for : } n \ge 2 \end{split}$$

**Averaging Example** This section shows how to use averaging and the compared display.

Performing Averaging. Average and display Trace 1.

- **1.** Display the spectrum of the measurement signal.
- 2. Pause data acquisition to simplify the operation. If you are acquiring the data in the continuous mode, press the **Run/Stop** key to stop the acquisition.
- 3. Press the Trace/Avg key on the front panel.
- 4. Press the Select Trace side key to select 1.
- 5. Press the Trace 1 Type side key and select Average for example.
- 6. Press the Number of Averages side key and set 64 for example.
- 7. Press the Run/Stop key to acquire data.
- 8. Press the **Reset Average** side key to perform averaging again.

The averaged waveform is displayed on screen (see Figure 3-125). The average count is indicated in the right top portion of the screen.

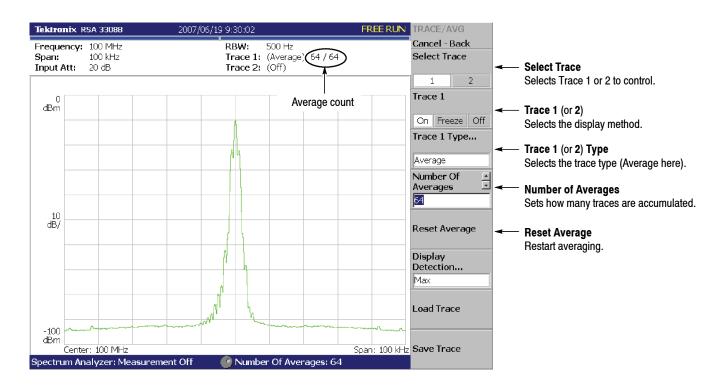


Figure 3-125: Displaying an averaged waveform

**Compared Display.** Display Trace 1 as a normal spectrum and Trace 2 as the averaged waveform concurrently, and compare them.

- 1. Display the spectrum of the measurement signal.
- 2. Press the Trace/Avg key on the front panel.
- **3.** Display Trace 1 as a normal spectrum.
  - a. Press the Select Trace side key to select 1.
  - b. Press the Trace 1 Type... side key and select Normal.
- 4. Display Trace 2 as the averaged waveform.
  - a. Press the Select Trace side key to select 2.
  - **b.** Press the **Trace 2 Type...** side key and select **Average**, **Max Hold** or **Min Hold**.

The waveform currently being acquired (trace 1 in yellow) and the averaged waveform (trace 2 in green) are displayed together.

Figure 3-126 shows an example of the concurrent display of the ordinal spectrum and its Max Hold waveform.

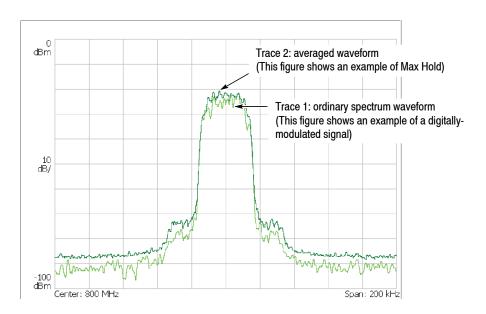


Figure 3-126: Compared display

### Saving/Loading Waveform Data

You can save the waveform data currently being acquired in a file, and load the saved data as trace 1 or 2.

- 1. Press the Trace/Avg key on the front panel.
- 2. Press the Select Trace side key to select Trace 1 or 2.
- **3.** Press the **Save Trace** side key and select the destination file. For information about file operations, refer to page 3-233 and subsequent pages.

**Loading a Trace** To load Trace 1 or 2 from a file:

- 1. Press the Trace/Avg key on the front panel.
- 2. Press the Select Trace side key to select Trace 1 or 2.
- **3.** Press the **Load Trace** side key and select the source file. (The trace is automatically frozen.) For information about file operations, refer to page 3–233 and subsequent pages.

#### **Trace Compression**

Although waveform data are acquired as 1024 points per frame, the acquired data are displayed after being compressed due to the limitation in the number of pixels on screen. The compression method and the procedure to select the method are described in the following section.

Relation between Frame,<br/>Bin, and PixelOne frame contains 1024 FFT output points. Part of the data in the frame is<br/>invalid for the calculation. The analyzer discards invalid data and displays only<br/>valid data. Valid data are placed into "bins". The number of bins depends on the<br/>span, as shown in Table 3-29.

Span	Number of bins
2 MHz or less	641
5 MHz	801
10 MHz	801
15 MHz	601
20 MHz (Baseband only)	801

Table 3-29: Number	of bins (	FFT P	points:	1024)
--------------------	-----------	-------	---------	-------


The number of bins is valid in all modes except Scalar mode. The number of bins is meaningless in Scalar mode because it uses several physical frames to display data. Generally, the following calculations are used:

Frequency bandwidth of one bin = (Sampling rate)/(Number of FFT points)

Number of bins = (Specified span)/(Frequency bandwidth of one bin) +1

The sampling rate varies depending on span. Refer to page A-7 in *Appendix A: Specifications* for more details.

**Compression Method** Generally, because the number of horizontal pixels on screen is smaller than the number of bins, bin data is compressed according to the number of pixels on screen when they are displayed (see Figure 3-127).

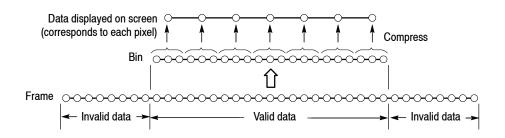

#### Figure 3-127: Relationships between frame, bin, and pixel

There are three types of compression methods: Max, Min, and Max-Min (see Figure 3-128 on page 3-185).

Selecting Compression<br/>MethodMax is the most commonly used compression method. Max-Min is used in the<br/>time domain waveform display in the Demod mode. (The compression method is<br/>always set to Max on the spectrogram view in three-dimensional modes.)


You can select the compression method using the following steps in S/A mode:

- 1. Press the Trace/Avg key on the front panel.
- 2. Press the Select Trace side key to select Trace 1 or 2.
- 3. Press the Display Detection... side key and select Max, Min, or Max-Min.

**NOTE**. Display Detection determines how to deal with data between spectrum display points. The analyzer performs RMS detection using the linear A/D converter regardless of the Display Detection setting.


#### Max


Acquires the maximum value of data points corresponding to each pixel.

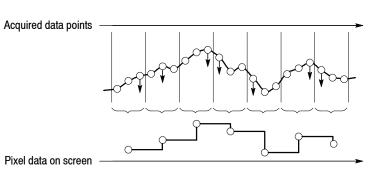



Acquires the minimum value of data points corresponding to each pixel.

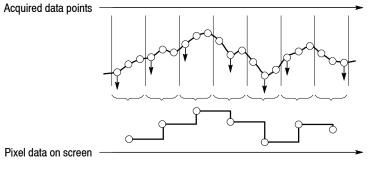

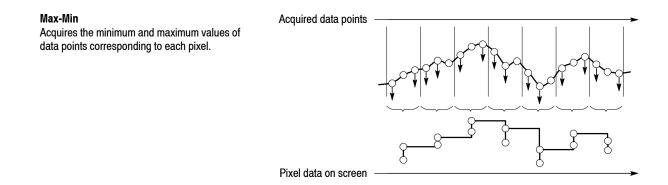



# **Setting Views**

This section describes scaling and formatting for the following view types:

- Spectrum view
- DPX Spectrum view
- Spectrogram view
- Time domain view
- CCDF view


For Option 21 only:

- Constellation view
- EVM view
- Symbol table
- Eye diagram
- AM/AM view
- AM/PM view
- PDF view
- Noisogram view

For information on views specific to the other optional analyses, refer to each user manual.

## **View Menu**

Use the following View keys on the front panel to scale and format views:

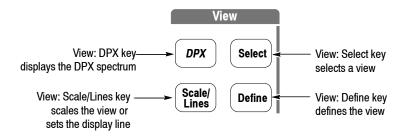

- **DPX** Displays the DPX spectrum view.
- **Select** Selects a view when displaying two or three views on screen. The view switches each time you press this key. The selected view is surrounded with a white box.
- **Define** Selects a view style, format, and content. The menu items depend on the measurement. Refer to each measurement description.
- **Scale/Lines** Press this key and then select one of the following items:
  - View Scale... Sets the horizontal and vertical axes for the view selected using the Select key.
  - View Lines... Sets the display line for the view selected using the Select key. Refer to *Display Line* on page 3-209 for the details.

## **Basic Procedure**

The following procedures set the scale or format the view using the View keys in single view or multiple views.

**Procedure for Single View.** When you are displaying one view on screen, press View: **Scale/Lines**  $\rightarrow$  **View Scale...** and set the scale. For the Scale menu of a specific view, refer to the following pages.

**Procedure for Multiple Views.** When you are displaying multiple views on the screen, follow these steps:

- **1.** Select a view by pressing the View: **Select** key. The selected view is surrounded by a white frame.
- 2. If necessary, change a multiple view display to a single view display:
  - **a.** Press the View: **Define** key.
  - b. Press the Show Views side key to select Single.

Only the selected view will be displayed.

- 3. Set the scale by pressing View: Scale/Lines  $\rightarrow$  View Scale....
- 4. If necessary, return to the multiple view display:
  - **a.** Press the View: **Define** key.
  - b. Press the Show Views side key to select Multi.

Refer to page 3-34 for information about changing the overview and the subview using the Define menu.

## **Spectrum View Setting**

The spectrum view is the frequency domain representation of a signal. It indicates frequency along the horizontal axis and power along the vertical axis.

**Scale Menu** The Scale menu for the spectrum view contains the following controls (see Figure 3-129):

**Auto Scale.** Sets the start value and the scale of the vertical axis automatically to fit the waveform to the screen.

Horizontal Scale. Sets the range of the horizontal axis.

Horizontal Start. Sets the minimum value (left edge) of the horizontal axis.

**Vertical Scale.** Sets the range of the vertical axis. Range: 100µ to 100 dB (default: 100 dB).

**Vertical Stop.** Sets the maximum value (top) of the vertical axis. Range: -100 to 100 dBm (default: 0 dBm).

Full Scale. Sets the scale of the vertical axis to the default full-scale value.

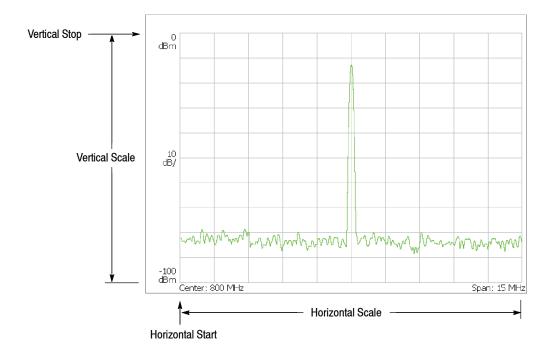


Figure 3-129: Setting the scale in spectrum view

## **DPX Spectrum View Setting**

The DPX Spectrum indicates frequency along the horizontal axis and power along the vertical axis. The Bitmap trace displays the density of acquired data by color grading.

# **Scale Menu** The Scale menu for the DPX spectrum view contains the following controls (see Figure 3-130):

**Vertical Scale.** Sets the full-scale of the vertical axis. Setting value: 20 or 100 dB (default: 100 dB).

**Vertical Stop.** Sets the maximum value (top) of the vertical axis. Range: -100 to 100 dBm (default: 0 dBm).

Bitmap Color. Selects the color scheme used for the Bitmap trace.

#### Table 3-30: Bitmap color scheme

Color scheme	Description
Temperature	Temperature color (low density in blue and high in red)
Spectral	Spectral color (low density in red and high in blue)
Gray	Gray scale (low density in black and high in white)
Binary Cyan	Binary cyan (lower density than the minimum value in black and higher in cyan)

**Maximum.** Sets the histogram value represented by the top color in the selected palette. Range: 1 to 100% in 1% steps (default: 100%) The Maximum value must be greater than Minimim.

**Minimum.** Sets the histogram value represented by the bottom color in the selected palette. Range: 0 to 99% in 1% steps (default: 0%) The Minimum value must be less than Maximum.

Full Scale. Sets Vertical Scale, Maximum and Minimum to default values.

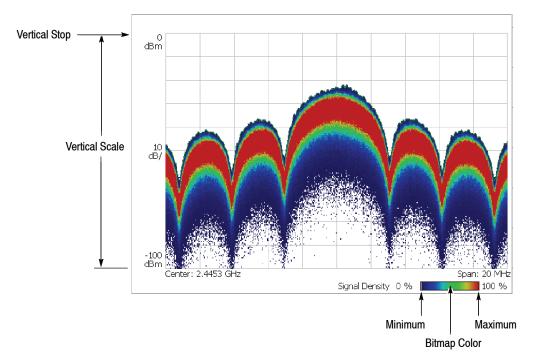


Figure 3-130: Setting the scale in spectrum view

### **Spectrogram View Setting**

The spectrogram view shows time-varying spectrum three-dimensionally. It indicates frequency along the horizontal axis, frame number along the vertical axis, and power level by using color. This view is used in the *S/A with Spectrogram* and *Real Time S/A* modes, or the Demod and Time modes as the overview content.

**NOTE**. When you select S/A with Spectrogram in the S/A mode, you cannot set the scale of the spectrogram.

Scale Menu The Scale menu for the spectrogram view contains the following controls (see Figure 3-131):

**Auto Scale.** Sets the start value and the scale of the color axis automatically to fit the spectrogram to the screen.

Horizontal Scale. Sets the range of the horizontal axis.

Horizontal Start. Sets the minimum value (left edge) of the horizontal axis.

Vertical Size. Sets the range of the vertical axis. Range: 87 to 89088 frames.

Vertical Start. Sets the start frame number for the vertical axis.

**Color Scale.** Sets the scale (value subtracting the minimum power value from the maximum power value) of the color axis. The level is represented in 100 steps (100 colors) from the minimum value (blue) to the maximum value (red) in the default state.

Color Stop. Sets the maximum value (top) of the color axis.

**Full Scale.** Sets the maximum value of the color axis to the reference level and the height to 100 dB.

The following menu items appear when the overview is a spectrogram in the Demod (modulation analysis) and the Time (time analysis) modes.

**Step to Spectrum Window.** Vertical Start is automatically set so that the frame displayed as a spectrum in the subview appears on the spectrogram.

**Step to Analysis Window.** Vertical Start is automatically set so that the frame analyzed and displayed in the main view appears on the spectrogram.

**Step to Trigger.** Vertical Start is automatically set so that the frame at the trigger position appears on the spectrogram.

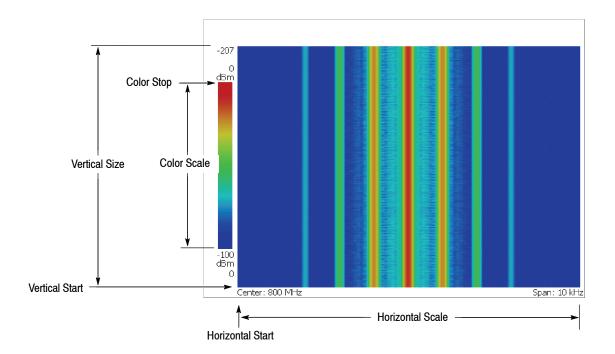


Figure 3-131: Scale and format settings in spectrogram view

#### **Time Domain View Setting**

The time domain view includes the following displays shown on the overview or the main view in the Demod and the Time modes:

- Power versus Time
- AM demodulation display (modulation factor versus time)
- FM demodulation display (frequency deviation versus time)
- PM demodulation display (phase deviation versus time)
- I/Q level versus Time

#### For Option 21 only:

- Random jitter versus Time
- Integrated phase noise versus Time
- C/N versus Time
- **Scale Menu** The Scale menu for the time domain view contains the following controls (see Figure 3-132):

**Auto Scale.** Sets the start value and the scale of the vertical axis automatically to fit the waveform to the screen.

Horizontal Scale. Sets the range of the horizontal axis.

Horizontal Start. Sets the minimum value (left edge) of the horizontal axis.

Vertical Scale. Sets the range of the vertical axis.

**Vertical Stop.** Sets the maximum value (top) of the vertical axis. Available when the vertical axis represents:


Power Random jitter Integrated phase noise C/N

**Vertical Offset.** Sets the center value ((maximum + minimum) / 2) of the vertical axis. Available when the vertical axis represents:

AM modulation factor FM frequency deviation PM phase deviation IQ voltage.

Full Scale. Sets the scale of the vertical axis to the default full-scale value.

Vertical axis: Power, Random jitter, Integrated phase noise, or C/N (Power shown)

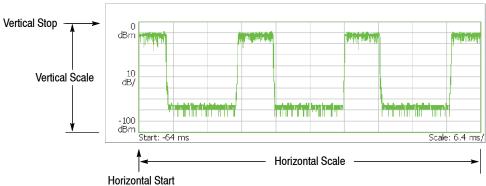

Vertical axis: AM modulation factor, FM frequency deviation, PM phase deviation, or IQ level (AM modulation factor shown)

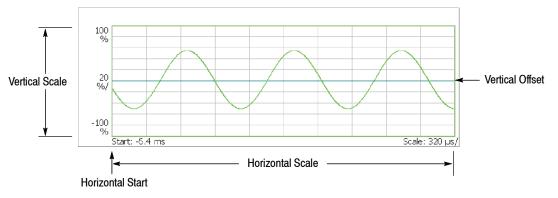


Figure 3-132: Setting the scale in time domain view

### **CCDF View Setting**

The CCDF view displays the CCDF measurement in the Digital Demod mode (digital modulation analysis, Option 21 only) and the Time mode (time analysis). The horizontal axis indicates amplitude and the vertical axis (logarithmic scale) indicates CCDF. Refer to page 3-88 for more information about the CCDF measurement.

# **Scale Menu** The Scale menu for the CCDF view contains the following controls (see Figure 3-133):

**Auto Scale.** Sets the starting value and the scale of the vertical axis automatically to fit the waveform to the screen.

Horizontal Scale. Sets the range of the horizontal axis.


Horizontal Start. Sets the minimum value (left edge) of the horizontal axis.

**Vertical Stop.** Sets the maximum value (top) of the vertical axis. Range: Twice of Vertical Start value to 100% in a 1-2-5 sequence.

**Vertical Start.** Sets the minimum value (bottom) of the vertical axis. Range:  $10^{-9}$  to 1/2 of Vertical Stop value in a 1-2-5 sequence.

Full Scale. Sets the scale of the vertical axis to the default full-scale value.

Sub Grid. *Time mode only*. Determines whether to display the sub-grid.

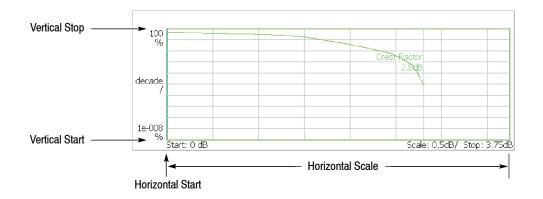



## **Constellation View Setting**

**Option 21 Only** 

The constellation view shows the signal represented by phase and amplitude in polar coordinate or IQ diagram. You can display it in the constellation analysis of the digital modulation analysis (refer to page 3–50).

The constellation displays I and Q components of the signal in two dimensions. As the overall signal level changes, both I and Q are automatically scaled to keep the size of the constellation relatively constant from one update to the next. This new scale is then normalized to a unitless scale.

**Scale Menu** The Scale menu for the constellation view contains the following controls (see Figure 3-134):

**Measurement Content...** Selects vector or constellation display. (see Figure 3-134).

- Vector. Selects the vector display. Signals represented by phase and amplitude (such as digitally modulated signals) are shown in a polar coordinate or IQ diagram. The red point indicates the symbol position of the measurement signal and the yellow trace indicates the locus of the signal between symbols. The magnitude of the error vector is evaluated by comparing the point where yellow traces are concentrated to the red point. The cross marks indicate symbol positions of the ideal signal.
- Constellation. Selects the constellation display. It is the same as the vector display, except that only symbols of the measurement signal are indicated in red, and the locus between symbols is not displayed.

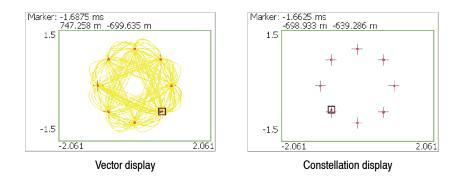


Figure 3-134: Vector and constellation displays


# **EVM View Setting**


**Option 21 Only** 

	The horizontal axis indicates time, and the vertical axis indicates EVM, amplitude, or phase in the EVM view. You can display this view in EVM analysis of the digital modulation analysis (refer to page 3-51).
Scale Menu	The Scale menu for the EVM view contains the following controls (see Figure 3-135):
	Auto Scale. Sets the starting value and the scale of the vertical axis automatically to size the waveform to the screen.
	Horizontal Scale. Sets the range of the horizontal axis.
	Horizontal Start. Sets the minimum value (left edge) of the horizontal axis.
	Vertical Scale. Sets the range of the vertical axis.
	<b>Vertical Start.</b> Available when the vertical axis represents EVM. Sets the minimum value (bottom) of the vertical axis.

**Vertical Offset.** Available when the vertical axis represents magnitude or phase error. Sets the center value ((maximum + minimum) / 2) of the vertical axis.


**Full Scale.** Sets the scale of the vertical axis to the default full-scale value.

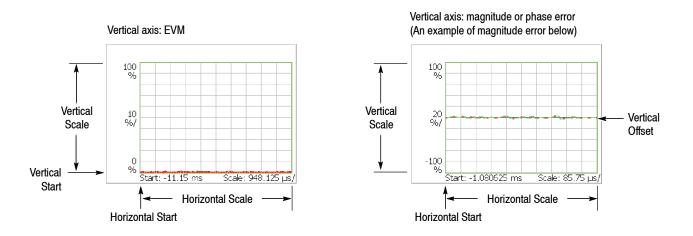

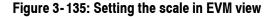

**Measurement Content...** Selects one of the following formats for view (see Figure 3-136).

- **EVM.** Displays changes of EVM (Error Vector Magnitude) in time series.
- Mag Error. Displays changes of magnitude error in time series.
- **Phase Error.** Displays changes of phase error in time series.

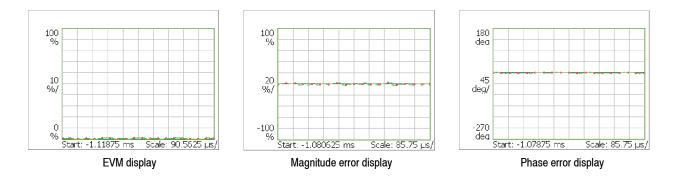

#### Figure 3-136: EVM, magnitude and phase error displays

You can select degree or radian for the angular unit by pressing System  $\rightarrow$  Instrument Setup... $\rightarrow$  Angular Units.

See Figure 3-137 for EVM, Mag Error, and Phase Error. The figure is an example of the constellation view in  $1/4\pi$ QPSK modulation. The cross marks are referred to as symbols and indicate phase positions of the ideal signal (amplitude is fixed in this case). In this modulation, the bit pattern is determined by movements from each position. For example, if the actual signal has shifted from the ideal symbol position to the  $\bullet$  position, you can evaluate vector signal quality as magnitudes of error in radius direction (amplitude), error in phase direction, and total error vector. These three types of errors correspond to the three types of views in the EVM view:

- EVM (%RMS): Root-mean-square value of EVM (Error Vector Magnitude)
- Mag Error (%RMS): Root-mean-square value of magnitude error
- Phase Error (%RMS): Root-mean-square value of phase error

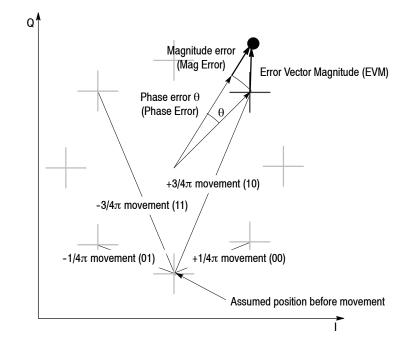


Figure 3-137: Constellation view and error vectors in 1/4 $\pi$  QPSK

## **Symbol Table Setting**

**Option 21 Only** 

The symbol table shows symbol change with time. You can display this table in the symbol table analysis in the digital modulation analysis (refer to page 3–54).

Marker: 199sym 01	
0: 01111010 01101110 01000101 00001010 11010011 11110110	×

Figure 3-138: Symbol table

**Scale Menu** The Scale menu for the symbol table contains the following controls (see Figure 3-138):

**Radix.** Selects the radix for displaying the table from hexadecimal digit (Hex), octal digit (Oct), and binary digit (Bin).

**Rotate.** Sets the value start position. The setting range is 0 to 3. This item is invalid in  $1/4\pi$  QPSK and GMSK modulations because absolute coordinates have no meaning.

**Decoding Start Position.** Selects the decoding start position for the ASK, FSK, and GFSK signals with the decoding format of Manchester or Miller.

- Auto. Determines the decoding start position automatically.
- 0. Starts decoding from the beginning of a symbol.
- +1. Delays the decoding start position by half a symbol.

**Sequence.** Selects how to determine the symbol value for the D8PSK signals. Valid when D8PSK is selected as Modulation Type in the Meas Setup menu.

- Code. Selects the code sequence that determines the symbol value by the phase difference between the current and the previous symbols using Gray code. For example, when the phase difference is  $\pi/2$ , the symbol value is 3.
- Phase. Selects the phase sequence that determines the symbol value by the phase angle directly. For example, when the phase angle is  $\pi/2$ , the symbol value is 2. See Figure 3-139.

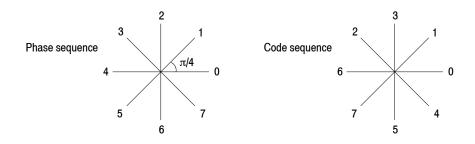


Figure 3-139: D8PSK symbol value

#### Eye Diagram Setting

**Option 21 Only** 

The eye diagram shows I or Q signal that is triggered on a clock signal aligned with the symbol rate, representing time along the horizontal axis and amplitude or phase along the vertical axis. You can display this diagram in the eye diagram analysis in the digital modulation analysis (refer to page 3–55).

**Scale Menu** The Scale menu for the eye diagram contains the following controls:

**Measurement Content...** Selects the vertical axis of the eye diagram (see Figure 3-140).

- I. Displays I signal level on the vertical axis (default).
- **Q.** Displays Q signal level on the vertical axis.
- **Trellis.** Indicates phase on the vertical axis.

**Eye Length.** Sets the number of symbols displayed on the horizontal axis. The time required for movement between symbols is defined as 1 (one).

Range: 1 to 16 (the default is 2)

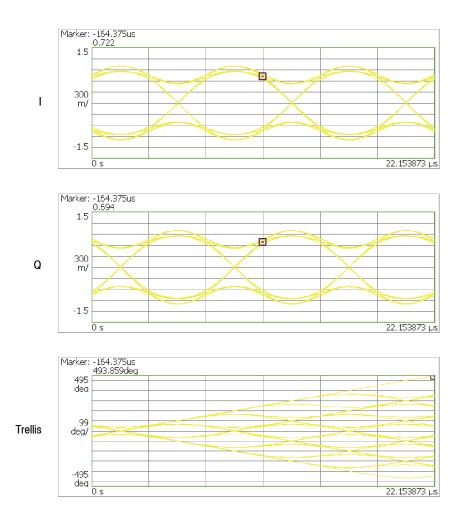


Figure 3-140: Eye diagram view

## **AM/AM View Setting**

#### **Option 21 Only**

The AM/AM view shows non-linearity characteristics of a device under test (DUT) such as an RF amplifier, representing reference amplitude along the horizontal axis and measured signal amplitude along the vertical axis. Refer to *AM/AM Measurement* on page 3–56.

**Scale Menu** The Scale menu for the AM/AM view contains the following controls:

Measurement Content... Selects a vector or dot display (see Figure 3-141).

- Vector. Displays yellow lines between the dots (default).
- **Dot.** Displays the calculated result as a series of red dots.

By default, the measurement results are displayed with the horizontal and vertical axes scaled automatically to fit the graph to the screen.

**Horizontal Start.** Sets the minimum value (left edge) of the horizontal axis to expand the graph.

Range: The initial minimum value (left edge) to the maximum value (right edge)

**Vertical Start.** Shows the minimum value (bottom) of the vertical axis (not settable). Equal to the Horizontal Start value.

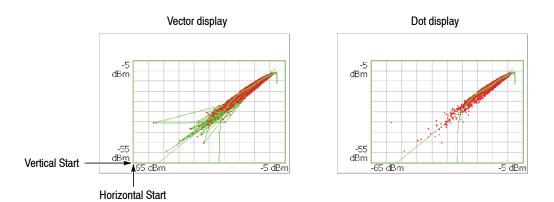


Figure 3-141: Vector and dot displays for AM/AM view

### **AM/PM View Setting**

#### **Option 21 Only**

The AM/PM view shows non-linearity characteristics of a device under test (DUT) such as RF amplifier, representing reference amplitude along the horizontal axis and phase error along the vertical axis. Refer to *AM/PM Measurement* on page 3–58.

**Scale Menu** The Scale menu for the AM/PM view contains the following controls:

Measurement Content... Selects a vector or dot display (see Figure 3-142).

- Vector. Displays yellow lines between the dots (default).
- **Dot.** Displays the calculated result as a series of red dots.

By default, the measurement results are displayed with the horizontal and vertical axes scaled automatically to fit the graph to the screen.

**Horizontal Start.** Sets the minimum value (left edge) of the horizontal axis to expand the graph horizontally.

Range: The initial minimum value (left edge) to the maximum value (right edge)

**Vertical Scale.** Sets the vertical full-scale value to expand the graph vertically. Range: The initial maximum value  $\times 0.05$  to the maximum value

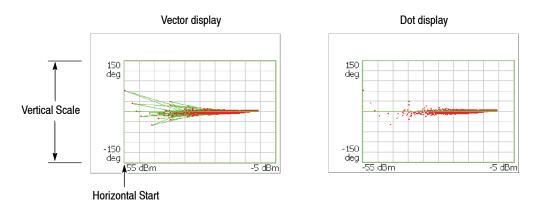


Figure 3-142: Vector and dot displays for AM/PM view

## **PDF View Setting**

**Option 21 Only** 

This view is displayed in the PDF (Probability Distribution Function) measurement in the digital modulation analysis (refer to page 3-60). The horizontal axis represents power level relative to the mean value at the center and the vertical axis represents occurrence probability.

**Scale Menu** The Scale menu for the PDF view contains the following controls (see Figure 3-143):

**Auto Scale.** Sets the starting value and the scale of the vertical axis automatically to fit the waveform to the screen.

Horizontal Scale. Sets the range of the horizontal axis.

Horizontal Start. Sets the minimum value (left edge) of the horizontal axis.

**NOTE**. By default, the horizontal axis is displayed in the range of  $\pm 12$  dB. If the amplitude exceeds 12 dB, the horizontal axis will be set automatically to fit the waveform to the screen.

Vertical Scale. Sets the range of the vertical axis.

Vertical Stop. Sets the maximum value (top) of the vertical axis.

**Full Scale.** Sets the scale of the horizontal and vertical axes to each full-scale value. You can change the scale of the horizontal and vertical axes within these full-scale values.

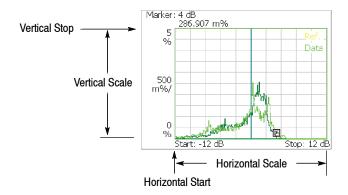


Figure 3-143: Setting the scale in PDF view

## **Noisogram View Setting**

**Option 21 Only** 

The noisogram provides a three dimensional representation of time-varying phase noise. It indicates frequency along the horizontal axis, frame number along the vertical axis, like a spectrogram, and C/N by using color. This view is used as the subview in the signal source analysis (refer to page 3-103).

## **Scale Menu** The Scale menu for the Noisogram contains the following controls (see Figure 3-144):

**Auto Scale.** Sets the start value and the scale of the color axis automatically to fit the noisogram to the screen.

**Horizontal Stop.** Sets the maximum value (right edge) of the horizontal axis. Range: (Horizontal Start)×2 to (Span)/2.

**Horizontal Start.** Sets the minimum value (left edge) of the horizontal axis. Range: (The minimum offset frequency of the analyzed data) to (Horizontal Stop)/2.

Vertical Size. Sets the range of the vertical axis. Range: 40 to 40960 frames.

**Vertical Start.** Sets the start frame number for the vertical axis. Range: - (the number of frames in the analysis range) to 0. Zero (0) represents the latest frame.

**Color Scale.** Sets the scale (value subtracting the minimum power value from the maximum power value) of the color axis. The level is represented in 100 steps (100 colors) from the minimum value (blue) to the maximum value (red) in the default state. Range: 10 to 100 dB in a 1-2-5 sequence.

**Color Stop.** Sets the maximum value (top) of the color axis. Range: [70 - (Color Scale)] to 70 dBc/Hz.

**Full Scale.** Sets the maximum value of the color axis to 0 and the height to 100 dBc/Hz.


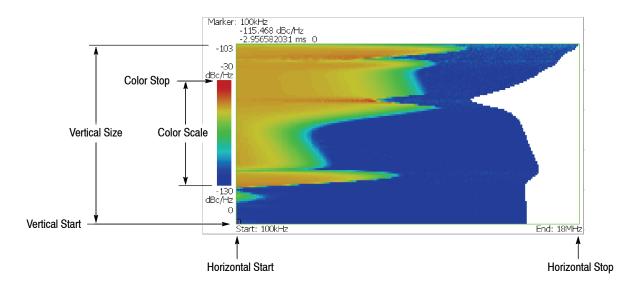



Figure 3-144: Setting the scale in noisogram view

## **Display Line**

The analyzer provides a convenient method to determine whether a signal peak is higher or lower than a specified level, or whether it falls within a specified range. Use the Lines menu to control the display lines. The Display Line feature displays horizontal and/or vertical line(s) at the position you specify.

There are two types of lines: horizontal line and vertical line as shown in Figure 3-145. In the S/A mode other than Real Time S/A, you can display one or two of each line. You can display the horizontal and vertical lines concurrently. The Real Time S/A mode has multi display lines function that shows multiple horizontal and vertical lines.

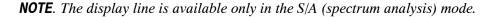

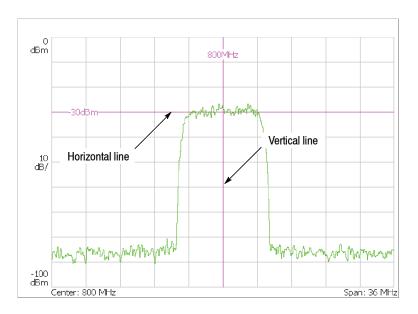

#### Figure 3-145: Display line

This section provides these two topics:

- Display Line Operation (Other than Real Time S/A)
- Multi Display Lines (Real Time S/A Only)

## **Display Line Operation (Other than Real Time S/A)**

Controls the display line. See Figure 3-146 for the Lines menu structure.

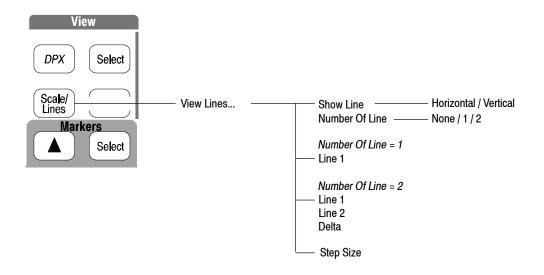


Figure 3-146: Lines menu structure (other than Real Time S/A)

**Lines Menu** The Lines menu in the S/A mode other than the Real Time S/A has the following controls:

Show Line. Selects the display line to be controlled:

- Horizontal
- Vertical

Number Of Line. Selects how many horizontal lines are displayed on the graph:

- None
- **1**
- **2**

Line 1. Sets the position of the first line.

Line 2. Sets the position of the second line.

**Delta.** Sets the difference between the first line and the second: (Value of Line 2) - (Value of Line 1)

#### **Horizontal Display Line**

- 1. Press the Scale/Lines key on the front panel, and then press the View Lines... side key.
- 2. Press the Show Line side key to select Horizontal.

**NOTE**. In the following steps, use the general purpose knob or the numeric input keypad to move the lines.

- 3. Perform one of the following options.
  - To display one horizontal line: Press the Number Of Line side key to select 1.

Press the Line 1 side key and move the line.

To display two horizontal lines:
 Press the Number Of Line side key to select 2.

Press the Line 1 side key and move line 1. Line 2 moves in parallel.

Press the Line 2 side key and move line 2.

Press the **Delta** side key and move line 2. The delta value indicates the following value: (value of Delta) = (value of line 2) - (value of line 1)

 To turn off the horizontal line: Press the Number Of Line side key to select None.

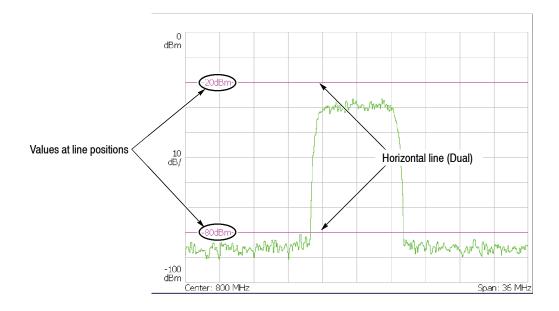


Figure 3-147: Two horizontal lines

#### Vertical Display Line

- 1. Press the Scale/Lines key on the front panel, and then press the View Lines... side key.
- 2. Press the Show Line side key to select Vertical.

**NOTE**. In the following steps, use the general purpose knob or the numeric input keypad to move the lines.

- **3.** Perform one of the following options.
  - To display one vertical line: Press the Number Of Line side key to select 1.

Press the Line 1 side key and move the line.

 To display two vertical lines: Press the Number Of Line side key to select 2.

Press the Line 1 side key and move line 1. Line 2 moves in parallel.

Press the Line 2 side key and move line 2.

Press the **Delta** side key and move line 2. The delta value indicates the following value: (value of Delta) = (value of line 2) - (value of line 1)

 To turn off the vertical line: Press the Number Of Line side key to select None.

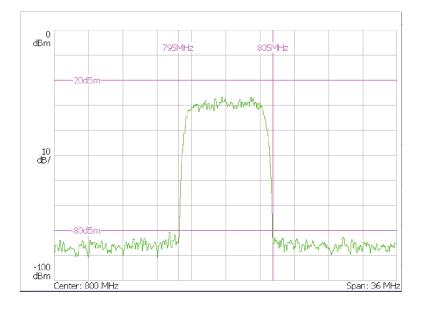


Figure 3-148: Two horizontal and two vertical lines

## Multi Display Lines (Real Time S/A Only)

The Real Time S/A mode has the multi display lines function that shows multiple horizontal and vertical lines as shown in Figure 3-149. The spectrum view uses the amplitude and frequency lines. The spectrogram uses the time and frequency lines. The display lines are placed at regular intervals from the reference line.

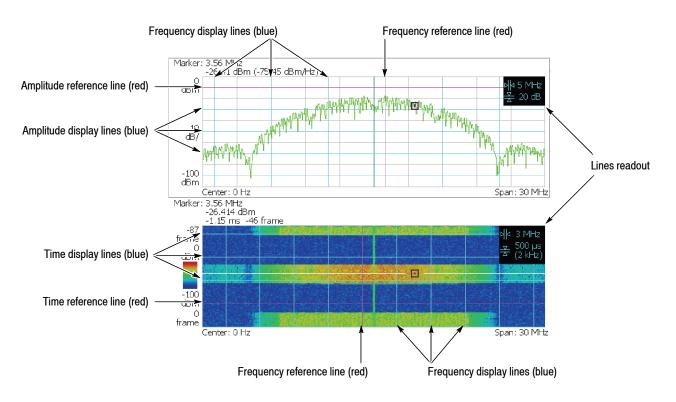


Figure 3-149: Multi display lines (Real Time S/A mode)

**Lines Menu** The Lines menu for the Real Time S/A mode has different controls between the spectrum and spectrogram view. The settings are shown in the lines readout at the upper right of the view (see Figure 3-149).

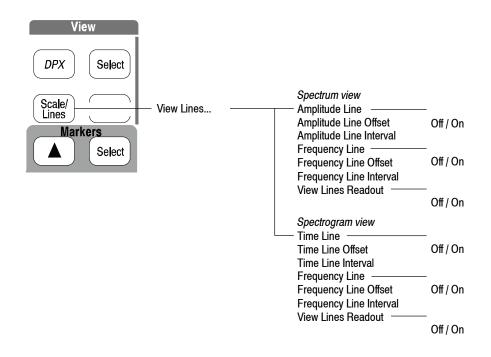


Figure 3-150: Lines menu structure (Real Time S/A)

#### Spectrum View.

**Amplitude Line.** Determines whether to turn on or off the amplitude display lines. The default is Off.

**Amplitude Line Offset.** Sets the offset of the amplitude reference line. Range: -100 to 0 dBm (the default is 0 dBm)

**Amplitude Line Interval.** Sets the interval of the amplitude display lines. Range: 0 to 100 dB (the default is 0 dB)

**Frequency Line.** Determines whether to turn on or off the frequency display lines. The default is Off.

**Frequency Line Offset.** Sets the offset of the frequency reference line. Range: Center frequency  $\pm$  Span/2 (Hz). The default value is the center frequency; the reference line is at center screen.

**Frequency Line Interval.** Sets the interval of the frequency display lines. Range: 0 to full span (the default is 0 Hz)

**View Lines Readout.** Determines whether to turn on or off the display lines readout. The default is On.

#### Spectrogram View.

**Time Line.** Determines whether to turn on or off the time display lines. The default is Off.

**Time Line Offset.** Sets the offset of the time reference line. Range: 0 second maximum (Zero represents the latest frame.) The minimum value depends on acquired data quantity.

**Time Line Interval.** Sets the interval of the time display lines. Range: 0 second minimum. The maximum value depends on acquired data quantity.


The other menu items (Frequency Line, Frequency Line Offset, Frequency Line Interval, and View Lines Readout) are the same as the spectrum view.

Display Line

## **Marker Operation and Peak Search**

The marker moves on a waveform to measure amplitude or frequency. One or two markers can be displayed on the screen. You can also use the reference cursor together with the marker. The marker is also used for peak search. For information on the band power marker used in the spectrum analysis, refer to *Spectrum Analysis (S/A mode)* on page 3-1.

One or two markers (active and fixed) and/or a reference cursor can be displayed on the screen as shown in Figure 3-151.

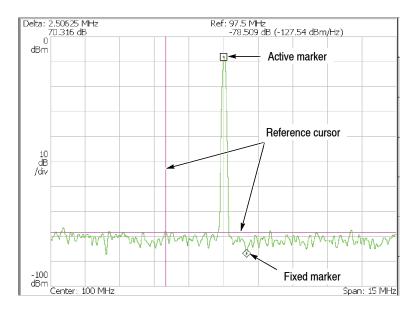

#### Figure 3-151: Marker display

#### Single marker mode


One marker (called Marker 1) indicated by " $\Box$ " moves on a waveform. It is used to measure absolute values.

#### Delta marker mode

Two markers (called Marker 1 and 2) indicated by " $\Box$ " and " $\diamondsuit$ " move on a waveform. " $\Box$ " and " $\diamondsuit$ " represent an active and a fixed marker, respectively. They are used to measure relative values.

## **Markers Menu**

The Markers menu controls the marker operation and peak search. Figure 3-152 shows the menu structure.

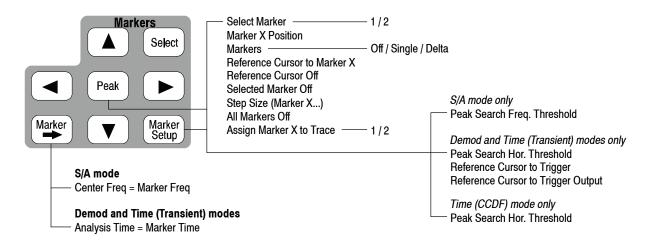

#### Figure 3-152: Markers menu structure

Select	t Selects which marker is being controlled in the delta marker mode. If markers are disabled, pressing this key enables Marker 1.	
Marker 🔿	Sets instrument parameters based on the marker position(s).	
<b>Center Freq = Marker Freq.</b> Changes the center frequency to match the marker position.		
	<b>Analysis Time = Marker Time.</b> <i>Demod and Time modes only.</i> Sets the origin of an analysis range with the marker in the overview. For the detail, refer to <i>Setting the Analysis Range</i> on page 3–30.	
Peak	Positions the marker to the highest peak signal.	
Marker Left (◀) / Right (►)	Moves the marker lower or higher in frequency to the next signal, respectively.	
Marker Up (▲) / Down (▼)	Moves the marker higher or lower in amplitude to the next signal, respectively.	
	The definition of "next" can be set in the Marker Setup menu, using Peak Search Freq. Threshold or Hor. Threshold. Refer to <i>Setting the Minimum Jump of the Marker</i> on page 3–226.	

**Marker Setup** Set up the behavior of the markers.

**Select Marker.** Selects which marker is being controlled in the delta marker mode. Same as Markers: **Select**.

Marker X Position. Sets the horizontal position of the selected marker.

Markers. Selects the marker mode:

- Off. Displays no marker.
- **Single.** Displays one marker (Marker 1).
- **Delta.** Displays two markers (Marker 1 and 2).

**Reference Cursor to Marker X.** Enables the Reference Cursor, at the same position as the selected marker. Only available when the selected view can contain a reference cursor.

**Reference Cursor Off.** Turns off the Reference Cursor. Only available when the selected view can contain a reference marker.

Selected Marker Off. Turns off the selected marker.

All Markers Off. Turns off both markers, the reference cursor, and all their readouts.

Assign Marker X to Trace. When two traces are being displayed, moves the selected marker to another trace with the same horizontal position. Available only if two traces are being displayed.

**Peak Search Freq. Threshold.** *S*/*A mode only*. Sets the minimum frequency jump that should be made when selecting the next signal to the left/right/up/ down.

**Peak Search Hor. Threshold.** *Demod and Time modes only.* Selects the minimum horizontal jump that should be made when selecting the next signal to the left/right/up/down.

**Marker X Vertical.** Spectrogram view in the Real Time S/A mode only. Specifies the number of the frame in the spectrogram to display the spectrum. The latest frame is number zero. Older frames have larger negative numbers. Same as Spectrum Offset in the Acquisition/Analysis menu.

**Reference Cursor to Trigger.** *Demod and Time modes only.* Enables the reference cursor at the same position as the trigger in the overview.

**Reference Cursor to Trigger Output.** *Demod and Time modes only.* Enables the reference cursor at the same position as the trigger output in the overview.

### **Using Markers**

This section shows some examples of marker operation.

**NOTE**. If multiple views are displayed on the screen, select the view to operate the markers using the View: **Select** key on the front panel. The selected view is shown with a white frame.

Measuring with a Single Marker

- To measure amplitude or frequency with a single marker, follow these steps:
  - 1. Press the Marker Setup key on the front panel.
  - 2. Press the Makers side key to select Single. Only Marker 1 is enabled.

Note that the Marker X Position menu item is already selected.

**3.** Rotate the general purpose knob (or enter the value using the numeric keypad) to move the marker to the measurement location.

The marker readout is displayed in the top left portion of the screen (see Figure 3-153).

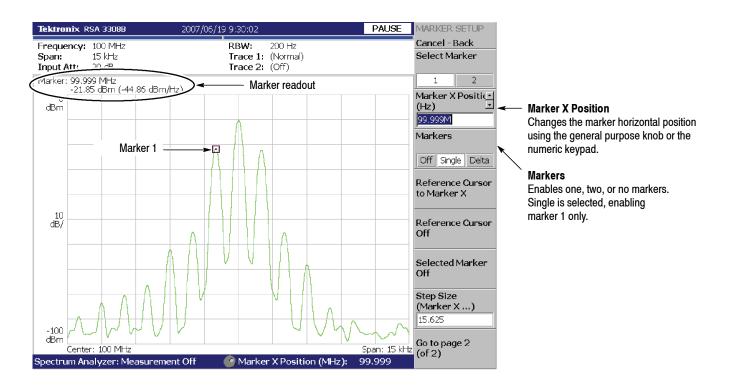


Figure 3-153: Measurement with a single marker

## Measuring with the<br/>Delta MarkerTo measure differences in amplitudes or frequencies using the delta marker,<br/>follow these steps:


- 1. Press the Marker Setup key on the front panel.
- 2. Press the Makers side key to select Delta.

Marker 1 and 2 are displayed, with the two markers overlapping each other. **Select Marker** is set to **1** (Marker 1) by default.

Note that the Marker X Position menu item is already selected.

- **3.** Rotate the general purpose knob (or enter the value using the keypad) to move the active marker to the reference point.
- 4. Change the active marker by pressing the **Select Marker** side key to select **2**.
- 5. Rotate the general purpose knob (or enter the value using the keypad) to move the active marker to the measurement point.

As shown in Figure 3-154, the marker readout is shown in the top left portion of the screen:

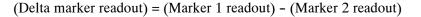

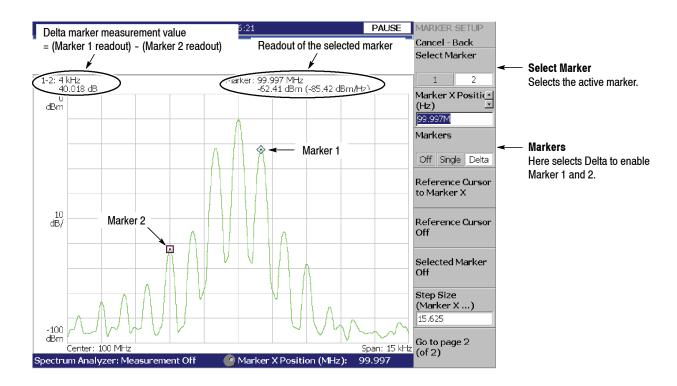


Figure 3-154: Using the delta marker to take measurements

# Measuring with the<br/>Reference CursorThe reference cursor is another way to measure relative amplitude or frequency.<br/>The reference cursor is positioned with the marker and fixed on the screen.

- 1. Press the Marker Setup key on the front panel.
- 2. Press the Makers side key to select Single or Delta.
- **3.** Press the **Marker X Position** side key and rotate the general purpose knob (or enter the value using the keypad) to move the marker to the reference point.
- **4.** Press the **Reference Cursor to Marker X** side key to display the cursor at the marker position.

#### Demod and Time modes only:

On the overview, you can also place the reference cursor at a trigger output position with the **Reference Cursor to Trigger** side key.

5. Press the Marker X Position side key and rotate the general purpose knob (or enter the value using the keypad) to move the marker to the measurement point.

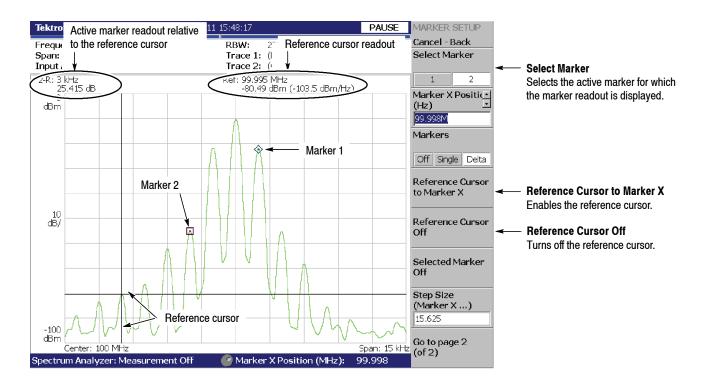


Figure 3-155: Measurement with the reference cursor

As shown in Figure 3-155, the marker readout relative to the reference cursor is shown in the top left portion of the screen for the active marker selected using the **Select Marker** side key.


To turn off the reference cursor, press the Reference Cursor Off side key.

**Changing the Trace** Use the following steps to change the trace on which you operate the marker when displaying two traces on a view. Trace 1 is displayed in yellow and Trace 2 in green.

In the S/A (spectrum analysis) mode, two traces are shown when Trace 2 is turned on using the Trace/Avg menu.

In Demod (modulation analysis) mode and Time (time analysis) mode, two traces are shown in IQ vs. Time view (I level vs. time and Q level vs. time).

- 1. Press the Marker Setup key on the front panel.
- 2. Press the Go to page 2 (of 2) (bottom) side key to show the next menu page.
- **3.** Press the **Assign Marker X to Trace** side key to select Trace **1** or **2**, as shown in Figure 3-156.

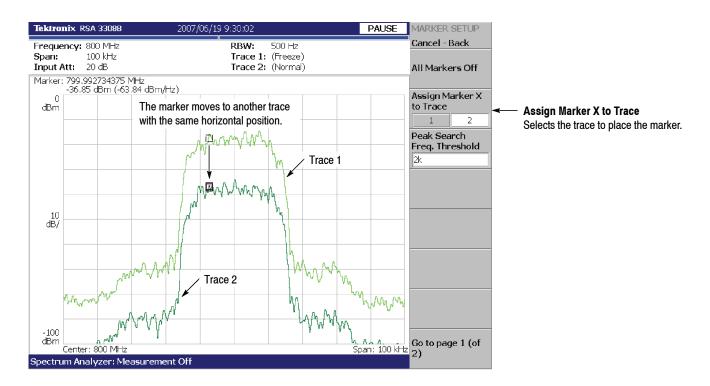

The marker moves to another trace with the same horizontal position.

Figure 3-156: Changing the Trace

#### **Markers in Multiple Views**

When you display multiple views on one screen, markers displayed in those views are locked together. Figure 3-157 shows a concurrent display of the spectrum and the spectrogram. If you move the marker to the left on the spectrum, the marker on the spectrogram moves to the left in accordance with that movement. Conversely, if you move the marker on the spectrogram to the left, the marker on the spectrum moves to the left.


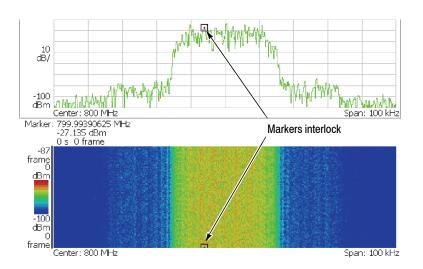


Figure 3-157: Interlocked markers

### Marker Operation in the DPX Spectrum


The marker operation in the DPX spectrum is the same as in the normal spectrum display. But note that when a marker is on the Bitmap trace, it is placed on the most popular amplitude value (see Figure 3-158).



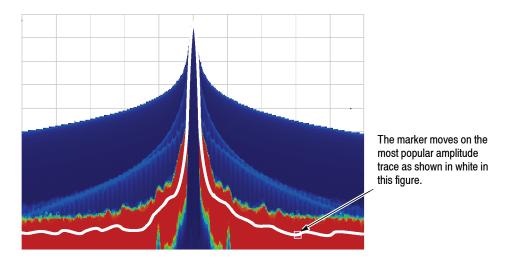

### **Peak Search**

The peak search function searches for the peak on the waveform and moves the marker to the peak position. Use the **Peak** and marker left/right/up/down  $(\blacktriangleleft \blacktriangleright \blacktriangle \bigtriangledown)$  keys on the front panel for peak search.

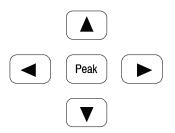

▼


#### Using the Peak Search keys


The peak search keys have the following functions (see Figure 3-159):

**Peak** Positions the marker to the highest peak signal.


Moves the marker lower in frequency to the next signal.

- Moves the marker higher in frequency to the next signal.
- Moves the marker higher in amplitude to the next signal.
- Moves the marker lower in amplitude to the next signal.

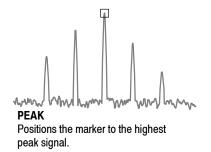


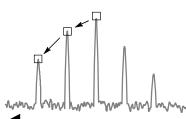


Moves the marker lower in frequency to the next signal.

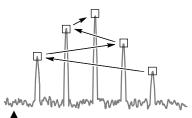

Moves the marker higher in amplitude to the next signal.

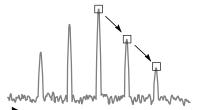

Moves the marker higher in frequency to the next signal.

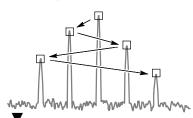

Moves the marker lower in amplitude to the next signal.

#### Figure 3-159: Functions of the peak search keys

#### Setting the Minimum Jump of the Marker


The definition of "next" can be set in the Marker Setup menu using the following parameters:

**Peak Search Freq. Threshold.** *S/A mode only.* Sets the minimum frequency jump that will be made when selecting the next signal to the left/right/up/down.

**Peak Search Hor. Threshold.** *Demod and Time modes only.* Sets the minimum horizontal jump that will be made when selecting the next signal to the left/right/up/down.

For example, when Peak Search Freq. Threshold is set to 1 kHz, the next peak is recognized only if it is 1 kHz or more distant from the first peak (see Figure 3-160).

Example: Peak Search Freq. Threshold = 1 kHz

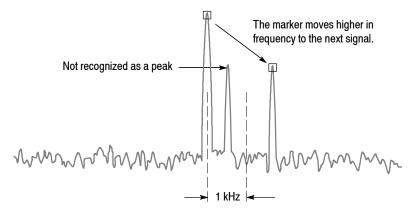

#### Figure 3-160: Setting the minimum frequency jump

Use the following steps to set the minimum jump:

- 1. Press the Marker Setup key on the front panel.
- 2. Press the Go to page 2 (of 2) (bottom) side key to show the next menu page.
- **3.** Press either of the following side keys depending on the measurement mode and then set the minimum jump:
  - Peak Search Freq. Threshold in the S/A mode
  - Peak Search Hor. Threshold in the Demod and Time modes

## **Using the Online Help**

This section describes the analyzer help system and how to access it.

The online help system is integrated as part of the analyzer user interface, using the Windows help system. You can use the front-panel arrow keys (as in the following procedure), a mouse, or a keyboard to navigate through the help system. See page 3-230 for information about using a mouse or keyboard to navigate through the online help.


The help system provides information about operating this analyzer, providing the following help resources online:

- Front-panel key descriptions
- Online user manual

### **Displaying the Online Help**

Use the following procedure to show the online help:

1. Press the **Help** key on the front panel.

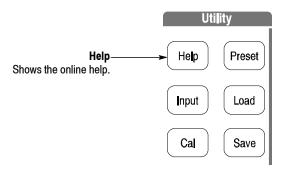

- 2. Select the type of online help using the side keys:
  - View Front Panel Button Help. Provides descriptions of the front-panel keys.
  - View Online User Manual. Displays the analyzer user manual, formatted for online use.

- **3.** Depending on the type of online manual you selected, use one of the following two procedures:
  - If you selected *View Front Panel Button Help:* Press the front-panel key about which you want to receive information. For example, press the **Measure** key to show information about that key (see Figure 3-161).

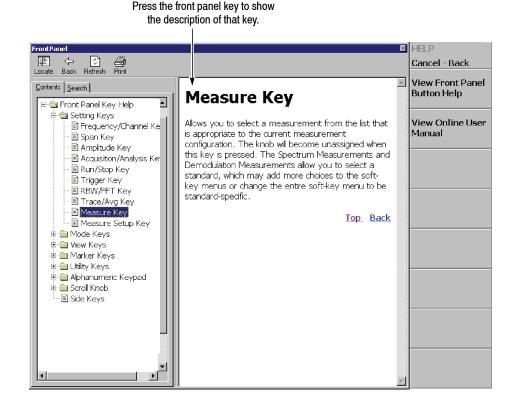


Figure 3-161: Online help for the front panel key

If you selected *View Online User Manual:* Using the front panel keys listed in Table 3-31, select the topic from the contents window on screen (see Figure 3-162). The key functions are different in the contents (left) window than in the descriptions (right) window.

	Function	
Front panel key	Contents window	Descriptions window
General purpose knob	Selects a topic.	Scrolls the display up or down.
Up key (🔺)	Selects a topic on top in the window.	Scrolls the display up.
Down key (▼)	Selects a topic on bottom in the window.	Scrolls the display down.
Markers:	Selects the next topic up.	Scrolls the display up.
Markers: 🔻	Selects the next topic down.	Scrolls the display down.
Markers:	Moves to the next level up.	Scrolls the display to the left.
Markers: Moves to the next level down.		Scrolls the display to the right.
Peak	Switch the window to operate the cursor.	
Marker 🔿	-	Selects a character string with hyperlink.
BKSP	Moves to the next level up.	Returns to the previous display.
Enter	Accepts the selection and displays the corresponding content.	

Table 3-31: Front panel keys for online help

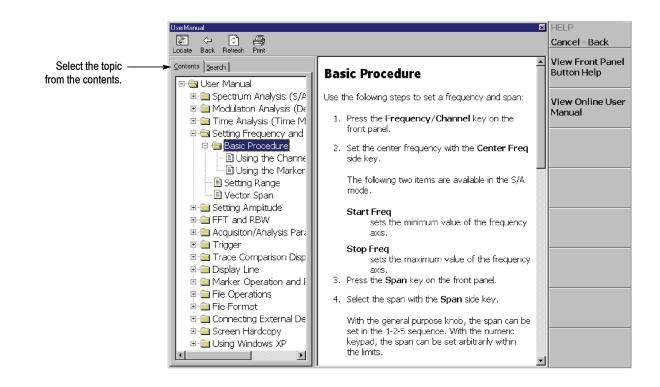


Figure 3-162: Online user manual

4. To exit the help system, press the Cancel-Back (top) side key.

## Using a Mouse and Keyboard

You can control the help system using a mouse and a keyboard instead of the front panel keys. Using a mouse, click the topic in the contents window to show the description. Using a keyboard, select a topic or enter a search word (see Figure 3-163). Use the arrow keys to move the cursor, and the enter key to accept the selection.

Refer to *Side-Panel Interface* on page 2-5 on connecting a USB mouse and keyboard to the analyzer.


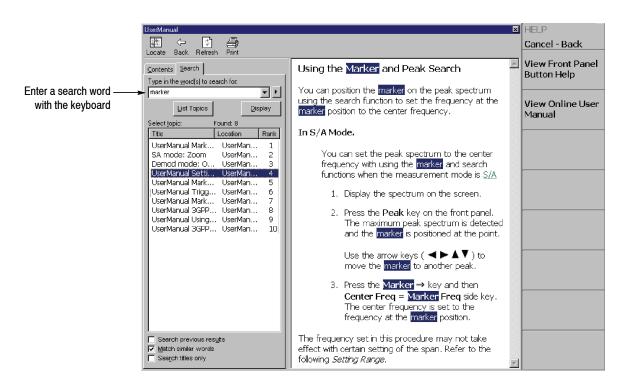


Figure 3-163: Word search using the keyboard

## **Selecting Input Source**

You can select from three input sources in the Input menu: RF, IQ, and Cal.



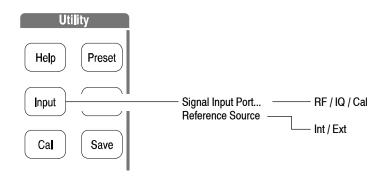

### **Input Menu**

The Input menu contains the following controls:

Signal Input Port... Selects which input is used for signal connection.

- **RF.** *Default.* Uses the signal from the INPUT connector on the front panel. See Figure 2-1 on page 2-2 for the front panel connectors.
- IQ. Option 03 only. Uses the signal from the I INPUT and Q INPUT connectors on the rear panel. See Figure 2-2 on page 2-4 for the rear panel connectors.

**NOTE**. *IQ INPUT* (*Option 03 only*). *The IQ input gain of the analyzer is set in 10 dB steps. To maximize the dynamic range, you may need to adjust the signal level externally, or insert attenuators in the I and Q signal paths.* 

When you calibrate the IQ offset, set the IQ input signal level to zero externally. (The IQ offset calibration is included in the center offset calibration. Refer to Calibrating Center Offset on page 1-26.)

 Cal. Uses the internal 50 MHz calibration signal as the signal source. Refer to *Functional Check* on page 1-16 for using the cal signal to check the system performance. **Reference Source.** Selects the source for reference frequency.

- Int. Uses the internal clock (10 MHz simulated sine wave).
- Ext. Uses the 10 MHz sine wave of -10 to +6 dBm from the REF IN connector on the rear panel when you want to synchronize the analyzer with other instrument(s).

The reference clock is available from the REF OUT connector on the rear panel.

## **File Operations**

The instrument settings and waveform data can be saved in or loaded from files on the hard disk or 3.5 inch disk. This section contains the following topics:

- File Type
- Load/Save Menu
- Saving and Loading Files
- Entering a File Name
- Deleting a File

Making/Deleting a Directory

### File Type

The analyzer can make files with the following extensions:

- **.STA** (Status file) A status file used to save all the current menu settings. You can save the settings you frequently use and reset the analyzer any time by loading the settings from this file.
- IQT (Data file)
   A file used to save time-domain waveform data acquired on the data memory in the Demod (Modulation Analysis) and the Time (Time Analysis) modes.
   I, Q, and T stand for in-phase, quadrature phase, and time domain, respectively.
- .TRC (Trace file)
   A file used to save Trace 1 or 2 in the S/A (Spectrum Analysis) mode. The data in this file is loaded as the reference waveform in comparison displays between Trace 1 and 2, for example.
- .DPT (Bitmap trace file) A file used to save the Bitmap trace in the DPX spectrum display.
- COR (Amplitude correction file)
   A file used in the S/A (Spectrum Analysis) mode to save an amplitude correction table. The data in this file is loaded when the amplitude correction is performed. Refer to page 3-132 for the amplitude correction.

## Load/Save Menu

Figure 3-165 shows the Load/Save menu structure.

Preset		
$\frown$	S/A mode except Real Time	
Load	Load State	── <del>≻</del> Load From File Menu
	Load Trace 1 —	Select File
$\frown$	Load Trace 2 —	Load File Now
Save	_ Load Correction →	Folder
		Load from XXXA
$\frown$	Real Time S/A, Demod and Time modes	Load from XXXB
Print	Load State	Load from XXXC
	Load Data ———	
	S/A mode except Real Time	
	Save State	→ Save To File Menu
	Save Trace 1	Filename
	Save Trace 2 —	Save File Now
	Save Correction —	Overwrite
		Folder
	Real Time S/A mode	Save to XXXA
	Save State	Save to XXXB
	Save Data ———	Save to XXXC
		Where XXX =
	Real Time S/A with Zoom, Demod and Time modes	State for a state file
	Save State	IQData for an IQ data file
	Save Data — All Blocks / Current Block / Current Area —	Trace for a trace file
	Save Data(.CSV) — Current Block / Current Area —	Corr for a correction file
	Save Data(.MAT) ———— Current Block / Current Area ——	Bitmap for a bitmap trace

#### Figure 3-165: Load/Save menu structure

**Load Menu** Loads waveform data or instrument settings from a file.

Load State. Loads the instrument settings.

**Load Data.** *Real Time S/A, Demod, and Time modes only.* Loads waveform data (IQ data in the time domain).

**NOTE**. The following menu items (Load Trace 1, Load Trace 2, and Load Correction) are available in the S/A mode except Real Time S/A.

Load Trace 1. Loads data for Trace 1 of a graph.

Load Trace 2. Loads data for Trace 2 of a graph.

**Load Correction.** Loads a correction table of Frequency/Amplitude pairs used for adjusting amplitude values.

file

**Save Menu** Saves waveform data or instrument settings to a file.

Save State. Saves the instrument settings.

**Save Data.** *Real Time S/A only.* Saves waveform data (IQ data in the time domain) to a data file (.iqt).

**Save Data...** *Demod and Time modes only.* Saves waveform data (IQ data in the time domain) to a data file (.iqt).

- All Blocks. Saves all the blocks.
- Current Block. Saves the block currently displayed.
- **Current Area.** Saves the data displayed in the main view.

**Save Data(.CSV)...** *Real Time S/A with Zoom, Demod and Time modes only.* Saves waveform data (IQ data in the time domain) to a file in the CSV (Comma Separated Values) format, allowing you to export the file into Microsoft Excel or other database systems.

- **Current Block.** Saves the block currently displayed.
- Current Area. Saves the data displayed in the main view.

**Save Data(.MAT)...** Real Time S/A with Zoom, Demod and Time modes only. Saves waveform data (IQ data in the time domain) to a file in the MATLAB format, allowing you to export the file into the MATLAB technical computing environment.

- **Current Block.** Saves the block currently displayed.
- **Current Area.** Saves the data displayed in the main view.

**NOTE**. The data capacity for the .CSV and .MAT files is 256 MB maximum. If exceeding 256 MB, the **Save Data** side key is greyed out to disable saving data.

**NOTE**. The following menu items (Save Trace 1, Save Trace 2, and Save Correction) are available in the S/A mode except Real Time S/A.

Save Trace 1. Saves data for Trace 1 of a graph.

**Save Trace 2.** Saves data for Trace 2 of a graph.

**Save Correction.** Saves a correction table of Frequency/Amplitude pairs used for adjusting amplitude values.

## **Saving and Loading Files**

The following procedures show how to save and load files.

- **Using a Disk** You can use a 2HD (1.44 M bytes) or 2DD (720 K bytes) 3.5-inch disk formatted for MS-DOS. The disk drive is on the side panel. See Figure 2-3 on page 2-5.
- **Saving a File** Use the following steps to save settings or waveforms to a file. There are two ways to save data to a file: using the preset file name and entering a new file name. To return to the previous menu, press the **Cancel Back** side key at any time.
  - 1. Press the Save key on the front panel.
  - 2. Press one of the following side keys, depending on the type of data to save.

Meas. mode Side key		Save content	Extension
S/A	Save State	Current instrument settings	.sta
(except Real Time S/A)	Save Trace 1	Trace 1 waveform	.trc, .dpt 1
. ,	Save Trace 2	Trace 2 waveform	.trc, .dpt 1
	Save Correction	Amplitude correction data	.cor
Real Time S/A	Save State	Current instrument settings	.sta
	Save Data	Acquired data in time domain	.iqt
Real Time S/A	Save State	Current instrument settings	.sta
with Zoom, Demod and Time	Save Data	Acquired data in time domain	.iqt
	All Blocks	Data of all blocks	
	Current Block	Data of the block currently displayed	
	Current Area	Data currently shown in the main view	
	Save Data(.CSV)	Acquired data in time domain	.CSV
	Current Block	Data of the block currently displayed	
	Current Area	Data currently shown in the main view	
	Save Data(.MAT)	Acquired data in time domain	.mat
	Current Block	Data of the block currently displayed	
	Current Area	Data currently shown in the main view	

#### Table 3-32: File save operation

<sup>1</sup> .dpt is only for the Bitmap trace in the DPX spectrum.

**3.** *To use a preset file name*: Press the **Save to** *file name* side key. (To enter a new file name, skip to step 4.)


There are three preset names for each file type as shown in Table 3-33.

Save file name
StateA, StateB, StateC
IQDataA, IQDataB, IQDataC
TraceA, TraceB, TraceC
BitmapA, BitmapB, BitmapC
CorrA, CorrB, CorrC
TableA, TableB, TableC

<sup>1</sup> DPX spectrum only.

The file extension is automatically added according to the file type.


**NOTE**. The preset files are always saved in the My Documents folder on the *C*: drive even if you select another folder.

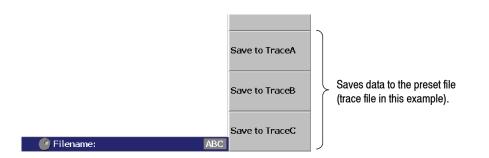

#### Figure 3-166: Save to the preset file (lower right of the screen)

- 4. To enter a new file name:
  - **a.** Press the **Folder...** side key to select the destination folder using the following menu items (see Figure 3-167):
    - +Open Folder opens the selected folder.
    - Select Folder selects a folder.

- -Close Folder closes the selected folder.
- **Done** accepts the selected file.

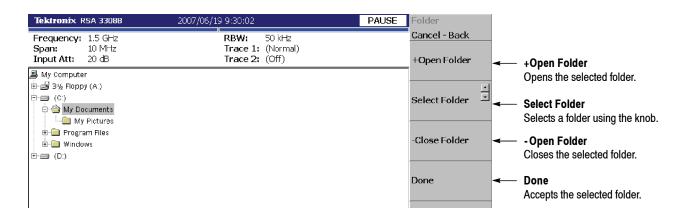

#### Figure 3-167: Selecting the folder

- **b.** After selecting the folder, press the **Done** side key to accept it. The list of files in the selected folder appears, displaying only files with the same extension as the file to be stored.
- **c.** Press the **Filename** side key (see Figure 3-168) and enter a file name using the front-panel keypad. The file extension is automatically added. Refer to *Entering a File Name* on page 3-242 for details about how to enter the name.

Here are instructions to input "TRACE1" for example, using the front-panel keypad.

- Press the TUV (number 2) key once to select the letter "T", and then press the Enter key.
- Press the PQRS (number 1) key three times to select the letter "R", and then press the Enter key.
- Press the ABC (number 8) key once to select the letter "A", and then press the Enter key.
- Press the ABC (number 8) key three times to select the letter "C", and then press the Enter key.
- Press the DEF (number 9) key twice to select the letter "E", and then press the Enter key.
- Press the PQRS (number 1) key five times to select the letter "1", and then press the Enter key.

Tektronix F	ISA 33088	2	007/06/19 9:30:02	PAUSE	Save To File		
Frequency: Span: Input Att:	1.5 GHz 5 MHz 15 dB		Acquisition Length: 320 µs		Cancel - Back Filename	-	- Filename
File name:	Size	Туре	Modified		TRACE1		Enters a file name.
⊠ a.iqt ॼ b.iqt ॼ IQDataA.iqt ॼ IQDataB.iqt	87KB 385KB	IQT ファイル IQT ファイル	07/07/13 09:43 07/07/13 09:43 07/07/13 09:43 07/07/13 09:43		Save File Now	◄	<ul> <li>Save File Now</li> <li>Stores the data to the specified file.</li> </ul>
					Overwrite	•	<ul> <li>Overwrite</li> <li>Overwrites an existing file.</li> </ul>
					Folder (C:¥Documents and)	•	- Folder Shows the selected directory.

#### Figure 3-168: Save to File menu

**d.** After entering the file name, press the **Save File Now** side key. The data is saved to the specified file.

*If the file already exists in the specified folder.* The message "File already exists. Do you want to overwrite it?" appears at the bottom of the screen. When you overwrite the file, press the **Overwrite** side key. Otherwise go back to step c to enter another name.

- **Loading a File** Use the following steps to load settings or waveform from a file. There are two ways to load data from a file: selecting the preset file and selecting a file by name. When you want to return to the previous menu, press the **Cancel Back** side key at any time.
  - 1. Press the Load key on the front panel.
  - 2. Press one of the following side keys according to the type of data to load:

### Table 3-34: File load operation

Meas. mode	Side key	Load content	Extension
S/A	Load State	Instrument settings	.sta
(except Real Time S/A)	Load Trace 1	Trace 1 waveform	.trc, .dpt 1
	Load Trace 2	Trace 2 waveform	.trc, .dpt <sup>1</sup>
	Load Correction	Amplitude correction data	.cor
Real Time S/A,	Load State	Instrument settings	.sta
Demod or Time	Load Data	Acquired data in time domain	.iqt


<sup>1</sup> .dpt is only for the Bitmap trace in the DPX spectrum.

**NOTE**. In the Real Time S/A, Demod, and Time modes, when you abort data acquisition before it is completed, such as in a trigger-armed state, the last data block is empty. Therefore, when you save all the blocks to a file and then load the file, the waveform does not appear at first. If you have acquired multiple blocks, select an old block to display the waveform.

**3.** *To select the preset file:* (To select a file by name, skip to step 4.) To select the preset file to which you have saved data in step 3 on page 3-237, press the **Load from** *file name* to load the data (see Figure 3-169).

**NOTE**. The preset file must be in the My Documents folder on the C: drive.

If you have not saved any data to the preset file that you select, the error message "File name not found" appears on the bottom left of the screen.

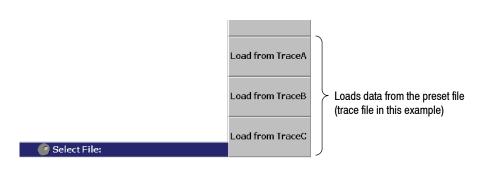




- 4. To select an existing file:
  - **a.** Press the **Folder...** side key. Select the folder containing the file to load. (For a file select procedure, see step 4 on page 3-237.)
  - **b.** After selecting the folder, press the **Done** side key to accept it.

The list of files in the selected folder appears, displaying only files with the same extension as the file to be loaded.

- **c.** Press the **Select File** side key and then select the file from the list using the general purpose knob.
- **d.** After selecting the file, press the **Load File Now** side key. See Figure 3-170.

The data is loaded from the specified file.

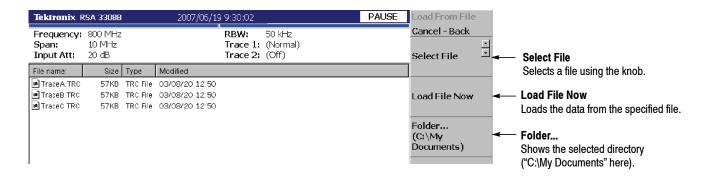

#### Figure 3-170: Load from File menu

**Displaying the File Name on the Load Data Side Key.** *RFID analysis (Option 21) only.* In the RFID analysis, when you load an .iqt file, the file name appears on the **Load Data** side key. Because of the limitation of display area, if the number of characters exceed the limit, only the part within the allowable range will be displayed with line feed inserted at space if any (see Figure 3-171).



File name: "Example.iqt" Displayed normally.



File name: "ExampleOfLongName withSpaces Inserted.iqt.iqt" Only the part within the allowable range is displayed with line feed inserted at space. The third and the subsequent lines are omitted.

Figure 3-171: Displaying the file name on the Load Data side key

# **Entering a File Name**

This section shows how to enter a file name in the **Filename** side key for saving the file using the Save menu.

**NOTE**. You can overwrite an existing file; after entering the file name, press Save File Now  $\rightarrow$  Overwrite.

Use the keypad on the front panel to input a file name (see Figure 3-172).

- 0 to 9 and "." keys are used to enter alphabetic characters, numbers, punctuation and special characters indicated above each key. For example, the 8 key inputs A, B, and C, as well as 8. Each time you press the key, the selected character moves through the list: A → B → C → 8 ...
- **Caps Lock** key toggles a character between uppercase and lowercase.
- **BKSP** (backspace) key deletes the character just before the cursor.
- Enter/Next key accepts the input character.

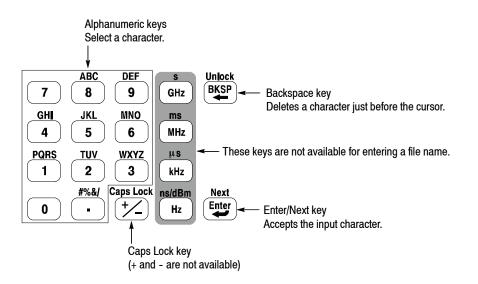


Figure 3-172: Alphanumeric keypad

**Procedure** To enter a new name, SAMPLE1, use the following steps:

1. Press the Filename side key.

The cursor appears on the left end in the text box.

File	name	

2. Press the PQRS key on the keypad four times to enter "S".

Filename	
S	

3. Press the ABC key on the keypad once to enter "A".

Filename
SA

4. Press the MNO key on the keypad once to enter "M".

Filename	
SAM	

5. Enter the remaining characters in the same way.

Filename	
SAMPLE1	

If you make a mistake, you can correct it any time before pressing the **Enter** side key by repeatedly pressing **BKSP** (the backspace key) until the incorrect character disappears and then typing the correct value. If you want to restart from the beginning, press the **Cancel - Back** (top) side key. You can also move the insertion cursor using the general purpose knob or up/down keys.

# **Deleting a File**

You can delete a file using Windows, but not using the analyzer menu. Refer to page 2-26 for information about using Windows on the analyzer. Refer to Windows documentation for additional information about operating Windows.

# Making/Deleting a Directory

You can make or delete a directory using Windows, but not using the analyzer menu. Refer to page 2-26 for information about using Windows on the analyzer. Refer to Windows documentation for additional information about operating Windows.

# **File Format**

This section describes the structure of the following files:


- Data file (\*.IQT)
- Trace file (\*.TRC)
- DPX Bitmap trace (\*.DPT)
- CSV (Comma Separated Values) file (\*.CSV)
- MAT (MATLAB) file (\*.MAT)

# **Data File Format**

File Structure

The data file normally consists of five blocks. The date and time may be added to the end for data logging (see Figure 3-173).

Data file (\*.IQT)

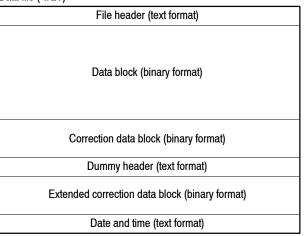

### Figure 3-173: Data file structure

The file is normally created after data acquisition is complete. When logging data is acquired continuously, a data block is added each time the data is acquired.

In data logging, since the internal program adds data blocks while acquiring data, when the program creates the file header, the date and time for the last frame is not known. The program, therefore, adds date and time data to the end of the file. If you examine the file size and find the date and time at the end, use them in place of DateTime in the file header. Refer to *DateTime* on page 3-249 for the format of date and time.

The number of frames ValidFrames (refer to page 3-249) is not known at the time the program creates the file header, so the internal program temporarily writes ValidFrames=0. If the value of ValidFrames of the file header is 0, examine the file size and calculate the true value of ValidFrames. In this case, the correction data block, the dummy header, and the extended correction data block are always added.

The subsections below provide details on each block.

**File Header** The following is an example of the file header. Although Type is always written at the beginning, the order of the other items is not fixed and some new items may be added.

40416Type=RSA3303B or RSA3308B IQT FrameReverse=Off FramePadding=Before Band=RF3 MemoryMode=Zoom FFTPoints=1024 Bins=801 MaxInputLevel=0 LevelOffset=0 CenterFrequency=7.9G FrequencyOffset=0 Span=5M BlockSize=40 ValidFrames=40 FramePeriod=160u UnitPeriod=160u FrameLength=160u DateTime=2005/01/10@13:21:16 GainOffset=-82.3326910626668 MultiFrames=1 MultiAddr=0 IOffset=-0.0475921630859375 QOffset=0.12628173828125

The first character "4" indicates that the number of bytes of the file header is expressed with four characters after the second character.

In the preceding example:

Number of bytes of the file header = 1 + 4 + 0416 = 421

The number of bytes is therefore 421. Data starts from the 422<sup>nd</sup> byte.

Descriptions of the file header items follow.

**Type.** Shows the type of data. In the RSA3303B or RSA3308B analyzer, only this type is used:

RSA3303B or RSA3308BIQT. The data block contains I and Q values in the time domain. It is upward compatible with the \*.IQT file of the WCA300 series.

**FrameReverse.** Shows the frame order. This item is always Off in the RSA3303B or RSA3308B analyzer. The following parameters are provided:

- Off. Frames are stored in the order of acquisition. The last frame in the data block is the latest acquired frame.
- **On.** Frames are stored in the reverse order of acquisition. The first frame in the data block is the latest acquired frame.

**FramePadding.** The analyzer adds dummy frames when acquired frames do not fill BlockSize (block size). Figure 3-174 shows the dummy frames. The FramePadding parameter is always set to Before in the RSA3303B or RSA3308B analyzer.

- **Before.** Adds dummy frames before valid frames, but not in the first block.
- After. Adds dummy frames after valid frames, but not in the last block.

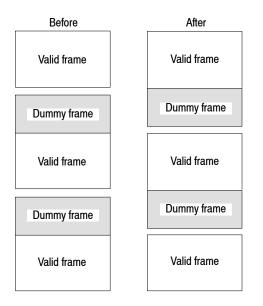


Figure 3-174: Adding dummy frames

**Band.** Shows the frequency band that was set when the analyzer acquired the data. It is necessary only when the analyzer reloads the data.

**MemoryMode.** Shows the memory mode that was set when the analyzer acquired the data. It is necessary only when the analyzer reloads the data.

**FFTPoints.** Shows the FFT point setting that was set when the analyzer acquired the data. This item is always 1024 in the RSA3303B or RSA3308B analyzer.

**Bins.** Shows the number of bins. The same value is also entered in *bins* of each frame header of the data block (refer to *Frame Header* on page 3–250).

**MaxInputLevel.** Shows the reference level setting in dBm when the analyzer acquired the data.

**LevelOffset.** Shows the level offset setting in dB when the analyzer acquired the data.

**CenterFrequency.** Shows the center frequency setting in Hz when the analyzer acquired the data.

**FrequencyOffset.** Shows the frequency offset setting in Hz when the analyzer acquired the data.

Span. Shows the span setting in Hz when the analyzer acquired the data.

BlockSize. Shows the block size setting when the analyzer acquired the data.

**ValidFrames.** Shows the number of frames in the data block. This value divided by MultiFrames represents the number of frames that are scanned and synthesized into one frame.

**FramePeriod.** Shows the frame period setting in seconds. The actual period is obtained by multiplying UnitPeriod by the difference of time stamp ticks of each frame in the data block.

**UnitPeriod.** Shows the unit time in seconds of time stamp Ticks of each frame in the data block.

FrameLength. Shows the time in second necessary to acquire one frame.

**DateTime.** Shows the time when the analyzer acquired the last frame in a data block. You can change "@" to "" (space) for display because files may have many "@" characters.

**GainOffset.** Shows the gain offset. It is used for calculating the amplitude (refer to page 3-252).

**MultiFrames.** Shows the number of frames in multi-frame mode. For example, when MultiFrames = 20, scanning 20 times with a span of 5 MHz makes a span of 100 MHz.

**MultiAddr.** Shows the last frame address in the multi-frame mode. The range is 0 to MultiFrames -1. MultiFrames -1 indicates that the data ends just at the end of scans.

**IOffset.** Shows the offset value of I data. It is used for calculating the data value (refer to page 3-252).

**QOffset.** Shows the offset value of Q data. It is used for calculating the data value (refer to page 3-252).

**Data Block** Each data block contains several pairs of frame header and frame data (see Figure 3-175). The number of pairs in a block is indicated by ValidFrames in the file header. The frame order is determined by FrameReverse in the file header.

Frame header	Frame data		
Frame header	Frame data		
		l	≻ Number indicated by ValidFrames
Frame header	Frame data		
Frame header	Frame data		

### Figure 3-175: Data block

**Frame Header** The frame header is defined by the following structure:

struct frameHeader\_st {
 short reserved1;
 short validA;
 short validP;
 short validI;
 short validQ;
 short validQ;
 short bins;
 short reserved2;
 short triggered;
 short overLoad;
 short lastFrame;
 unsigned long ticks;
};

Descriptions of the frame header items follow.

short reserved1. Internal use only.

**short validA, validP, validI, validQ.** These parameters indicate whether data of amplitude, phase, I or Q (value converted into a two-byte signed integer) is written.

- 0. Indicates that data is not written in the file.
- **-1.** Indicates that data is written in the file.

There are seven combinations as listed in Table 3-35.

validA	validP	validl	validQ
0	0	0	0
-1	0	0	0
0	-1	0	0
-1	-1	0	0
0	0	-1	0
0	0	0	-1
0	0	-1	-1

Table 3-35: Combinations of valid A, P, I and Q

short bins. Shows the number of bins. It is the same as Bins in the file header.
short reserved2. Internal use only.

short triggered. Indicates whether the frame is before or after the trigger.

- 0. Indicates that the frame is before the trigger (pre-trigger).
- Indicates that the frame is after the trigger (post-trigger).

short overLoad. Indicates whether an input overload occurred.

- 0. Indicates that the MaxInputLevel value in the file header was correct.
- -1. Indicates that the MaxInputLevel value in the file header was too low.

short lastFrame. The analyzer can divide its memory into  $(100 \text{ frames}) \times (40 \text{ blocks})$ , for example. LastFrame indicates the last frame in a block.

- 0. Indicates that the frame is not the last in the block.
- **-1.** Indicates that the frame is the last in the block.

**unsigned long ticks.** Shows a time stamp using UnitPeriod in the file header (not FramePeriod) as a unit time.

```
Frame Data
                        A frame contains 1024 pairs of I and Q data in the time domain in order of data
                        acquisition.
                        Definition of Bin. The bin is defined by the following structure:
                            struct iqBin_st {
                                 short q;
                                 short i;
                            };
                        Definition of Frame. The frame is defined by the following structure:
                            struct iqFrame1024 st {
                                 struct iqBin st iq[1024];
                            };
Calculation of Data
                        All the data of amplitude, phase, I and Q are transformed to 2-byte signed
                        integers, and then written in the file.
                        Amplitude. The amplitude is calculated using i and q in the IQT file using these
                        formulas:
                            i = i - IOffset
                            q = q - QOffset
                            Amplitude = 10 * \ln(i * i + q * q) / \ln(10)
                                             + GainOffset + MaxInputLevel + LevelOffset [dBm]
                        Phase. The phase is calculated using i and q in the IQT file using these formulas:
                            i = i - IOffset
                            q = q - QOffset
                            Phase = atan2(q, i) * (180 / Pi) [degree]
                        I, Q. I and Q are calculated using i and q in the IQT file using these formulas.
                            i = i - IOffset
                            q = q - QOffset
                            IQScale = Sqrt(Power(10, (GainOffset + MaxInputLevel
                                             + LevelOffset) / 10) / 20 * 2)
                            I = i * IQScale [V]
                            Q = q * IQScale [V]
```

### Correction and Extended Correction Data Blocks

The correction data block contains gain and phase correction data as one frame of data block in the frequency domain. The extended correction data block contains the next lower byte of the correction data.

When these blocks are added, the amplitude and phase are calculated using the following formulas. Be careful about the sign for phase correction.

Amplitude [dBm] = (Original data) - (((Gain correction data << 8) | (Extended gain correction data & 0x000000ff)) / (128 \* 256))

Phase [degrees] = (Original data) + (((Phase correction data << 8) | (Extended phase correction data & 0x00000ff)) / (128 \* 256))

**Definition of Bin.** The bin is defined by the following structure:

```
struct apBin_st {
    short a;
    short p;
};
```

**Definition of Frame.** The frame is defined by the following structure:

```
struct apFrame1024_st {
    struct apBin_st ap[1024];
};
```


**Definition of Extended Correction Data Block.** The extended correction data block is defined by the following structure:

```
struct extendedCorrectionData_st {
    unsigned char a[1024];
    unsigned char p[1024];
};
```

**Dummy Header** The dummy header contains the ASCII codes of "40000" with no line feed code.

# **Trace File Format**

File Structure The trace file have two blocks in the text format as shown in Figure 3-176.

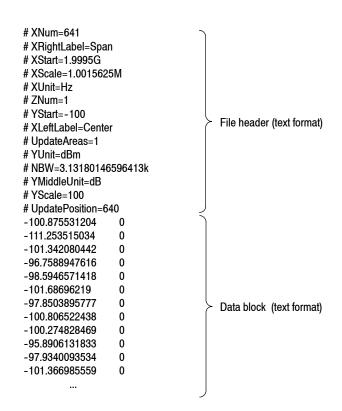



**File Header** The file header contains the following information. *"Required"* means that the item is prerequisite for displaying waveform.

XNum. Required. Represents the number of data in the data block.

**XRightLabel.** Represents the label displayed on the right of the horizontal axis: Span (span) when XLeftLabel is Center (center frequency), or Stop (stop frequency) when XLeftLabel is Start (start frequency).

XStart. Required. Represents the minimum (left) edge of the horizontal axis.

**XScale.** *Required.* Represents the scale of the horizontal axis.

XUnit. Represents the unit of the horizontal axis: fixed to Hz.

**ZNum.** *Required.* Internal use only: fixed to 1.

YStart. Required. Represents the minimum (bottom) edge of the vertical axis.

**XLeftLabel.** Represents the label on the left of the horizontal axis: Center (center frequency) when XLeftLabel is Span (span), or Start (start frequency) when XLeftLabel is Stop (stop frequency).

UpdateAreas. Internal use only.

YUnit. Represents the unit of the vertical axis: dBm, dBµV, V, or W.

**NBW.** Specifies NBW (Noise Bandwidth) when setting FFT parameters, or RBW (Resolution Bandwidth) when setting RBW parameters.

**YMiddleUnit.** Represents the unit of the vertical axis scale. Specify dB when YUnit is dBm. For all other values of YUnit, the YUnit value is used.

**YScale.** *Required.* Represents the scale of the vertical axis.

UpdatePosition. Internal use only.

**Data Block** Pairs of a power value and a mask value are written in order, with a tab between the power and the mask (see the figure below). The number of lines is indicated by XNum in the file header.

Power 1	(tab)	Mask 1	)
Power 2	(tab)	Mask 2	
Power 3	(tab)	Mask 3	
			$\succ$ Number indicated by XNum
	•		
	•		
Power N	(tab)	Mask N	J

Part of the data block might look like this:

-100.875531204	0
-111.253515034	0
-101.342080442	0
-96.7588947616	0
-98.5946571418	0

For example, the first line indicates that the power is -100.875531204 dBm and the mask value is 0.

Mask Value. The mask value represents whether to display the data.

- 0. Indicates that the data is displayed.
- -1. Indicates that the data is not displayed.

# **Bitmap Trace File Format**

**File Structure** The Bitmap trace file have two blocks in the text format as shown in Figure 3-177.

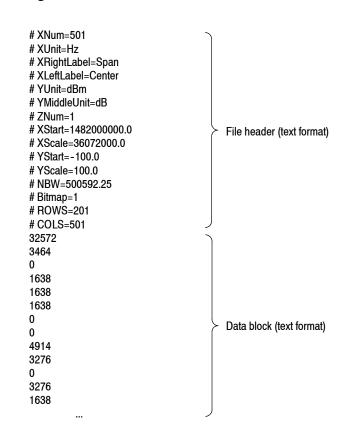


Figure 3-177: Bitmap trace file (.DPT) structure

**File Header** The file header contains the following information. *"Required"* means that the item is prerequisite for displaying waveform.

XNum. Required. Fixed to 501 (= COLS).

XUnit. Represents the unit of the horizontal axis: fixed to Hz.

**XRightLabel.** Represents the label displayed on the right of the horizontal axis: Span (span) when XLeftLabel is Center (center frequency), or Stop (stop frequency) when XLeftLabel is Start (start frequency).

**XLeftLabel.** Represents the label on the left of the horizontal axis: Center (center frequency) when XLeftLabel is Span (span), or Start (start frequency) when XLeftLabel is Stop (stop frequency). YUnit. Represents the unit of the vertical axis: fixed to dBm.

YMiddleUnit. Represents the unit of the vertical axis scale: fixed to dB.

**ZNum.** *Required.* Internal use only: fixed to 1.

XStart. Required. Represents the minimum (left) edge of the horizontal axis.

**XScale.** *Required.* Represents the scale of the horizontal axis.

YStart. Required. Represents the minimum (bottom) edge of the vertical axis.

YScale. Required. Represents the scale of the vertical axis.


NBW. Specifies RBW (Resolution Bandwidth).

Bitmap. Required. Indicates that the trace is the bitmap: fixed to 1.

**ROWS.** Represents the number of data points along the vertical axis on the Bitmap display: fixed to 201.

**COLS.** Represents the number of data points along the horizontal axis on the Bitmap display: fixed to 501.

**Data Block** The normalized number of occurrences is written in integer values for each data point in the order shown in Figure 3-178: from bottom up for each column from left to right on the Bitmap display.

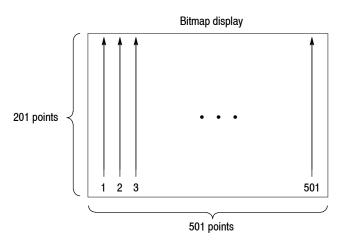



# **CSV File Format**

The acquisition data files saved in CSV format have two sections. At the top of the file is the following header information:

SamplingFrequency. The sampling frequency in Hertz.

NumberSamples. The number of IQ samples in the file.

DateTime. When the data was acquired.

Frequency. The center frequency in Hertz.

AcquisitionBandwidth. The acquisition bandwidth in Hertz.

Following the AcquisitionBandwidth are the data. The data are IQ pairs, in Volts.

# **MAT File Format**

The acquisition data files saved in MatLab format contain the following MatLab variables:

InputCenter. The center frequency in Hertz.

**XDelta.** The sample period in seconds.

Y. A complex array containing IQ pairs.

**InputZoom.** Always 1, and indicates complex data.


# **Screen Copy**

This section explains how to send a screen copy to a printer or a file. Use a printer that complies with USB specifications. For files, data in the bitmap format (.BMP) are created. The following topics are described in this section:

- Print Menu
- Printing a Screen Copy
- Creating a File

# **Print Menu**

The Print menu contains the following controls:

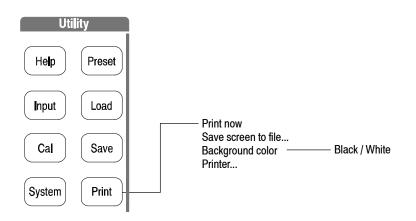

### Figure 3-179: Print menu structure

Print now. Starts printing a copy of the analyzer screen to an attached printer.

**Save screen to file...** Opens the Save menu to save a bitmap file. For file operations, refer to page 3–233.

Background color. Selects the background color for printing.

- **Black.** Prints the background of screen in black.
- White. Reverses the background of screen to white.

Printer... Selects a destination printer.

# **Printing a Screen Copy**

You must connect the USB printer to the analyzer and install the printer driver before you can print.

**Connecting a Printer** Connect the USB cable of the printer to the USB port on the side panel of the analyzer. You can use either of two ports. Refer to *Side-Panel Interface* on page 2–5 for connecting the USB devices.

If the analyzer is connected to a network through a LAN interface, you can also use a network printer.

# **Installing Printer Driver** Install the printer driver using the Windows XP printer wizard. Refer to the printer manual for the installation method. Refer to page 2–26 for using Windows XP.

- **Print** To print a copy of the analyzer screen, follow these steps:
  - 1. Press the **Print** key on the front panel, and select **Black** or **White** as the background color with the **Background color** side key.
  - 2. Display the measurement to be printed and stop data acquisition.
  - **3.** Press the **Print** key again. The screen image data is captured on the internal memory with the background color set in step 1.
  - 4. Press the **Printer...** side key and select a destination printer.
  - 5. Press the **Print now** side key to start printing.

If the printer driver displays a message during printing, follow the instructions.

# **Creating a File**

You can save a screen image to a file as bitmap data and read the data into a word processor to create a report, for example.

To save the file to a floppy disk, insert it to the floppy disk drive.

- 1. Press the **Print** key on the front panel, and select **Black** or **White** of the background color with the **Background color** side key.
- 2. Display the measurement to be copied and stop data acquisition.
- **3.** Press the **Print** key again. The screen image data is captured on the internal memory with the background color set in step 1.
- 4. Press the Save screen to file... side key. The file selection screen appears.
- 5. Do one of the following to output the screen image to a specified file:
  - Using a preset file name: Press the Save to BitmapA, Save to BitmapB, or Save to BitmapC side key. The image data is stored to the file BitmapA.bmp, BitmapB.bmp, or BitmapC.bmp in the C:\My Documents folder, respectively.
  - Entering a new file name: Enter a file name with the Filename side key, referring to Entering a File Name on page 3-242. The file extension
     .BMP is automatically assigned. When you have finished entering the name, press the Save File Now side key to store the copy in the file.

After the screen copy is saved in the file, the screen returns to the original waveform display.

Screen Copy

# **Appendices**

# **Appendix A: Specifications**

This appendix lists the electrical, physical, and environmental characteristics of the analyzer, and specifies the performance requirements for those characteristics. The specifications are common to RSA3303B and RSA3308B, unless otherwise noted.

Unless otherwise stated, the following tables of electrical characteristics and features apply to the analyzer after a 20 minute warm-up period (within the environmental limits) and after all calibration procedures have been carried out.

# **Electrical Characteristics**

### Table A-1: Frequency

Characteristics	Description			
Measurement frequency				
Frequency range	Baseband: DC to 20 MHz           RF:         15 MHz to 3 GHz         (RSA3303B)           RF1:         15 MHz to 3.5 GHz         (RSA3308B)           RF2:         3.5 GHz to 6.5 GHz         (RSA3308B)           RF3:         5 GHz to 8 GHz         (RSA3308B)			
Center frequency setting resolution	0.1 Hz			
Frequency marker readout accuracy	$ \begin{array}{ll} \text{Baseband:} & \pm(\text{RE}\times\text{MF}+0.001\times\text{Span}+0.2)\ \text{Hz}\\ \text{RF, RF1 to 3:} & \pm(\text{RE}\times\text{MF}+0.001\times\text{Span}+2)\ \text{Hz}\\ \text{RE: Reference frequency error; MF: Marker frequency [Hz]} \end{array} $			
At specified frequency	Baseband, Frequency = 10 MHz, Span = 1 MHz Marker: ±1 kHz; Carrier frequency measurement: ±1.2 Hz RF/RF1, Frequency = 2 GHz, Span = 1 MHz Marker: ±1.2 kHz; Carrier frequency measurement: ±210 Hz RF2, Frequency = 5 GHz, Span = 1 MHz (RSA3308B only) Marker: ±1.5 kHz; Carrier frequency measurement: ±510 Hz RF3, Frequency = 7 GHz, Span = 1 MHz (RSA3308B only) Marker: ±1.7 kHz; Carrier frequency measurement: ±710 Hz			
Residual FM (Typical)	2 Нz р-р			
Span accuracy	±1 bin			
RBW filter bandwidth	0.1%			
Reference frequency				
Aging per day	$1 \times 10^{-9}$ (after 30 days of operation)			
Aging per year	$1 \times 10^{-7}$ (after 30 days of operation)			
Temperature drift	$1 \times 10^{-7}$ (10 to 40 °C)			
Total frequency error	$2 \times 10^{-7}$ (within one year after calibration)			
Reference output level	>0 dBm			
External reference input	10 MHz, -10 to +6 dBm, Spurious level <-80 dBc within 100 kHz offset			

### Table A-2: Spectrum purity

Characteristics	Description
Spectrum purity (Frequency = 1500 MHz)	
Carrier offset = 10 kHz, Span = 100 kHz	-100 dBc/Hz
Carrier offset = 100 kHz, Span = 1 MHz	-105 dBc/Hz
Carrier offset = 1 MHz, Span = 5 MHz	-125 dBc/Hz

### Table A-3: Noise sideband

Characteristics	Description		
	Noise sideband	Offset	
Frequency = 1000 MHz	-100 dBc/Hz -105 dBc/Hz -105 dBc/Hz -105 dBc/Hz -112 dBc/Hz -132 dBc/Hz -135 dBc/Hz -135 dBc/Hz	1 kHz 10 kHz 20 kHz 30 kHz 100 kHz 1 MHz 5 MHz 7 MHz	
Frequency = 2000 MHz	-96 dBc/Hz -104 dBc/Hz -105 dBc/Hz -105 dBc/Hz -112 dBc/Hz -132 dBc/Hz -135 dBc/Hz -135 dBc/Hz	1 kHz 10 kHz 20 kHz 30 kHz 100 kHz 1 MHz 5 MHz 7 MHz	
Frequency = 6000 MHz (RSA3308B only)	-87 dBc/Hz -104 dBc/Hz -105 dBc/Hz -105 dBc/Hz -112 dBc/Hz -128 dBc/Hz -130 dBc/Hz -130 dBc/Hz	1 kHz 10 kHz 20 kHz 30 kHz 100 kHz 1 MHz 5 MHz 7 MHz	
Typical	Noise sideband	Offset	
Frequency = 1000 MHz	-103 dBc/Hz -108 dBc/Hz -108 dBc/Hz -108 dBc/Hz -108 dBc/Hz -115 dBc/Hz -135 dBc/Hz -138 dBc/Hz -138 dBc/Hz	1 kHz 10 kHz 20 kHz 30 kHz 100 kHz 1 MHz 5 MHz 7 MHz	

## Table A-3: Noise sideband (Cont.)

Characteristics	Description		
Frequency = 2000 MHz	-99 dBc/Hz -107 dBc/Hz -108 dBc/Hz -108 dBc/Hz -115 dBc/Hz -135 dBc/Hz -138 dBc/Hz -138 dBc/Hz	1 kHz 10 kHz 20 kHz 30 kHz 100 kHz 1 MHz 5 MHz 7 MHz	
Frequency = 6000 MHz (RSA3308B only)	-90 dBc/Hz -107 dBc/Hz -108 dBc/Hz -108 dBc/Hz -115 dBc/Hz -131 dBc/Hz -133 dBc/Hz -133 dBc/Hz	1 kHz 10 kHz 20 kHz 30 kHz 100 kHz 1 MHz 5 MHz 7 MHz	

# Table A-4: Input

Characteristics	Description	
Signal input		
Input connector	N type (RF and baseband input); BNC type (Option 03 IQ input)	
Input impedance	50 Ω	
VSWR	<1.4 (2.5 GHz, RF attenuation≥10 dB) <1.8 (7.5 GHz, RF attenuation≥10 dB, RSA3308B only)	
Typical	<1.4 (300 kHz to 10 MHz, RF attenuation≥10 dB) <1.3 (10 MHz to 3 GHz, RF attenuation≥10 dB) <1.4 (3 GHz to 8 GHz, RF attenuation≥10 dB, RSA3308B only)	
Maximum input level		
Maximum DC voltage	±0.2 V (RF, RF1 to 3) ±5 V (Baseband) ±5 V (Option 03 IQ input)	
Maximum input power	+30 dBm (RF, RF1 to 3, RF attenuation≥10 dB)	
nput attenuator		
RF/Baseband attenuator	0 to 50 dB (2 dB step in RF/RF1/Baseband; 10 dB step in RF2/RF3)	
I/Q attenuator	0 to 30 dB (10 dB step)	

## Table A-5: Amplitude

Characteristics	Description	
Reference level		
Reference level setting range	-30 to +20 dBm (2 dB step, Baseband) -51 to +30 dBm (1 dB step, RF/RF1) -50 to +30 dBm (1 dB step, RF2/RF3, RSA3308B only) -10 to +20 dBm (10 dB step, Option 03 IQ input)	
Accuracy (-10 to -50 dBm)	$\pm$ 0.2 dB (at 50 MHz, 10 dB attenuation, 20 to 30 °C)	
Frequency response (RF attenuation≥10 dB)		
at 20 to 30 °C	±0.5 dB (Baseband) ±1.2 dB (RF/RF1) ±1.7 dB (RF2, RSA3308B only) ±1.7 dB (RF3, RSA3308B only)	
Typical	±0.3 dB (100 kHz to 20 MHz) ±0.5 dB (10 MHz to 3 GHz) ±0.5 dB (10 MHz to 3.5 GHz, RSA3308B only) ±1.0 dB (3.5 GHz to 6.5 GHz, RSA3308B only) ±1.0 dB (5 GHz to 8 GHz, RSA3308B only)	
at 10 to 40 °C	±0.7 dB (Baseband) ±1.5 dB (RF/RF1) ±2.0 dB (RF2, RSA3308B only) ±2.0 dB (RF3, RSA3308B only)	
Absolute amplitude accuracy at calibration point (0 dB attenuation, 20 to 30 °C)	±0.3 dB (Baseband at 10 MHz, -10 dBm signal) ±0.5 dB (at 50 MHz, -20 dBm signal)	
Input attenuator setting uncertainty	±0.5 dB (at 50 MHz)	
Level linearity in display range	±0.2 dB (0 to -40 dBfs) ±0.2 dB (0 to -50 dBfs, <i>Typical</i> ) ±0.12 dB (0 to -50 dBfs, <i>Typical</i> )	
Channel power measurement accuracy for W-CDMA signal at 20 to 30 °C ( <i>Typical</i> )	$\pm 0.6$ dB (Signal frequency: 1900 to 2200 MHz; Signal power: +10 to -30 dBm; RF attenuation: 0 to 20 dB; After Auto Level operation at 10 MHz span)	

# Table A-6: Dynamic range

Characteristics	Description
1 dB compression input	+2 dBm (RF attenuation = 0 dB, 2 GHz)
Second harmonic distortion (-30 dBm tone at input mixer)	-56 dBc (10 MHz to 1500 MHz) -56 dBc (10 MHz to 1750 MHz, RSA3308B only)
3rd order intermodulation distortion (Reference	e level = +5 dBm, RF attenuation = Adjusted, Total signal power = -7 dBm)
Center frequency = 2 GHz	-74 dBc
100 MHz to 3 GHz	-74 dBc
3 GHz to 8 GHz (RSA3308B only)	-72 dBc

## Table A-6: Dynamic range (Cont.)

Characteristics	Description
Displayed average noise level	-144 dBm/Hz (1 kHz to 10 kHz) -151 dBm/Hz (10 kHz to 10 MHz) -151 dBm/Hz (10 MHz to 100 MHz) -151 dBm/Hz (100 MHz to 1 GHz) -150 dBm/Hz (1 GHz to 2 GHz) -150 dBm/Hz (2 GHz to 3 GHz) -142 dBm/Hz (3 GHz to 5 GHz) -142 dBm/Hz (5 GHz to 8 GHz)
ACLR (3GPP down link, test model 1, 16 ch)	60 dB (5 MHz offset) 63 dB (10 MHz offset)
Typical	66 dB (5 MHz offset) 70 dB (10 MHz offset)
Local feed-through to input connecter (Typical)	-40 dBm (local frequency 4.2 to 5 GHz) -55 dBm (local frequency 5 to 6 GHz) -60 dBm (local frequency 6 to 7 GHz) -60 dBm (local frequency 7 to 7.7 GHz, RSA3308B only)

# Table A-7: Spurious response

Characteristics	Description		
Image Suppression (Typical)	mage Suppression (Typical)		
1 <sup>st</sup> IF	75 dB (RF/RF1) 70 dB (RF2/RF3, RSA3308B only)		
2 <sup>nd</sup> and 3 <sup>rd</sup> IF	80 dB (RF/RF1 ) 75 dB (RF2/RF3, RSA3308B only)		
Alias suppression (Typical)	65 dB (Baseband)		
Residual response (Reference level = -30 dBm,	RBW = 100 kHz)		
Baseband, 1 to 20 MHz	-93 dBm (Span = 20 MHz)		
RF, 0.5 to 3 GHz (RSA3303B)	-90 dBm (Span = 2.5 GHz)		
RF1, 0.5 to 3.5 GHz (RSA3308B)	-90 dBm (Span = 3 GHz)		
RF2, 3.5 to 6.5 GHz (RSA3308B)	-85 dBm (Span = 3 GHz)		
RF3, 5 to 8 GHz (RSA3308B)	-85 dBm (Span = 3 GHz)		
Spurious response (Span = 10 MHz, Reference	level = 0 dBm, RBW = 50 kHz)		
Baseband, 10 MHz	-73 dBc (Signal frequency = 10 MHz, Signal level = -5 dBm)		
RF/RF1, 2 GHz	-73 dBc (Signal frequency = 2 GHz, Signal level = -5 dBm)		
RF2, 5 GHz (RSA3308B)	-70 dBc (Signal frequency = 5 GHz, Signal level = -5 dBm)		
RF3, 7 GHz (RSA3308B)	-70 dBc (Signal frequency = 7 GHz, Signal level = -5 dBm)		

### Table A-8: Acquisition

Characteristics	Description	
Acquisition mode	Single and Continuous	
Acquisition memory size	64 MB (Standard); 256 MB (Option 02)	
Block size	1 to 16,000 frames (Standard); 1 to 64,000 frames (Option 02)	
Number of data samples in one frame	1024 (Vector mode)	
A/D converter	14 bits, 51.2 Msps	
Vector span	20 MHz (Baseband), 15 MHz (RF, RF1 to 3), 20 MHz (Option 03 IQ input)	
Real-time capture bandwidth	20 MHz (Baseband), 15 MHz (RF, RF1 to 3), 10 MHz (Option 03 IQ input)	

## Table A-9: Sampling rate

Characteristics	Description
Sampling rate (Real Time S/A, Demod, and Tim	e modes)
20 MHz span (Baseband)	25.6 Msps
15 MHz span (RF, RF1 to 3)	25.6 Msps
10 MHz span	12.8 Msps
5 MHz span	6.4 Msps
2 MHz span	3.2 Msps
1 MHz span	1.6 Msps
500 kHz span	800 ksps
200 kHz span	320 ksps
100 kHz span	160 ksps
50 kHz span	80 ksps
20 kHz span	32 ksps
10 kHz span	16 ksps
5 kHz span	8 ksps
2 kHz span	3.2 ksps
1 kHz span	1.6 ksps
500 Hz span	800 sps
200 Hz span	320 sps
100 Hz span	160 sps

### Table A-10: Frame time

Characteristics	Description
rame time (Real Time S/A, Demod, and	d Time modes)
20 MHz span (Baseband)	40 µs
15 MHz span (RF, RF1 to 3)	40 µs
10 MHz span	80 µs
5 MHz span	160 μs
2 MHz span	320 μs
1 MHz span	640 μs
500 kHz span	1.28 ms
200 kHz span	3.2 ms
100 kHz span	6.4 ms
50 kHz span	12.8 ms
20 kHz span	32 ms
10 kHz span	64 ms
5 kHz span	128 ms
2 kHz span	320 ms
1 kHz span	640 ms
500 Hz span	1.28 s
200 Hz span	3.2 s
100 Hz span	6.4 s

# Table A-11: Analog demodulation

Characteristics	Description
Accuracy (Typical)	
AM demodulation	±2% (-10 dBfs input at center, 10 to 60% modulation depth)
PM demodulation	±3° (-10 dBfs input at center)
FM demodulation	±1% of span (-10 dBfs input at center)

## Table A-12: RBW (Resolution Bandwidth)

Characteristics	Description
Filter shape	Gaussian, Rectangle, Root Nyquist
Setting range	1 Hz to 10 MHz
Minimum Resolution Bandwidth (S/A m	ode)
>2 GHz span	100 kHz
1 to 1.99 GHz span	50 kHz
500 to 990 MHz span	20 kHz
200 to 490 MHz span	10 kHz
100 to 190 MHz span	10 kHz
50 to 90 MHz span	10 kHz
20 to 40 MHz span	10 kHz
10 MHz span	1 kHz
5 MHz span	1 kHz
2 MHz span	1 kHz
1 MHz span	1 kHz
500 kHz span	500 Hz
200 kHz span	200 Hz
100 kHz span	100 Hz
50 kHz span	50 Hz
20 kHz span	20 Hz
10 kHz span	10 Hz
5 kHz span	5 Hz
2 kHz span	2 Hz
1 kHz span	1 Hz
500 Hz span	1 Hz
200 Hz span	1 Hz
100 Hz span	1 Hz

## Table A-13: Trigger

Characteristics	Description
Trigger mode	Free run, Triggered
Trigger event source	Power (Span BW), Frequency mask (Option 02), External
Pre/Post trigger setting	Trigger position can be set from 0 to 100% of total data length.
Power trigger	0 to -40 dBfs
Frequency mask trigger (Option 02)	
Mask resolution	1 bin
Level range	0 to -60 dBfs
Real-time event detection bandwidth	Baseband: 20 MHz (1024 point FFT, 50% overlapping) RF: 15 MHz (1024 point FFT, 50% overlapping)
External trigger	
Threshold voltage	-1.5 V to +1.5 V adjustable
Setting resolution	0.1 V
Input impedance	>2 kΩ
Trigger output voltage	High: >2.0 V, Low: <0.4 V (output current <1 mA)
Trigger marker position timing uncertainty	±2 sample points (Power/External trigger)

### Table A-14: Measurement function

Characteristics	Description
S/A mode	Noise power, Channel power, Adjacent channel power ratio, Occupied bandwidth, Emission bandwidth, Carrier to Noise ratio, Carrier frequency, Spurious
Demod mode	Analog demodulation (AM, FM, PM), Digital demodulation (refer to Table A-15), I/Q vs. Time, Power vs. Time, Frequency vs. Time
Time mode	I/Q vs. Time, Power vs. Time, Frequency vs. Time, CCDF, Pulse measurement
Pulse measurement	Pulse width, Pulse peak power, On/Off ratio, Pulse ripple, Pulse repetition interval, Duty cycle, Pulse-Pulse phase, Channel power, OBW, EBW, Frequency deviation
Pulse length	Minimum 20 samples; Maximum 260,000 samples

Characteristics	Description		
Demodulator			
Carrier type	Continuous and B	urst	
Modulation format	BPSK, QPSK, π/4 GMSK, GFSK, AS	I shift DQPSK, 8PSK, 16QAI SK and FSK	M, 64QAM, 128QAM, 256QAM,
Measurement filter	Root cosine		
Reference filter	Cosine, Gauss		
Filter parameter	α/BT: 0.0001 to 1,	, 0.0001 step	
Maximum symbol rate	12.8 Msps		
Standard setup	PDC, PHS, NADO	C, TETRA, GSM, CDPD, and	Bluetooth
Display format			
Vector diagram	Symbol/Locus dis	play, Frequency error measu	rement, and Origin offset measurement
Constellation diagram	Symbol display, F	requency error measurement	, and Origin offset measurement
Eye diagram	I, Q, and Trellis di	splays ( 1 to 16 symbols )	
Error vector diagram		error, Phase error, Waveform neasurement, and Origin offs	
Symbol table	Binary, Octal, Hexadecimal		
AM/AM	Measured amplitu	de vs. Reference amplitude,	1 dB compression measurement
AM/PM	Phase error vs. Re	eference amplitude	
CCDF	Probability of excess vs. Power level, Crest factor measurement		
PDF	Probability of occurrence vs. Power level		
Accuracy			
PDC (100 kHz span)	EVM ≤1.2%, Mag	nitude error ≤1.0%, Phase e	rror ≤0.8°
PHS (1 MHz span)	EVM $\leq$ 1.4%, Magnitude error $\leq$ 1.2%, Phase error $\leq$ 0.8°		rror ≤0.8°
GSM (1 MHz span)	EVM ≤1.8%, Mag	nitude error ≤1.2%, Phase e	rror ≤1.0°
64QAM, 5.3 Msps, 1 GHz carrier (15 MHz span, <i>Typical</i> )	EVM ≤2.5%		
QPSK, 3.84 Msps, 2 GHz carrier (15 MHz span, <i>Typical</i> )	EVM ≤2.5%		
QPSK	EVM (Typical)	Symbol rate	
Center frequency = 1 GHz	0.5% 0.5% 1.2% 2.7%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 2 GHz	0.5% 0.5% 1.2% 2.7%	100 kHz 1 MHz 4 MHz 10 MHz	

### Table A-15: Digital demodulation (Option 21 only)

Table A-15: Digital demodulation (Option 21 only) (Cont.)

aracteristics	Description		
Center frequency = 3 GHz	0.7% 0.7% 1.5% 2.9%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 5 GHz (RSA3308B only)	0.7% 0.7% 1.5% 3.0%	100 kHz 1 MHz 4 MHz 10 MHz	
π/4 DQPSK	EVM (Typical)	Symbol rate	
Center frequency = 1 GHz	0.6% 0.6% 1.2% 2.7%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 2 GHz	0.6% 0.6% 1.2% 2.7%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 3 GHz	0.7% 0.7% 1.5% 2.9%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 5 GHz (RSA3308B only)	0.7% 0.7% 1.5% 3.0%	100 kHz 1 MHz 4 MHz 10 MHz	
16QAM	EVM (Typical)	Symbol rate	
Center frequency = 1 GHz	0.9% 0.5% 1.2% 2.2%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 2 GHz	0.9% 0.5% 1.2% 2.2%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 3 GHz	0.9% 0.5% 1.2% 2.5%	100 kHz 1 MHz 4 MHz 10 MHz	
Center frequency = 5 GHz (RSA3308B only)	0.9% 0.5% 1.2% 2.5%	100 kHz 1 MHz 4 MHz 10 MHz	

### Table A-16: ACLR measurement (Option 30 only)

Characteristics	Description	
ACLR (3GPP down link, Test model 1, 1 DPCH ch, 1800 to 2200 MHz carrier frequency)		
Real-Time measurement	66 dB (5 MHz offset) 68 dB (10 MHz offset)	
Sweep measurement with noise correction (Typical)	70 dB (5 MHz offset) 72 dB (10 MHz offset)	

### Table A- 17: Pulse measurement

Characteristics	Description
Measurement functions	Pulse Width, Pulse Peak Power, On/Off Ratio, Pulse Ripple, Pulse Repitition Interval, Duty Cycle, Pulse-Pulse Phase, Channel Power, OBW, EBW, Frequency Deviation
Minimum pulse length	20 samples
Maximum pulse length	260,000 samples

### Table A- 18: Digital phosphor spectrum processing (DPX)

Characteristics	Description
Spectrum processing rate (Typical)	≥10,000 per second (span independent)
Minimum signal duration for 100% probability of intercept ( <i>Typical</i> )	≤31 μs
Span range	100 Hz to 20 MHz (1-2-5 sequence), 36 MHz
Minimum RBW settings	300 kHz (36 MHz span) 200 kHz (20 MHz span) 100 kHz (10 MHz span) 50 kHz (5 MHz span) 20 kHz (2 MHz span) 20 kHz (2 MHz span) 10 kHz (1 MHz span) 5 kHz (500 kHz span) 2 kHz (200 kHz span) 1 kHz (100 kHz span) 500 Hz (50 kHz span) 200 Hz (20 kHz span) 100 Hz (10 kHz span) 50 Hz (5 kHz span) 20 Hz (2 kHz span) 10 Hz (1 kHz span) 5 Hz (500 Hz span) 2 Hz (200 Hz span) 2 Hz (200 Hz span) 1 Hz (100 Hz span)

### Table A-18: Digital phosphor spectrum processing (DPX) (Cont.)

Characteristics	Description
RBW accuracy	+1%, -7%
RBW range	1 Hz to 5 MHz (1,2,3,5 sequence)
Span accuracy	±0.3%
Amplitude accuracy	±0.5 dB

### Table A-19: Display

Characteristics	Description	
View		
Number of views	1, 2, 3, or 4	
Number of display traces	2	
LCD		
Size	213 mm (8.4 in)	
Resolution	$800 \times 600$ pixels	
Color	Maximum 256 colors	
Display detector	Positive peak, Negative peak, and Positive-Negative peak	

### Table A-20: Marker, trace, and display line

Characteristics	Description
Marker type	Normal, Reference, and Band power
Search function	Peak right, Peak left, Maximum, Larger peak, and Smaller peak
Trace content	Active, Average, Max Hold, Min Hold, View, and Off
Display line	Horizontal line 1 and 2, Vertical line 1 and 2

### Table A-21: Controller and interface

Characteristics	Description	
Controller		
CPU	Intel Celeron 1.2 GHz	
DRAM	512 MB DIMM	
OS	Windows XP	
System bus	PCI, ISA	

Characteristics	Description
Storage medium	
Hard disk	≥20 GB 2.5 inch IDE
Floppy disk	1.44 MB 3.5 inch
Interface	
Printer port	USB
GPIB	IEEE 488.1
LAN	10/100 BASE-T (IEEE 802.3)
Mouse	USB
Keyboard	USB
Monitor out	VGA (D-SUB 15 pins)

### Table A-21: Controller and interface (Cont.)

### Table A-22: Power requirements

Characteristics	Description
Rating voltage	100 to 240 VAC
Voltage range	90 to 250 VAC
Line frequency	47 to 63 Hz
Mains Fuse Data	Densei-Lambda supplies: 5 A, Time-delayed, 250 V (not operator replaceable) Cosel supply: 2 A, Time-delayed, 250 V (not operator replaceable)
Heat dissipation	· · · · ·
Maximum power	350 VA
Maximum line current	5 A rms at 50 Hz (90 V line with 5% clipping)
Surge current	Maximum 52 A peak (25 °C) for ≤5 line cycles after the product has been turned off for at least 30 s.

### Table A-23: Power connector

Characteristics	Description	
Preamp power connector		
Connector type	LEMO 6 poles	
Pin assignment	Pin 1: NC, Pin 2: ID1, Pin 3: ID2, Pin 4: -12 V, Pin 5: GND, Pin 6: +12 V	

## **Physical Characteristics**

### Table A-24: Physical characteristics

Characteristics	Description		
Dimensions			
Width	425 mm (16.7 in) without belts		
Height	215 mm (8.5 in) without feet		
Depth	425 mm (16.7 in) without cover and feet		
Net weight	19 kg		

### **Environmental Characteristics**

### Table A-25: Environmental characteristics

Characteristics	Description
Temperature	
Operating	+10 to +40 °C
Nonoperating	-20 to +60 °C
Relative humidity	
Operating and nonoperating	20 to 80% (no condensation), maximum wet-bulb temperature 29 $^\circ\text{C}$
Altitude	
Operating	Up to 3000 m (10000 ft)
Nonoperating	Up to 12000 m (40000 ft)
Vibration	
Operating	2.16 m/s <sup>2</sup> rms (0.22 G rms), 5 to 500 Hz (class 8)
Nonoperating	22.3 m/s <sup>2</sup> rms (2.28 G rms), 5 to 500 Hz (class 5)
Shock	
Nonoperating	196 m/s <sup>2</sup> (20 G), half-sine, 11 ms duration Three shocks in each direction along each major axis, total of 18 shocks
Cooling clearance	
Bottom	20 mm (0.79 in)
Both sides	50 mm (1.97 in)
Rear	50 mm (1.97 in)

Characteristic	Description		
EC Declaration of Conformity - EMC	Meets intent of Directive 2004/108/EEC for Electromagnetic Compatibility. Compliance was demonstrated to the following specifications as listed in the Official Journal of the European Communities:		
	EN 61326:1997.	EMC requirements for Class A electrical equipment for measurement, control and laboratory use. Annex D <sup>1,2,3</sup>	
	IEC 61000-4-2:1999.Electrostatic discharge immunityIEC 61000-4-3:2002.RF electromagnetic field immunityIEC 61000-4-4:2004.Electrical fast transient/burst immunitIEC 61000-4-5:2005.Power line surge immunityIEC 61000-4-6:2003.Conducted RF immunityIEC 61000-4-11:2004.Voltage dips and interruptions immunity		
	EN 61000-3-2:2002. AC power line harmonic emissions EN 61000-3-3:1995. Voltage changes, fluctuations, and flicker		
Australia/New Zealand Declaration of Conformity - EMC	Complies with EMC provision of Radiocommunications Act per the following standar in accordance with ACMA:         EN 61326:1997.       EMC requirements for electrical equipment for measurement, control, and laboratory use.		
EC Declaration of Conformity - Low Voltage	Compliance was demonstrated to the following specification as listed in the Official Journal of the European Union:		
	Low Voltage Directive 2006	/95/EC	
	EN 61010-1:2001.	Safety requirements for electrical equipment for measurement control and laboratory use.	
U.S. Nationally Recognized Testing Laboratory Listing	UL61010-1.	Standard for electrical measuring and test equipment.	
Canadian Certification	CAN/CSA C22.2 No. 61010	GA C22.2 No. 61010-1-04. Safety requirements for electrical equipment for measurement, control, and laboratory use. Part 1.	
Additional Compliance	ANSI/ISA 82.02.01:2004.	Safety standard for electrical and electronic test, measuring, controlling, and related equipment.	
	IEC61010-1.	Safety requirements for electrical equipment for measurement, control, and laboratory use.	

### Table A-26: Certifications and compliances

<sup>1</sup> This product is intended for use in nonresidential areas only. Use in residential areas may cause electromagnetic interference.

<sup>2</sup> Emissions which exceed the levels required by this standard may occur when this equipment is connected to a test object.

- <sup>3</sup> To ensure compliance with the EMC standards listed here, high quality shielded interface cables should be used.
- <sup>4</sup> Performance Criteria with disturbance signal frequencies within 50 MHz of EUT Center Frequency: Residual spurious signals can increase to -60 dBm with exposure to the disturbance levels of the test.

Characteristic	Description	Description			
Overvoltage Category Description		Terminals on this product may have different installation (overvoltage) category designations. The installation categories are:			
		pution-level mains (usually permanently connected). Equipment at this is typically in a fixed industrial location.			
	applia	level mains (wall sockets). Equipment at this level includes inces, portable tools, and similar products. Equipment is usually connected.			
	CAT I Seco	ndary (signal level) or battery operated circuits of electronic equipment.			
Overvoltage Category	Overvoltage Ca	tegory II (as defined in IEC61010-1)			
Pollution Degree Description	a product. Typic	ne contaminates that could occur in the environment around and within cally the internal environment inside a product is considered to be the ternal. Products should be used only in the environment for which they			
	Pollution Degre	e 1 No pollution or only dry, nonconductive pollution occurs. Products in this category are generally encapsulated, hermetically sealed, or located in clean rooms.			
	Pollution Degre	e 2 Normally only dry, nonconductive pollution occurs. Occasionally a temporary conductivity that is caused by condensation must be expected. This location is a typical office/home environment. Temporary condensation occurs only when the product is out of service.			
	Pollution Degre	e 3 Conductive pollution, or dry, nonconductive pollution that becomes conductive due to condensation. These are sheltered locations where neither temperature nor humidity is controlled. The area is protected from direct sunshine, rain, or direct wind.			
	Pollution Degre	e 4 Pollution that generates persistent conductivity through conductive dust, rain, or snow. Typical outdoor locations.			
Pollution Degree	Pollution Degre	e 2 (as defined in IEC61010-1). Note: Rated for indoor use only.			
Safety Certification Compliance					
Equipment Type	Test and measu	iring equipment			
Safety Class	Class I (as defin	Class I (as defined in IEC61010-1) - grounded product			
Operating Temperature Range	+5 to +40 °C				

### Table A-26: Certifications and compliances (Cont.)

## **Appendix B: Menu Tree**

This section shows the structure of the menus and submenus displayed by pressing the menu keys as illustrated in Figure B-1. The menu keys are divided into the following groups:

- Settings
- Mode
- Markers
- View
- Utility

Some of the submenus are used for programming or servicing; information about these submenus is in the programmer or service manual (optional). See *Appendix E* for manual part numbers.



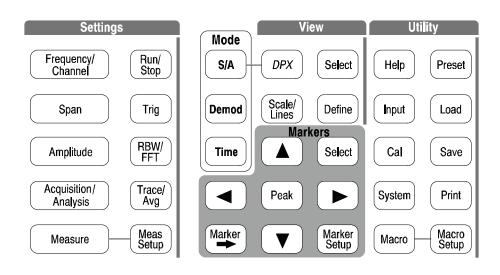


Figure B-1: Menu keys

### Settings Menu (1 of 3)


### S/A Mode except Real Time S/A

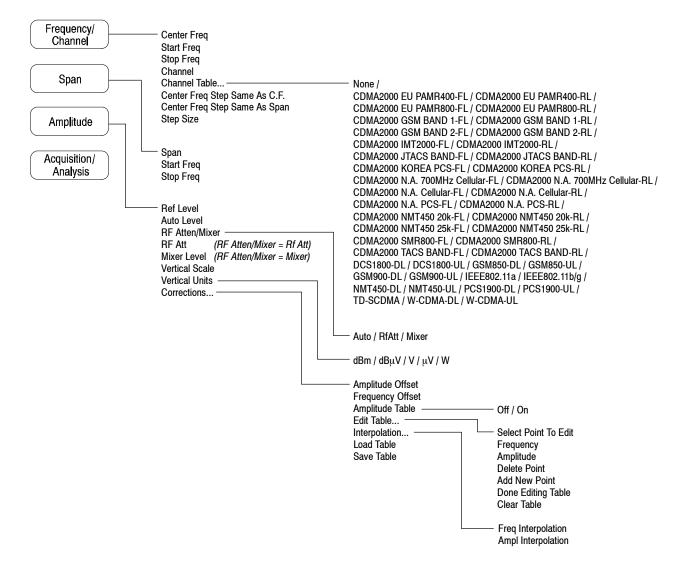

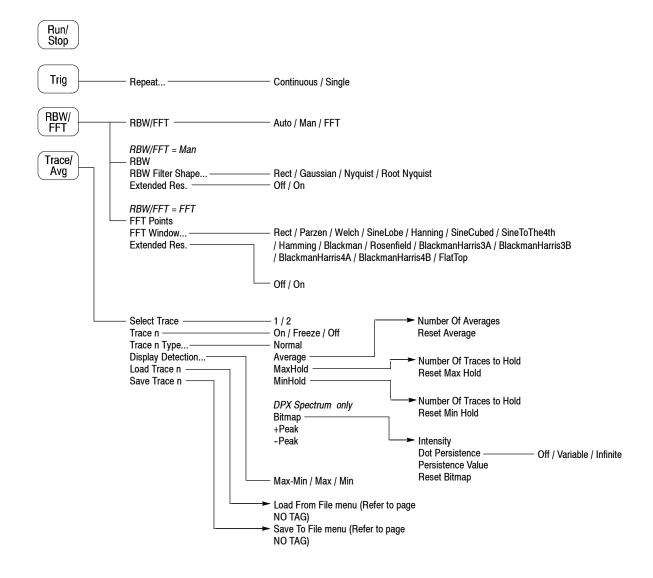

### Settings Menu (2 of 3)


### Real Time S/A, Demod and Time Modes

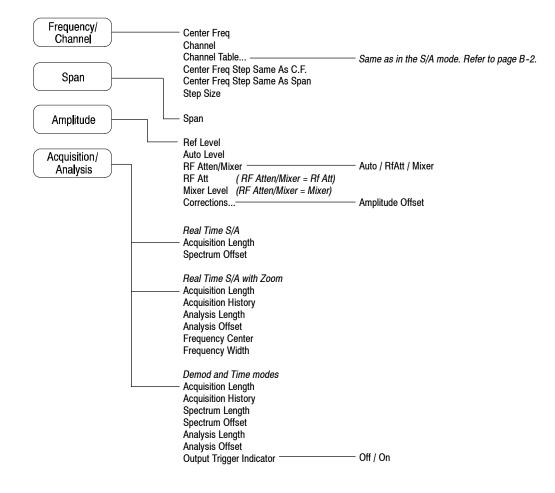


Run/ Stop		
Trig	<ul> <li>Mode</li> <li>Repeat</li> <li>Stop and Show Results</li> </ul>	<sup>—</sup> Free Run / Triggered — Continuous / Single
FFT Trace/ Avg	Mode = Triggered Source Save on Trigger	– Power / Freq Mask * / External  * <i>Option 02 only.</i> – On / Off – On / Off
	Source = Power, External Level Slope Position	- Rise / Fall / Rise and Fall / Fall and Rise
	Source = Freq Mask Define Mask Slope Position	<ul> <li>Select Next Point</li> <li>Set Selected Point X</li> <li>Set Selected Point Y</li> <li>Delete Selected Point</li> <li>Insert New Point</li> <li>Set All Points to Maximum</li> <li>Set All Points to Minimum</li> <li>Reset Mask to Default</li> </ul>
		- In / Out / In and Out / Out and In
	Real Time S/A Mode - FFT Start Point FFT Overlap FFT Window	- Same as in the S/A mode. Refer to page B-3.
	Pulse Measurements in the Time M - FFT Window Rolloff Ratio	
	<i>Demod and Time Modes</i> – Average – Average Count	- Off / On
	Average Term Control	<sup>—</sup> Expo / Repeat

### Settings Menu (3 of 3)

### S/A Mode

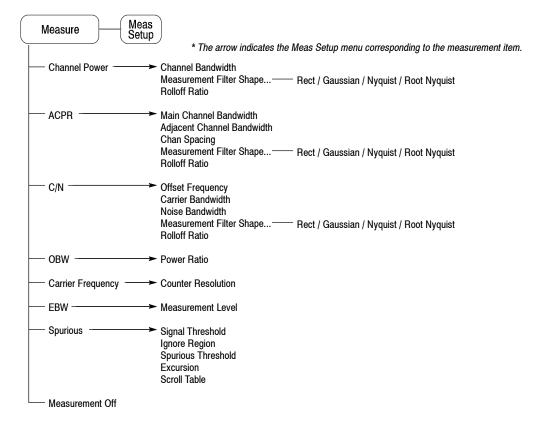


### **Demod Mode**

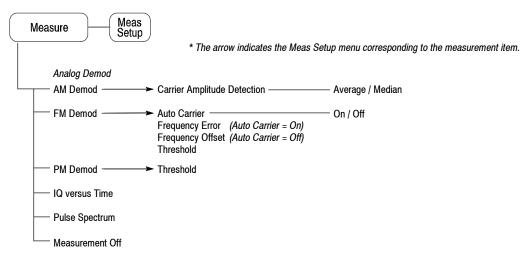


### Demod Mode

Measure Meas	* The arrow indicates the Meas Setup	menu corresponding to the measurement ite	m.
Digital Demod (Option 2	21 only)		
	Parameter Presets		A / GSM / CDPD / Bluetooth
EVM	Modulation Type		
IQ/Frequency	Modulation Parameters	1/4PI_QPSK / BPSK / QPSK / 8P	
versus Time	Auto Carrier On / Off	/ 16QAM / 32QAM / 64QAM / 128	
Power versus Time – Symbol Table –	Frequency Error <i>(Auto Carrier = On)</i> Frequency Offset <i>(Auto Carrier = Off)</i>	/ GMSK / GFSK / ASK / FSK / OC	JL 2V
Eye Diagram		Symbol Rate	
AM/AM	AM/AM, AM/PM, CCDF, PDF	Measurement Filter	— None / Root Raised Cosine
	Linear Signal Region Unit —— dB / dBm		— None / Raised Cosine / Gaussian
PDF	Linear Signal Region	Filter Parameter	
Measurement Off	AM/AM, AM/PM	Modulation Type = ASK	
	Scroll Table	Auto Modulation Depth	— On / Off
	Maximum Coefficients	Modulation Depth	
	Display Reference Line On / Off Display Best-Fit Line On / Off		
	Display Best-Fit Line — On / Off	Modulation Type = FSK, GFSK Auto Frequency Deviation	
	AM/AM only	Frequency Deviation	
	Linear Signal Region Mask — On / Off		
		Modulation Type = ASK, FSK, GF	FSK
	CCDF, PDF	Decoding Format	— NRZ / Manchester / Miller
	Horizontal Division		
	CCDF only		
	Display Gaussian Line ——— On / Off		
RFID Measurements (C	Dption 21 only)		
└─── Carrier Analyz	e ACPF	R Main Channel Bandwidth	
	er Resolution	Adjacent Channel Bandwidth	
	Ratio for the OBW — % / dB	Chan Spacing	
	Ratio for the OBW ude Offset	Measurement Filter Shape — Rolloff Ratio	
	idth Ch Power	Scroll Table	Root Nyquist
Spurious> Same	as Spurious in the S/A mode.		
Power On/Down	➤ Analyze		
RF Envelope		8000-4-1 / 18000-6-A / 18000-6-B / 18000-6-	
Constellation —— Eye Diagram ——		18092(494k) / 15693-2 / C0G1 / C1G1 / Mar	nual
Symbol Table	Modulation Type In In In In	iterrogator / Tag	
Measurement Off	Auto Bit Rate	nerrogator / Tag	
		SK / DSB-ASK / SSB-ASK / PR-ASK / FSK /	OOK / Subcarrier OOK
	Tari / S	Subcarrier BPSK	
	Settling Error Width +-		
		lanchester / Miller / Miller (M_2, M_4, or M_8 PIE (Type A or C) / NRZ / NRZ-L (2, 4, or 8	
		1 out of 4 / 1 out of 256 / SSC (High or Low)	
	- /	· · · · · · · · · · · · · · · · · · ·	
	L 0	n / Off	

### **Time Mode**

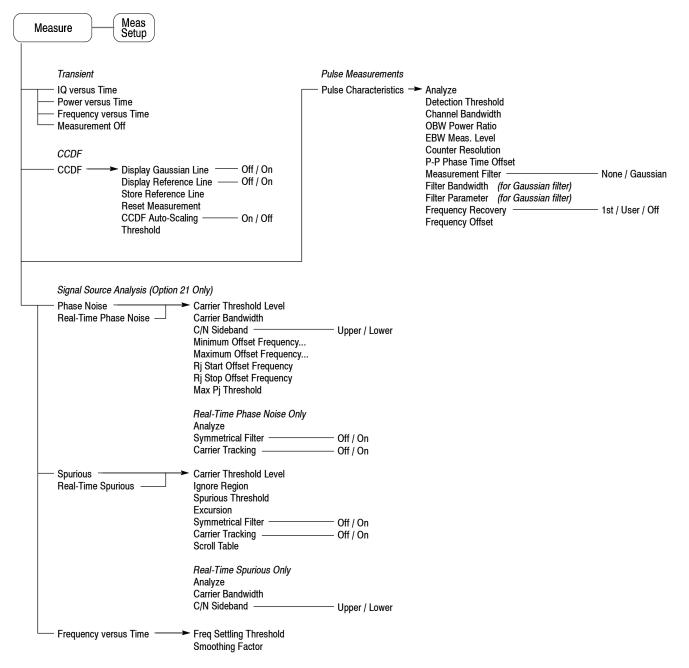


### Mode Menu

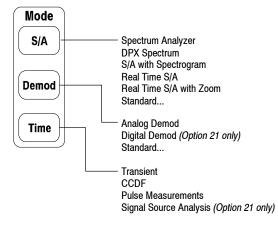


\* For details about the Standard... menu, refer to the optional software user manual.

### **Markers Menu**

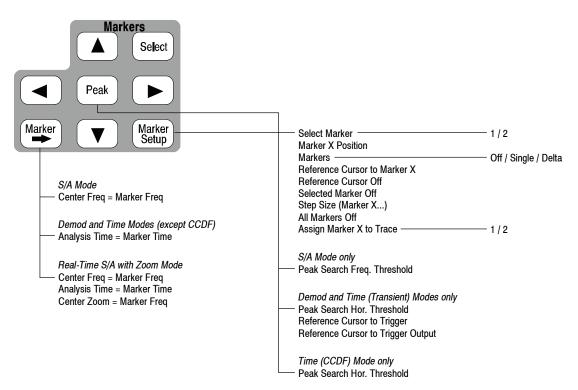


### Display: View Menu (1 of 2)

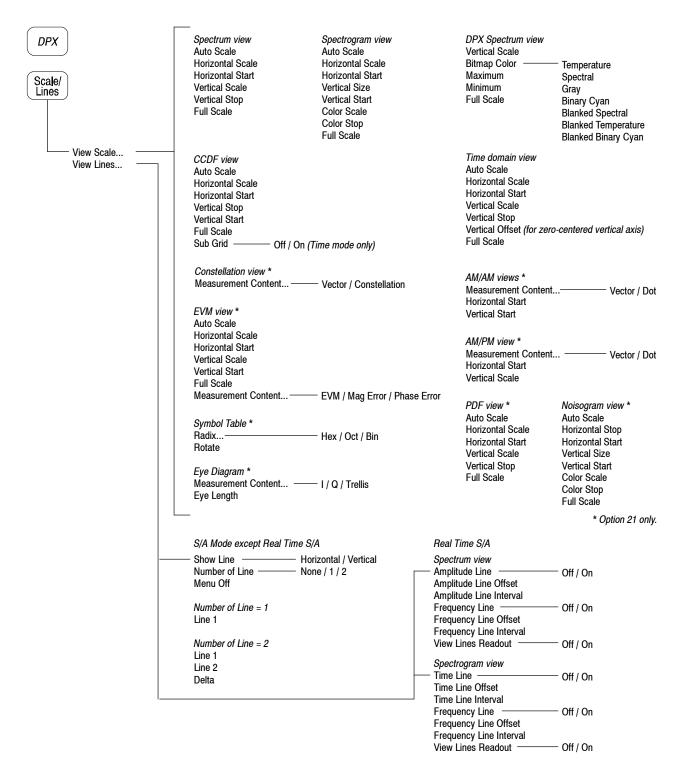


### Display: View Menu (2 of 2)

### S/A Mode

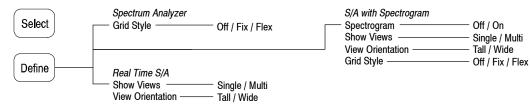


### **Demod and Time Modes**

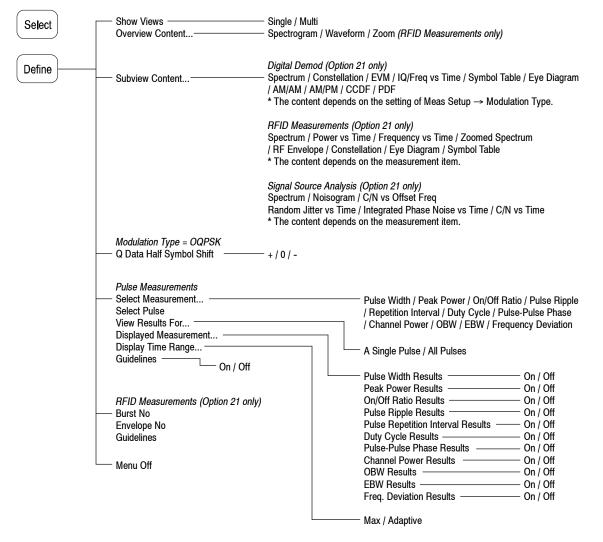

### **Utility Menu**


Help	<ul> <li>View Front Panel Button Help</li> <li>View Online User Manual</li> <li>View Online Programmer Manual</li> </ul>				
	- Signal Input Port Reference Source	RF / IQ * / Cal100 Int / Ext	0M / Cal25M	* Option	03 only.
Cal	Calibrate All     Calibrate Gain     Calibrate Center Offset     Calibrate DC Offset     Calibrate IF Flatness     Auto Calibration     Service (Refer to the Service Manual)	— Yes / No			
	Display Brightness     Reset All to Factory Defaults     Remote Setup (Refer to the Programmer Manual)     Version and Installed Options Instrument Setup	<ul> <li>— Scroll</li> <li>Option Key</li> <li>Change Option Key</li> <li>— Angular Units (Ph Digital IQ Output</li> </ul>	ey ase) —————	Degrees On / Off	/ Radians
Preset Load Save	S/A mode except Real Time  Load State Load Trace 1 Load Trace 2 Load Correction  Real Time S/A, Demod and Time modes Load State Load State		Load From File I Select File Load File Now Folder Load from XXXA Load from XXXC	Menu	
Print	Load State         Load Data         Save State         Save Trace 1         Save Trace 2         Save Correction         Real Time S/A mode         Save State         Save Data         Real Time S/A with Zoom, Demod and Time modes         Save State         Save Data         Real Time S/A with Zoom, Demod and Time modes         Save Data         Save Data         Save Data         Save Data         Current Block / Current Block / Current Area         Save Data (.MAT)         Current Block / Current Area         Save Data (.MAT)		Save To File Mer Filename Save File Now Overwrite Folder Save to XXXA Save to XXXB Save to XXXC		here XXX = State for a state file IQData for an IQ data file Trace for a trace file Corr for a correction file Bitmap for a bitmap file
	Save screen to file Background color ——— Black / White Printer				

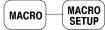

\* Contact your local Tektronix distributor or sales office for installing user-specific macros.

## **Appendix C: Digital Demodulation Symbol Mapping**

This section shows the symbol mapping for each digital modulation technique.

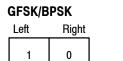

QPSK	
Left	Right
1	3
0	2

8	BPSK Left				Right
			3		
		2		1	
	6				0
		7		4	
			5		

<b>QAM</b> .eft			Right
3	2	1	0
7	6	5	4
В	А	9	8
F	E	D	С

3	20	<u>ک</u>	М

Left					Right
	3	2	1	0	
9	8	7	6	5	4
F	E	D	С	В	Α
15	14	13	12	11	10
1B	1A	19	18	17	16
	1F	1E	1D	1C	

64QAM

Left							Right
7	6	5	4	3	2	1	0
F	E	D	С	В	А	9	8
17	16	15	14	13	12	11	10
1F	1E	1D	1C	1B	1A	19	18
27	26	25	24	23	22	21	20
2F	2E	2D	2C	2B	2A	29	28
37	36	35	34	33	32	31	30
3F	3E	3D	3C	3B	3A	39	38

Figure C-1: Symbol mapping: GFSK/BPSK, QPSK, 8PSK, 16QAM, 32QAM, and 64QAM

Appendix	C: Digital	Demodulation	Symbol	Mapping

Left															Right
EF	FD	EB	F9	E7	F5	E3	F1	0F	3F	4F	7F	8F	BF	CF	FF
CE	DC	CA	D8	C6	D4	C2	D0	0C	3C	4C	7C	8C	BC	СС	FC
AF	BD	AB	B9	A7	B5	A3	B1	0B	3B	4B	7B	8B	BB	СВ	FB
8E	9C	8A	98	86	94	82	90	08	38	48	78	88	B8	C8	F8
6F	7D	6B	79	67	75	63	71	07	37	47	77	87	B7	C7	F7
4E	5C	4A	58	46	54	42	50	04	34	44	74	84	B4	C4	F4
2F	3D	2B	39	27	35	23	31	03	33	43	73	83	B3	C3	F3
0E	1C	0A	18	06	14	02	10	00	30	40	70	80	B0	C0	F0
E1	D1	A1	91	61	51	21	11	01	13	05	17	09	1B	0D	1F
E2	D2	A2	92	62	52	22	12	20	32	24	36	28	ЗA	2C	3E
E5	D5	A5	95	65	55	25	15	41	53	45	57	49	5B	4D	5F
E6	D6	A6	96	66	56	26	16	60	72	64	76	68	7A	6C	7E
E9	D9	A9	99	69	59	29	19	81	93	85	97	89	9B	8D	9F
EA	DA	AA	9A	6A	5A	2A	1A	A0	B2	A4	B6	A8	BA	AC	BE
ED	DD	AD	9D	6D	5D	2D	1D	C1	D3	C5	D7	C9	DB	CD	DF
EE	DE	AE	9E	6E	5E	2E	1E	E0	F2	E4	F6	E8	FA	EC	FE

Figure C-2: Symbol mapping: 256QAM

### Table C-1: $\pi/4$ shift DQPSK

Rotation angle (degree)	Bit
+45	0
+135	1
-45	2
-135	3

### Table C-2: GMSK

Rotation direction	Bit
Left	1
Right	0

## **Appendix D: Restoring the Instrument Operating System**

You can restore the instrument operating system from a file on the instrument hard disk drive.



**CAUTION.** Using the restore process reformats the hard drive and reinstalls the operating system. All saved data will be lost. If possible, save important files to external media before performing a system restore.

### **Restoring the Operating System from the Instrument Hard Disk**

Perform the following steps to restore the operating system:

1. Restart the instrument. During the boot-up process you will see the following message at the bottom of the screen:

Starting Acronis Loader... press F5 for Acronis Startup Recovery Manager

2. Repeatedly press the F5 key until the Acronis True Image Tool opens.

There is a three-second time period from when the message appears until the instrument proceeds with the normal instrument startup. If the instrument does not open the Acronis application, power off the instrument, then power on the instrument and try again.

- 3. Click Restore.
- 4. In the Confirmation dialog box, click **Yes** to restore the instrument operating system, or No to exit the restore process.

The restore process takes approximately 30 minutes; the actual time depends on the instrument configuration.

## **Appendix E: Inspection and Cleaning**

Inspect and clean the exterior of the instrument as often as operating conditions require. Dirt acts as an insulating blanket, preventing efficient heat dissipation. Regular cleaning may prevent instrument malfunction and enhance reliability.



**WARNING.** To avoid personal injury, unplug the power cord from line voltage before cleaning the instrument. To avoid getting moisture inside the instrument during external cleaning, use only enough liquid to dampen the applicator.

### Inspecting the Exterior

Inspect the outside of the instrument for damage, wear, and missing parts, using Table E-1 as a guide. Instruments that appear to have been dropped or otherwise abused should be checked thoroughly to verify correct operation and performance. Immediately repair defects that could cause personal injury or lead to further damage to the instrument.

ltem	Inspect for	Repair action
Cabinet, front panel, and cover	Cracks, scratches, deformations, damaged hardware or gaskets. Contact your loca Tektronix distribu	
Front-panel knobs	Missing, damaged, or loose knobs.	sales office.
Connectors	Broken shells, cracked insulation, and deformed contacts. Dirt in connectors.	
Carrying handle	Correct operation.	
Accessories	Missing items or parts of items, bent pins, broken or frayed cables, and damaged connectors.	



**CAUTION.** Avoid the use of chemical cleaning agents that might damage the plastics used in this instrument. Organic solvents such as benzene and acetone must never be used.

### **Cleaning the Exterior**

Remove the power cord before cleaning the instrument. Clean the exterior surfaces of the instrument with a dry, lint-free cloth or a soft-bristle brush. If dirt remains, use a cloth or swab dampened with a 75% isopropyl alcohol solution. A swab is useful for cleaning in narrow spaces around the controls and connectors. Do not use abrasive compounds on any part of the instrument.



**CAUTION.** Do not allow moisture inside the instrument. During exterior cleaning, use only enough solution to dampen the cloth or swab.

### Lubrication

There is no periodic lubrication required for this instrument.

### **Cleaning the Interior**

Consult your Tektronix Service Center or representative for cleaning the analyzer interior.

## **Appendix F: Options and Accessories**

This appendix describes the options and standard and optional accessories that are available for the analyzer.

### Options

The following options can be ordered for the analyzer:

- Option 1A: External preamplifier (20 dB gain to 3 GHz)
- Option 1R: Rack Mount Kit (includes: hardware and instructions for converting to 19-inch rackmount configuration)
- Option L5: Provides Japanese language user manual and programmer manual instead of English manuals
- Option 02: 256 MB data memory with frequency mask trigger
- Option 03: IQ input function
- Option 21: Advanced measurement suite
- Option 24: GSM/EDGE analysis software
- Option 25: cdma2000 analysis software
- Option 26: cdma2000 1xEV-DO analysis software
- Option 28: TD-SCDMA analysis software
- Option 30: 3GPP Release 99 and Release 5 uplink and downlink analysis software
- Option 40: 3GPP Release 6 (HSUPA) analysis software
- Option C3: Calibration service, 3 years
- Option C5: Calibration service, 5 years
- Option D1: Calibration data report
- Option D3: Calibration data report, 3 years (requires Option C3)
- Option D5: Calibration data report, 5 years (requires Option C5)
- Option R3: Repair service, 3 years
- Option R5: Repair service, 5 years

■ Power cord options listed in Table F-1

### Table F-1: Power cord identification

Plug configuration	Normal usage	Option number
A CAR	North America 120 V	A0
The second secon	Universal Euro 230 V	A1
	United Kingdom 230 V	A2
T CB	Australia 240 V	A3
- Contraction of the second se	North America 240 V	A4
	Switzerland 220 V	A5
	Japan 100 V	A6
Tool Section 1997	China 240 V	A10
	No power cord supplied.	A99

### **Standard Accessories**

All analyzers in this series are shipped with the accessories listed in Table F-2.

Table F-2: Standard accessories

Accessory	Part number
English language user manual	071-2363-XX
English language programmer manual <sup>1</sup>	071-2382-XX
U.S. power cord	161-0066-XX
USB mouse	119-6936-XX
USB keyboard	119-B146-XX
BNC-N adapter	103-0045-XX
Front cover	200-A524-50
Documents CD	063-4089-XX

<sup>1</sup> PDF documents are stored on the Documents CD.

### **Optional Accessories**

You can order the accessories listed in Table F-3 to use with your analyzer.

### Table F-3: Optional accessories

Accessory	Part number
Japanese language user manual	071-2362-XX
Japanese language programmer manual	071-2383-XX
English language service manual	071-2367-XX
Accessory bag	016-A330-00
Preamplifier (20 dB gain to 3 GHz)	650-A900-XX
RTPA2A real-time probe adapter (P7000 series probes recommended)	
Rack mount kit (for field conversion)	

### Optional Software User Manuals

Optional software includes the user manual listed in Table F-4.

### Table F-4: Accessories for optional software

Accessory	Part number
Option 24 GSM/EDGE analysis software user manual	071-2398-XX
Option 25 cdma2000 analysis software user manual	071-2400-XX
Option 26 1xEV-DO analysis software user manual	071-2402-XX
Option 28 TD-SCDMA analysis software user manual	071-2404-XX
Option 30 3GPP Release 99 and Release 5 uplink and downlink analysis software user manual	071-2396-XX
Option 40 3GPP Release 6 (HSUPA) analysis software user manual	071-2407-XX

# Glossary

## Glossary

#### Accuracy

The closeness of the indicated value to the true value.

#### Acquisition

A series of time-contiguous frames. This is also called a Block.

### **Acquisition Time**

The length of time contained in one acquisition.

#### Amplitude

The magnitude of an electrical signal.

#### **Amplitude Modulation (AM)**

The process, or result of a process, in which the amplitude of a sine wave (the carrier) is varied in accordance with the instantaneous voltage of a second electrical signal (the modulating signal).

#### **Analysis Time**

A subset of time-contiguous samples from one Block Time, used as input to an Analysis View.

#### **Analysis View**

Any view which produces measurement results except for Spectrum, Power versus Time, and Spectrogram.

#### ASK

Acronym for Amplitude Shift Keying. The process, or result of a process, in which the amplitude of the carrier is varied in accordance with the state of a digital input signal.

### Block

See Acquisition.

### Calibrator

A signal generator producing a specified output used for calibration purposes.

### **Carrier, Carrier Signal**

The electrical signal, typically a sine wave, upon which modulation is impressed.

### **Carrier Frequency**

The frequency of the carrier signal.

### Carrier-to-Noise Ratio (C/N)

The ratio of carrier signal power to average noise power in a given bandwidth surrounding the carrier; usually expressed in decibels.

#### **Center Frequency**

That frequency which corresponds to the center of a frequency span, expressed in hertz.

### Clear (Erase)

Presets memory to a prescribed state, usually that denoting zero.

#### Controls

User-accessible mechanisms to modify setup parameters.

#### dBm

A unit of expressed power level in decibels referenced to 1 milliwatt.

#### dBmV

A unit to express voltage levels in decibels referenced to 1 millivolt.

### dBµV

A unit to express voltage levels in decibels referenced to 1 microvolt.

#### **Decibel (dB)**

Ten times the logarithm of the ratio of one electrical power to another.

### **Depth of Modulation**

The difference between the maximum and the minimum of the RF envelope amplitude expressed as a percentage of the maximum RF envelope level.

#### **Display Law**

The mathematical law that defines the input-output function of the instrument. The following cases apply:

#### Linear

A display in which the scale divisions are a linear function of the input signal voltage.

### Square Law (Power)

A display in which the scale divisions are a linear function of the input signal power.

### Logarithmic

A display in which the scale divisions are a logarithmic function of the input signal voltage.

#### **Display Line**

A horizontal or vertical line on a waveform display, used as a reference for visual (or automatic) comparison with a given level, time, or frequency.

#### **Display Reference Level**

A designated vertical position representing a specified input level. The level may be expressed in dBm, volts, or any other units.

#### Distortion

Degradation of a signal, often a result of nonlinear operations, resulting in unwanted signal components. Harmonic and intermodulation distortion are common types.

#### **Dynamic Range**

The maximum ratio of the levels of two signals simultaneously present at the input which can be measured to a specified accuracy.

#### Export

Save data to a file in a format other than application-native.

#### **FFT Time**

The length of time covered by one frame, whether or not an RBW filter is applied to the FFT. This was previously called a "multi-frame".

#### Filter

A circuit which separates electrical signals or signal components based on their frequencies.

#### Frame

A series of time-contiguous samples, long enough in duration and at a sufficient sample rate to produce a spectrum view of a specified span and RBW.

#### Frequency

The rate at which a signal oscillates, or changes polarity, expressed as hertz or number of cycles per second.

#### **Frequency Band**

The continuous range of frequencies extending between two limiting frequencies, expressed in hertz.

#### **Frequency Domain Representation**

The portrayal of a signal in the frequency domain; representing a signal by displaying its sine wave components; the signal spectrum.

#### **Frequency Drift**

Gradual shift or change in displayed frequency over the specified time due to internal changes in the spectrum analyzer, where other conditions remain constant. Expressed in hertz per second.

#### **Frequency Modulation (FM)**

The process, or result of a process, in which the frequency of an electrical signal (the carrier) is varied in accordance with some characteristic of a second electrical signal (the modulating signal or modulation).

#### **Frequency Range**

The range of frequencies over which the performance of the instrument is specified.

#### **Frequency Span (Dispersion)**

The magnitude of the frequency band displayed; expressed in hertz or hertz per division.

#### FSK

Acronym for Frequency Shift Keying. The process, or result of a process, in which the frequency of the carrier is varied in accordance with the state of a digital input signal.

#### **GPIB**

Acronym for General Purpose Interface Bus, the common name for the communications interface system defined in IEEE Std 488.

#### Graticule

The calibrated grid overlaying the display screen of spectrum analyzers, oscilloscopes, and other test instruments.

#### **Grayed Out**

An on-screen control is "grayed out" if it is not adjustable.

#### Import

Bring data into the application from a file of some format other than application-native.

#### **Input Impedance**

The impedance at the desired input terminal. Usually expressed in terms of VSWR, return loss, or other related terms for low impedance devices and resistance-capacitance parameters for high impedance devices.

#### **Intermediate Frequency (IF)**

In a heterodyne process, the sum or difference frequency at the output of a mixer stage which will be used for further signal processing.

#### **Key Contacts**

A pattern of interlaced fingers on the front-panel circuit board that form a contact closure when mated with a conductive pad under a mechanical front panel button.

#### Link

Connect to a file containing settings or other data to be used by the application. Maintain the connection.

#### Live Trace

Any combination of the A trace and/or the B trace when SAVE A is off.

#### Load

Copy settings or other data from a file, but don't maintain any connection with that file.

#### Local Oscillator (LO)

An oscillator which produces the internal signal that is mixed with an incoming signal to produce the IF signal.

#### Marker

A visually identifiable point on a waveform trace, used to extract a readout of domain and range values represented by that point.

#### MAX HOLD

Digitally stored display mode which, at each frequency address, compares the incoming signal level to the stored level and retains the greater level. In this mode, the display indicates the peak level at each frequency after several successive sweeps.

#### MIN HOLD

A spectrum analyzer feature which captures the minimum signal amplitude at all displayed frequencies over a series of sweeps.

#### Modulate

To regulate or vary a characteristic of a signal.

#### **Modulating Signal**

The signal which modulates a carrier. The signal which varies or regulates some characteristic of another signal.

#### Modulation

The process of varying some characteristic of a signal with a second signal.

#### Noise

Unwanted random disturbances superimposed on a signal which tend to obscure it.

#### **Noise Floor**

The self-noise of an instrument or system that represents the minimum limit at which input signals can be observed. The spectrum analyzer noise floor appears as a "grassy" baseline in the display, even when no signal is present.

#### **Open (Recall, Restore)**

Bring data into the application from a file of application-native format.

#### **OQPSK**

Acronym for Offset QPSK (Quadrature Phase Shift Keying).

#### PDF

Acronym for Probability Distribution Function.

#### **Peak Detection**

A detection scheme wherein the peak amplitude of a signal is measured and displayed. In spectrum analysis, 20 log (peak) is often displayed.

#### **Platform Utility**

User-available software applications provided with the instrument.

#### **Primary Marker**

The marker displayed in the Single Marker mode whose frequency and/or position is changed when tuning with the general purpose knob.

#### PSK

Acronym for Phase Shift Keying. The process, or result of a process, in which the carrier phase is varied discretely in accordance with a digital code.

#### QAM

Acronym for Quadrature Amplitude Modulation. The process, or result of a process, in which the amplitude and phase of the carrier are varied concurrently by synthesizing two orthogonal ASK waves (see *ASK*).

#### **Reference Level**

The signal level required to deflect the CRT display to the top graticule line.

#### **Residual FM (Incidental FM)**

Short term displayed frequency instability or jitter due to instability in the spectrum analyzer local oscillators. Given in terms of peak-to-peak frequency deviation and expressed in hertz or percent of the displayed frequency.

#### **Residual Response**

A spurious response in the absence of an input signal. (Noise and zero pip are excluded.)

#### **Resolution Bandwidth (RBW)**

The width of the narrowest filter in the IF stages of a spectrum analyzer. The RBW determines how well the analyzer can resolve or separate two or more closely spaced signal components.

#### Save

Commit data to memory in application-native format.

#### **Secondary Marker**

The "second" marker displayed only in the Delta Marker mode.

#### Sensitivity

Measure of a spectrum analyzer's ability to display minimum level signals, expressed in volts or decibels. Internal frequency (IF) bandwidth, display mode, and any other influencing factors must be given.

#### **Shape Factor (Skirt Selectivity)**

The ratio of the frequency separation of the two (60 dB/6 dB) down points on the response curve to the static resolution bandwidth.

#### Single Sweep

Operating mode in which the sweep generator must be reset for each sweep. Especially useful for obtaining single examples of a signal spectrum.

#### Span Per Division, Span/Div

Frequency difference represented by each major horizontal division of the graticule.

#### Spectrum

The frequency domain representation of a signal wherein it is represented by displaying its frequency distribution.

#### **Spectrum Analysis**

The technique or process of determining the frequency distribution of a signal.

#### Spectrum Analyzer

A device for determining the frequency components of a signal.

#### Spectrum Length

The amount of time required to acquire one frame.

#### **Spurious Response**

A response to a spectrum analyzer wherein the displayed frequency is not related to the input frequency.

#### Vertical Scale Factor, Vertical Display Factor

The number of units (such as dB or V), represented by one vertical division of a spectrum analyzer display screen.

#### View

An area or window of the display screen which contains information all of the same type (such as a Spectrum View or a Power versus Time View). Glossary

# Index

# Index

# A

A/D overflow, 3-131 AC line frequency, 1-12 Accessories, F-1 optional, F-4 standard, F-3 ACPR measurement, 3-6 RFID analysis (Option 21 only), 3-68 Acquisition history, 3-142 length, 3-142 Acquisition & Analysis menu, zoom, 3-23 Acquisition/Analysis menu, 3-141 AM signal measurement, 3-39 AM/AM measurement, 3-56 view setting, 3-204 AM/PM measurement, 3-58 view setting, 3-205 Amplitude correction, 3-132 menu, 3-127, 3-128 setting, 3-127 Analog modulation analysis, 3-37 basic procedure, 3-37 Analysis length, 3-142 offset, 3-142 Analysis range, setting, 3-30 Angular units, System menu, 2-23 Architecture, 1-7 Arrow keys, Markers menu, 3-218 ASK Modulation Type, 3-45 note for measurements, 3-49 Auto level, 3-128 Average, 3-175 Demod mode, 3-179 DPX spectrum, 3-177 menu, 3-176 Time mode, 3-179 type, 3-180

### B

Backing up user files, 1-22 Bin, relation between frame, bin, and pixel, 3-183 Block diagram, 1-7 Bluetooth, 3-45 Brightness adjusting, 1-28 uneven, 1-21

### С

C/N measurement, 3-7 C4FM, 3-45 Cal menu, 1-24 Calibration, 1-23 Carrier frequency measurement, 3-10 Carrier measurement, RFID analysis (Option 21 only), 3-64 Carrier tracking, Signal source analysis (Option 21 only), 3-109 Caution, RF INPUT, 1-14 CCDF Demod mode, 3-59 Time mode, 3-88 view setting, 3-197 Center offset calibration, 1-26 Changing overview, 3-34 subview, 3-34 Channel power measurement, 3-5 Channel table, using, 3-122 Characteristics electrical, A-2 environmental, A-16 physical, A-16 Cleaning, E-1 Compression, trace data, 3-183 Constellation measurement, 3-50 RFID analysis (Option 21 only), 3-80 view setting, 3-198 Continuous, trigger, 3-149 Conventions, xxi Copy, screen, 3-259 Correction, amplitude description, 3-132 menu, 3-133 CSV file format, 3-258

### D

DC offset calibration, 1-27 Decoding format, RFID analysis (Option 21 only), 3-71 Default settings, restoring, 2-22 Demod mode, 3-27 Difference between RSA3303B and RSA3308B, 1-2 Digital modulation analysis, 3-44 basic procedure, 3-49 symbol mapping, C-1 Dimensions, A-16 Display, dead pixel, 1-21 Display line, 3-209 multi display lines (Real Time S/A only), 3-213 operation, 3-210 Documentation list, xx PDF, xx DPX spectrum average, 3-177 displaying, 3-14 marker operation, 3-224 view setting, 3-191

### Ε

EBW measurement, 3-10 Electrical characteristics, A-2 Environmental characteristics, A-16 Ethernet connector, side panel, 2-5 EVM measurement, 3-51 view setting, 3-199 Extended Res, 3-165 Eye diagram analysis, 3-55 RFID analysis (Option 21 only), 3-80 setting, 3-202

### F

Features, 1-1 FFT description, 3-163 Extended Res., 3-165 menu, 3-164 overlap, 3-167 points, 3-168 window, 3-169 File Bitmap trace file format, 3-256 data file format, 3-246

entering a name, 3-242 format, 3-245 loading, 3-236 operation, 3-233 overwrite, 3-239 saving, 3-236 trace file format, 3-254 type, 3-233 Floppy disk drive, side panel, 2-5 using, 3-236 FM signal measurement, 3-40 Frame, relation between frame, bin, and pixel, 3-183 Free Run, trigger mode, 3-149 Frequency menu, 3-120 setting, 3-119 Frequency versus Time measurement, 3-87 Signal source analysis (Option 21 only), 3-117 Frequency/Channel, 3-120 FSK Modulation Type, 3-45 note for measurements, 3-49 Functional check, 1-16

### G

Gain calibration, 1-25 General purpose knob, 2-3

### Η

Hard copy, screen, 3-259 Help, menu, 3-227 How to adjust center offset, 1-26 DC offset, 1-27 display brightness, 1-28 gain, 1-25 How to create amplitude correction files, 3-135 directory, 3-244 files, 3-261 trigger mask, 3-155 How to delete directories, 3-244 files, 3-244 How to display averaged waveforms, 3-181 data by trigger and repeat modes, 3-153 DPX spectrum, 2-46 horizontal display line, 3-211

multiple views, 3-189 online help, 3-227 options, 2-24 trace 1 and 2, 3-179 trigger point, 3-161 versions, 2-24 vertical display line, 3-212 Windows desktop, 2-26 How to enter file name, 3-242 numeric data, 2-17 How to install keyboard, 2-6 mouse, 2-6 printer driver, 3-260 USB printer connection, 3-260 How to load files, 3-236 menu settings, 3-236 trace, 3-183 How to save amplitude correction table, 3-132, 3-233 data to a PC via a network, 2-5 files, 3-236 menu settings, 3-236 preset files, 3-237 screen copy to a disk, 3-261 trace, 3-183 How to set, zoomed area, 3-24

# I

```
Input
menu, 3-231
over-voltage, 3-131
source, 3-231
Inspection, E-1
Installation, 1-11
Integrated phase noise, description, 3-107
Interface, maps, 2-1
IQ
input, note, 3-231
IQ versus Time measurement, 3-42, 3-85
IQ/Frequency versus Time measurement, 3-52
```

### K

Key, lock, 2-11 Keyboard, operation, 2-7 Keypad, 2-3 entering a file name, 3-242 using, 2-19

### L

LAN Ethernet connector, side panel, 2-5 Lines menu, 3-210 menu (Real Time S/A), 3-214 Load, menu, 3-234

### Μ

Main view, modulation analysis screen, 3-29 Manuals list, xx PDF, xx Map, interface, 2-1 Marker Bitmap trace, 3-224 menu, 3-218 operation, 3-217 Marker Setup, 3-219 MAT file format, 3-258 Max Hold, description, 3-180 Meas Setup menu digital modulation analysis, 3-45 FM Demod measurement, 3-39, 3-40 PM Demod measurement, 3-41 pulse measurement, 3-100 Measurement, mode, 2-21 Measurements ACPR, 3-6 AM demod, 3-39 AM/AM. 3-56 AM/PM. 3-58 analog modulation analysis, 3-37 C/N, 3-7 carrier frequency, 3-10 CCDF, 3-59, 3-88 channel power, 3-5 constellation, 3-50 digital modulation analysis (Option 21 only), 3-44 EBW, 3-10 EVM, 3-51 eye diagram, 3-55 FM demod, 3-40 frequency versus time, 3-87 IQ versus time, 3-42, 3-85 IQ/frequency versus time, 3-52 modulation analysis (Demod mode), 3-27 OBW, 3-8 PDF, 3-60 PM demod, 3-41 power versus time, 3-53, 3-86 pulse, 3-93

pulse spectrum, 3-43 real-time analysis, 3-19 spectrum analysis (S/A mode), 3-1 spurious signal, 3-12 symbol table, 3-54 time analysis (Time mode), 3-83 Menu Acquisition/Analysis, 3-141 Amplitude, 3-127 amplitude correction, 3-133 Average, 3-176 basic operation, 2-16 FFT, 3-164 Frequency, 3-120 front panel keys, 2-3 Help, 3-227 Input, 3-231 Lines, 3-210 Lines (Real Time S/A), 3-214 Load, 3-234 Markers, 3-218 Print, 3-259 RBW, 3-164 Save, 3-234 Span, 3-120 System, 2-23 Trace, 3-176 tree, B-1 Trigger, 3-148 View, 3-188 Min Hold, description, 3-180 Minimum jump, setting, 3-226 Mode, selecting, 2-21 Modulation analysis (Demod mode), 3-27 Mouse, operation, 2-7 Mult display lines (Real Time S/A only), 3-213

### Ν

NBW, setup display, 2-12 Noisogram Signal source analysis (Option 21 only), 3-111 view setting, 3-207 Numeric input, 2-17

### 0

OBW measurement, 3-8 Online help, using, 3-227 Option, displaying, 2-24 Option key, 2-25 Optional accessories, F-4 Options, F-1 OQPSK, shifting Q data, 3-47 Overlap, FFT, 3-167 Overrange, 3-131 Overview changing, 3-34 modulation analysis screen, 3-29 Over-voltage input, 3-131

### Ρ

P25 C4FM, 3-45 PDF measurement, 3-60 view setting, 3-206 PDF manual, xx Peak, 3-218 Peak search, 3-217 operation, 3-225 setting minimum jump, 3-226 Phase noise measurement, Signal source analysis (Option 21 only), 3-105 Physical characteristics, A-16 Pixel, relation between frame, bin, and pixel, 3-183 PM signal measurement, 3-41 Power AC requirements, 1-12 connecting the cord, 1-12Power on/down measurement, RFID analysis (Option 21 only), 3-69 Power switch, 1-13 Power versus Time measurement, 3-86 Demod mode, 3-53 Powering off, 1-20 Powering on, 1-12, 1-13 Preset, 2-22 Principal power switch, 1-13 Print menu, 3-259 screen copy, 3-259 Product description, 1-1 Pulse measurements, 3-93 Pulse spectrum, 3-43

### Q

Q Data Half Symbol Shift, 3-47

### R

Random jitter, description, 3-107

RSA3303B & RSA3308B Analyzers User Manual

RBW description, 3-163 menu, 3-164 Real Time S/A, operation, 3-19 Real-time phase noise measurement, Signal source analysis (Option 21 only), 3-111 Real-time spurious measurement, Signal source analysis (Option 21 only), 3-115 Real-time analysis, 1-3 S/A mode, 3-19 Rear panel, 2-4 Repeat, trigger, 3-149 Restart, 1-21 Restoring the operating system, D-1 RF envelope measurement, RFID analysis (Option 21 only), 3-76 RF INPUT, Caution, 1-14 RFID analysis (Option 21 only), 3-61 ACPR measurement, 3-68 carrier measurement, 3-64 constellation, 3-80 eye diagram, 3-80 Power on/down measurement, 3-69 RF envelope measurement, 3-76 spurious measurement, 3-66 symbol table, 3-80 RMS, description, 3-180 RMSExpo, description, 3-180 RUN/STOP key, 2-22

### S

S/A mode, 3-1 S/A with Spectrogram, operation, 3-17 Save, menu, 3-234 Save on trigger, 3-154 Scalar mode, span, 3-126 Scale menu AM/AM, 3-204 AM/PM, 3-205 CCDF, 3-197 Constellation, 3-198 DPX spectrum, 3-191 EVM, 3-199 eye diagram, 3-202 noisogram, 3-207 PDF, 3-206 spectrogram, 3-193 spectrum, 3-190 symbol table, 3-201 time domain view, 3-195 Scan Disk, 1-21

Screen, elements, 2-8 Screen copy, 3-259 Seamless acquisition, 3-145 Select, Markers, 3-218 Setting, analysis range, 3-30 Setup, display, 2-12 Shifting Q data (OQPSK), 3-47 Side keys, 2-3 Side panel, 2-5 Signal amplitude limit, 1-14 Signal processing, 1-7 Signal source analysis (Option 21 only), 3-103 frequency versus time measurement, 3-117 phase noise measurement, 3-105 real-time phase noise measurement, 3-111 real-time spurious measurement, 3-115 spurious measurement, 3-109 Single, trigger, 3-149 Single view, 3-36 Span, 3-121 menu, 3-120 setting, 3-119 vector span, 3-126 Specifications, A-1 Spectrogram S/A with Spectrogram, 3-17 view setting, 3-193 Spectrum analysis, S/A mode, 3-1 length, 3-142 offset, 3-143 view setting, 3-190 Spurious measurement RFID analysis (Option 21 only), 3-66 Signal source analysis (Option 21 only), 3-109 Spurious signal measurement, 3-12 Stand, setting up, 1-15 Standard accessories, F-3 Status, display, 2-10 Step size center frequency, 2-20 setting, 2-20 Stop and Show Results, 3-150 Subview changing, 3-34 modulation analysis screen, 3-29 Symbol mapping, C-1 Symbol table analysis, 3-54 RFID analysis (Option 21 only), 3-80 setting, 3-201

Symmetrical filter, Signal source analysis (Option 21 only), 3-109
Synchronous operation, external trigger, 3-162
System menu, 2-23

### Т

Tari, RFID analysis (Option 21 only), 3-70 Time analysis (Time mode), 3-83 Time domain, view setting, 3-195 Time variation measurement, 3-84 basic procedure, 3-84 Trace comparison, 3-175 displaying, 3-179 menu, 3-176 saving/loading, 3-183 Trigger, 3-147 creating a mask, 3-155 level, 3-151 menu, 3-148 mode, 3-149 position, 3-150 save on trigger, 3-154 slope, 3-152 source, 3-150 specification, A-10 synchronous operation, 3-162 trigger output, indicating, 3-161 trigger point indicator, 3-161 Triggered, trigger mode, 3-149

Troubleshooting, 1-21 Tutorial, 2-29

### U

UNCAL, 1-14 Unpacking, 1-11 USB connectors, side panel, 2-6

### V

Vector mode, span, 3-126 Vector span, 3-126 Version, displaying, 2-24 VGA output connector, side panel, 2-6 Video filter, RFID analysis (Option 21 only), 3-67 View menu, 3-188 setting, 3-187

### W

Windows XP, using, 2-26

### Ζ

Zoom Acquisition & Analysis menu, 3-23 function, 3-23 setting a zoomed area, 3-24